

Manuals+

[Q & A](#) | [Deep Search](#) | [Upload](#)

manuals.plus /

› [AERLANG](#) /

› [AERLANG Cordless Neck and Back Massager BK02 User Manual](#)

AERLANG BK02

AERLANG Cordless Neck and Back Massager with Heat (Model: BK02)

User Instruction Manual

INTRODUCTION

The AERLANG Cordless Neck and Back Massager with Heat (Model BK02) is designed to provide targeted relief for muscle soreness, stiffness, and fatigue in the neck, shoulders, and back. Featuring 3D Shiatsu kneading nodes, adjustable intensity levels, bi-directional rotation, and soothing heat, this device offers a customizable massage experience. Its cordless design enhances portability and convenience, allowing for use in various settings.

IMPORTANT SAFETY INFORMATION

Please read all safety instructions carefully before using this product. Failure to follow these instructions may result in injury or damage to the device.

- Always unplug the massager from the power source before use.
- Do not use the massager while sleeping or taking a bath.
- Keep the massager away from water and heat sources.
- Do not use the massager if you have any medical conditions, are pregnant, or have any injuries or skin conditions in the area to be massaged. Consult a physician before use if you have any health concerns.
- Always follow the instructions provided with the massager.

SETUP

Charging the Massager

Before initial use, fully charge the massager for 5 hours. For subsequent uses, charge for 2-3 hours. A flashing red light on the remote indicates that charging is needed. The red light will flash during charging and remain solid red when fully charged. This product cannot be operated during charging.

Multipurpose Use

Can not be operated while it is charging



Neck



Abdomen



Waist



Arm



Thigh



Calf

Image: Remote control with charging indicator and function buttons. Note: The massager cannot be operated while charging.

Placement and Securing

The massager is designed to be versatile. Use the adjustable strap on the back of the unit to secure it to a chair or position it comfortably on a sofa or bed. Ensure the massager is stable before use.

Your browser does not support the video tag.

Video: Demonstrates the AERLANG massager being placed on a sofa for neck massage and its cordless operation.

OPERATING INSTRUCTIONS

The AERLANG massager features a multi-function remote control for easy operation.


1. **Power On/Off:** Press the **ON/OFF** button to turn the massager on or off.
2. **Rotation:** The massager features bi-directional rotation. Press the **ROTATION** button to change the direction of the kneading nodes. The rotation direction automatically changes every few minutes during operation.
3. **Intensity:** Select from three adjustable intensity settings (Low/Medium/High) by pressing the **INTENSITY** button.
4. **Heat Function:** Activate the soothing heat function by pressing the **HEAT** button. The massager offers two adjustable heat levels

and dual heating zones for both neck and back.

Recommended Usage


- Place the massager on the desired body area: neck, shoulders, back, abdomen, waist, arm, thigh, or calf.
- For safety, do not use the massager for more than 15 minutes at a time. The device will automatically shut down after 15 minutes of continuous use.

Perfect Fit the Neck and Back for More Comfort





Other Massager



Our Massager


Can not be operated while it is charging

Image: The massager's ergonomic design provides a perfect fit for the neck and back, enhancing comfort during use.

Heated Massager

Promotes blood circulation
and pain relief

2 Levels | 2 Heating zones

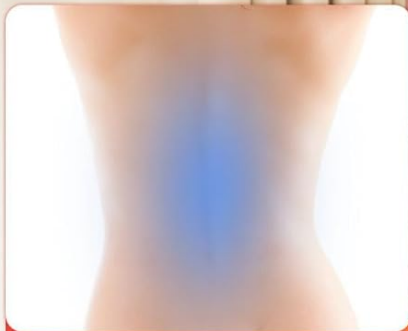

Before



After

Can not be operated while it is charging

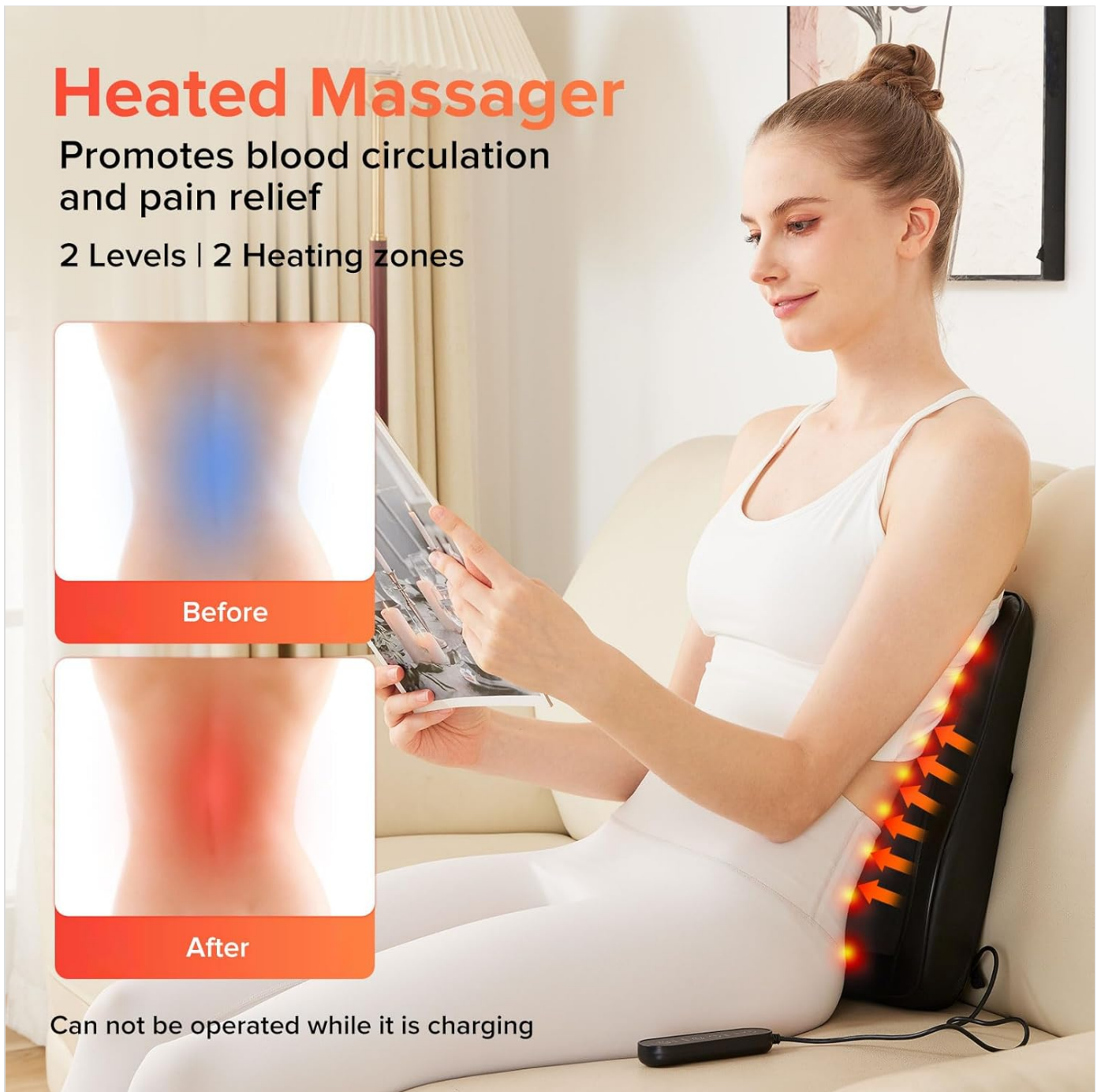


Image: The heated massager feature, illustrating its effect on blood circulation and pain relief with two heat levels and two heating zones.

Remote Control is Easy to Operate

Can not be operated while it is charging



Wireless Back Massager

Rechargeable + remote control is easy to operate



INDICATOR LIGHT

ON/OFF

ON/OFF BUTTON

ROTATION

ROTATION BUTTON

INTENSITY

INTENSITY BUTTON

HEAT

HEAT BUTTON

Built-in rechargeable battery

1. Charging time is 5 hours.
A flashing red light indicates charging is needed.

2. The red light flashes during charging and remains solid red when fully charged.

3. Please note that this product cannot be operated during charging. (Automatically shut down within 15Min.)

Image: Demonstrates the multipurpose use of the massager on various body areas including the neck, abdomen, waist, arms, thighs, and calves.

MAINTENANCE

Cleaning

- Ensure the massager is turned off and unplugged before cleaning.
- Wipe the surface of the massager with a soft, damp cloth. Do not use abrasive cleaners or immerse the device in water.
- Allow the massager to air dry completely before storing or next use.

Storage

- Store the massager in a cool, dry place away from direct sunlight and extreme temperatures.
- Avoid placing heavy objects on the massager during storage.

TROUBLESHOOTING

- **Massager does not turn on:** Ensure the device is fully charged. Check if the charging cable is securely connected if attempting to

charge. Remember, the massager cannot be operated while charging.

- **Heat function is not working:** Verify that the heat button has been pressed. Allow a few minutes for the heat to become noticeable.
- **Massage intensity feels too weak/strong:** Adjust the intensity level using the intensity button on the remote control.
- **Unusual noise during operation:** Discontinue use and contact customer support if you experience any abnormal sounds.

SPECIFICATIONS

Feature	Detail
Model Number	BK02
Product Dimensions	12.8 x 4 x 15 inches
Item Weight	3.7 Pounds
Power Source	Battery Powered (1 Lithium Ion battery included)
Material	Leather, Nylon
Massage Nodes	4 x 3D Shiatsu Kneading Nodes
Intensity Levels	3 (Low/Medium/High)
Heat Levels	2
Automatic Shut-off	15 minutes

WARRANTY AND SUPPORT

The AERLANG Cordless Neck and Back Massager is backed by a 1-year warranty. For any product inquiries, technical assistance, or warranty claims, please contact AERLANG professional customer support. Refer to the product packaging or the official AERLANG website for contact details.

Legal Disclaimer: Statements regarding dietary supplements have not been evaluated by the FDA and are not intended to diagnose, treat, cure, or prevent any disease or health condition.

