

HX-601/603

Wireless Knee Massager HX-601/603 Instruction Manual

Model: HX-601/603 | Brand: Generic

1. IMPORTANT SAFETY INFORMATION

Please read all instructions carefully before using the device. Keep this manual for future reference.

- Do not use if you have a pacemaker, artificial joints, or electronic medical devices.
- Consult a physician before use if you have any pre-existing medical conditions, skin sensitivities, or are pregnant.
- Do not use on open wounds, swollen areas, or areas with skin irritation.
- Keep out of reach of children.
- Do not immerse the device in water or use it in excessively humid environments.
- Use only the provided charging cable and power adapter.
- Discontinue use immediately if you experience any discomfort or pain.

2. PRODUCT OVERVIEW

The HX-601/603 Wireless Knee Massager is designed to provide therapeutic massage and integrated heating for your knees. Its ergonomic and portable design allows for comfortable use at home, in the office, or during travel.

2.1 Components

- Main Massager Unit with Control Panel
- Adjustable Straps
- Magnetic Battery Compartment
- USB Charging Cable

2.2 Features

- **Wireless Operation:** Enjoy massage without tangled cords.
- **Integrated Heating:** Therapeutic heat to improve circulation and relax muscles.
- **3 Intensity Levels:** Customize your massage experience.
- **15-Minute Timer:** Automatic shut-off for safety and optimal session duration.

- **Digital Control Panel:** Easy adjustment of settings with LED display.
- **Ergonomic Design:** Conforms to the knee for maximum comfort and effectiveness.
- **Comfortable Fabrics:** Made with breathable materials for extended use.



Forte atração magnética


Difícil romper

Tecidos confortáveis

Confortável e bonito,
respirável e não abafado



The image displays the wireless knee massager worn on a knee, illustrating the integrated heating function with red arrows indicating heat flow.

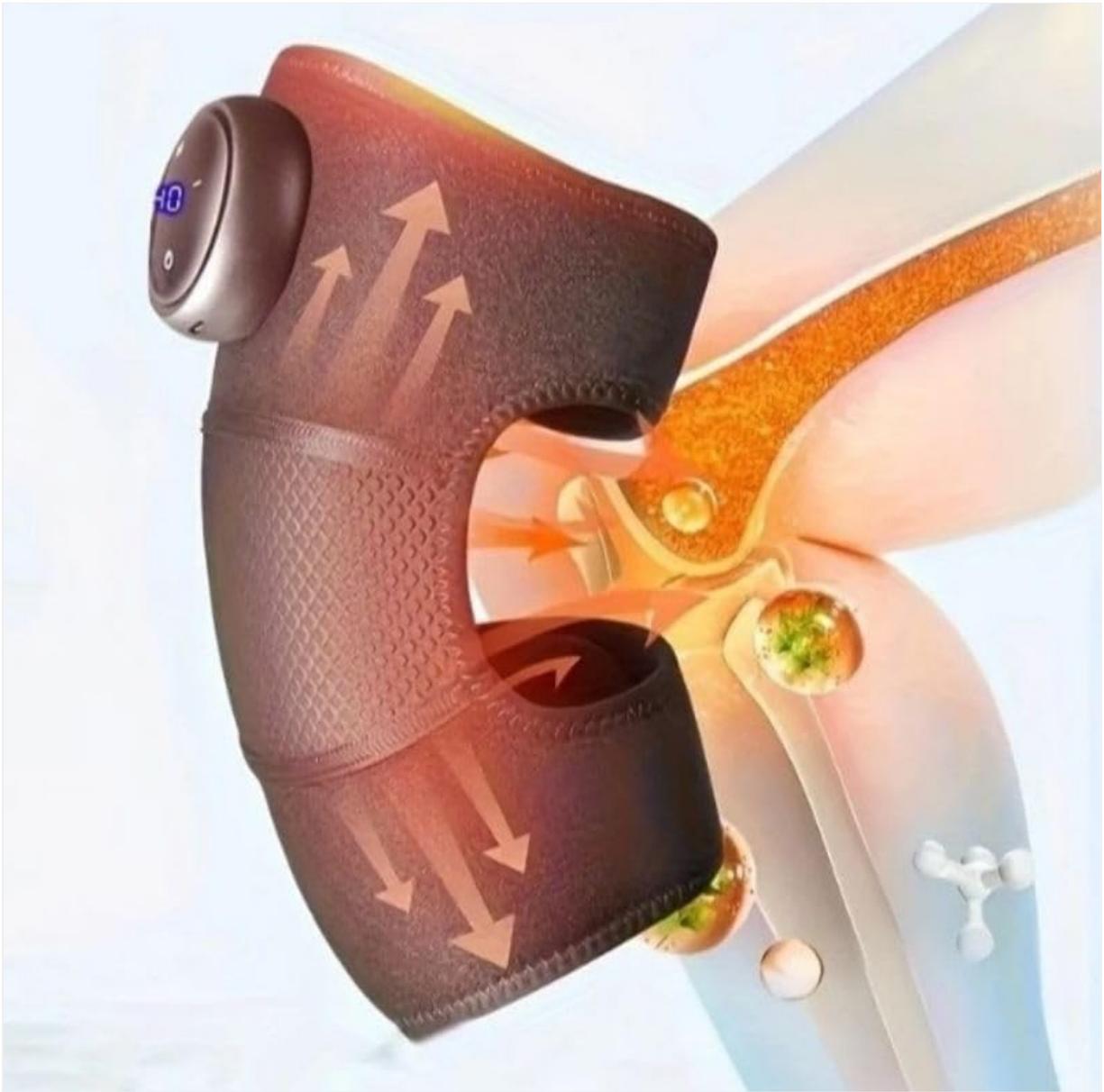

This image highlights the device's control panel with a digital display and buttons, along with the magnetic battery compartment for wireless operation.

Compartimento magnético da bateria

Desfrute de uma massagem sem fio



A detailed view of the massager's fabric, emphasizing its comfortable, breathable, and non-stuffy texture.

3. SETUP

3.1 Charging the Device

1. Locate the charging port on the massager unit.
2. Connect the USB charging cable to the device and to a compatible USB power adapter (not included) or computer USB port.
3. The LED indicator on the control panel will show charging status. A full charge typically takes several hours.
4. Ensure the device is fully charged before first use for optimal battery life.

3.2 Attaching to the Knee

1. Unfold the knee massager.
2. Position the main unit over your knee, ensuring the massage elements are centered.
3. Secure the massager using the adjustable straps, ensuring a snug but comfortable fit. Do not overtighten.

4. OPERATING INSTRUCTIONS

4.1 Power On/Off

- To turn on: Press and hold the Power button (usually marked with a circle or power symbol) on the control panel for 2-3 seconds. The LED display will illuminate.
- To turn off: Press and hold the Power button again for 2-3 seconds. The display will turn off.

4.2 Adjusting Massage Intensity

- Once powered on, press the Intensity button (usually marked with '+' or '-') to cycle through the 3 available intensity levels.
- The current intensity level will be indicated on the LED display.

4.3 Activating Heating Function

- Press the Heat button (usually marked with a wave or heat symbol) to activate or deactivate the integrated heating function.
- The heating indicator on the display will light up when active.

4.4 Automatic Timer

The device is equipped with an automatic 15-minute timer. After 15 minutes of continuous operation, the massager will automatically shut off. You can restart the device for another session if desired, but it is recommended to take breaks between sessions.

5. MAINTENANCE

5.1 Cleaning

- Ensure the device is turned off and unplugged before cleaning.
- Wipe the surface of the massager with a soft, damp cloth.
- Do not use abrasive cleaners, solvents, or harsh chemicals.
- Do not immerse the device in water.

5.2 Storage

- Store the device in a cool, dry place, away from direct sunlight and extreme temperatures.
- Keep it away from sharp objects that could damage the fabric or internal components.

6. TROUBLESHOOTING

Problem	Possible Cause	Solution
Device does not turn on.	Low battery or device not charged.	Charge the device fully. Ensure the power button is pressed and held for 2-3 seconds.
No heat or weak heat.	Heating function not activated or low battery.	Ensure the heat button is pressed. Charge the device.
Massage intensity is too weak/strong.	Incorrect intensity setting.	Adjust the intensity using the '+' or '-' buttons.

Problem	Possible Cause	Solution
Device stops unexpectedly.	Automatic 15-minute timer activated or low battery.	This is normal operation due to the timer. Recharge the device if the battery is low.

7. SPECIFICATIONS

Model Number	HX-601/603
Brand	Generic
Power Source	Battery Powered
Voltage	5V
Power	7W
Product Dimensions	25 x 14 x 8 cm
Product Weight	100 Grams
Material	Electronic Components, Plastic, Fabric
Timer	15 minutes (automatic shut-off)
Intensity Levels	3
Country of Origin	Brazil



This image provides the physical dimensions of the knee massager, indicating a length of 50cm and a width of 28cm when laid flat.

8. WARRANTY AND SUPPORT

This product is manufactured by Knee. For warranty information or technical support, please contact your retailer or the manufacturer directly. Please retain your proof of purchase.

Note: Specific warranty terms may vary by region and retailer.