

KINGJOY KJF1

KINGJOY KJF1 Foot Massager Instruction Manual

Model: KJF1

1. IMPORTANT SAFETY INFORMATION

Please read all instructions carefully before using the KINGJOY KJF1 Foot Massager. Retain this manual for future reference.

- **Medical Consultation:** Consult a physician before use if you have heart disease, neurological conditions, diabetes, or vascular issues.
- **Power Supply:** Always use the provided original adapter and a power supply matching the specified voltage and power rating (Input: 12 V, 3 A; Rated Power: 36 W).
- **Water and Moisture:** Do not use the massager in wet or damp environments. Keep away from water.
- **Children and Pets:** Keep out of reach of children and pets. This appliance is not intended for use by persons with reduced physical, sensory, or mental capabilities unless supervised.
- **Damaged Cord/Plug:** Do not operate the massager if the power cord or plug is damaged. Contact customer support for assistance.
- **Overheating:** If the massager feels excessively hot, discontinue use immediately.
- **Intended Use:** This device is designed for personal foot massage and relaxation. Do not use for any other purpose.

2. PRODUCT OVERVIEW

The KINGJOY KJF1 Foot Massager is designed to provide comprehensive foot pain relief through a combination of deep shiatsu kneading, rolling, and U-shaped air compression, along with dual-zone heat

therapy. It is suitable for various foot discomforts, including those associated with plantar fasciitis and neuropathy.



Image: The KINGJOY KJF1 Foot Massager, a gray unit with foot inserts and a digital display, accompanied by its remote control.

Key Features:

- **Full-Coverage Foot Pain Relief:** Utilizes deep shiatsu kneading, rolling, and U-shaped air squeeze techniques.
- **Dual-Zone Heat:** Adjustable heat (104-131°F) applied to both the top and bottom of the feet.
- **3D Shiatsu Rollers:** Ergonomic design with rollers that conform to the foot's arch, heel, and forefoot.
- **Easy Control:** Features a large LCD touchscreen and a wireless remote control for convenient operation.
- **Customizable Settings:** 3 kneading intensities, 3 air squeeze levels, 3 heat levels, and 3 timer settings (15/20/30 minutes).
- **Easy Maintenance:** Removable, machine-washable foot covers.

- **Compact Design:** Lightweight and portable for easy storage and use.



3-IN-1 MULTIFUNCTION
Get The Value of 3 Products for 1 Price!



Foot Boots
(Air Compression)

Foot Spa
(Rolling Massager)

Foot Warmers
(Heating)

Image: An illustration highlighting the 3-in-1 multifunctionality of the massager, comparing it to foot boots (air compression), a foot spa (rolling massager), and foot warmers (heating).

3. SETUP

1. **Unpacking:** Carefully remove the foot massager from its packaging. Inspect for any signs of damage.
2. **Placement:** Place the massager on a flat, stable surface. Ensure there is enough space around the unit for proper ventilation.

3. **Power Connection:** Connect the provided power adapter to the massager's power input port, then plug the adapter into a suitable electrical outlet.
4. **Remote Control Battery:** The remote control requires 1 nonstandard battery (included). Ensure the battery is correctly inserted.
5. **Foot Covers:** Ensure the removable foot covers are properly installed before use.

4. OPERATING INSTRUCTIONS

The KINGJOY KJF1 Foot Massager can be operated using either the integrated LCD touchscreen or the wireless remote control.

A ONE-OF-A-KIND DESIGN
The protective film on the panel can be removed.

Ambient Light *Design*

Large LCD Touch Panel

Elegant & Smooth Body Lines

Image: Close-up view of the KINGJOY KJF1 Foot Massager's LCD touchscreen control panel and the accompanying wireless remote control.

Getting Started:

1. **Power On:** Press the 'Power' button on the control panel or remote to turn on the massager. The LCD screen will illuminate.
2. **Insert Feet:** Place your feet comfortably into the foot inserts. The massager is designed to fit up to men's size 14.

Adjusting Settings:

Use the corresponding buttons on the control panel or remote to adjust the massage settings:

- **Kneading (Knead):** Press the 'Knead' button to cycle through 3 intensity levels of shiatsu kneading massage.
- **Air Compression (Air):** Press the 'Air' button to cycle through 3 intensity levels of U-shaped air squeeze compression.
- **Heat (Heat):** Press the 'Heat' button to cycle through 3 heat levels (104-131°F). Heat is applied to both soles and tops of feet. It is recommended to start at the lowest setting.
- **Timer (Timer):** Press the 'Timer' button to select a massage duration of 15, 20, or 30 minutes. The default setting is 15 minutes.

YOUR FEET IN WARMTH



Image: A person using the foot massager, illustrating the dual-zone heating feature with visual indicators for heat on both the top and bottom of the feet.




Image: A diagram showing the 360-degree U-shaped air compression system of the foot massager, highlighting how it hugs the heels and provides compression.

3D LIFELIKE SHIATSU MASSAGE

3-Section Comprehensive
Massage for Ultimate Relaxation

INSTEP

Air Compression

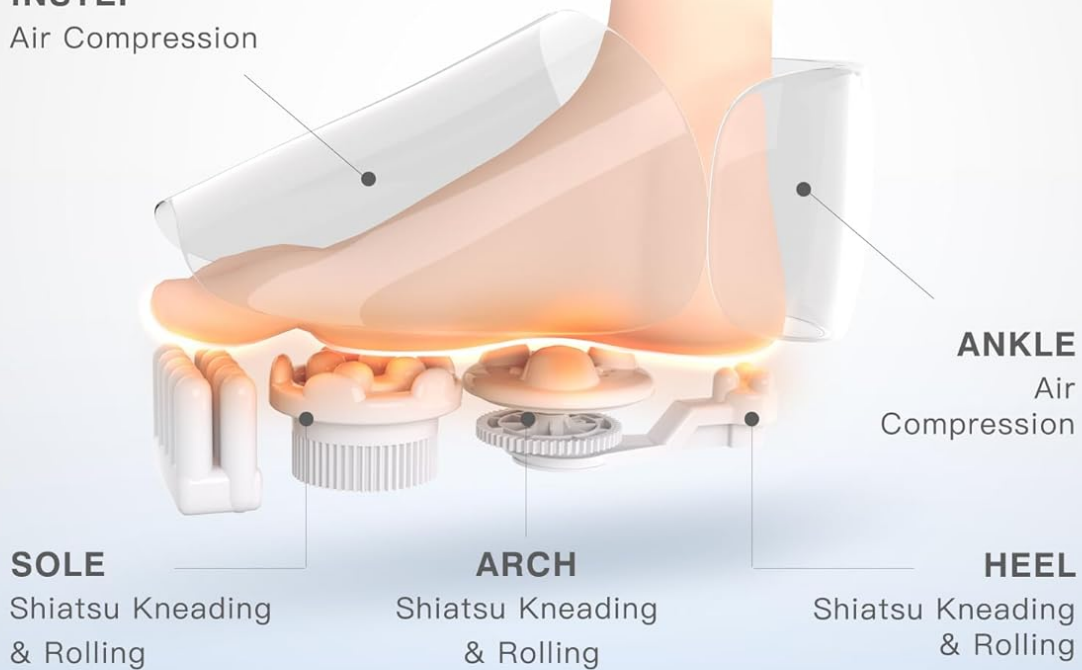

SOLE

Shiatsu Kneading
& Rolling

ARCH

Shiatsu Kneading
& Rolling

HEEL

Shiatsu Kneading
& Rolling



Image: A detailed diagram illustrating the 3D shiatsu rollers and their placement to target the sole, arch, and heel, along with instep and ankle air compression.

Recommendation:

For initial use and muscle relaxation, it is recommended to start with the lowest intensity settings for kneading and air compression, and a 15-minute session.

Power Off:

Press the 'Power' button again to turn off the massager. The unit will also automatically shut off once the timer concludes.

5. MAINTENANCE

Cleaning:

- The foot massager features removable, machine-washable foot covers for easy hygiene.
- To remove, unzip the covers from the foot inserts.
- Wash the covers according to the care instructions on their label.
- Ensure the covers are completely dry before reattaching them to the massager.
- Wipe the exterior of the massager with a soft, damp cloth. Do not use abrasive cleaners or immerse the unit in water.

Storage:

- Store the massager in a cool, dry place when not in use.
- Its compact and lightweight design makes it easy to store in a closet or under furniture.

6. TROUBLESHOOTING

If you encounter any issues with your KINGJOY KJF1 Foot Massager, please refer to the following common solutions:

- **Massager does not power on:**
 - Ensure the power adapter is securely plugged into both the massager and a working electrical outlet.
 - Check if the power outlet is functional by plugging in another device.
- **No heat or weak heat:**
 - Verify that the heat function is activated and set to the desired level.
 - Allow a few minutes for the heating elements to warm up.
- **Massage feels too intense or too weak:**
 - Adjust the kneading and air compression intensity levels using the respective buttons.
 - Ensure your feet are positioned correctly within the massager.
- **Remote control not working:**
 - Check and replace the battery in the remote control if necessary.
 - Ensure there are no obstructions between the remote and the massager.

If the issue persists after attempting these steps, please contact KINGJOY customer support for further assistance.

7. SPECIFICATIONS

Feature	Detail
Brand	KINGJOY
Model Number	KJF1
Product Dimensions	16.22 x 13.31 x 7.01 inches
Item Weight	6.61 Pounds
Power Source	Corded Electric
Rated Power	36 W
Input Voltage	12 V, 3 A
Material	Plastic
Fits Foot Size	Up to Men's Size 14
Heat Range	104-131 °F
Batteries (Remote)	1 Nonstandard Battery (included)

**IDEAL GIFT
FOR ANY OCCASION**

One Size Fits Whole Family

16.2"

13.3"

Image: The KINGJOY KJF1 Foot Massager with its dimensions (16.2 inches length, 13.3 inches width) indicated, alongside a group of happy people, suggesting it as an ideal gift.

8. WARRANTY AND SUPPORT

Warranty Information:

The KINGJOY KJF1 Foot Massager comes with a **1-year warranty**. This warranty covers free replacement for any damage that occurs within the first year from the date of purchase.

Customer Support:

For any questions, concerns, or warranty claims, please contact KINGJOY customer support through the retailer where the product was purchased or visit the official KINGJOY website for contact information.