

[Manuals.plus](#) /

› [FOREVIVE](#) /

› FOREVIVE Baby Finger Toothbrush Instruction Manual (Model B0FNVR5QKR)

FOREVIVE B0FNVR5QKR

FOREVIVE Baby Finger Toothbrush Instruction Manual

Model: B0FNVR5QKR

Brand: FOREVIVE

1. INTRODUCTION AND OVERVIEW

The FOREVIVE Baby Finger Toothbrush is designed to support early oral development for infants and toddlers aged 0 to 3 years. This ultra-soft, dual-function silicone brush serves as both a gentle gum massager for teething babies and an effective tongue cleaner and toothbrush for emerging teeth. It promotes healthy oral hygiene habits from an early age.



Image 1.1: The FOREVIVE Baby Finger Toothbrush, designed for young children.

2. SAFETY INFORMATION

- Always use this product under adult supervision.
- Inspect the toothbrush before each use. Discard immediately if any signs of wear, tear, or damage are observed.
- This product is made from 100% medical-grade soft silicone, which is BPA-free and latex-free.
- Do not allow the child to chew on the toothbrush unsupervised.
- Keep out of reach of children when not in use.

3. PRODUCT FEATURES

The FOREVIVE Baby Finger Toothbrush offers several key features:

- **Gentle and Safe:** Made from medical-grade soft silicone, it is ultra-gentle on delicate gums, emerging teeth, and sensitive mouths.
- **Dual Functionality:** Acts as a teething massager for soothing swollen gums (0-12 months) and a tongue cleaner/toothbrush for cleaning milk/formula residue and introducing brushing motions (1-3

years).

- **Ergonomic Design:** Slides comfortably onto an adult's index finger, providing full control during cleaning. The open-back design prevents hand sweating.
- **Promotes Healthy Habits:** Helps transition from gum wiping to regular brushing, fostering good dental hygiene from infancy through toddlerhood.
- **Hygienic and Durable:** Easy to clean, dishwasher safe (top rack), and designed for repeated use.



Professional Oral Cleaning Toothbrush

- ✓ One-piece molding No hair loss
- ✓ Liquid silicone, environmentally friendly and non-toxic
- ✓ Highly transparent and elastic
- ✓ Innovative shape is highly accepted by babies

Gauze cleaning

- ✗ Rough and lint free
- ✗ Chemical glue
- ✗ Easy to come off with water
- ✗ Lint in the mouth

Image 3.1: The silicone toothbrush offers a professional oral cleaning alternative to gauze.

Innovative Silicone Designed For Young Children

BPA-free

One-piece molding



Image 3.2: Innovative silicone design with gum-friendly, soft massaging bristles.

4. SETUP

Before first use, and after each subsequent use, thoroughly clean the finger toothbrush. Refer to the 'Maintenance and Care' section for detailed cleaning instructions.

5. OPERATING INSTRUCTIONS

To use the FOREVIVE Baby Finger Toothbrush:

1. **Preparation:** Ensure the toothbrush is clean. Wash your hands thoroughly.
2. **Placement:** Gently slide the silicone toothbrush onto your index finger. The flexible material is designed to fit snugly.
3. **For Infants (0-12 months):**
 - Gently massage your baby's gums to soothe teething discomfort.
 - Use the soft bristles to carefully clean milk or formula residue from the tongue and inner cheeks.
4. **For Toddlers (1-3 years):**

- Introduce brushing motions by gently cleaning emerging teeth.
- Encourage your child to become familiar with the sensation of brushing.

5. **After Use:** Remove the toothbrush from your finger and clean it immediately.



Image 5.1: The finger toothbrush effectively cleans the tongue, as shown in before and after images.

Your browser does not support the video tag.

Video 5.2: Demonstration of the FOREVIVE Baby Finger Toothbrush in use, showing its flexibility and how it cleans gums and teeth.

6. MAINTENANCE AND CARE

- **Cleaning:** After each use, rinse the toothbrush thoroughly under warm running water. It is also dishwasher safe (top rack recommended).
- **Drying:** Allow the toothbrush to air dry completely before storing.
- **Storage:** Store in a clean, dry place away from direct sunlight and extreme temperatures.
- **Inspection:** Regularly inspect the toothbrush for any signs of wear, such as cracks, tears, or discoloration. Replace the toothbrush if any damage is detected to ensure safety and hygiene.

7. TROUBLESHOOTING

- **Issue: The toothbrush does not stay securely on the finger.**

Solution: Ensure your finger is dry before sliding on the toothbrush. The silicone is designed to stretch and fit various finger sizes. If it consistently feels loose, it may not be the right fit for your finger.

- **Issue: Child resists brushing.**

Solution: Make brushing a positive experience. Start with short, gentle sessions. Allow the child to hold the toothbrush or play with it (under supervision) to get accustomed to it. Sing songs or use a mirror to make it fun. Consistency is key, even if it's just for a few seconds initially.

8. SPECIFICATIONS

- **Product Dimensions:** 0.5 x 0.5 x 1.9 inches
- **Item Weight:** 1.76 ounces
- **Material:** 100% Medical-Grade Soft Silicone
- **Age Range:** 0 to 3 Years
- **Special Feature:** Lightweight
- **Item Firmness:** Extra Soft
- **Safety:** BPA-free, Latex-free

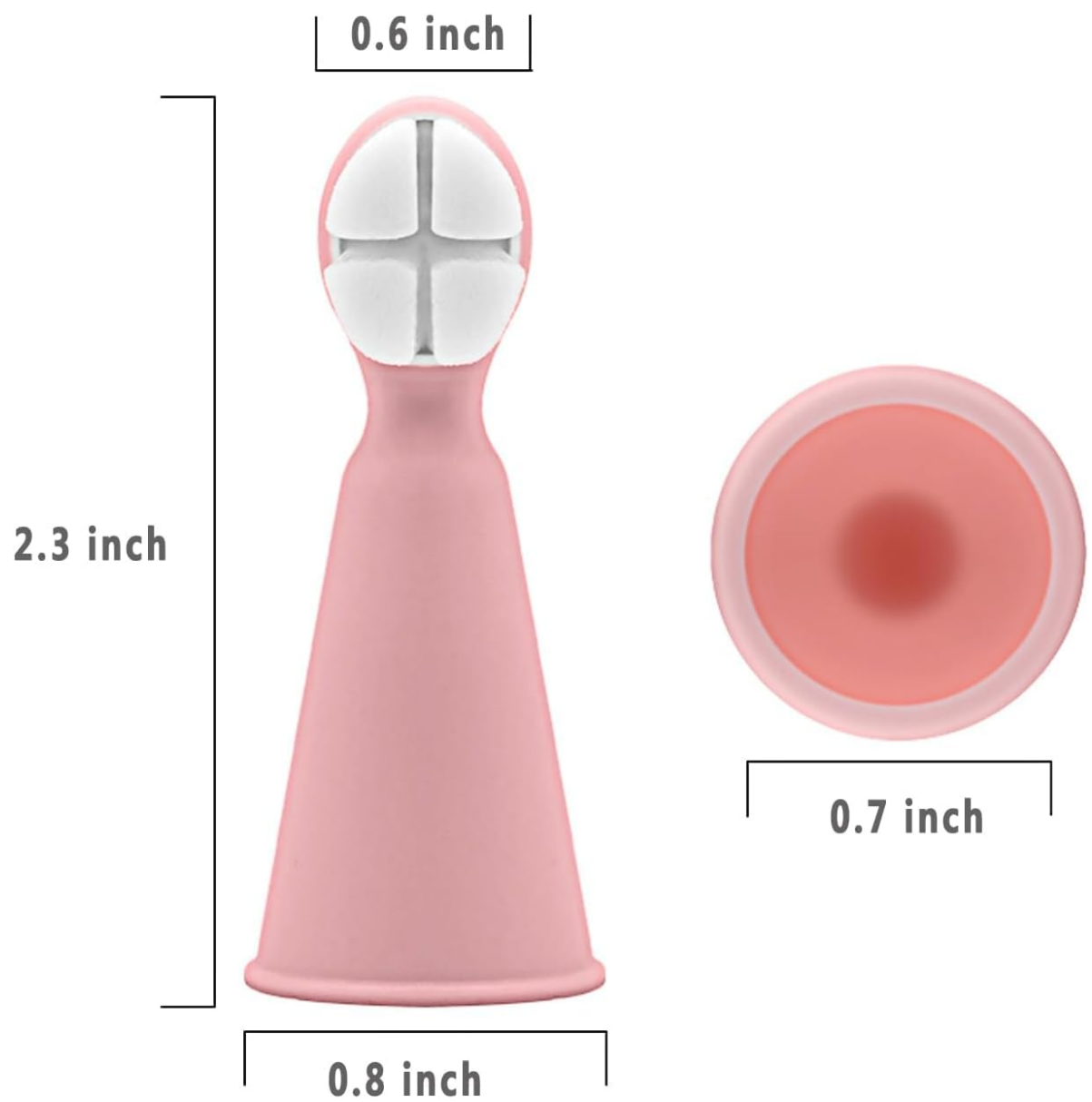


Image 8.1: Product dimensions for the FOREVIVE Baby Finger Toothbrush.

9. WARRANTY AND SUPPORT

For any questions, concerns, or support regarding your FOREVIVE Baby Finger Toothbrush, please contact FOREVIVE customer service. Specific warranty details may vary; please refer to the product packaging or the official FOREVIVE website for the most current information.
