

COMFIER B0FP2BZWBP

COMFIER Shiatsu Foot and Back Massager (Model: B0FP2BZWBP) User Manual

Your guide to safe and effective use of your COMFIER massager.

IMPORTANT SAFETY INSTRUCTIONS

Please read all instructions carefully before using this COMFIER massager. Failure to follow these instructions may result in electric shock, fire, or serious injury.

- Always unplug the appliance from the electrical outlet immediately after using and before cleaning.
- Do not use while bathing or in a shower.
- Do not place or store appliance where it can fall or be pulled into a tub or sink.
- Do not place in or drop into water or other liquid.
- Never use pins or other metallic fasteners with this appliance.
- Close supervision is necessary when this appliance is used by, on, or near children, invalids, or disabled persons.
- Use this appliance only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.
- Never operate this appliance if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water. Return the appliance to a service center for examination and repair.
- Keep the cord away from heated surfaces.
- Never use while sleeping or drowsy.
- Never drop or insert any object into any opening.
- Do not operate where aerosol (spray) products are being used or where oxygen is being administered.
- Do not use on bare skin.
- Do not use on infants, invalids, or on a sleeping or unconscious person.
- Do not use on insensitive skin or on a person with poor blood circulation.
- Do not crush. Avoid sharp folds.
- This appliance is equipped with a heated surface. Persons insensitive to heat must exercise caution when

using the appliance.

- Consult your doctor before using this product if you are pregnant, have a pacemaker, or have any concerns about your health.

Legal Disclaimer: Statements regarding dietary supplements have not been evaluated by the FDA and are not intended to diagnose, treat, cure, or prevent any disease or health condition. This product is not a medical device.

PRODUCT OVERVIEW

The COMFIER Shiatsu Foot and Back Massager is designed to provide a soothing massage experience with optional heat. It features multiple shiatsu massage nodes that rotate to deliver a deep kneading massage, suitable for both feet and the lower back. The device includes a washable cover for hygiene and comfort.



Image 1: Top view of the COMFIER Shiatsu Foot and Back Massager. The device is black with a mesh surface, featuring six sets of red illuminated shiatsu massage nodes arranged in a grid pattern. Each set of nodes has three points. Below the massage area, there is a control panel with three buttons: a power button, a heat button, and a massage intensity/mode button. The COMFIER logo is visible at the top center.

Key Features:

- **2-in-1 Design:** Functions as both a foot massager and a back massager.
- **Shiatsu Kneading Massage:** Multiple rotating nodes provide a deep tissue massage.
- **Soothing Heat Function:** Optional heat therapy to enhance relaxation.
- **Washable Cover:** Removable and washable fabric cover for easy maintenance.
- **Portable Design:** Compact and lightweight for easy storage and transport.

SETUP

1. **Unpack:** Carefully remove the massager from its packaging.
2. **Inspect:** Check the device and power cord for any signs of damage. Do not use if damaged.
3. **Placement:**
 - For foot massage: Place the massager on a flat, stable surface in front of a chair or sofa.
 - For back massage: Place the massager on a chair or sofa, ensuring it is stable and will not slip during use.

4. **Connect Power:** Plug the power adapter into the massager's power input, then plug the other end into a standard electrical outlet (AC 100-240V, 50/60Hz).
5. **Ensure Comfort:** Position yourself comfortably before turning on the device.

OPERATING INSTRUCTIONS

Familiarize yourself with the control panel before operation. The control panel typically includes a Power button, a Heat button, and a Massage Mode/Intensity button.



Image 2: A detailed view of the control panel on the COMFIER Shiatsu Foot and Back Massager. From left to right, there is a button with three vertical lines (indicating massage intensity/mode), a circular power button with a universal power symbol, and a button with a wavy heat symbol. All buttons are illuminated when active.

1. Power On/Off:

- Press the **Power Button** (circular button with universal power symbol) once to turn on the massager. The massage nodes will begin to rotate.
- Press the **Power Button** again to turn off the massager.

2. Heat Function:

- With the massager on, press the **Heat Button** (button with wavy heat symbol) to activate the soothing heat function. The heat indicator light will illuminate.
- Press the **Heat Button** again to turn off the heat function.

3. Massage Mode/Intensity:

- Press the **Massage Mode/Intensity Button** (button with three vertical lines) to cycle through different massage modes or intensity levels. Refer to the product's specific settings for details (e.g., low, medium, high intensity, or different rotation patterns).
- The corresponding indicator light will change to reflect the selected mode/intensity.

4. **Usage Duration:** It is recommended to use the massager for sessions of no more than 15-20 minutes at a time. The device may have an automatic shut-off feature after a certain period for safety.

MAINTENANCE

Cleaning:

- Always unplug the massager before cleaning.
- Allow the massager to cool down completely if the heat function was used.
- The fabric cover is removable and washable. Follow the care instructions on the label for washing.
- Wipe the main unit with a soft, slightly damp cloth. Do not use abrasive cleaners, brushes, gasoline, kerosene, glass/furniture polish, or paint thinner to clean the massager.
- Do not immerse the massager in water or any other liquid.

Storage:

- Store the massager in a cool, dry place.
- Avoid contact with sharp edges or pointed objects which might cut or puncture the fabric.
- Do not wrap the power cord around the appliance.

TROUBLESHOOTING

Problem	Possible Cause	Solution
Massager does not turn on.	<ul style="list-style-type: none"> • Not properly plugged in. • Power outlet not working. • Power button not pressed. 	<ul style="list-style-type: none"> • Ensure the power cord is securely plugged into both the massager and the electrical outlet. • Test the outlet with another device. • Press the Power button firmly.
Heat function is not working.	<ul style="list-style-type: none"> • Heat function not activated. • Device needs time to warm up. 	<ul style="list-style-type: none"> • Ensure the Heat button is pressed and its indicator light is on. • Allow 5-10 minutes for the heat to become noticeable.
Massage feels too intense or too weak.	<ul style="list-style-type: none"> • Incorrect massage mode/intensity selected. • Positioning of body part. 	<ul style="list-style-type: none"> • Press the Massage Mode/Intensity button to adjust the setting. • Adjust your feet or back position on the massager to find a comfortable intensity.
Massager stops unexpectedly.	<ul style="list-style-type: none"> • Automatic shut-off timer activated. • Overheating protection. 	<ul style="list-style-type: none"> • This is normal. Wait a few minutes and restart if desired. • Unplug the device and let it cool down for at least 30 minutes before restarting.

SPECIFICATIONS

- **Model:** B0FP2BZWBP
- **Brand:** COMFIER

- **Power Source:** Corded Electric
- **Material:** Polyester
- **Color:** Black
- **Product Dimensions:** 14.4 x 3.55 x 13.4 inches
- **Item Weight:** 5 Pounds
- **Manufacturer:** Comfier


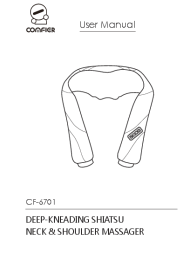

WARRANTY AND SUPPORT


COMFIER products are designed for reliability and performance. For warranty information, product support, or to inquire about replacement parts, please contact COMFIER customer service directly.

Please refer to the official COMFIER website or the packaging for the most up-to-date contact information and warranty details.

Online Support: [Visit the COMFIER Store on Amazon](#)

Related Documents - B0FP2BZWBP

 <p>COMFIER Shiatsu Neck and Back Massager</p> <p>User Manual CF-2309A</p>	<p>Comfier Shiatsu Neck and Back Massager CF-2309A User Manual</p> <p>User manual for the Comfier Shiatsu Neck and Back Massager, model CF-2309A. Includes safety instructions, features, technical data, setup and operation, controller instructions, care and cleaning, troubleshooting, and warranty information.</p>
 <p>COMFIER User Manual</p> <p>CF-6701 DEEP-KNEADING SHIATSU NECK & SHOULDER MASSAGER</p>	<p>Comfier CF-6701 Deep-Kneading Shiatsu Neck & Shoulder Massager User Manual</p> <p>User manual for the Comfier CF-6701 Deep-Kneading Shiatsu Neck & Shoulder Massager, covering product features, technical data, setup, operation, safety instructions, care, troubleshooting, and warranty information.</p>
 <p>COMFIER Leg and Foot Massager</p> <p>User Manual CF-F25H</p>	<p>Comfier CF-F25H Leg and Foot Massager User Manual</p> <p>Comprehensive user manual for the Comfier CF-F25H Leg and Foot Massager, detailing safety instructions, product features, setup and operation, maintenance, troubleshooting, warranty information, and FCC compliance.</p>

<div><div><div>COMFIER</div><div>Shiatsu massage pillow with heat</div></div><div><div><div>CE</div><div>User Manual</div><div>CF-1228</div></div></div></div>	<p>Comfier CF-1228 Shiatsu Massage Pillow with Heat User Manual</p> <p>User manual for the Comfier CF-1228 Shiatsu massage pillow with heat, detailing safety instructions, operating procedures, features, maintenance, disposal, technical specifications, and warranty information.</p>
<div><div><div>COMFIER</div><div>Shiatsu Neck and Back Massager with heat</div></div><div><div><div>CE</div><div>User Manual</div><div>CF-2303-KR</div></div></div></div>	<p>COMFIER CF-2303-KR Shiatsu Neck and Back Massager with Heat User Manual</p> <p>Comprehensive user manual for the COMFIER CF-2303-KR Shiatsu Neck and Back Massager with Heat, covering safety instructions, features, technical data, setup, operation, controller functions, care, cleaning, troubleshooting, and warranty information.</p>
<div><div><div>COMFIER</div><div>Shiatsu massagecushion with heat Shiatsu massage cushion with heat Coussin de massage Shiatsu avec chaleur Cushion per massage Shiatsu riscaldato Cajín masajeador Shiatsu con calefacción</div></div><div><div><div>CE</div><div>User Manual</div><div>CF-2606MC</div></div></div></div>	<p>COMFIER CF-2606MC Shiatsu Massage Cushion with Heat User Manual</p> <p>Comprehensive user manual for the COMFIER CF-2606MC Shiatsu massage cushion with heat. Includes safety instructions, operating guide, features, maintenance, and warranty information.</p>