

Manuals+

[Q & A](#) | [Deep Search](#) | [Upload](#)

manuals.plus /

- › [pedete](#) /
- › [Pedete Red Light Therapy Eye Mask User Manual](#)

pedete Red Light Therapy Eye Mask

Pedete Red Light Therapy Eye Mask User Manual

Model: Red Light Therapy Eye Mask

INTRODUCTION

The Pedete Red Light Therapy Eye Mask is designed to provide targeted light therapy for the delicate skin around the eyes. This device utilizes advanced light technology to address various eye area concerns, promoting a refreshed appearance. This manual provides essential information for the safe and effective use of your eye mask.



Image: The Pedete Red Light Therapy Eye Mask in use, demonstrating its ergonomic design.

SAFETY INFORMATION

Please read all safety instructions before using the device.

- **Disclaimer:** Statements regarding this device have not been evaluated by the FDA and are not intended to diagnose, treat, cure, or prevent any disease or health condition.
- **Discontinue Use:** If you experience any discomfort or observe an abnormal skin reaction during use, immediately stop using the device.
- **Eye Protection:** The unique design of the eye mask incorporates features to block light from entering your eyes, allowing for comfortable use without direct eye strain.
- **Cleaning:** Always clean the silicone part of the mask with a soft towel or a wet wipe after each use.

PACKAGE CONTENTS

- Pedete Red Light Therapy Eye Mask

- Charging Cable
- Storage Bag
- Adjustable Head Strap Extender

SETUP

1. **Charge the Device:** Before first use, ensure the eye mask is fully charged using the provided charging cable. The device contains a 1000 mAh Lithium Ion battery.
2. **Prepare Your Skin:** Cleanse and dry your face thoroughly before applying the mask.
3. **Position the Mask:** Place the eye mask over your eyes.
4. **Adjust the Strap:** Adjust the Velcro head strap to a comfortable and secure position. An extender is included if needed for a larger fit.

OPERATING INSTRUCTIONS

The Pedete Red Light Therapy Eye Mask features three distinct light modes, each designed for specific benefits.

Power On/Off and Mode Selection

- **To Power On:** Press and hold the On/Off button located on the side of the mask.
- **To Select Mode:** Once powered on, click the On/Off button to cycle through the available light modes.
- **Automatic Shut-off:** The system defaults to a 10-minute treatment session and will automatically shut off after this period.
- **Continuous Use:** To continue a session beyond 10 minutes, press the On/Off button again to restart the timer.

Light Modes and Wavelengths

3 Light Modes for Multi Eyecare

Fight dark circles, puffiness, and wrinkles in Just 10 Minutes a day

1 Red Light Mode

Lifting and firming, tighten the skin around eyes



2 Blue Light Mode

Calming and soothing, relieve eye fatigue



3 Yellow Light Mode

Increase skin elasticity and gloss, enhance skin cell vitality



Image: Visual representation of the three light modes: Red, Blue, and Yellow, and their associated benefits.

- **Red Light + Infrared Mode (620-850nm):** This mode is designed for lifting and firming the skin, helping to tighten the area around the eyes.
- **Blue Light Mode (460-465nm):** This mode provides a calming and soothing effect, intended to relieve eye fatigue and reduce swelling.
- **Yellow Light Mode (600-605nm):** This mode aims to increase skin elasticity and gloss, enhancing skin cell vitality.

4 Light Source Wavelengths

Boost collagen production for firmer, smoother skin

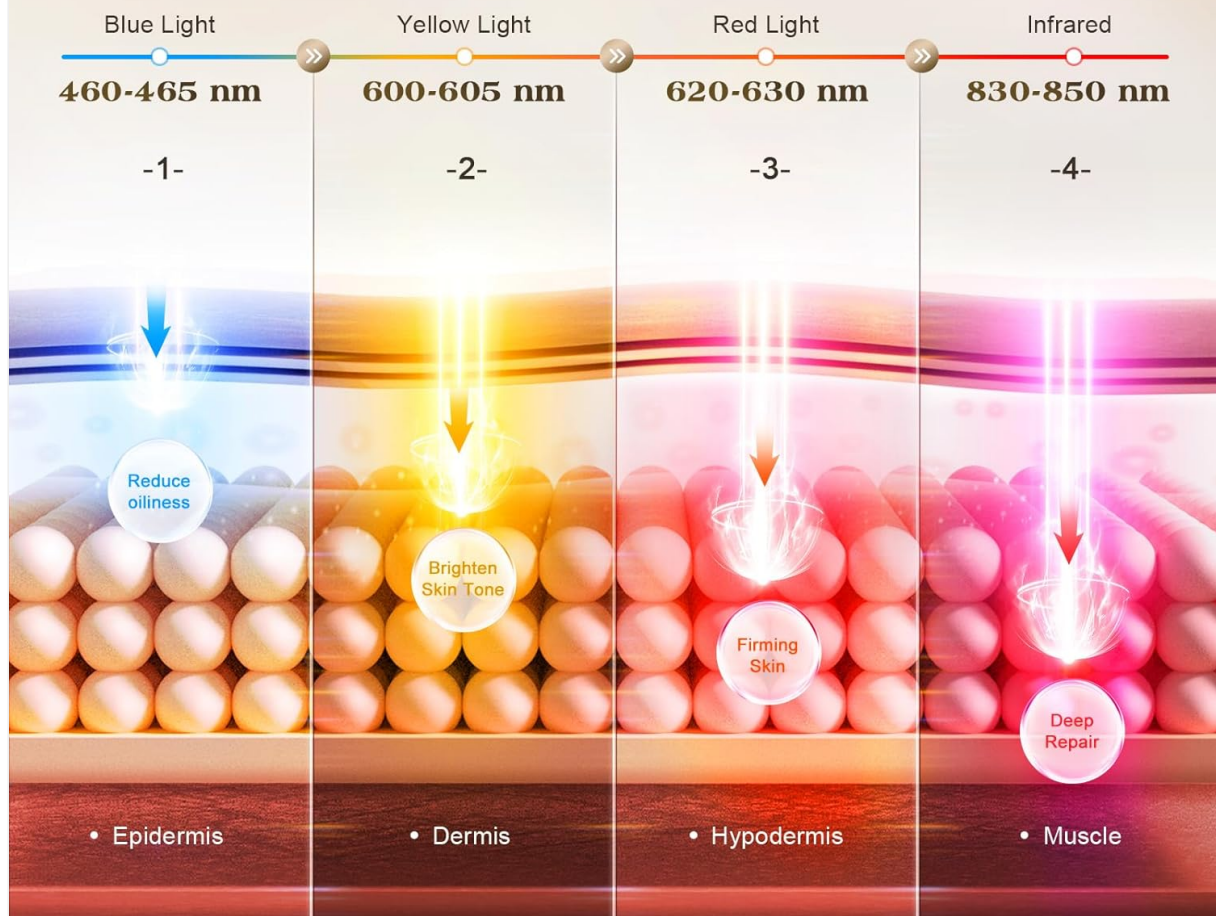



Image: Illustration of different light source wavelengths and their penetration depths into the skin. Each light wavelength penetrates to different depths of the skin, targeting specific concerns.

Design and Eye Protection

The eye mask features a full-coverage design that targets the entire eye area. It includes specially designed eye protection elements to block light from entering your eyes, ensuring comfort and safety during use.

Full Coverage Eye Mask Design

Targets on the full eye area with high-density LED eyecare

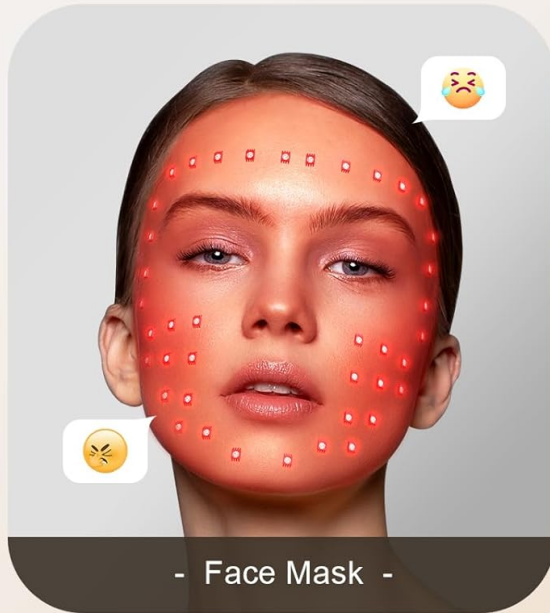

Specially Designed Eye Mask

Blocks light from entering your eyes



Image: Details of the eye mask's full coverage design and integrated eye protection.

Recommended Usage

Use the eye mask for 10 minutes daily. The lightweight and comfortable design allows for use during various activities such as resting, working, reading, or during entertainment time.

Enjoy Anytime Anywhere

Enjoy eye massager anytime anywhere



Rest Time

Promote eye cream absorption



Entertainment Time

Relieve eye fatigue



Working Time

Refresh your eyes



Reading Time

Enhance eye relaxation



Image: Examples of using the eye mask during daily activities.

MAINTENANCE

- **Cleaning:** After each use, gently wipe the silicone part of the eye mask with a soft, damp cloth or a wet wipe. Do not submerge the device in water.
- **Storage:** Store the eye mask in the provided storage bag to protect it from dust and damage when not in use.
- **Charging:** Recharge the device as needed. Avoid overcharging for extended periods.

TROUBLESHOOTING

- **Device Not Turning On:** Ensure the device is fully charged. If the battery is low, connect it to the charging cable.
- **Discomfort During Use:** If you experience any irritation, redness, or discomfort, discontinue use immediately. Consult a healthcare professional if symptoms persist.
- **Light Not Emitting:** Check if the device is powered on and if a mode has been selected. If the issue persists after charging, contact customer support.

SPECIFICATIONS

Brand	pedete
Model Name	Red Light Therapy Eye Mask
ASIN	B0FP2DRZLG
UPC	850077067389
Battery Type	1 Lithium Ion (included)
Battery Capacity	1000 mAh
Red Light Wavelength	620-850 nm (Red + Infrared)
Blue Light Wavelength	460-465 nm
Yellow Light Wavelength	600-605 nm
Default Session Time	10 minutes
Weight	Approximately 5 ounces (141.75 grams)
Manufacturer	pedete

CONTACT INFORMATION

For any questions or support regarding your Pedete Red Light Therapy Eye Mask, please contact the manufacturer or seller through your purchase platform.