

[Manuals.plus](#) /

› [COSLUS](#) /

› COSLUS C50 Portable Water Flosser Instruction Manual

COSLUS C50

COSLUS C50 Portable Water Flosser Instruction Manual

Model: C50

1. IMPORTANT SAFETY INFORMATION

Read all instructions before using this appliance. When using electrical products, especially when children are present, basic safety precautions should always be followed, including the following:

- Do not immerse the unit in water or other liquids.
- Do not use if the power cord or plug is damaged.
- Keep out of reach of children.
- Do not direct water under the tongue, into the ear, nose, or other delicate areas.
- Use this product only for its intended use as described in this manual.

Legal Disclaimer

Statements regarding dietary supplements have not been evaluated by the FDA and are not intended to diagnose, treat, cure, or prevent any disease or health condition.

2. PRODUCT OVERVIEW

The COSLUS C50 Portable Water Flosser is designed for effective oral hygiene. It features 4 cleaning modes and 5 intensity levels, providing 20 unique cleaning combinations for various oral care needs. This device is ADA Accepted and suitable for sensitive gums, braces, and daily oral care. Its cordless and waterproof design makes it convenient for both home and travel use.



Image: COSLUS C50 Portable Water Flosser in black.

Key Features:

- **4 Cleaning Modes:** Clean, Soft, Orthodontic, Jet.
- **5 Intensity Levels:** Adjustable water pressure for personalized comfort.
- **ADA Accepted:** Recognized for effectiveness in removing plaque and preventing gingivitis.
- **IPX7 Waterproof:** Safe for use in the shower and easy to clean.
- **4 Functional Nozzles:** Includes 2 standard, 1 orthodontic, and 1 periodontal tip.
- **Portable Design:** Cordless and compact for travel and on-the-go use.
- **Large Water Tank:** 270ml capacity for uninterrupted cleaning.

3. WHAT'S INCLUDED

The COSLUS C50 Water Flosser package typically includes:

- COSLUS C50 Water Flosser Unit
- 2 Standard Nozzle Tips
- 1 Orthodontic Nozzle Tip
- 1 Periodontal Nozzle Tip
- Charging Cable (USB-A or Type-C, randomly shipped)

Port: Type-C or USB-A (Randomly Shipped)

Please confirm your power port before purchase to avoid inconvenience



- ① If the cable provided does not match your needs, simply contact us and we will send the suitable version to you promptly.
② The adapter is not included to avoid wasting resources, but you can charge it using any 5V adapter on hand.

Image: The COSLUS C50 Water Flosser unit with its various nozzle tips and charging cable.

4. CHARGING THE DEVICE

Before first use, fully charge the water flosser. The device uses either a USB-A or Type-C charging cable. The charging adapter is not included; use any 5V adapter.

1. Locate the charging port on the device.
2. Connect the provided charging cable to the device.
3. Plug the other end of the cable into a compatible 5V USB power adapter (not included).
4. The battery indicator lights will show charging progress. A full charge typically takes a few hours.

Note: Due to a product upgrade, the charging cable type (USB-A or Type-C) may vary. If the provided cable does not match your needs, contact customer support for assistance.



Image: Illustration of the two possible charging port types (Type-C or USB-A) for the COSLUS C50 Water Flosser.

5. SETUP

Follow these steps to prepare your water flosser for use:

1. **Remove the irrigator tip:** Gently pull the tip from the handle if it's already attached.
2. **Fill the water tank:** Open the water tank cap and fill the tank with warm water. Do not exceed 40°C (104°F) for enhanced teeth protection.
3. **Reattach the tip:** Insert your desired nozzle tip firmly into the handle until it clicks into place.
4. **Adjust mode and intensity:** Use the control buttons to select your preferred cleaning mode and intensity level.

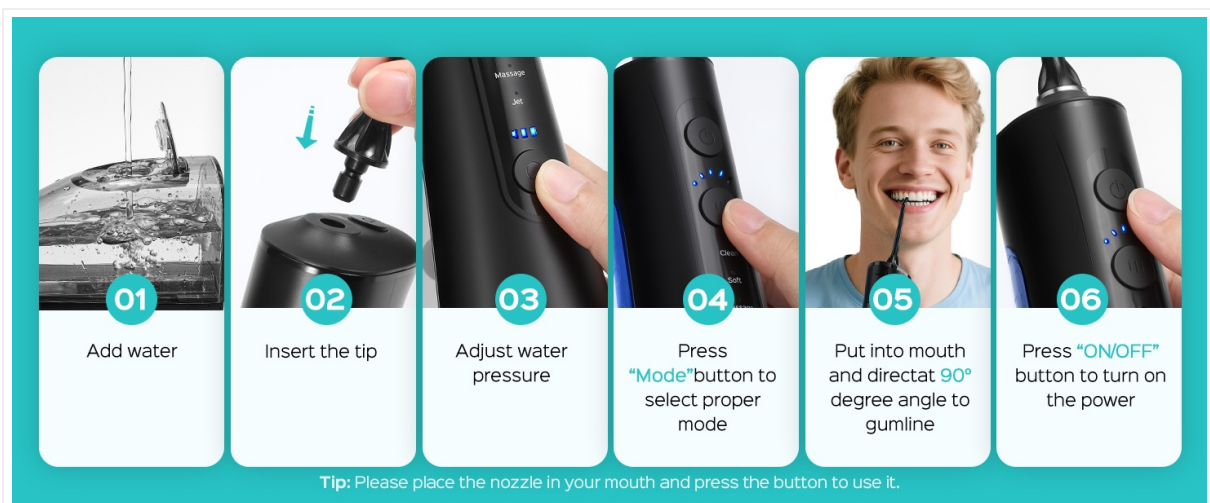


Image: Visual guide demonstrating how to add water, insert the tip, adjust pressure, select mode, position the nozzle, and turn on the device.

6. OPERATING INSTRUCTIONS

The COSLUS C50 offers 4 modes and 5 intensity levels. Select the combination that best suits your oral health needs.

Modes:

- **Clean:** Effectively removes deposits between teeth.
- **Soft:** Gently cleans interdental spaces and removes deposits. Suitable for first-time users, children, and individuals with sensitive gums.
- **Orthodontic:** Cleans and massages gums, ideal for braces users.
- **Jet:** This mode prevents water splashing. To activate, hold the mode button pressed for 2 seconds. The device will auto-off after 2 minutes of inactivity.

Usage Steps:

1. Lean over the sink and place the tip in your mouth, pointing it at your teeth.
2. Turn on the device using the power button.
3. Guide the tip along your gum line, moving it slowly and carefully. Ensure the water stream is directed at a 90-degree angle to your gumline.
4. Keep your mouth partially closed to prevent splashing while allowing water to drain into the sink.
5. Floss all areas of your mouth, including between teeth and along the gumline.
6. Turn off the device when finished.

4 x 5 = 20

Modes Intensities Pressures

Solve the issue of inappropriate water pressure for users



5 Intensity Levels

Can be adjusted as needed



Clean

Effectively removes deposits between teeth



Soft

Gently cleans interdental spaces and removes deposits
Suitable for first-time users, children, and people with sensitive gums



Orthodontic

Cleans and massages gums



Jet

This mode prevents water splashing
To activate, hold the button pressed



Image: Close-up of the COSLUS C50 control panel, highlighting the power button, mode selection, and intensity indicators.

Upgraded Jet Mode

No More Splash



Jet Mode

Minimizes Water Spray Effect



Other

Tip: This mode requires long-pressing the power button and automatically switches off after several seconds of inactivity



Image: A user demonstrating the 'Jet Mode' feature, designed to reduce water splashing during use.

7. NOZZLE TYPES AND USAGE

The COSLUS C50 comes with four specialized nozzle tips to cater to different oral care needs:

- **Standard Nozzle (x2):** For daily deep cleaning and plaque removal.
- **Orthodontic Nozzle (x1):** Designed for cleaning around braces wires and orthodontic appliances.
- **Periodontal Nozzle (x1):** Provides gentle care for gum pockets and sensitive areas.

4 Multifunctional Nozzle Tips



x2 Standard nozzle
For daily dental cleaning



x1 Orthodontic nozzle
For braces users



x1 Periodontal nozzle
Deeply cleans between teeth, ideal for periodontitis care

360°
Rotation



Image: A visual representation of the four multifunctional nozzle tips, detailing their specific uses.

8. MAINTENANCE AND CLEANING

Proper maintenance ensures the longevity and hygiene of your water flosser.

1. **After Each Use:** Remove the nozzle tip and rinse it under clean water.
2. **Clean the Water Tank:** Empty any remaining water from the tank. Rinse the tank thoroughly with clean water.
3. **Air Dry:** Allow both the tip and the water tank to air dry completely before reassembling or storing.
4. **Wipe Down Unit:** Wipe the main unit with a soft, damp cloth. Do not use abrasive cleaners or harsh chemicals.
5. **Storage:** Store the device with the tip and tank attached in a clean, dry place.



Larger Removable Water Tank

20% Larger Capacity Than Comparable Devices



Ergonomically Designed
Convenient Grip



BPA Accepted
Safe and protected

Tip : Water temperature should not exceed 40°C (104°F) for enhanced teeth protection

Image: The COSLUS C50 Water Flosser highlighting its larger, removable water tank for easy filling and cleaning.

9. TROUBLESHOOTING

If you encounter issues with your COSLUS C50 Water Flosser, consider the following general troubleshooting steps:

- **Device not turning on:** Ensure the device is fully charged. If the battery is low, connect it to the charger.
- **Weak water pressure:** Check if the water tank is sufficiently filled. Ensure the nozzle tip is securely attached. Clean the nozzle tip for any blockages.
- **Water leakage:** Verify that the water tank is correctly seated and the cap is tightly closed. Ensure the nozzle tip is properly inserted.
- **Unusual noise:** If the device makes unusual noises, turn it off immediately and contact customer support.

For persistent issues, refer to the warranty and support section for contact information.

10. SPECIFICATIONS

Feature	Specification
Model Number	C50
Product Dimensions	8.97 x 3.5 x 2.55 inches

Item Weight	0.94 Pounds (15.04 ounces)
Water Tank Volume	270 Milliliters
Power Source	Battery Powered (1 Nonstandard Battery included)
Waterproof Rating	IPX7
Modes	4 (Clean, Soft, Orthodontic, Jet)
Intensity Levels	5
Material Type	Plastic
Color	Black
Certification	ADA Accepted

11. WARRANTY AND SUPPORT

The COSLUS C50 Portable Water Flosser comes with a **1-year manufacturer's warranty**. Additionally, COSLUS offers **24-hour after-sales service** and a **free 30-day refund/replacement** policy.

For warranty claims, technical support, or any product-related inquiries, please contact COSLUS customer service through the retailer where the product was purchased or visit the official COSLUS website for contact details.