

## Manuals+

[Q & A](#) | [Deep Search](#) | [Upload](#)

[manuals.plus](#) /

> [mornwell](#) /

> [Mornwell Water Dental Flosser F5106 User Manual](#)

## mornwell F5106

# Mornwell Water Dental Flosser F5106 User Manual

Model: F5106

[Overview](#)   [Setup](#)   [Operation](#)   [Introduction](#)   [Maintenance](#)   [Safety Information](#)   [Troubleshooting](#)   [Product](#)   [Specifications](#)   [Warranty & Support](#)

## 1. INTRODUCTION

The Mornwell Water Dental Flosser F5106 is designed to provide effective oral hygiene by delivering a pulsating water stream to clean between teeth and below the gumline. This device features 2 modes and 10 adjustable pressure settings, making it suitable for various users, including those with braces, implants, or sensitive gums. The 600ML water tank ensures a complete cleaning session without frequent refills.

## 2. IMPORTANT SAFETY INFORMATION

### 2.1 General Safety Precautions

- Read all instructions before using this appliance.
- Do not immerse the unit in water or other liquids.
- Do not use this product if you have an open wound on your tongue or in your mouth.
- Keep out of reach of children.
- Use only as directed in this manual.

### 2.2 Electrical Safety

- Ensure the power adapter is compatible with your local power supply.
- Unplug the device before cleaning or when not in use for an extended period.
- Do not handle the power plug with wet hands.

## 3. PRODUCT OVERVIEW AND COMPONENTS

The Mornwell Water Dental Flosser F5106 consists of several key components designed for ease of use and effective cleaning.

# Product Details



Figure 1: Exploded view of the Mornwell Water Dental Flosser F5106, highlighting the 600ML water tank, anti-leak design, magnetic attraction for the handle, 10 adjustable pressure button, ON/OFF switch, and non-slip design.

## 3.1 Main Components

- **Main Unit:** Houses the motor, pump, and control panel.
- **600ML Water Tank:** Detachable for easy refilling and cleaning. Features measurement markings.
- **Handle with ON/OFF Switch:** Ergonomically designed for comfortable grip, includes a switch to control water flow.
- **Pressure Adjustment Dial:** Allows selection of 10 pressure levels and modes (Clean/Massage).
- **Nozzles Storage Box:** Integrated into the lid for hygienic storage of tips.
- **Power Adapter:** For connecting the device to an electrical outlet.

# WHAT YOU GET

- 1 1 x Box
- 2 1 x Water Dental Flosser
- 3 1 x Adapter
- 4 1 x User Manual
- 5 8 x Nozzles



Figure 2: Contents of the Mornwell Water Dental Flosser F5106 package, including the main unit, adapter, user manual, and 8 nozzles.

## 3.2 Included Nozzles (8 Multi-Functional Tips)



# 8 NOZZLES

with Different Functions

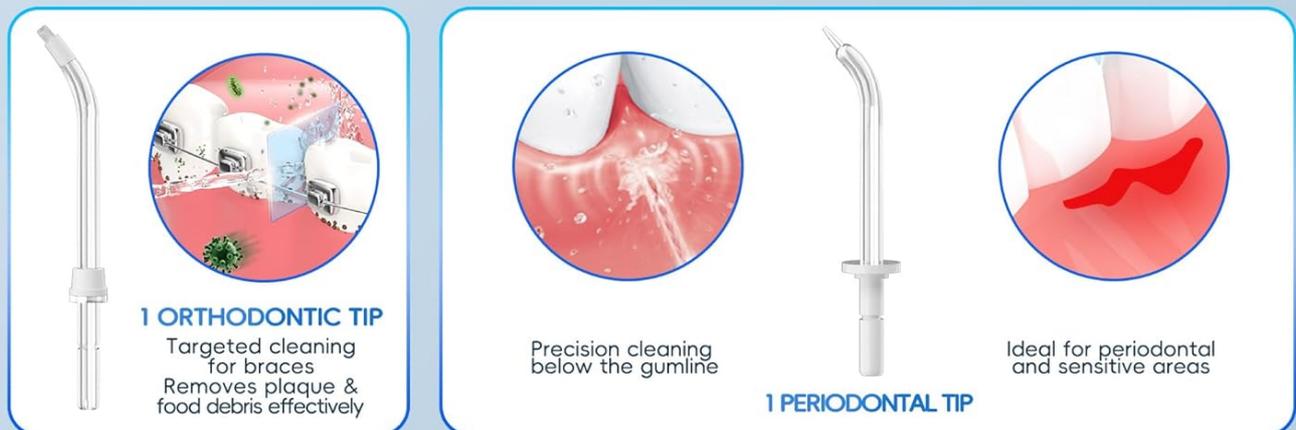


Figure 3: Illustration of the 8 different nozzles included with the Mornwell Water Dental Flosser F5106, detailing their specific functions.

- **Standard Tips (4):** For daily cleaning of teeth and gumline with an ultra-fine jet stream.
- **Tongue Cleaner Tip (1):** Designed to clean the tongue effectively, reducing bad breath.
- **Toothbrush Tip (1):** Combines brushing and water flossing for comprehensive daily oral care.
- **Orthodontic Tip (1):** For targeted cleaning around braces, effectively removing plaque and food debris.
- **Periodontal Tip (1):** For precision cleaning below the gumline and in periodontal pockets, ideal for sensitive areas.

## 4. SETUP INSTRUCTIONS

Follow these steps for the initial setup of your Mornwell Water Dental Flosser F5106:

1. **Connect Power:** Plug the power adapter into the port at the base of the unit, then plug the adapter into a suitable electrical outlet.
2. **Fill Water Tank:** Remove the water tank from the base. Fill the 600ML tank with lukewarm water. You may also use mouthwash as recommended by your dentist. Ensure the tank is securely placed back on the base.

3. **Select and Attach Tip:** Choose the desired jet tip from the nozzles storage box. Insert the tip firmly into the handle until it clicks into place.
4. **Prime the Unit (First Use):** Direct the tip into the sink. Turn the unit ON using the power button on the base. Set the pressure to the lowest setting (1). Turn the handle switch to ON to allow water to flow until a steady stream is achieved. Then turn the handle switch OFF.

Your browser does not support the video tag.

Video 1: Demonstration of Mornwell Water Dental Flosser F5106 setup and initial use, including connecting power, filling the water tank, and attaching a tip.

## 5. OPERATING INSTRUCTIONS

To use your Mornwell Water Dental Flosser F5106 effectively, follow these steps:



Figure 4: Step-by-step visual guide on how to use the Mornwell Water Dental Flosser F5106, from adding water to cleaning teeth.

1. **Power On:** Press the power button on the main unit to turn on the device. The indicator light for the selected mode (Clean or Massage) will illuminate.
2. **Select Mode and Pressure:** Rotate the pressure adjustment dial to select your desired pressure level (1-10) and mode (Clean or Massage). Start with a lower pressure setting, especially if you have sensitive gums or are a new user.
3. **Position Tip:** Lean over the sink and place the tip in your mouth. Close your lips enough to prevent splashing, but allow water to flow from your mouth into the sink.
4. **Begin Flossing:** Turn the switch on the handle to the ON position to start the water flow. Aim the water stream at a 90-degree angle to your gumline. Move the tip slowly along your teeth, pausing briefly between teeth and cleaning both the front and back surfaces.
5. **Complete Cleaning:** Continue until you have cleaned all areas of your mouth. The 600ML tank typically provides enough water for a thorough cleaning.
6. **Power Off:** Once finished, turn the switch on the handle to OFF, then press the power button on the main unit to turn off the device.

### 5.1 Modes and Pressure Levels

The device offers two primary modes and 10 adjustable pressure levels:

# 10 Adjustable Levels

## Fine-Tuned Water Pressure Control

Rapid startup and ultra-smooth knob design deliver effortless control with no lag or stiffness.



Figure 5: Guide to the 10 adjustable water pressure levels and their recommended uses for different users and cleaning needs.

- **Clean Mode:** Provides a steady, powerful pulse for deep cleaning and plaque removal.
- **Massage Mode:** Offers a gentle pulse to stimulate and massage gums, promoting gum health.

### Pressure Recommendations:

- **Levels 1-2:** Ideal for new users or children, providing a very gentle clean.
- **Levels 3-5:** Suitable for sensitive gums.
- **Levels 6-8:** Recommended for daily cleaning.
- **Levels 9-10:** For deep cleaning and reaching blind areas, suitable for experienced users.

## 6. MAINTENANCE AND CLEANING

Regular cleaning and maintenance ensure the longevity and hygienic operation of your water flosser.

### 6.1 Daily Cleaning

## Daily Cleaning

- **Empty Water Tank:** After each use, empty any remaining water from the tank.
- **Rinse Tank:** Rinse the water tank with clean water.
- **Wipe Unit:** Wipe the exterior of the main unit and handle with a soft, damp cloth.
- **Remove and Rinse Tip:** Press the tip eject button on the handle to remove the tip. Rinse the tip under running water.

## 6.2 Weekly Cleaning

- **Deep Clean Tank:** The 600ML tank is detachable for thorough cleaning. Wash it with mild soap and water, then rinse completely.
- **Clean Internal Components:** To prevent mineral deposits, fill the tank with 2 parts water and 1 part white vinegar. Run the solution through the unit until the tank is empty. Rinse by running a full tank of clean water through the unit.

**10X More Effective Plaque Removal**  
Healthier Gums Start Here

**Up to 1400**  
Pulsations per minute

**0.6 mm micro-bubble**  
Water column

**30-110 psi**  
Strong water pressure

**mornwell**

EASY TO UP TO **10X**

**1200**  
Pulsations per minute

**0.66 mm**  
Water column

**40-90 psi**  
Strong water pressure

**others**

EASY TO UP TO **5X**

Figure 6: The detachable 600ML water tank, designed for easy refilling and cleaning.

## 6.3 Tip Replacement

## 6.3 Tip Replacement

It is recommended to replace the standard jet tips every 3 months, and specialized tips (orthodontic, periodontal, tongue cleaner, toothbrush) every 6 months, or sooner if they show signs of wear.

## 7. TROUBLESHOOTING

If you encounter issues with your Mornwell Water Dental Flosser F5106, refer to the following common problems and solutions:

Problem	Possible Cause	Solution
Unit does not turn on	Not plugged in; Power outlet issue; Unit malfunction	Ensure power adapter is securely plugged into the unit and a working electrical outlet. Test the outlet with another device. If problem persists, contact customer support.
No water flow	Water tank empty; Tip not properly inserted; Handle switch OFF; Unit not primed	Fill the water tank. Ensure the tip is firmly clicked into place. Turn the handle switch to ON. Prime the unit as described in the Setup section.
Weak water pressure	Low pressure setting; Mineral buildup; Clogged tip	Increase the pressure setting using the dial. Perform a deep clean with vinegar solution. Replace the tip if it appears clogged or damaged.
Water leaks from unit	Water tank not seated correctly; Damaged O-ring; Loose connection	Ensure the water tank is securely placed on the base. Check for any visible damage to the O-ring or connections. Contact customer support if leaks persist.

If you experience any issues not listed here or if the suggested solutions do not resolve the problem, please contact Mornwell customer support for assistance.

## 8. SPECIFICATIONS

Feature	Detail
Model Number	F5106
Water Tank Capacity	600ML
Pressure Settings	10 Adjustable Levels
Modes	Clean, Massage
Pulse Rate	1,250–1,700 pulses per minute
Pressure Range	10–120 PSI
Included Tips	8 Multi-Functional Tips
Power Source	Corded Electric
Package Dimensions	8.9 x 6.81 x 5.83 inches
Item Weight	2.07 Pounds

Feature	Detail
Manufacturer	Prooral Health Technology (Shenzhen) Ltd

## 9. WARRANTY AND SUPPORT

---

### 9.1 Product Warranty

Mornwell products are manufactured to the highest quality standards. This product is covered by a limited warranty against defects in materials and workmanship from the date of purchase. Please retain your proof of purchase for warranty claims.

The warranty does not cover damage caused by misuse, accident, unauthorized modification, or normal wear and tear. For specific warranty terms and conditions, please refer to the warranty card included in your product packaging or visit the official Mornwell website.

### 9.2 Customer Support

If you have any questions about your Mornwell Water Dental Flosser F5106, require technical assistance, or need to make a warranty claim, please contact our customer support team:

- **Email:** [support@mornwell.com](mailto:support@mornwell.com) (Example)
- **Phone:** 1-800-XXX-XXXX (Example)
- **Website:** [www.mornwell.com/support](http://www.mornwell.com/support) (Example)

Please have your model number (F5106) and purchase date ready when contacting support.