

## FAQ 501

# FAQ 501 Red LED Face Light Instruction Manual

Model: 501

## 1. PRODUCT OVERVIEW

The FAQ 501 Red LED Face Light is an at-home beauty device designed to enhance skin vitality and appearance through concentrated light innovation. This device supports a refreshed, glowing, and smoother-looking complexion. Key features include:

- **Visible Skin Boost:** Utilizes concentrated light to enhance natural skin vitality.
- **Smoother Looking Finish:** Helps refine surface appearance and minimize texture irregularities.
- **Balanced Tone & Radiance:** Supports a healthy-looking glow by improving tone and reducing dullness.
- **Calming Experience:** Contributes to a balanced look and feel with each use.
- **Intuitive At-Home Design:** Lightweight, USB-rechargeable, and app-connected for guided sessions.



Image: The FAQ 501 Red LED Face Light device, showcasing its sleek design.

## 2. SAFETY INFORMATION

**Please read all instructions carefully before using the FAQ 501 device.**

- If you have a skin condition or any medical concerns, consult a dermatologist before use.
- FAQ treatments should be comfortable. If you experience discomfort, discontinue use and consult your physician.
- For hygiene reasons, do not share your device with others.
- Avoid leaving your device in direct sunlight and never expose it to extreme heat or boiling water.
- Discontinue use if the product appears damaged in any way.
- Use this device for its intended purpose only.
- Always wear the provided protective glasses during treatment to shield your eyes from the LED light.

## 3. WHAT'S IN THE BOX

Your FAQ 501 package includes the following items:

## CONTENU DE LA BOÎTE

Guide de  
démarrage  
rapide




Lunettes de protection



Pochette



Câble  
de charge USB

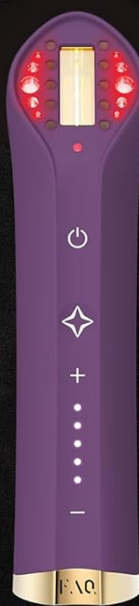

FAQ 501

Image: All items included in the FAQ 501 product box.

- FAQ 501 Device
- USB Charging Cable
- Protective Glasses
- Protective Pouch
- Quick Start Guide

## 4. SETUP AND CHARGING

### 4.1 Initial Charge

Before first use, fully charge your FAQ 501 device. Connect the provided USB charging cable to the device and a suitable USB power adapter. The indicator lights on the device will show charging progress.

## MODO DE USO

1



### PREPARACIÓN

Comience con la piel limpia y seca, y utilice gafas protectoras. Aplique el sérum FAQ Red Light Peptide para un efecto adicional.

2



### SELECCIONE INTENSIDAD

Encienda con el botón de encendido y ajuste con + / -. Explore sesiones guiadas a través de la aplicación FAQ Swiss.

3



### BRILLO & RENOVACIÓN

Coloque sobre la piel. Mantenga presionado el flash para áreas grandes, presione para áreas pequeñas.

Image: The FAQ 501 device being charged via its USB cable.

## 4.2 App Connection

For guided sessions and full functionality, download the FAQ Swiss app. The device is designed to connect with the app for personalized treatment plans and intensity adjustments.

- Download the FAQ Swiss app from your device's app store.
- Follow the in-app instructions to pair your FAQ 501 device.

## BENEFICIOS ANTI-EDAD



**La piel se ve más suave y firme, con un aspecto más refinado y juvenil.**

Image: The FAQ 501 device next to a smartphone running the FAQ Swiss app.

## 5. OPERATING INSTRUCTIONS

Follow these steps for an effective treatment with your FAQ 501 device:

### 1. 5.1 Preparation

Begin with clean, dry skin. It is recommended to apply the [FAQ Red Light Peptide Serum](#) (sold separately) for enhanced light absorption and optimal results. Always wear the provided protective glasses before starting your treatment.



Image: A woman touching her face, demonstrating skin preparation before using the device.

## 2. **5.2 Select Intensity**

Turn on the device using the power button. Adjust the light intensity using the '+' and '-' buttons. You can explore guided sessions and recommended intensity levels through the FAQ Swiss app.



Image: A hand holding the FAQ 501 device, illustrating the power and intensity controls.

## 3. **5.3 Treatment Application**

Place the device gently on your skin. For larger areas, press and hold the flash button to activate continuous pulsed light. For smaller, targeted areas, press the flash button once for a single pulse.



Image: A woman using the FAQ 501 device on her cheek, wearing protective eyewear.

### LUZ DE ESPECTRO COMPLETO

Combina luz ancha y estrecha para mejorar la claridad general de la piel.

### LED ROJO

Favorece una apariencia más energizada y radiante.

### BOTÓN FLASH

Activa la luz pulsada: presione una vez para modo de sellado, mantenga presionado para modo deslizante.

### IMPULSO DE LUZ PULSADA

Estimula suavemente la piel para un aspecto revitalizado y fresco.



Image: Diagram illustrating the full-spectrum light, red LED, flash button, and pulsed light impulse features of the FAQ 501.

## FAQ 501

La luz roja de espectro completo actúa en armonía con el sérum para refrescar y energizar su piel.

## FAQ Red Light Peptide Serum

La fórmula conductora mejora la absorción de la luz para una piel visiblemente más suave, equilibrada y radiante.



\*SE VENDE POR SEPARADO

Image: The FAQ 501 device shown with the FAQ Red Light Peptide Serum, highlighting their complementary use.

## 6. MAINTENANCE

### 6.1 Cleaning

After each use, clean your FAQ 501 device. Use a soft, damp cloth with a mild, non-abrasive cleanser. Ensure the device is completely dry before storing.

### 6.2 Storage

Store the device in its protective pouch in a cool, dry place, away from direct sunlight. Do not expose the device to extreme temperatures or boiling water.

## 7. TROUBLESHOOTING

- **Device Not Turning On:** Ensure the device is fully charged. If the issue persists, contact customer support.
- **Discomfort During Use:** If you experience any discomfort, immediately discontinue use and consult a

physician.

- **Visible Damage:** If the device appears damaged, discontinue use and do not attempt to repair it yourself. Contact customer support.
- **App Connectivity Issues:** Ensure your device's Bluetooth is enabled and the FAQ Swiss app is up to date. Try restarting both the device and the app.

## 8. SPECIFICATIONS

| Feature            | Detail                             |
|--------------------|------------------------------------|
| Model              | FAQ 501                            |
| Color              | Berry                              |
| Power Source       | 1 Lithium Ion battery (included)   |
| Item Weight        | 0.55 Kilograms (1.21 Pounds)       |
| Manufacturer       | VIS (SHANGHAI) TECHNOLOGY CO., LTD |
| Package Dimensions | 7.87 x 5.55 x 3.31 inches          |

## 9. SUPPORT

For further assistance, product registration, or to access the full product manual, please visit the official FAQ Swiss support page or download the app:

- **Support Website:** [faqswiss.com/support](https://faqswiss.com/support)
- **FAQ Swiss App:** [faqswiss.com/app](https://faqswiss.com/app)