

AUVON TU1080

AUVON Rechargeable TENS Unit Muscle Stimulator (Model TU1080) - Instruction Manual

Your guide to safe and effective pain relief.

1. IMPORTANT SAFETY INFORMATION

WARNING: This product is not suitable for individuals with certain medical conditions or circumstances. Consult your doctor before use if you have any medical conditions.

- Do not use if you have an implanted defibrillator, cardiac pacemaker, or any other implanted metallic or electronic device.
- Do not use if you have heart rhythm problems, epilepsy, brain diseases, or neurological disorders.
- Do not use if you have undiagnosed chronic pain.
- Do not use if you are pregnant or in labor.
- Do not use if you have hypertension (high blood pressure) or thrombosis (blood clots).
- Do not use if you have cancer.
- Do not use if you are sensitive to electrotherapy.
- Do not use excessively.
- Do not use above the neck, especially on the head or neck area.
- Do not use if you are under 22 years of age.

Individual symptoms, situations, circumstances, and results may vary. The information provided is not intended for medical diagnosis or treatment or as a substitute for professional medical advice. Please consult your physician or a qualified healthcare provider regarding your condition and appropriate medical treatment.



Live More with Less Pain

AUVON user-friendly TENS unit offers a drug-free pain relief method with an easy operation, helping you live a healthy, pain-free life.

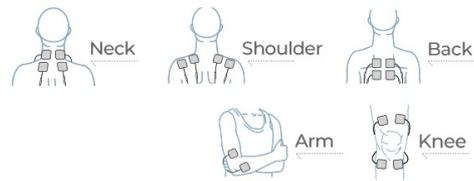


Image: Visual representation of areas where the TENS unit should not be applied, emphasizing safety precautions.

2. PACKAGE CONTENTS

Your AUVON TENS Unit package includes the following items:

- 1 x AUVON TENS Muscle Relaxer (Main Device)
- 8 x Square Electrode Pads (2"x2" Premium)
- 2 x Lead Wires (1/2 Wire)
- 1 x Pads Holder
- 1 x Type-C Cable
- 1 x Pouch
- 1 x Professional User Manual

What You Receive

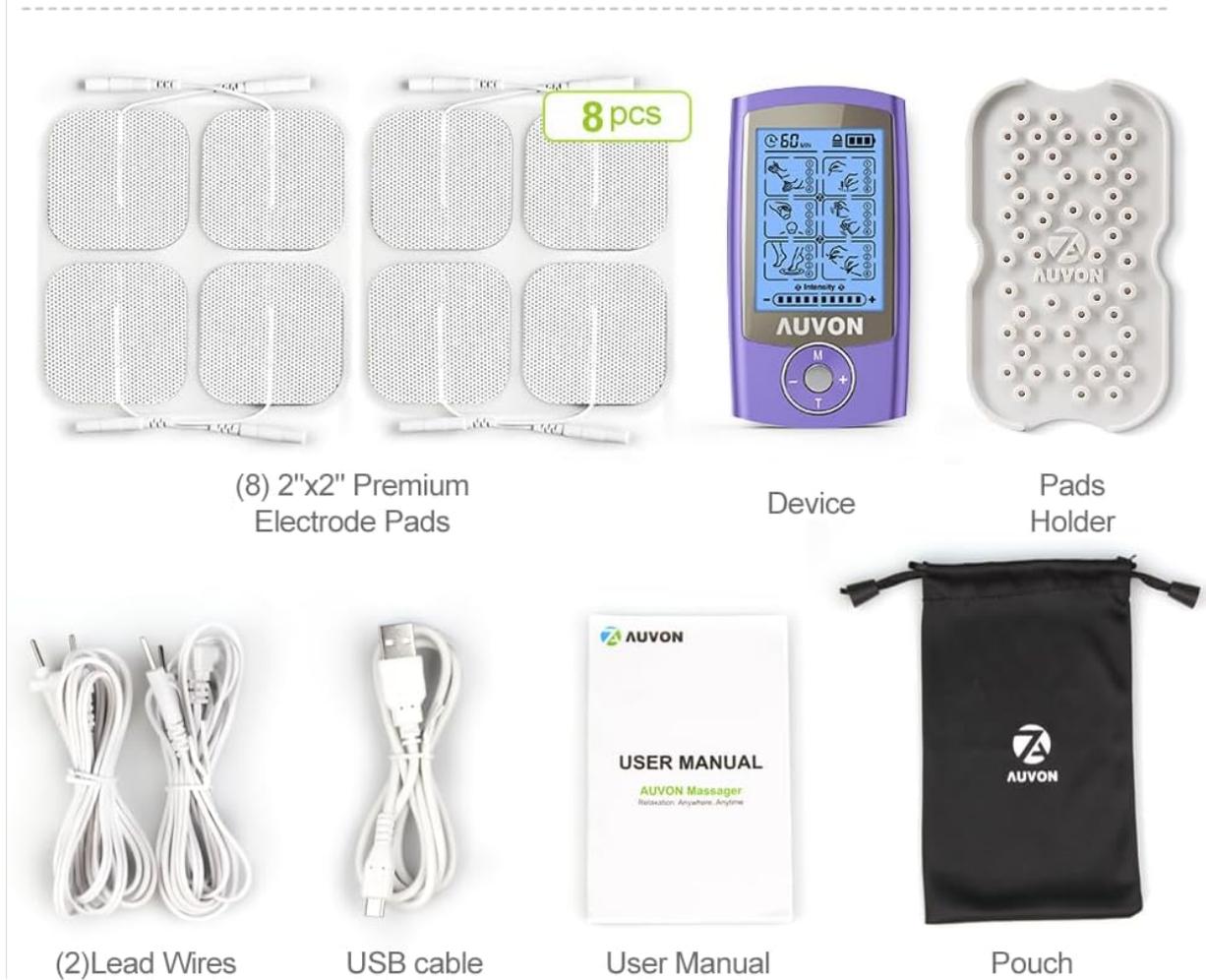


Image: All components included in the AUVON TENS Unit package.

3. PRODUCT OVERVIEW

The AUVON Rechargeable TENS Unit Muscle Stimulator is a professional-grade electrotherapy device designed for pain relief and muscle relaxation. It features 24 user-friendly pre-programmed massage modes and adjustable intensity levels.

Key Features:

- **24 Pre-programmed Modes:** Offers a variety of massage techniques for different needs.
- **20 Intensity Levels:** Allows for precise adjustment of stimulation strength.
- **Rechargeable Lithium Battery:** Provides up to 10 hours of continuous use on a single charge.
- **Power Indicator:** Keeps you informed about remaining battery life.
- **Upgraded Gel Pads:** Low-impedance solid gel for premium self-stick performance and extended lifespan (up to 45 uses).

Up to 24 Mode Choices

24 mode options and controls to give you a relaxing massage with little fuss.



24 Modes



20 Intensity



Pocket Size




10h Usage



Steady Current



Rechargeable

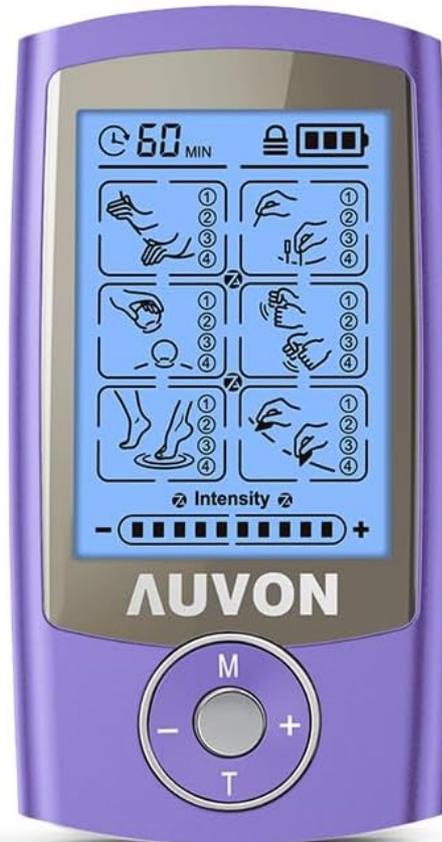

* The 24 sub-mode choices selection is under the 6 massage modes.

Image: Overview of the AUVON TENS Unit's features, highlighting its 24 modes and rechargeable battery.

4. SETUP

4.1 Connecting Electrode Pads

1. Ensure the TENS unit is turned OFF.
2. Connect the lead wires to the electrode pads. Each lead wire typically has two connectors for two pads.
3. Carefully remove the protective film from the electrode pads.
4. Apply the pads to clean, dry skin on the desired treatment area. Ensure pads are placed at least one inch apart.
5. Plug the lead wires into the output ports on the top of the TENS unit. The device has two output ports (A and B) allowing for up to 4 pads to be used simultaneously.

4.2 Charging the Device

The AUVON TENS Unit is equipped with a rechargeable lithium battery. To charge:

1. Turn the unit ON first.
2. Connect the provided Type-C USB cable to the charging port on the device.
3. Plug the other end of the USB cable into a compatible wall charger (not included) or a USB port.

- The power indicator on the screen will show charging status. A full charge typically takes about 1 hour and provides up to 10 hours of continuous use.

Set Your TENS Easily

Step 1



Main Mode: Press "M"
Sub-mode: Press round button

Step 2



Intensity: Press +/-

Step 3



Timer: Press T
Lock: Long press the round button

***Important: Do NOT use if you have a pacemaker, are pregnant, or have epilepsy.**

Image: The AUVON TENS Unit being charged via its Type-C USB cable.

5. OPERATING INSTRUCTIONS

The AUVON TENS Unit is designed for ease of use with pre-programmed settings.

5.1 Turning On/Off

- To turn ON the device, slide the power switch located on the top of the unit. The screen will illuminate.
- To turn OFF, slide the power switch back to its original position.

5.2 Selecting Modes

- Press the 'M' button (Mode) to cycle through the 6 main massage modes.
- Within each main mode, press the round center button to select one of the 4 sub-modes (totaling 24 sub-modes). The selected mode will blink on the screen.

5.3 Adjusting Intensity

- Use the '+' button to increase the intensity of the electrical stimulation.
- Use the '-' button to decrease the intensity.
- Adjust the intensity to a comfortable level where you feel a strong but non-painful sensation.

5.4 Setting the Timer

- Press the 'T' button (Timer) to adjust the treatment duration. Available times typically range from 10 to 60 minutes.
- The device will automatically shut off once the set time is complete.

Rechargeable TENS Unit

The long-lasting battery frees you from the worry of the frequently dead battery and gives you uninterrupted therapy.



Image: Simplified guide to setting up the AUVON TENS Unit's modes, intensity, and timer.

6. ELECTRODE PAD CARE

Proper care of your electrode pads ensures optimal performance and longevity.

- The upgraded low-impedance solid gel pads are designed for premium self-stick performance and can be reused up to 45 times with proper care.
- After each use, place the pads back onto the provided pad holder or the plastic film to maintain their adhesive quality.
- If the pads lose their stickiness, gently wipe the adhesive surface with a damp cloth or a few drops of water. Allow them to air dry completely before reapplying to the holder.
- Store pads in a cool, dry place away from direct sunlight.
- Replacement pads (2mm pigtail TENS Unit pads) are widely available.



Good Stickiness

Upgrade Gel, More Reuses

Lower Impedance, Better Conductivity.

	AUVON GEL	OTHERS GEL
Number of reuse	45 times	30 times
Resistance	45 Ω	50 Ω
Biocompatibility	✓	×
Production Tool	Automated Machine	Hands
Certificate Report	510K Biocompatibility Report	×

Image: Details on the durability and quality of AUVON's upgraded gel pads.

7. APPLICABLE USERS

The AUVON TENS Unit is suitable for a wide range of individuals experiencing muscle soreness, tension, and everyday aches. This includes office workers, individuals who are overworked, homemakers, fitness enthusiasts, and middle-aged or elderly people.

Applicable User

Office Worker / Overworked /
Housewife / Fitness /
Middle aged and
elderly people etc.

* Important: Do NOT use if you have a
pacemaker, are pregnant, or have epilepsy.

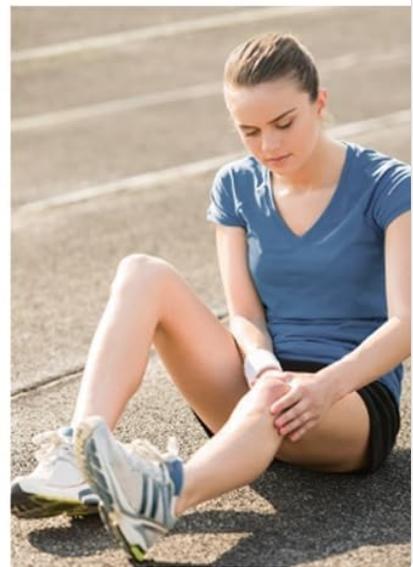
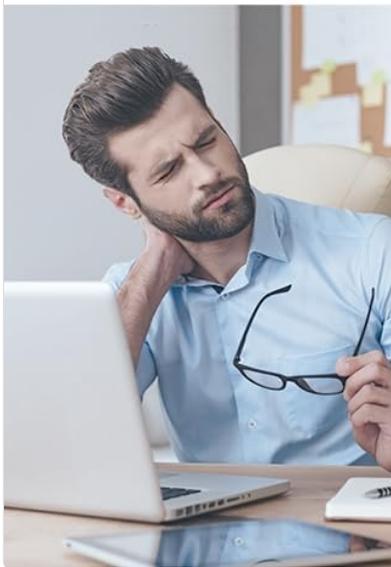


Image: Examples of individuals who may find the AUVON TENS Unit beneficial for pain relief.

8. SPECIFICATIONS

Model Number	TU1080
Package Dimensions	5.71 x 3.31 x 2.52 inches
Item Weight	8.47 ounces
Batteries	1 Lithium Ion battery (included)
Manufacturer	AUVON
Date First Available	September 15, 2025

9. WARRANTY AND SUPPORT

AUVON provides a 12-month warranty for the TENS Unit Muscle Stimulator. For any technical support or customer service inquiries, please refer to the contact information provided in your Professional User Manual or visit the AUVON official website.

