

[Manuals.plus](#) /

› [Myoglow](#) /

› MyoGlow LED Sculpting Device & Serum Kit User Manual

Myoglow MyoGlow

MyoGlow LED Sculpting Device & Serum Kit User Manual

Brand: Myoglow

Model: MyoGlow

INTRODUCTION

The MyoGlow LED Sculpting Device, paired with Niacinamide and NoTox Serums, offers an at-home solution for skin care. This manual provides instructions for safe and effective use of your device and accompanying serums.

SAFETY INFORMATION

WARNING:


- Do not use on throat or esophagus.
- For home use only.

PRODUCT OVERVIEW

The MyoGlow device is designed for use on the arm, neck, and face. It incorporates LED light therapy, gentle vibrations, and thermal warmth to support skin appearance.

Key Features:

- **LED Technology:** Three different LED light modes: Amber, Blue, and Red.
- **Gentle Vibrations:** The device vibrates only when pressed against the skin.
- **Convenient & Ergonomic:** Designed to fit comfortably in hand and is compact for portability.
- **Lymphatic Drainage Massager Tool:** Functions to help reduce facial puffiness.

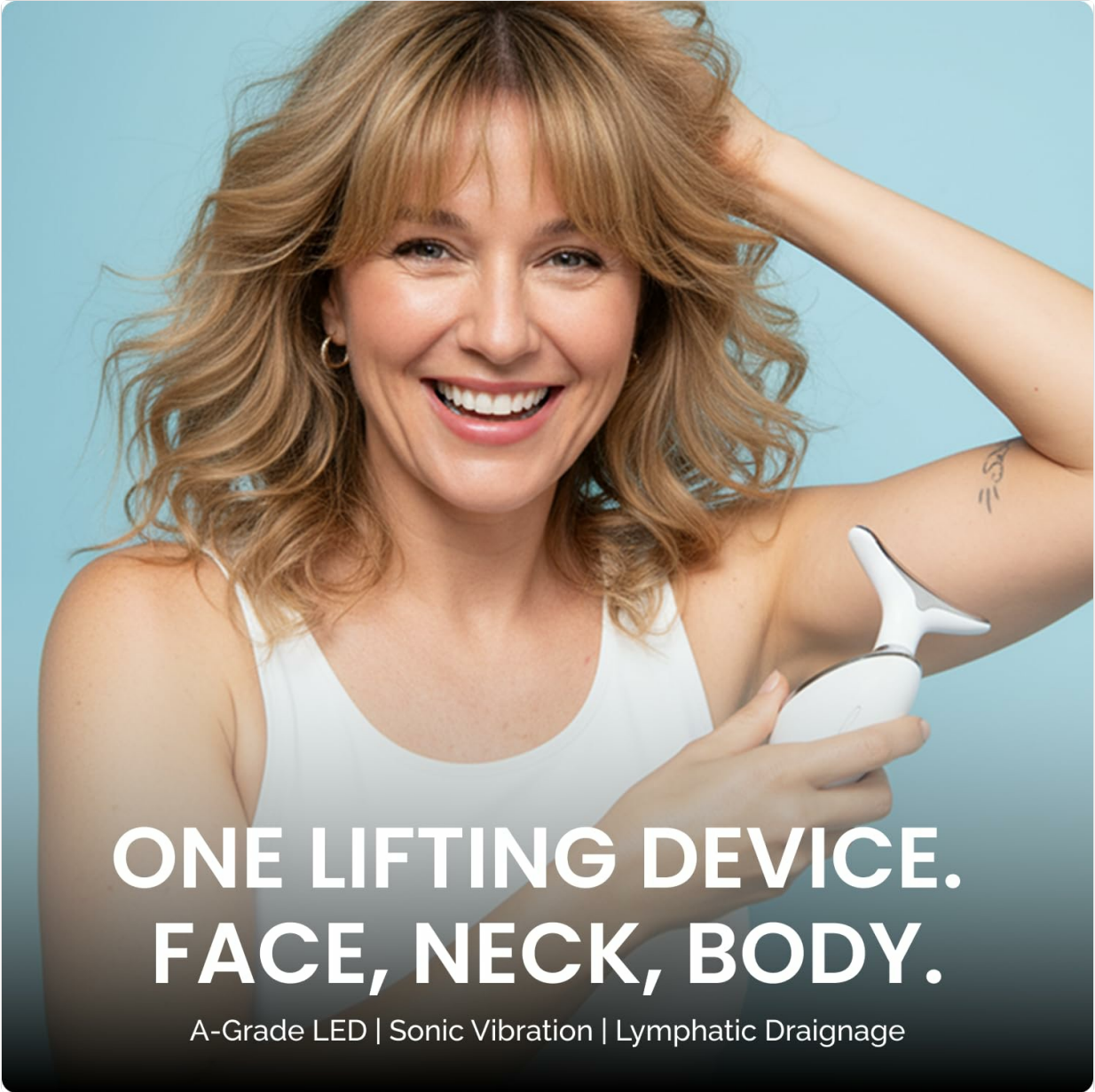

ONE LIFTING DEVICE. FACE, NECK, BODY.

A-Grade LED | Sonic Vibration | Lymphatic Drainage

The MyoGlow device features a sleek, ergonomic design for comfortable handling during use.



FACE, NECK AND ARMS

- ✓ LED Light Therapy
- ✓ Sonic Vibration
- ✓ Thermal Therapy
- ✓ Myofascial Release

The MyoGlow device is versatile and can be applied to various areas including arms, face, legs, and neck for comprehensive skin care.

INCLUDED SERUMS

The MyoGlow device is designed to be paired with MyDermaDream serums for optimal results.

Niacinamide Serum:

- Moisturizes, brightens, and revitalizes skin.
- Contains 5% Niacinamide.



3

Turn on MyoGlow

And choose treatment mode of choice

[◀Click Through to Explore ▶](#)

The Niacinamide Serum helps to moisturize, brighten, and revitalize the skin.

NoTox Serum:

- Formulated with Matrixyl peptides, Hyaluronic Acid, and Vitamin C.
- Aims to provide a firming effect.



2

Apply Activator Serum

We recommend No-Tox Serum

[◀Click Through to Explore ▶](#)

The NoTox Serum is formulated with key ingredients to support skin firmness.

SETUP

Charging the Device:

The MyoGlow device requires 2 Lithium Ion batteries (included). Ensure the device is fully charged before first use. Use the provided charging cable.

Preparation:

1. Cleanse your skin thoroughly.
2. Apply a small amount of your chosen MyDermaDream serum (Niacinamide or NoTox) to the desired treatment area, allowing it to absorb.



1

Take a Shower

To cleanse your skin properly

[◀ Click Through to Explore ▶](#)

Apply serum evenly to clean skin before using the MyoGlow device.

OPERATING INSTRUCTIONS

Follow these steps for effective use of your MyoGlow device:

1. **Power On:** Press and hold the power button on the device to turn it on.
2. **Select Mode:** Briefly press the power button to cycle through the three LED light modes: Red, Blue, and Amber.
3. **Apply to Skin:** Gently press the device against your skin. The device will begin to vibrate when it makes contact.
4. **Massage:** Use the device in a gliding, uplifting motion across the desired areas (face, neck, arms, legs).
5. **Power Off:** Press and hold the power button to turn off the device after use.

LED Light Modes:

- **Red Light Therapy:** Improves fine lines and wrinkles, promotes collagen production.
- **Blue Light Therapy:** Kills acne-causing bacteria, soothes inflammation, improves tone and texture.
- **Amber Light Therapy:** Promotes lymphatic drainage, reduces redness, restores radiance.

Experience Real Results

99% Saw improvement in their flabby arms and legs

95% Saw improvement in their skin texture and tone

88% Received compliments on their youthful skin

Visible Results Reported from 230 Sample-Size test*



The MyoGlow device offers three distinct LED light modes, each with specific benefits for skin health.

Unveil a healthy, beautifully balanced complexion.

Effective and gentle skincare for all ages and skin types



The red light mode is used for improving skin texture and promoting collagen.

SIMPLE 4-STEP PROCESS



① Take a Shower

To cleanse your skin properly



② Apply Activator Serum

We recommend No-Tox Serum



③ Turn on MyoGlow

And choose treatment mode of choice



④ Apply on desired area

Use in a gliding, uplifting motion

The blue light mode targets acne-causing bacteria and helps soothe inflammation.

This video demonstrates the MyoGlow device and its use with the ultimate skincare bundle.

This video shows how to achieve a spa-like skincare treatment at home using the MyoGlow device.

This video highlights the benefits of the MyoGlow skincare bundle for achieving optimal skin health.

This video provides a review of the MyoGlow ultimate skincare bundle, including the device and serum.

This video demonstrates a routine for achieving glowing skin at home using the MyoGlow device.

This video offers a relaxing demonstration of the MyoGlow device in use.

This video showcases the simple, convenient, and effective use of the MyoGlow device.

This video presents the MyoGlow device as an ultimate skincare tool.

This video demonstrates a spa-like skincare treatment at home using the MyoGlow device.

This video illustrates the multi-purpose use of the MyoGlow device and its visible effects.

MAINTENANCE

- After each use, gently wipe the device head with a soft, damp cloth.
- Do not submerge the device in water.
- Store the device in a cool, dry place away from direct sunlight.

TROUBLESHOOTING

Device Not Turning On:

- Ensure the device is fully charged.
- Press and hold the power button for a few seconds to activate.

No Vibration:

- The device is designed to vibrate only when in contact with skin. Ensure firm contact.

Reduced Performance:

- Ensure the device is clean and free from serum residue.
- Recharge the device if battery life is low.

SPECIFICATIONS

- **Package Dimensions:** 6.85 x 5.43 x 4.06 inches; 1.12 Pounds
- **Batteries:** 2 Lithium Ion batteries required (included)
- **UPC:** 850063274654
- **Manufacturer:** My DermaDream
- **ASIN:** B0FRB151L3

WARRANTY AND SUPPORT

- **Customer Support:** Responsive customer support is available through Amazon 7-days a week.
- **Return Policy:** 30-day refund/replacement policy.