

[Manuals.plus](#) /

› [Waterpulse](#) /

› Waterpulse WP62 Portable Oral Irrigator Instruction Manual

Waterpulse WP62

Waterpulse WP62 Portable Oral Irrigator Instruction Manual

Model: WP62

IMPORTANT SAFETY INFORMATION

ATTENTION: Read this entire manual before using the device.

- This appliance is designed for personal oral hygiene. Do not use it for any other purpose.
- Do not immerse the main unit in water. The unit has an IPX4 rating, meaning it is resistant to splashing water, but not full immersion.
- Ensure the power cord is not damaged before use. If damaged, contact customer support.
- Keep out of reach of children. Children should be supervised to ensure they do not play with the appliance.
- Use only with water or dentist-recommended mouthwash.
- Do not direct water under the tongue, into the nose, ear, or other sensitive areas, as this may cause harm.
- Unplug the device when not in use and before cleaning.

PRODUCT OVERVIEW

The Waterpulse WP62 Portable Oral Irrigator is designed for deep, safe, and effortless oral cleaning. It is ideal for daily oral health, especially for users with orthodontic braces, implants, or crowns. Featuring 10 adjustable pressure levels and a precise, powerful water jet, it effectively removes food debris and bacterial plaque, reaching areas that traditional brushing or flossing cannot. Its 600ml reservoir ensures complete hygiene without interruptions, and the 1.2-meter hose provides freedom of movement during use. With a modern and durable design, this irrigator is an essential tool for daily oral hygiene.



The main product image of the Waterpulse WP62 Oral Irrigator, showcasing its compact design and highlighting features such as 6 exclusive nozzle types, 600ml water capacity, 3-minute smart timer, adjustable pressure, 1-year warranty, 360-degree rotating nozzle, and flexible hose.

Key Features:

- Deep Oral Cleaning: Removes bacteria from hard-to-reach areas.
- 10 Adjustable Pressure Levels (21-145 PSI): Choose the ideal intensity for your gums.
- 8 Interchangeable Nozzles: Includes orthodontic, periodontal, tongue cleaner, and standard tips for family use.
- 360° Rotating Nozzle System: Reaches any area of the mouth with greater efficiency.
- Ideal for Braces and Implants: Removes debris around brackets, crowns, and prostheses.
- 600ml Reservoir: Provides longer autonomy for extended cleaning sessions.
- Automatic Bivolt Power Supply (110V/220V): Compatible with any outlet.
- Quiet Operation: Powerful cleaning without disturbance.
- Easy to Refill and Clean: Removable lid and moisture-resistant finish.

Package Contents:

- 01 Main Unit
- 01 Instruction Manual

- 03 Standard Nozzles
- 01 Periodontal Nozzle
- 01 Orthodontic Nozzle
- 01 Dental Plaque Nozzle
- 01 Tongue Cleaner
- 01 Rotating Brush Nozzle

PRODUCT COMPONENTS

Familiarize yourself with the parts of your Waterpulse WP62 Oral Irrigator:

- **Main Unit:** Houses the motor, pump, and control panel.
- **Water Reservoir:** 600ml capacity for water or mouthwash.
- **Handle with Nozzle:** Ergonomic handle with ON/OFF switch and pressure control.
- **Nozzle Storage Compartment:** Conveniently located on top of the unit.
- **Power Cord:** 1.2-meter cord for flexible placement.
- **Pressure Adjustment Dial:** Located on the main unit, offering 10 intensity levels.

Interchangeable Nozzles:

The WP62 comes with 8 specialized nozzles for comprehensive oral care:

- **Standard Nozzles (x3):** For general daily cleaning.
- **Periodontal Nozzle (x1):** Designed for deep cleaning of periodontal pockets.
- **Orthodontic Nozzle (x1):** Specifically for cleaning around braces and orthodontic appliances.
- **Dental Plaque Nozzle (x1):** For effective removal of plaque from tooth surfaces.
- **Tongue Cleaner (x1):** Helps remove bacteria and freshen breath.
- **Rotating Brush Nozzle (x1):** Provides additional cleaning action.

waterpulse


IRRIGADOR ORAL



*6 TIPOS DE BICOS
EXCLUSIVOS*



DESIGN COMPACTO

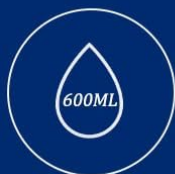

*MAIOR CAPACIDADE DE
ÁGUA, COM MENOS
NECESSIDADE DE
RECARGAS*



*TEMPO INTELIGENTE
DE 3 MINUTOS*



*AJUSTE DE PRESSÃO
COM MÚLTIPLAS
CONFIGURAÇÕES*



GARANTIA DE UM ANO



*BICO GIRATÓRIO
DE 360 GRAUS*



AJUSTE A PRESSÃO



*MANGUEIRA SANFONADA
DOBRÁVEL*

This image displays the complete set of 8 specialized nozzles included with the Waterpulse WP62 Oral Irrigator, designed for various oral hygiene needs, including standard, orthodontic, periodontal, and tongue cleaning tips. It also shows the nozzle release button and power switch on the handle.

SETUP

1. **Unpack and Place:** Carefully remove all components from the packaging. Place the main unit on a flat, stable surface near an electrical outlet and a sink.
2. **Fill the Reservoir:** Lift the reservoir lid and fill the 600ml reservoir with lukewarm water. You may also add a small amount of mouthwash as recommended by your dentist. Close the lid securely.
3. **Attach a Nozzle:** Select the desired nozzle. Insert the nozzle into the hole at the top of the handle until it clicks into place. To remove, press the release button on the handle and pull the nozzle out.
4. **Connect Power:** Plug the power cord into a suitable electrical outlet. The device is Bivolt (110V/220V) and will automatically adjust.

**ESCOLHA SEU
NÍVEL DE
PRESSÃO**



10 REGULAÇÃO DE
PRESSÃO

MENOR



MAIOR



The image displays the 600ml water reservoir of the Waterpulse WP62, designed for extended cleaning sessions. It also points out the 1-meter flexible hose, allowing for comfortable movement during use.

OPERATION

1. **Adjust Pressure:** Rotate the pressure adjustment dial on the main unit to your desired intensity level (1-10). Start with a lower setting and gradually increase as you get accustomed to the water jet.
2. **Position Nozzle:** Lean over the sink and place the nozzle tip in your mouth, aiming it at your gumline at a 90-degree angle. Close your lips slightly to prevent splashing, but allow water to flow into the sink.
3. **Start Irrigation:** Slide the ON/OFF switch on the handle to the "ON" position to start the water flow.
4. **Guide Water Jet:** Guide the water jet along your gumline, pausing briefly between teeth. Clean both the inside and outside surfaces of your teeth, as well as the chewing surfaces. For braces, implants, or crowns, carefully direct the jet around these areas. The 360° rotating nozzle allows for easy access to all areas.
5. **Finish:** When finished, slide the ON/OFF switch on the handle to the "OFF" position.
6. **Empty Reservoir:** It is recommended to empty any remaining water from the reservoir after each use to prevent bacterial growth.

VEM COM 8 BICOS

- 03 BICO PADRÃO
- 01 BICO PERIODONTAL
- 01 BICO ORTODÔNTICO
- 01 BICO DE PLACVAS DENTÁRIAS
- 01 LIMPADOR DE LÍNGUA
- 01 ESCOVA GIRATÓRIA



This view highlights the anti-slip rubber bases on the bottom of the irrigator, ensuring stability on various surfaces. It also shows the pressure adjustment dial with 10 levels, from gentle to intense.

MAINTENANCE AND CLEANING

Regular cleaning ensures optimal performance and hygiene of your Waterpulse WP62 Oral Irrigator.

1. After Each Use:

- Turn off the device and unplug it from the power outlet.
- Empty any remaining water from the reservoir.
- Remove the reservoir and rinse it thoroughly with warm water. Allow it to air dry.
- Remove the nozzle from the handle and rinse it under running water.

2. Cleaning the Main Unit:

- Wipe the exterior of the main unit with a soft, damp cloth. Do not use abrasive cleaners or immerse the unit in water.

3. Deep Cleaning (Monthly):

- To prevent mineral deposits, fill the reservoir with 2 parts water and 1 part white vinegar.
- Run the solution through the unit until the reservoir is empty.

- Rinse by running a full reservoir of clean water through the unit.

TROUBLESHOOTING

Problem	Possible Cause	Solution
Device does not turn on.	Not plugged in; power outlet issue; damaged cord.	Ensure the power cord is securely plugged into a working outlet. Check the outlet with another appliance. If the cord is damaged, contact customer support.
Weak or no water pressure.	Low water level in reservoir; nozzle clogged; pressure setting too low.	Refill the reservoir. Remove and clean the nozzle. Increase the pressure setting using the dial. Perform a deep cleaning if mineral buildup is suspected.
Water leakage.	Reservoir not properly seated; nozzle not securely attached.	Ensure the water reservoir is firmly placed on the main unit. Make sure the nozzle is clicked securely into the handle.
Excessive noise during operation.	Air in the pump; unit not on a stable surface.	Ensure the reservoir is full. Place the unit on a flat, stable surface. If noise persists, contact customer support.

SPECIFICATIONS

Feature	Detail
Brand	Waterpulse
Model	WP62
Input Voltage	AC100 – 240V (Bivolt)
Frequency	50/60Hz
Power	18W
Water Pressure	21 to 145 PSI
Pressure Levels	10 adjustable levels
RPM	1450±20%
Noise Level	Less than 75 Db
Startup Time	± 15 Seconds
Reservoir Capacity	600ml
Water Flow Rate	> 290ml/min
Number of Nozzles	8 specialized nozzles
Hose Length	1.2m

Feature	Detail
Product Dimensions (LxWxH)	20cm x 10cm x 9cm
Product Weight	765g
Material	Plastic
Special Feature	Lightweight
IP Rating	IPX4 (Splash-proof)

IRRIGADOR ORAL waterpulse





waterpulse
IRRIGADOR ORAL

Para uma vida saudável, dentes sempre limpos e sorriso mais bonito.

- LIMPEZA PROFUNDA Remove resíduos e placa bacteriana
- JATOS DE ÁGUA PULSANTES
- 10 NÍVEIS DE PRESSÃO
- BICO GIRATÓRIO 360°

Capacidade do tanque: 600ml / 19oz

IPX4
À PROVA D'ÁGUA

BENEFÍCIOS DO USO

- Ajuda a prevenir cáries, gengivite e mau hálito
- Melhora a saúde das gengivas
- Remove resíduos em áreas de difícil acesso
- Ideal para uso com aparelhos ortodonticos e implantes

ESPECIFICAÇÕES

- Tensão: 110V-240V - 50/60Hz
- Potência: 18W
- Pressão: 20-120 PSI
- Frequência: 1400 - 16000 puls/min

RoHS CE

PRESSIONE O BOTÃO PARA RETIRAR O BICO

Permite remover o bico com facilidade para troca ou limpeza.

BOTÃO DE LIGAR/DESLIGAR

Liga e desliga o aparelho e permite ajustar o nível de pressão.

- LIMPEZA PROFUNDA**
Remove resíduos e placa bacteriana
- JATOS DE ÁGUA PULSANTES**
Mais eficácia na limpeza bucal
- 10 NÍVEIS DE PRESSÃO**
Ajuste ideal para cada necessidade
- BICO GIRATÓRIO 360°**
Mais conforto e alcance total



This image illustrates the compact dimensions of the Waterpulse WP62 Oral Irrigator: 20cm height, 10cm width, and 9cm depth. It also highlights the 1.2-meter power cable, providing ample reach for convenient use.

WARRANTY INFORMATION

The Waterpulse WP62 Oral Irrigator comes with a **90-day warranty** from the date of purchase. This warranty covers manufacturing defects under normal use. Please retain your proof of purchase for warranty claims. The warranty does not cover damage caused by misuse, accidents, unauthorized repairs, or normal wear and tear.

CUSTOMER SUPPORT

If you have any questions, concerns, or require assistance with your Waterpulse WP62 Oral Irrigator, please contact Waterpulse customer support. Refer to the product packaging or the official Waterpulse website for the most current contact information.