

## FOLOKE B0FX3STD3X

# FOLOKE LED Light Therapy Mask Skin Care Mask

MODEL: B0FX3STD3X

## 1. Introduction

---

The FOLOKE LED Light Therapy Mask is designed to enhance your daily skincare and self-care routine. This cordless device utilizes red and near-infrared (NIR) light wavelengths to improve skin appearance, reduce wrinkles, and promote a healthier complexion. Its ergonomic design, flexible silicone construction, and adjustable straps ensure maximum comfort and portability, allowing for convenient use at home, during travel, or while engaging in various activities.

## 2. Product Features

---

- **Red and Near-Infrared (NIR) Light:** The device produces red and near-infrared light in the visible spectrum, specifically intended to improve the appearance of wrinkles.
- **Flexible Silicone Design:** Constructed from high-quality, flexible silicone, the mask contains LEDs and a controller, conforming perfectly to your face and neck for maximum comfort.
- **Integrated Neck Attachment:** The system includes a neck attachment for comprehensive treatment, extending benefits to the neck and shoulders.
- **Cordless and Portable:** Enjoy hands-free skincare anytime, anywhere, whether relaxing, working, traveling, or practicing yoga.
- **USB-C Rechargeable:** Equipped with a modern Type-C port, a full charge takes only 3-4 hours, providing extended use.

## 3. Safety Information

---

This device is intended for cosmetic use only. Consult a healthcare professional before use if you have any pre-existing skin conditions, are pregnant, or have any medical concerns. Avoid direct eye contact with the LED lights. Do not use the device if it is damaged. Keep out of reach of children.

## 4. Setup Guide

1. **Charge the Device:** Before first use, fully charge the mask using the provided USB-C cable. A full charge typically takes 3-4 hours.
2. **Attach Neck Piece:** If desired, connect the neck attachment to the main face mask. Ensure it is securely fastened.
3. **Adjust Straps:** Place the mask on your face and secure it with the adjustable Velcro straps for a comfortable and snug fit.



Figure 1: Components of the FOLOKE LED Light Therapy Mask, including the main mask, detachable neck piece, adjustable straps, and USB-C charging cable.

## 5. Operating Instructions

---

For best results, use the mask 3-4 times a week for 15-20 minutes per session.

1. **Cleanse Skin:** Ensure your face and neck are clean and dry before applying the mask.
2. **Wear the Mask:** Place the mask on your face and secure the straps. Ensure the neck piece is also comfortably positioned if using.
3. **Power On:** Press and hold the power button on the mask's controller to turn it on.
4. **Select Light Mode:** Choose your desired light mode. The mask offers various light wavelengths for different skin needs.
5. **Relax:** Relax for 15-20 minutes while the mask is in operation.
6. **Power Off:** After your session, press and hold the power button to turn off the device.
7. **Apply Skincare:** Remove the mask and apply your preferred skincare products.



Figure 2: The cordless design allows for flexible use during daily activities such as working, reading, yoga, or simply relaxing.

Your browser does not support the video tag.

Video 1: A seller-provided video demonstrating the use of the FOLOKE LED Light Therapy Mask, highlighting its features and ease of integration into a skincare routine.

Your browser does not support the video tag.

Video 2: A seller-provided video showcasing the FOLOKE LED Light Therapy Mask in action, emphasizing its benefits and user experience.

Your browser does not support the video tag.

Video 3: A seller-provided video introducing the FOLOKE LED Light Therapy Mask as a skincare partner, highlighting its ease of use and effectiveness.

## 6. Light Modes and Benefits

The mask features various light modes, each targeting specific skin concerns:

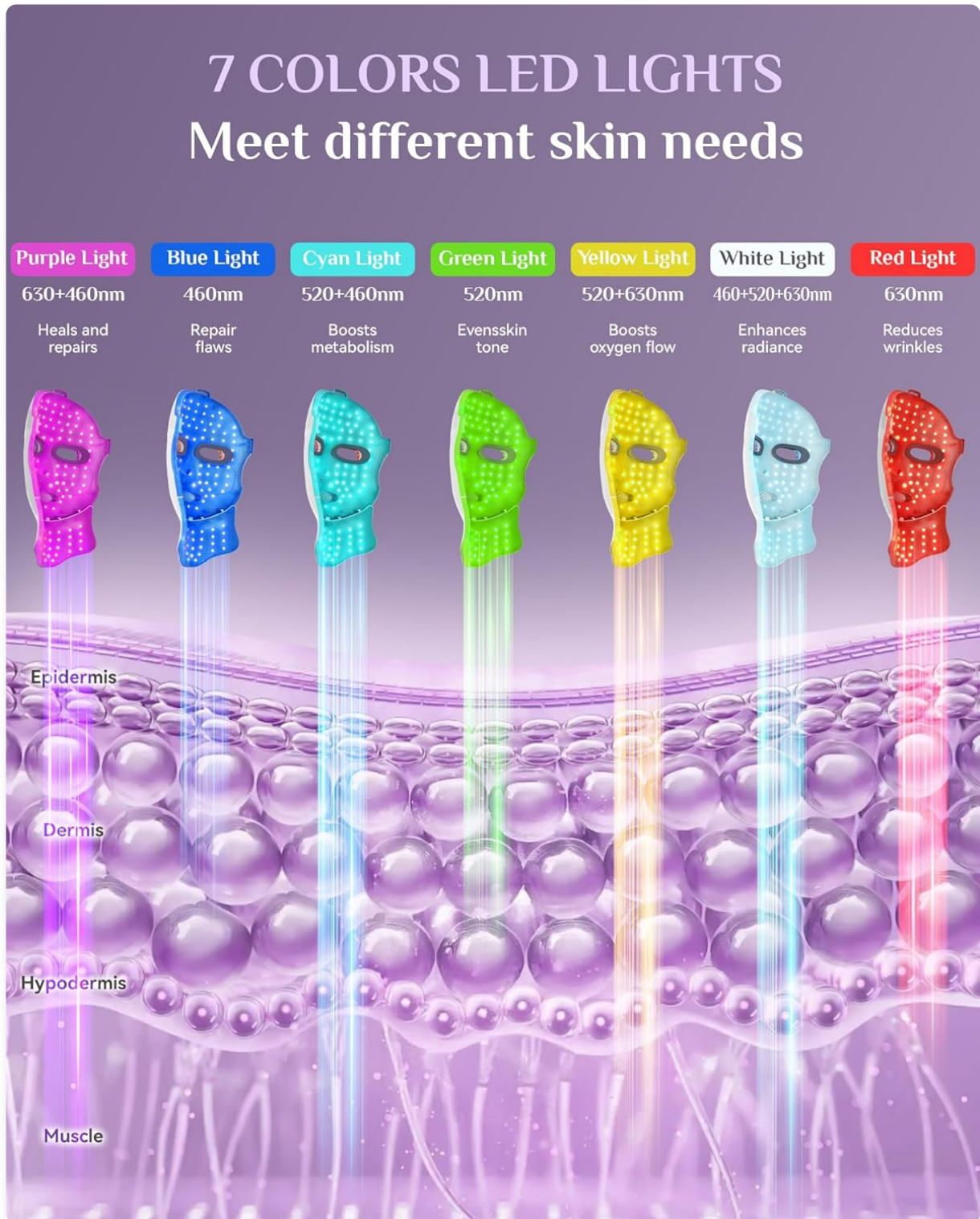


Figure 3: Detailed breakdown of the 7 LED light colors and their specific benefits for different skin layers.

- **Red Light (630nm):** Reduces wrinkles and stimulates collagen.
- **Blue Light (460nm):** Repairs flaws and helps with acne.
- **Green Light (520nm):** Evens skin tone.
- **Yellow Light (520+630nm):** Boosts oxygen flow.
- **Cyan Light (520+460nm):** Boosts metabolism.

- **Purple Light (630+460nm):** Heals and repairs skin.
- **White Light (460+520+630nm):** Enhances radiance.

## 7. Maintenance

---

- **Cleaning:** Wipe the mask and neck attachment with a soft, damp cloth after each use. Do not use abrasive cleaners or immerse the device in water.
- **Storage:** Store the device in a cool, dry place away from direct sunlight.
- **Charging:** Ensure the device is fully charged before storing for extended periods to maintain battery health.

## 8. Troubleshooting

---

- **Device not turning on:** Ensure the device is fully charged. Check the USB-C cable and adapter for damage.
- **Lights not working:** Verify that the mask is properly connected to the controller and fully charged.
- **Discomfort during use:** Adjust the straps for a more comfortable fit. If discomfort persists, discontinue use and consult a professional.

## 9. Specifications

---

Feature	Detail
Model	B0FX3STD3X
Light Wavelengths	Red (630nm), Blue (460nm), Green (520nm), Yellow (520+630nm), Cyan (520+460nm), Purple (630+460nm), White (460+520+630nm)
Material	Flexible Silicone
Power Source	USB-C Rechargeable (1 Lithium Polymer battery included)
Charging Time	3-4 hours
Usage Frequency	3-4 times a week
Session Duration	15-20 minutes
Package Dimensions	9.49 x 8.39 x 4.61 inches
Weight	1.7 Pounds

## 10. Warranty and Support

---

FOLOKE is committed to providing high-quality products and customer satisfaction. For any inquiries or support needs, please utilize our 24-hour responsive service. Contact information can typically be found on the product packaging or the official FOLOKE website.

