

## Manuals+

[Q & A](#) | [Deep Search](#) | [Upload](#)

manuals.plus /

› [MEGELIN](#) /

› [MEGELIN Duo-Lux LED Light Therapy Face and Neck Mask User Manual](#)

## MEGELIN Duo-Lux

# MEGELIN Duo-Lux LED Light Therapy Face and Neck Mask User Manual

Model: Duo-Lux

## 1. INTRODUCTION

---

The MEGELIN Duo-Lux LED Light Therapy Face and Neck Mask is designed for advanced skin care. This device combines 660nm (Red Light) and 1064nm (Near-Infrared Laser) wavelengths to support skin repair and promote a healthy complexion. It features a 3D contoured fit for comfort and optimized irradiance for effective treatment.



*Image 1: The Megelin Duo-Lux LED Light Therapy Face Mask and Neck Mask, showing both the face and neck components.*

## **2. IMPORTANT SAFETY INFORMATION**

Please read all instructions carefully before using the MEGELIN Duo-Lux LED Light Therapy Mask. Failure to follow these instructions may result in injury or damage to the device.

- Consult a healthcare professional before use if you have any pre-existing medical conditions, are pregnant, or are taking photosensitizing medications.
- Do not look directly into the LED or laser lights during operation.
- Keep the device out of reach of children.
- Do not use the device if it is damaged or appears to be malfunctioning.
- Use only the provided power adapter and accessories.
- Clean the device according to the maintenance instructions to prevent skin irritation or infection.

### 3. PACKAGE CONTENTS

---

Verify that all items are present in the package:

- MEGELIN Duo-Lux LED Light Therapy Face Mask
- MEGELIN Duo-Lux LED Light Therapy Neck Mask
- Controller Unit
- Charging Cable
- User Manual (this document)



*Image 2: The Megelin Duo-Lux LED Light Therapy Face Mask and Neck Mask, as packaged, showing the main components.*

### 4. SETUP AND CHARGING

---

1. **Charge the Controller:** Connect the charging cable to the controller unit and a suitable power source. The

controller typically takes approximately three hours for a full charge.

2. **Connect the Mask:** Ensure the controller is fully charged. Connect the mask's cable to the designated port on the controller unit.
3. **Prepare for Use:** Ensure your skin is clean and dry before placing the mask on your face and neck.

Your browser does not support the video tag.

*Video 1: A user unboxing and preparing the Megelin Duo-Lux LED Light Therapy Mask for use, demonstrating how to connect the mask to the controller.*

## 5. OPERATING INSTRUCTIONS

---

The MEGELIN Duo-Lux mask offers four specialized modes for various skin concerns. Follow these steps for operation:

1. **Wear the Mask:** Carefully place the face mask and neck mask onto your face and neck, ensuring a comfortable and secure fit. The 3D contoured design is intended to adapt perfectly to your facial structure.
2. **Power On:** Press and hold the power button on the controller unit to turn on the device.
3. **Select Mode:** Use the mode selection buttons on the controller to choose your desired treatment mode. The available modes are:
  - **Wrinkles:** Targets fine lines and wrinkles.
  - **Inflammation:** Addresses skin inflammation and redness.
  - **Skin Toning:** Improves overall skin tone and firmness.
  - **Daily Repair:** For general skin health and maintenance.
4. **Set Duration:** Adjust the treatment duration using the timer buttons. Recommended session times may vary by mode; refer to the specific mode guidelines below.
5. **Begin Treatment:** Once the mode and duration are set, the light therapy will begin. Relax during your session.
6. **Power Off:** The device will automatically turn off after the set duration. You can also manually power off by pressing and holding the power button.

# A New Era of Radiance

Engineered for professional facial anti-aging

Reduces wrinkles

---

Fades pigmentation

---

Increases skin tone

---

Calms redness

---

Relieves inflammation

---



*Image 3: A user wearing the Megelin Duo-Lux LED Light Therapy Face Mask, demonstrating its comfortable fit and ease of use.*

Your browser does not support the video tag.

*Video 2: This video illustrates the application of the Megelin Duo-Lux mask and highlights its advanced Laser & LED Fusion technology for targeted skin treatment.*

## 6. MODES AND TARGETED BENEFITS


The Duo-Lux mask utilizes a combination of 660nm (Red Light) and 1064nm (Near-Infrared Laser) wavelengths to address various skin concerns effectively.

### 6.1. Laser-Powered Facial Zone Mapping

The device employs laser technology for precise targeting of specific areas prone to wrinkles, such as forehead lines, crow's feet, and nasolabial folds, while LED light provides broader coverage for overall skin care.

# Laser-powered Facial Zone Mapping

Laser for wrinkles, LED for overall care.

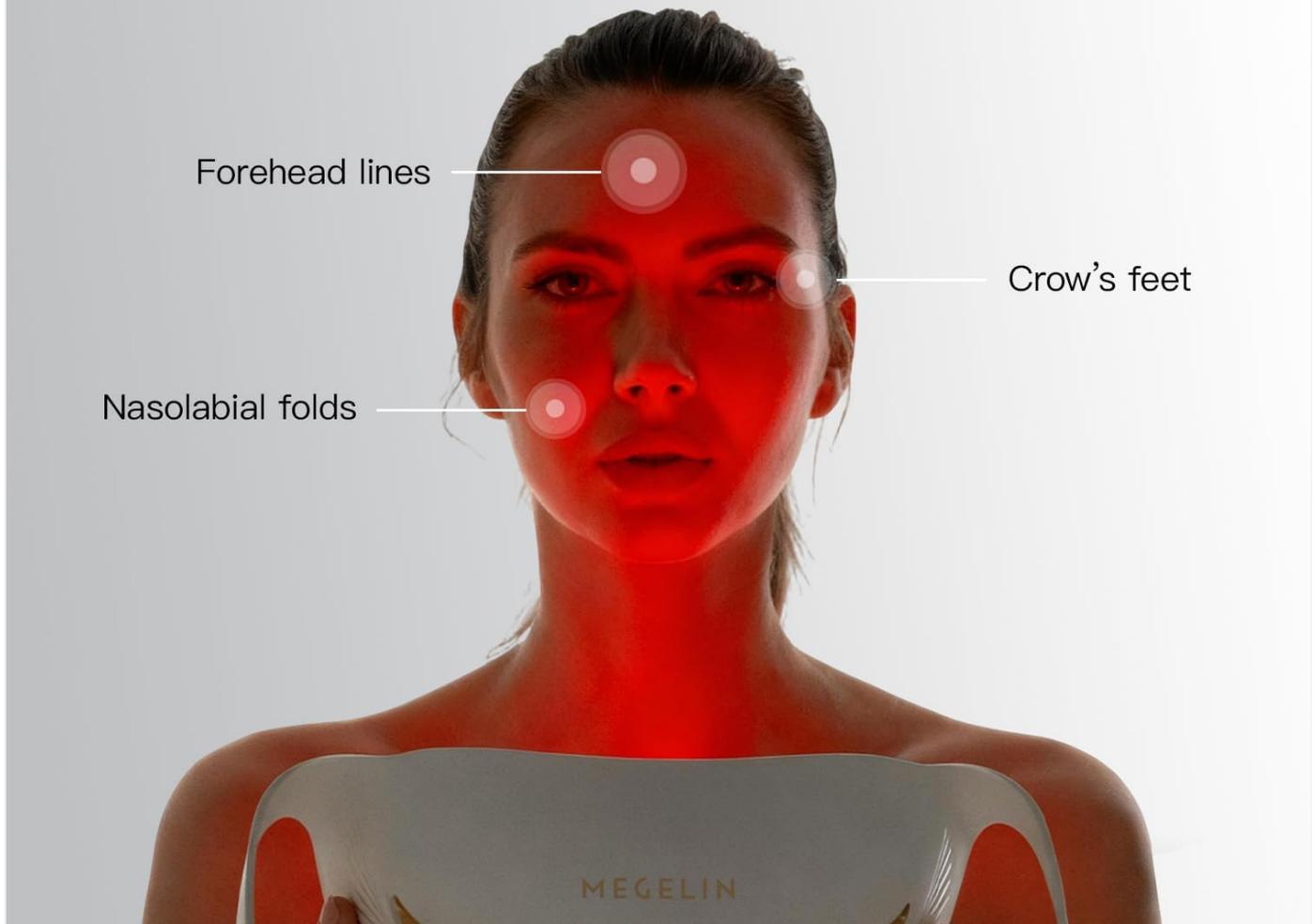


*Image 4: Diagram illustrating the laser-powered facial zone mapping feature, highlighting target areas for wrinkle treatment like forehead lines, crow's feet, and nasolabial folds.*

## 6.2. Duo-Lux™ Laser-LED Fusion

This technology integrates 660nm and 1064nm lasers with 660nm LEDs to deliver powerful and stable light energy at an optimized irradiance of 45mW/cm<sup>2</sup>. This fusion supports collagen production and cellular repair at both surface and deeper skin layers.

# Duo-Lux™ Laser&LED Fusion

Featuring 660+1064nm lasers  
paired with 660nm LED.

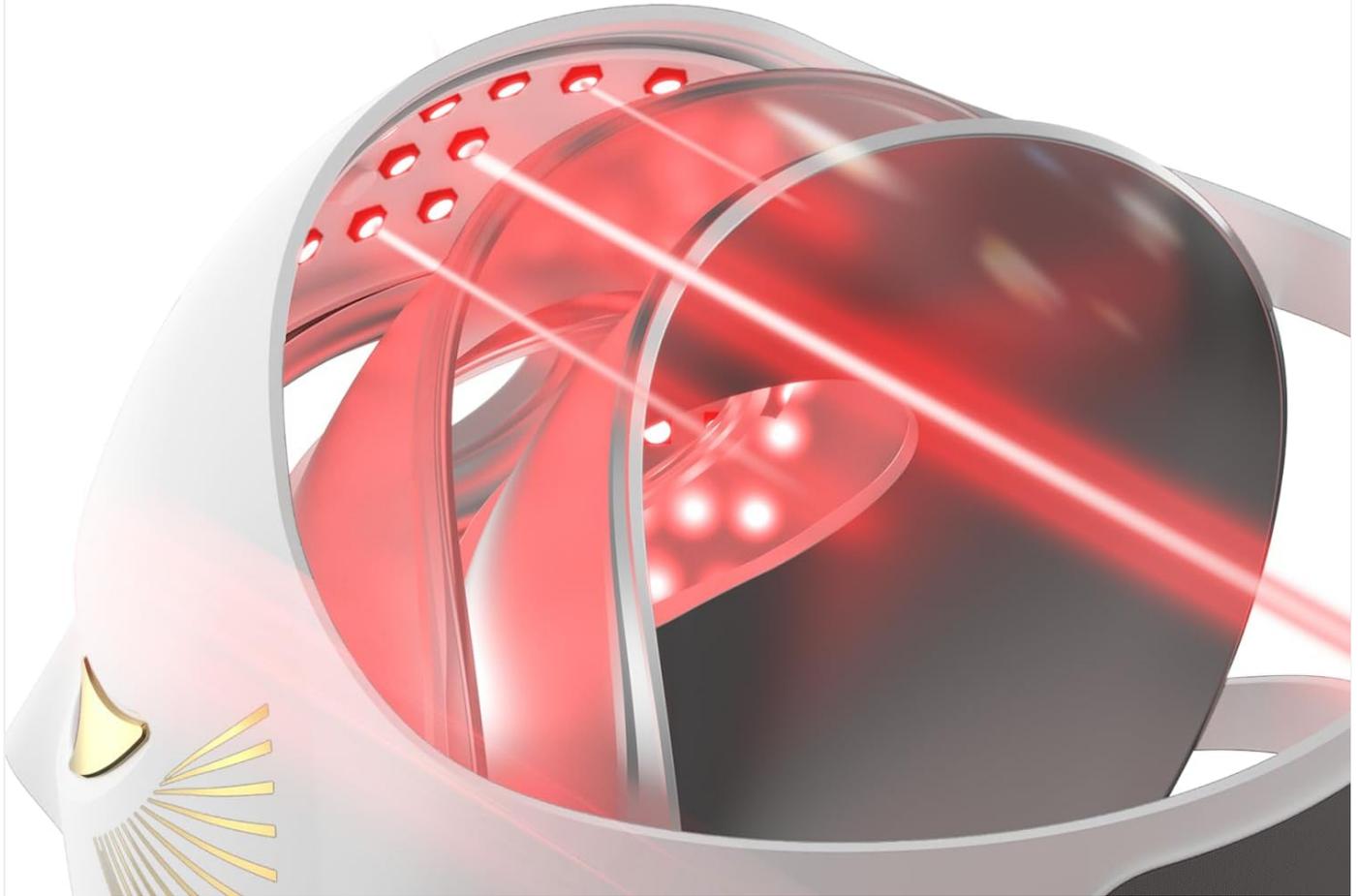

*Image 5: A detailed view of the Megelin Duo-Lux mask's internal structure, showcasing the integration of laser and LED components for comprehensive light therapy.*

Your browser does not support the video tag.

*Video 3: This video demonstrates the Megelin Duo-Lux LED Light Therapy Mask in use, highlighting its comfortable fit and the visible red light during operation.*

## 7. CARE AND MAINTENANCE

- **Cleaning:** After each use, gently wipe the interior and exterior surfaces of the mask with a soft, damp cloth. Do not use abrasive cleaners or immerse the device in water.
- **Storage:** Store the mask and controller in a cool, dry place, away from direct sunlight and extreme temperatures. Use the original packaging or a protective case to prevent damage.
- **Battery Care:** To prolong battery life, avoid fully discharging the controller frequently. Charge it regularly, even if not in daily use.

## 8. TROUBLESHOOTING

Problem	Possible Cause	Solution
Device does not turn on.	Low battery or power issue.	Ensure the controller is fully charged. Check the charging cable and adapter for damage.
Lights are not illuminating.	Incorrect mode selection or connection issue.	Verify the mask is securely connected to the controller. Re-select the desired treatment mode.
Discomfort during use.	Improper fit or sensitive skin.	Adjust the mask for a comfortable fit. If discomfort persists, discontinue use and consult a healthcare professional.

## 9. PRODUCT SPECIFICATIONS

---

- **Model:** Duo-Lux
- **Light Wavelengths:** 660nm (Red LED), 1064nm (Near-Infrared Laser)
- **Irradiance:** 45mW/cm<sup>2</sup> (optimized)
- **Fit:** 3D Contoured Fit Design
- **Power Source:** Rechargeable Controller Unit
- **Manufacturer:** Megelin
- **Package Dimensions:** 11.9 x 11.75 x 11.25 inches
- **Weight:** 6.37 Pounds

## 10. WARRANTY AND SUPPORT

---

For detailed warranty information, product support, or service inquiries, please refer to the official documentation included with your product or visit the MEGELIN official website. Contact information for customer service can typically be found in the product packaging or on the manufacturer's website.

