

Manuals+

[Q & A](#) | [Deep Search](#) | [Upload](#)

manuals.plus /

> [MYTREX](#) /

> [MYTREX Head Spa Deep X: Facial and Scalp Massager Instruction Manual](#)

MYTREX MTHSDX24B

MYTREX Head Spa Deep X: Facial and Scalp Massager Instruction Manual

Model: MTHSDX24B

1. INTRODUCTION

Thank you for choosing the MYTREX Head Spa Deep X. This device is designed to provide a comprehensive facial and scalp massage experience, incorporating EMS technology, red light therapy, and ion infusion for enhanced hair and skin care. Please read this manual thoroughly before use to ensure proper operation and maintenance.

2. SAFETY INFORMATION

- Do not use this device while sleeping.
- Avoid using the device in hot water baths. The device is IPX7 waterproof, suitable for shower use, but not for submersion in hot baths.
- Keep out of reach of children.
- Consult a physician before use if you have any pre-existing medical conditions, especially those related to skin, scalp, or circulatory issues.
- Discontinue use immediately if you experience any discomfort or adverse reactions.
- Use only the provided charging cable and adapter.
- Do not attempt to disassemble or repair the device yourself. Contact customer support for assistance.

3. PACKAGE CONTENTS

Verify that all items are present in your package:

- MYTREX Head Spa Deep X Device
- Charging Base
- Charging Cable
- Scalp Massage Heads (set)
- Facial Massage Head (optional, if included)

- User Manual



Image: MYTREX Head Spa Deep X device with its various accessories.

4. PRODUCT OVERVIEW

The MYTREX Head Spa Deep X features an ergonomic design for comfortable handling. It includes interchangeable massage heads and intuitive controls for various functions.



Image: Front view of the MYTREX Head Spa Deep X device.

Controls:

- **Power Button:** Press and hold to turn on/off.
- **Mode Button:** Press to cycle through massage speeds and EMS intensity levels.
- **LED Indicators:** Display current speed/intensity level and battery status.

Key Features:

- **EMS Technology:** Electrical Muscle Stimulation for scalp and facial muscles.
- **Red Light Therapy (630 nm):** Promotes scalp health and hair regeneration.
- **Ion Infusion & Cleansing:** Enhances absorption of serums and deep cleanses pores.
- **Interchangeable Massage Heads:** For scalp and facial applications.
- **Waterproof Design (IPX7):** Suitable for use in the shower.

Image: Illustration of the device's core technologies: EMS, Red Light, Ion Infusion, and Ion Cleansing.

5. SETUP

5.1 Initial Charging

Before first use, fully charge the device for approximately 3 hours. Place the device on its charging base, ensuring proper contact. The LED indicators will show charging progress and full charge status.

5.2 Attaching Massage Heads

The device comes with interchangeable massage heads. To attach or remove a head, align it with the corresponding slot on the device and gently push or pull until it clicks into place. Ensure the heads are securely attached before use.

Your browser does not support the video tag.

Video: Demonstration of turning on the device, cycling through modes, and changing massage heads. (Note: This video URL is a placeholder as no m3u8 URL was provided in the source data.)

6. OPERATING INSTRUCTIONS

6.1 Power On/Off

Press and hold the Power button for 2 seconds to turn the device on. The device will start in the default mode. To turn off, press and hold the Power button again for 2 seconds.

6.2 Selecting Mode and Speed

After powering on, press the Mode button to cycle through the three available speed/intensity levels. The LED indicators will illuminate to show the selected level. Choose the level that best suits your comfort and needs.

6.3 Application

Apply the device to the desired area (scalp or face). For scalp use, gently move the device across your scalp in small circular motions. For facial use, use the dedicated facial massage head and glide it over your skin. The device has a 15-minute safety timer and will automatically shut off after this period.

From Clogged to Clean with Ion Cleansing



Before: Buildup



After: Purified Pores

Image: Woman using the device for ion cleansing on her scalp, illustrating its effectiveness in purifying pores.

630nm Red Light for Radiant Skin



Image: Woman using the device on her face, demonstrating the red light therapy feature.

7. MODES AND FEATURES EXPLAINED

7.1 EMS Technology

The Electrical Muscle Stimulation (EMS) feature gently stimulates muscles, promoting blood circulation and enhancing the effectiveness of treatments. This can be beneficial for both scalp and facial areas.

7.2 Red Light Therapy (630 nm)

The integrated 630 nm red light therapy helps to revitalize the scalp, strengthen hair roots, and accelerate hair regeneration. For facial use, it supports skin health and collagen production.

7.3 Ion Infusion and Cleansing

This innovative technology uses ions to deeply cleanse pores by removing impurities and significantly enhances the absorption of hair and skin serums. Studies indicate up to a 206% increase in absorption for certain treatments like

Deeper Minoxidil Absorption, Stronger Results.



Image: Chart illustrating the enhanced absorption of minoxidil when used with the Deep X infusion feature.

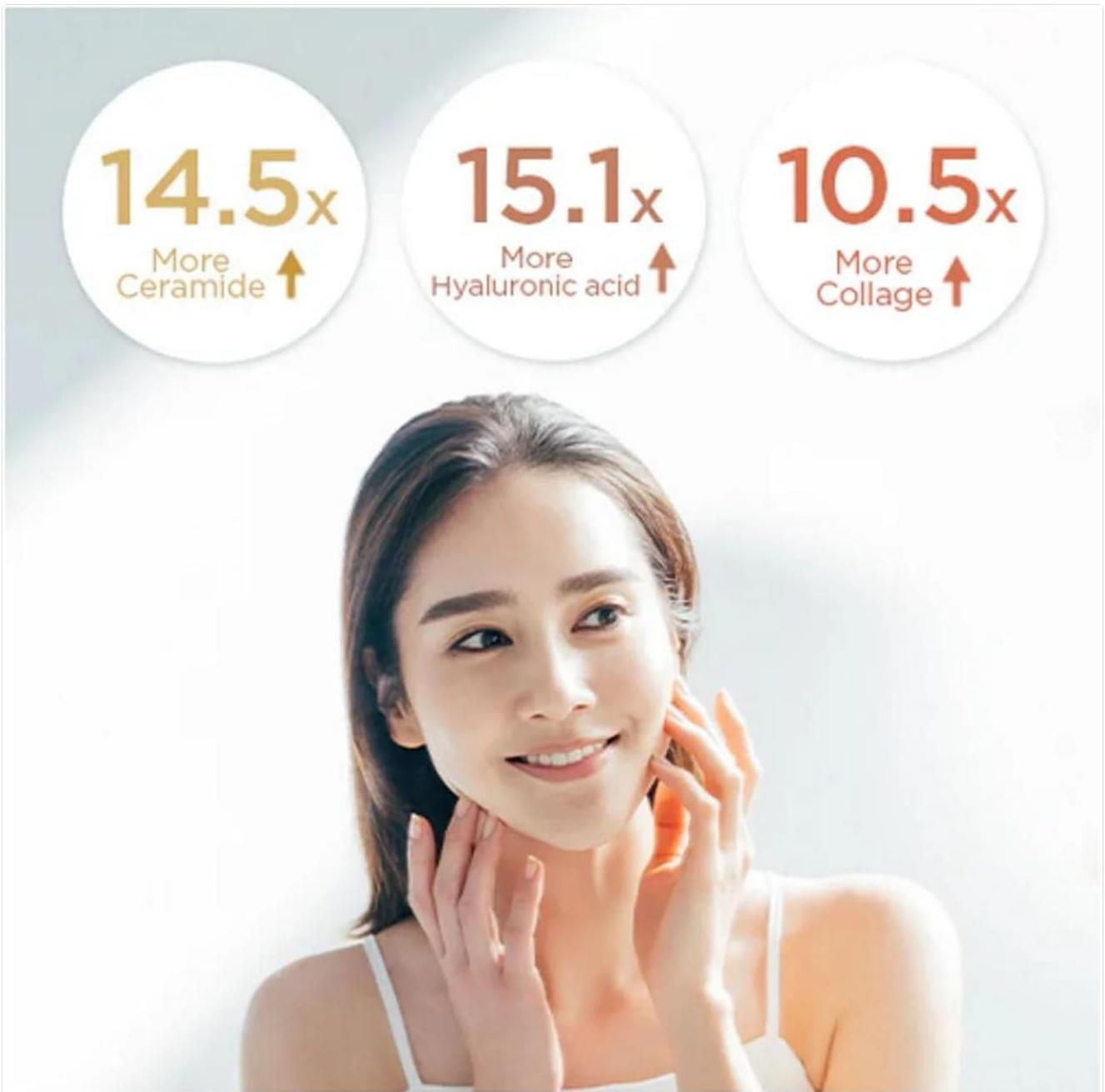


Image: Visual representation of increased absorption rates for key skincare ingredients like Ceramide, Hyaluronic acid, and Collagen.

8. MAINTENANCE AND CLEANING

Regular cleaning ensures optimal performance and hygiene of your device.

- After each use, especially if used with hair products or serums, remove the massage heads.
- Rinse the massage heads under running water. You may use a mild soap if necessary.
- Wipe the main body of the device with a soft, damp cloth. Do not submerge the main body in water, only the massage heads are designed for direct rinsing.
- Ensure all parts are completely dry before reattaching the heads or storing the device.
- Store the device in a cool, dry place away from direct sunlight.

9. TROUBLESHOOTING

Problem	Possible Cause	Solution
Device does not turn on.	Low battery or not charged.	Charge the device fully (approx. 3 hours).
Weak massage intensity.	Low battery or incorrect mode selected.	Charge the device. Press the Mode button to increase intensity.
Red light not visible.	Red light mode not activated or device malfunction.	Ensure the device is on and operating. If still not visible, contact support.
Massage heads detach easily.	Improper attachment.	Ensure massage heads are securely clicked into place.

If the problem persists after attempting these solutions, please contact MYTREX customer support.

10. SPECIFICATIONS

Feature	Detail
Model Number	MTHSDX24B
Dimensions	17.5 x 13.2 x 4 cm
Weight	313 g
Power Source	Rechargeable Battery (1 non-standard battery included)
Material	Polycarbonate
Color	Black
Water Resistance	IPX7 (suitable for shower use)
Red Light Wavelength	630 nm
Automatic Shut-off	15 minutes
Manufacturer	Shenzhen Chuangtong Electronic Instruments Co., Ltd.
Country of Origin	China

11. WARRANTY AND SUPPORT

The MYTREX Head Spa Deep X comes with a standard manufacturer's warranty. Please refer to the warranty card included in your package for specific terms and conditions. For technical support, troubleshooting assistance, or warranty claims, please contact MYTREX customer service through the official website or your retailer.

MYTREX Official Website: [Visit MYTREX Store on Amazon](#)