

[Manuals.plus](#) /

› [Philips Sonicare](#) /

› Philips Sonicare 7300 Series Electric Toothbrush HX7420/07 User Manual

Philips Sonicare HX7420/07

Philips Sonicare 7300 Series Electric Toothbrush User Manual

Model: HX7420/07

INTRODUCTION

This manual provides essential information for the safe and effective use of your Philips Sonicare 7300 Series Electric Toothbrush. Please read these instructions thoroughly before first use and retain them for future reference.

WHAT'S INCLUDED

Your Philips Sonicare 7300 Series Electric Toothbrush package includes the following components:

- 1x Philips Sonicare 7300 Series electric toothbrush handle
- 3x A3 Premium All-in-One brush heads
- 1x Charging travel case
- 1x Charging stand (power adapter not included)

What's in the box

3x A3 Premium All-in-One
brush head

7300 Series
Electric toothbrush

Charging travel case

USB Charging stand
Power adaptor not included



Image: The Philips Sonicare 7300 Series Electric Toothbrush, three A3 Premium All-in-One brush heads, a charging travel case, and a USB charging stand.

SETUP

Charging the Toothbrush

Before first use, fully charge your toothbrush. A full charge provides up to 21 days of regular brushing. Use the provided charging stand or the charging travel case.

1. Connect the charging stand to a USB power adapter (not included) and plug it into a power outlet.
2. Place the toothbrush handle onto the charging stand. The battery indicator lights will illuminate to show charging progress.
3. For travel, use the charging travel case by connecting it to its dedicated power adapter.

Attaching the Brush Head

Align the brush head with the metal shaft on the toothbrush handle and push it down firmly until it clicks into place.



A3 Premium All-in-One brush head

Delivers optimal cleaning

Image: A close-up view of the A3 Premium All-in-One brush head, designed for optimal cleaning.

OPERATION

Brushing Modes and Intensity Settings

Your toothbrush offers 4 brushing modes and 3 intensity levels for a personalized cleaning experience.

- **Clean Mode:** Standard mode for daily cleaning.
- **Sensitive Mode:** Gentle cleaning for sensitive teeth and gums.
- **Gum Health Mode:** Focuses on improving gum health.
- **White Mode:** Designed for surface stain removal.

To change modes or intensity, press the mode/intensity button on the handle. The indicator lights will show the selected setting.



Image: The toothbrush handle displaying the power button, mode/intensity selection buttons, and indicators for 4 modes and 3 intensity levels, along with the battery life indicator.

Next-Generation Sonicare Technology

The toothbrush utilizes advanced sonic technology, providing 62,000 bristle movements per minute. The motor automatically adjusts power to maintain consistent performance, even in hard-to-reach areas. Sonicare Fluid Action supports the bristles by driving fluid deep between teeth and along the gumline for effective cleaning.

Next-Generation Sonicare Technology

Gentle yet effective
clean even in hard-
to-reach areas

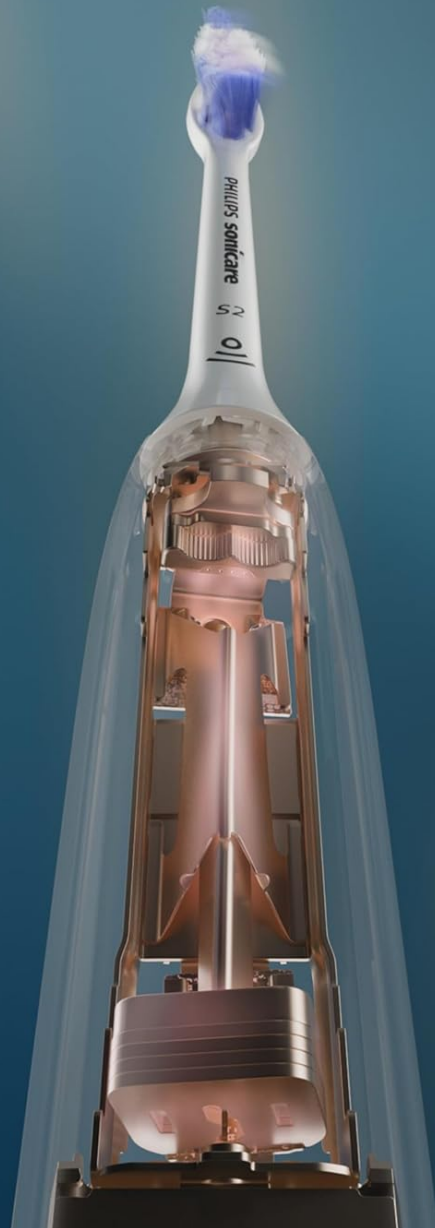


Image: An illustration showing the internal mechanism of the Philips Sonicare toothbrush, highlighting its next-generation sonic technology for gentle yet effective cleaning.

Sonicare Fluid Action

Supports bristles cleaning between teeth and along the gumline



Image: A close-up of toothbrush bristles in action, illustrating how Sonicare Fluid Action supports cleaning between teeth and along the gumline.

SmarTimer and BrushPacer

The built-in SmarTimer ensures you brush for the dentist-recommended 2 minutes. The BrushPacer prompts you to move to a new section of your mouth every 20 seconds, ensuring thorough cleaning of all areas.



SmarTimer

Complete your clean
after 2 minutes

BrushPacer

Prompts to clean
a new area every
20 seconds

Image: A person using the Philips Sonicare toothbrush, demonstrating the SmarTimer for a 2-minute clean and the BrushPacer for 20-second interval prompts.

Visual Pressure Sensor

To protect your gums, a visual pressure sensor at the base of the toothbrush handle will light up if you apply too much pressure. Reduce pressure when the light activates.

Visual pressure sensor

Helps to protect
your gums



Image: A hand holding the Philips Sonicare toothbrush, with the visual pressure sensor at the base glowing purple, indicating excessive brushing pressure.

Sonicare App Integration

Connect your toothbrush to the Sonicare app for guided brushing, tips, and personalized content to help achieve your oral health goals.

Sonicare App

Brushing insights
in the palm of your
hand



Image: A hand holding a smartphone displaying the Sonicare App interface, showing brushing insights and a 'Reduce Pressure' prompt.

Official Product Video

Your browser does not support the video tag.

Video: An official overview of the Philips Sonicare 7300 Series Electric Toothbrush, demonstrating its features and benefits.

MAINTENANCE

Brush Head Replacement

Dental professionals recommend replacing your brush head every 3 months. Your Philips Sonicare toothbrush includes a brush head replacement reminder that tracks how often and how hard you brush, then reminds you when it's time for a new brush head.

Cleaning the Toothbrush

After each use, rinse the brush head and bristles thoroughly. Wipe the handle with a damp cloth. Do not immerse the handle in

water.

TROUBLESHOOTING

If you encounter issues with your Philips Sonicare 7300 Series Electric Toothbrush, consider the following:

- **Toothbrush not turning on:** Ensure the toothbrush is fully charged.
- **Pressure sensor light activates frequently:** Reduce the pressure applied during brushing. The sensor is designed to help protect your gums from excessive force.
- **Reduced brushing effectiveness:** Check if the brush head needs replacement. The brush head replacement reminder will indicate when a new one is due.

SPECIFICATIONS

Model Number	HX7420/07
Brand	Philips Sonicare
Power Source	Battery Powered
Battery Life	Up to 21 days
Number of Settings	12 (4 modes, 3 intensity levels)
Special Features	Brushing Timer, Lightweight, Rechargeable, Sonic, Pressure Sensor
Item Weight	0.96 Pounds
Handle Material	Plastic
Bristle Material	Nylon

SAFETY INFORMATION

For your safety and optimal performance of the device, please adhere to the following guidelines:

- Replace brush head every 3 months.
- Do not use for purposes other than oral hygiene.

WARRANTY AND SUPPORT

Your Philips Sonicare 7300 Series Electric Toothbrush is covered by a **2-Year Manufacturer Warranty**. For warranty claims or technical support, please refer to the contact information provided in your product packaging or visit the official Philips Sonicare website.

ASK A QUESTION ABOUT THIS MANUAL

Ask about setup, troubleshooting, compatibility, parts, safety, or missing instructions. Manuals+ will review the question and use this page’s manual context to help answer it.

Name

Anonymous

Email

you@example.com

After a successful submission, we securely hash the supplied account or form email before sending it to Microsoft Advertising for conversion attribution. [Privacy details](#).

Question

Example: How do I reset this device or fix this error?

Details

Model number, symptoms, what you tried, or the section of the manual you are using.

Submit question