

XASAWHHSV-EC9C1380

Generic Replacement Brush Heads User Manual

For MES610, MES609, MES613, MES611 Electric Toothbrushes

1. INTRODUCTION

This manual provides instructions for the proper use and maintenance of your Generic Replacement Brush Heads. These brush heads are designed for compatibility with MES610, MES609, MES613, and MES611 electric toothbrush models, ensuring effective and comfortable oral hygiene.

2. PRODUCT FEATURES

- **Compatibility:** Designed to fit MES610, MES609, MES613, and MES611 toothbrush models.
- **Soft Nozzles:** Equipped with soft bristles for gentle contact, suitable for daily use and comfort.
- **Precision Fit:** Molded for a perfect match with MES series models, allowing easy installation and stable performance.
- **Convenient Packaging:** Lightweight and compact for easy storage and portability.

Efficient Cleaning Without The Hassle

Three-dimensional cut and design three brush heads to meet different cleaning needs.



Figure 1: Close-up view of the brush head bristles, illustrating the three-dimensional cut and design for varied cleaning needs.

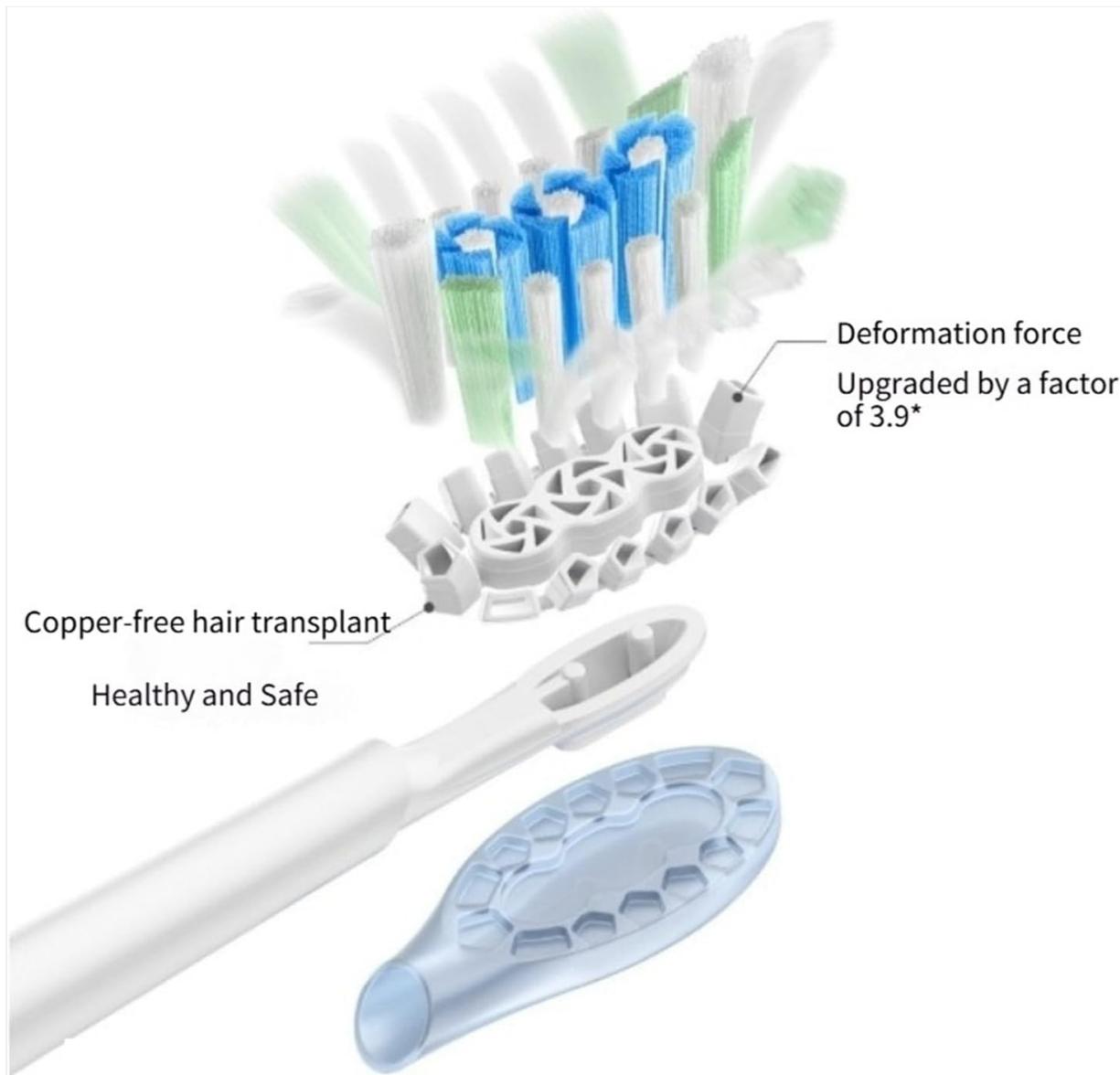


Figure 2: Exploded diagram of the brush head, highlighting the internal structure, copper-free hair transplant, and enhanced deformation force for effective cleaning.

3. SETUP AND INSTALLATION

1. **Remove Old Brush Head:** Gently pull the existing brush head straight off the toothbrush handle. Do not twist.
2. **Attach New Brush Head:** Align the new replacement brush head with the shaft of the toothbrush handle. Push the brush head down firmly until it clicks into place. Ensure it is securely seated.
3. **Verify Fit:** Lightly tug on the brush head to confirm it is properly attached and does not wobble.

4. OPERATING INSTRUCTIONS

Once the replacement brush head is securely installed, operate your electric toothbrush as per its original manufacturer's instructions. Apply a small amount of toothpaste to the bristles and guide the brush head gently across all surfaces of your teeth and gums. Avoid excessive pressure.

5. MAINTENANCE

- **Rinse After Use:** After each brushing session, rinse the brush head thoroughly under running water to remove any toothpaste residue and debris.

- **Air Dry:** Allow the brush head to air dry in an upright position.
- **Replacement:** Dentists recommend replacing brush heads every three months, or sooner if the bristles become frayed or worn.

6. TROUBLESHOOTING

- **Brush Head Does Not Fit:** Ensure you have the correct replacement brush head for your specific MES series toothbrush model (MES610, MES609, MES613, MES611).
- **Brush Head Feels Loose:** Re-attach the brush head, ensuring it clicks firmly into place. If it continues to feel loose, inspect both the brush head and the toothbrush handle for any damage.
- **Bristles Fraying Quickly:** This may indicate brushing too hard. Reduce pressure during brushing. Ensure proper maintenance as described in Section 5.

7. SPECIFICATIONS

Brand	Generic
Compatible Models	MES610, MES609, MES613, MES611
Material	Plastic (Soft Nozzles)
Age Group	Adults
Item Model Number	XASAWHHSV-EC9C1380
Country of Origin	China
Color	Purple (for this specific variant)

8. WARRANTY AND SUPPORT

For warranty information or product support, please refer to the documentation provided with your original electric toothbrush handle or contact the seller directly. This manual pertains specifically to the replacement brush heads.