



[Manuals.plus](#) /

› [Retekess](#) /

› Retekess TH017 Automatic Wi-Fi Pill Dispenser Instruction Manual

Retekess TH017

Retekess TH017 Automatic Wi-Fi Pill Dispenser Instruction Manual

Model: TH017

INTRODUCTION

Thank you for choosing the Retekess TH017 Automatic Wi-Fi Pill Dispenser. This device is designed to simplify medication management with its automatic dispensing, Wi-Fi connectivity, and customizable reminders. This manual provides essential information for safe and effective use of your pill dispenser.

IMPORTANT SAFETY INFORMATION

The product packaging may contain different information than what is shown on our website. Please always read the label, warnings, and instructions provided on the product before using or consuming it. Keep out of reach of children and pets. Do not expose the device to extreme temperatures or moisture. Use only the specified power sources.

PACKAGE CONTENTS

Please verify that all items are present in the package:

- Retekess TH017 Wi-Fi Smart Pill Box x1
- Pill Compartment (21 sections) x1 (pre-installed in the pill box)
- AA Batteries x4
- User Manual x1
- Quick Start Guide x1
- Key x2
- Dose Ring x9
- Charging Cable (Type-C) x1
- Adapter x1

Package Includes




TH017 WIFI Smart Pill Box *1



Pill Compartment *21
*Placed inside
the pill box



AA Batteries *4

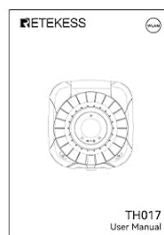


User Manual *1



Quick Start Guide *1



Key *2

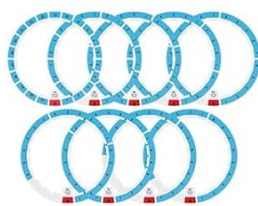

Dose Ring *9



Charging Cable *1



Adapter *1

Image: All components included in the Reteless TH017 Automatic Pill Dispenser package.

SETUP

1. Powering the Device

The TH017 offers dual power options: AA batteries or Type-C power.

- **Battery Power:** Insert the 4 included AA batteries into the battery compartment located on the underside of the device. Ensure correct polarity.
- **Type-C Power:** Connect the provided Type-C charging cable to the device and a power adapter, then plug it into a wall outlet.

The device features power-off memory, ensuring your medication schedule is not interrupted during power transitions or outages.



Image: The pill dispenser highlighting its Type-C power input and AA battery compartment.


2. Loading Medication

The dispenser has 21 independent compartments for a week's worth of medication. Use the provided dose rings to customize the dispensing frequency.

1. Unlock the top lid using one of the provided keys.
2. Carefully open the transparent lid to access the pill compartments.
3. Place your medication into the designated compartments. Each compartment is large enough for various capsules, tablets, and softgels.
4. Install the appropriate dose ring (e.g., 1x/day, 2x/day) to match your medication schedule.
5. Close the transparent lid and lock the top lid securely.

21 Doses, All in One

Effortlessly store your medication in one convenient box

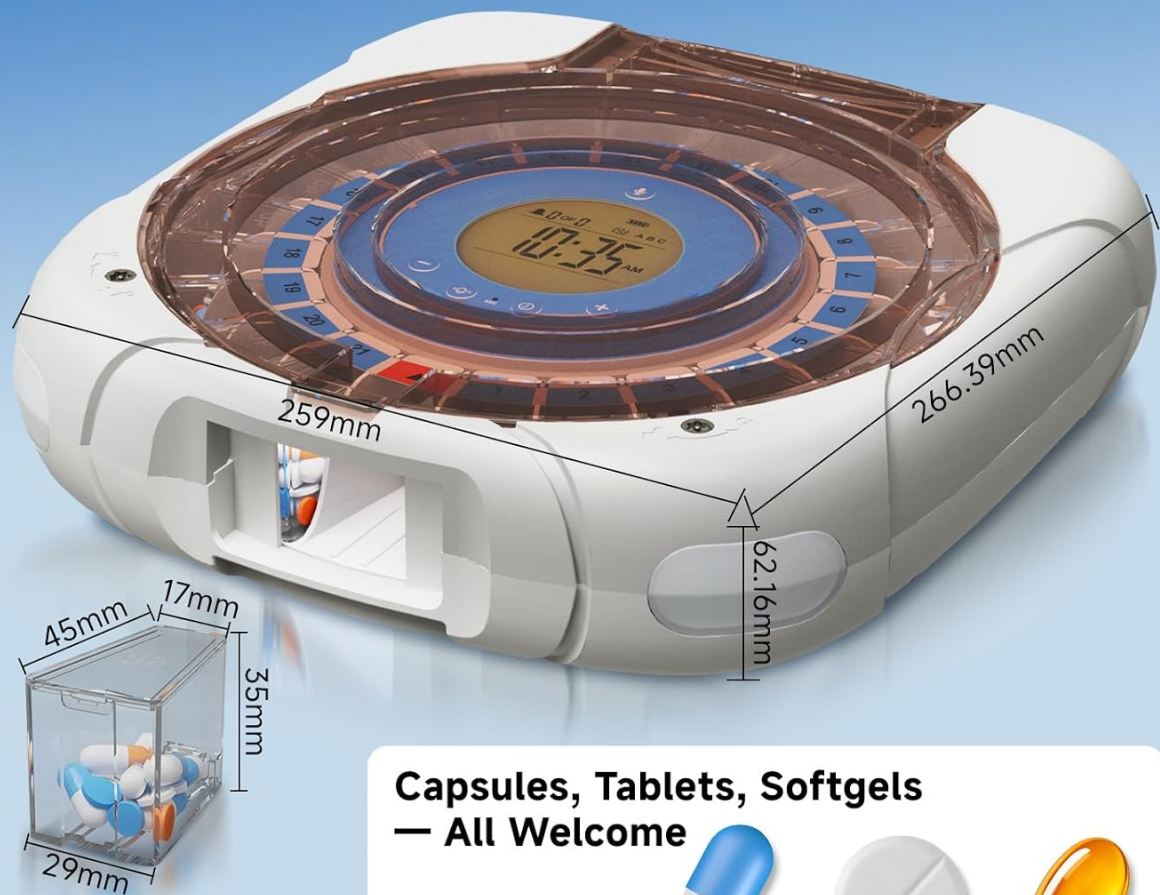


Capsules, Tablets, Softgels
— All Welcome



Image: The pill dispenser with an illustration of its 21 compartments and capacity for different pill types.

Personalized Medication Plan

Simplify medication management with
a customizable, automated system

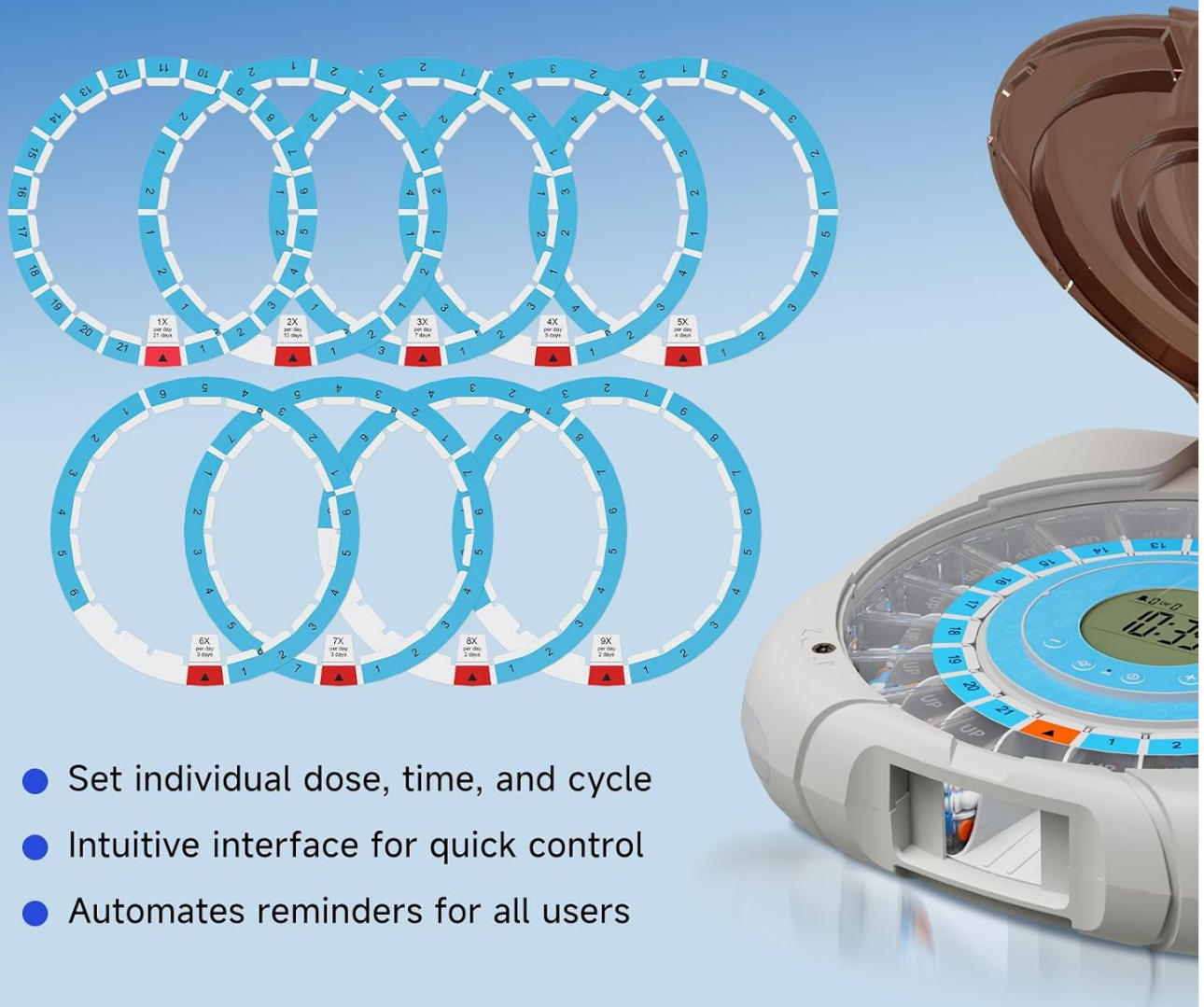

- Set individual dose, time, and cycle
- Intuitive interface for quick control
- Automates reminders for all users

Image: Visual representation of dose rings used for personalizing medication plans.

3. App Connectivity and Schedule Setup

The device connects to a mobile app via Wi-Fi for seamless management.

1. Download the Reteless companion app from your mobile device's app store.
2. Follow the in-app instructions to connect the TH017 to your home Wi-Fi network.
3. Once connected, use the app to set your personalized medication schedule, including dose times, frequency, and specific days.
4. The app provides multi-device alerts and real-time updates, allowing you to manage medication settings remotely.

Wireless Wi-Fi Communication

Manage your medication settings anytime, anywhere

Multi-Device Alerts

Sends app notifications to all bound phones



Remote Management


Set/modify schedules anytime via app



Image: The pill dispenser demonstrating wireless Wi-Fi communication and remote management via a mobile app.

Configurazione comoda tramite app

La connettività Wi-Fi wireless garantisce un funzionamento senza problemi

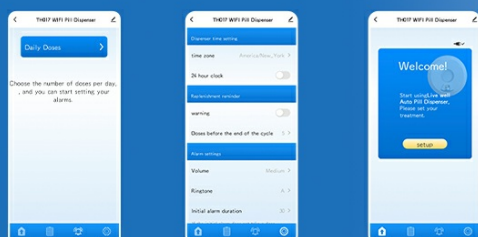

Regolazione dei promemoria

Personalizzazione delle impostazioni

Panoramica chiara



Image: The mobile app interface for convenient setup and customization of the medication schedule.

1. Automatic Pill Dispensing

When a scheduled medication time arrives, the dispenser will sound an alarm and automatically open the compartment containing the dose. This allows for easy access to your medication.

- When the alarm sounds, retrieve the pills from the opened compartment.
- To silence the alarm, simply tap the corresponding notification in the mobile app.



Image: The automatic pill dispensing mechanism in action, showing pills being ejected.

2. Customizable Alerts and Volume

The TH017 allows for personalized alerts to ensure you never miss a dose.

- **Ringtone Customization:** You can record your own voice or use preset ringtones for a more familiar and personalized reminder. This feature is managed through the mobile app.
- **Volume Levels:** The dispenser offers four adjustable volume levels (High, Medium, Low, Off) to suit your environment, whether quiet or noisy. Adjust this setting via the app or directly on the device if available.

Customizable Ringtone

Use family or personal voice for a warm, familiar alert



Image: The pill dispenser illustrating customizable ringtone options, including voice recording.



Image: Diagram illustrating the four adjustable volume levels for the device's alerts.

3. LCD Display Information

The clear LCD screen provides essential information at a glance:

- Current Time (12-hour / 24-hour format)

- Next Medication Time
- Total Number of Alarms Set
- Remaining Alarms for the Day
- Battery Level Indicator
- Wi-Fi Connection Status
- Volume Setting
- Ringtone Setting

Large LCD Display

Easy-to-read LCD display and easy setup for everyone

- Battery
- Volume
- Wi-Fi Connection
- Total Number of Alarms
- Ringtone
- Remaining Alarms for the Day
- Time (12-hour / 24-hour)
- Next Medication Time

Screen Size (Round)

60mm



Image: The large LCD display showing various status indicators and medication details.

Pannello LCD da 60 mm

Informazioni di facile lettura a colpo d'occhio



Image: Detailed view of the 60mm LCD panel and its displayed information.

MAINTENANCE

1. Cleaning the Dispenser

To ensure optimal performance and hygiene, clean your pill dispenser regularly.

- Turn off the device and disconnect it from any power source.
- Wipe the exterior with a soft, damp cloth. Do not use abrasive cleaners or immerse the device in water.
- For the pill compartments, remove any residue with a dry brush or cloth. Ensure all compartments are dry before reloading medication.

2. Secure Storage

The TH017 is designed with double-layer protection to keep your medication safe and effective.

- The UV-resistant top lid protects medication from light exposure, which can degrade certain drugs.
- The bottom layer features a double mechanical lock to prevent accidental openings, ensuring safety, especially around children and pets.

Secure Storage

Dependable protection for your medications



Blocks Light & Moisture
For full potency



Dual Mechanical Locks
For mis-touch prevention

Image: The pill dispenser highlighting its secure storage features, including UV protection and mechanical locks.

TROUBLESHOOTING

Problem	Possible Cause	Solution
Device not powering on	Low or dead batteries; Type-C cable not connected properly.	Replace AA batteries; ensure Type-C cable is securely connected to both device and power source.
Pills not dispensing automatically	Compartment jammed; schedule not set correctly; low power.	Check for obstructions in the compartment; verify schedule in the app; ensure adequate power supply.
Wi-Fi connection issues	Incorrect Wi-Fi password; device too far from router; router issues.	Re-enter Wi-Fi password in the app; move device closer to router; restart your Wi-Fi router.

Problem	Possible Cause	Solution
Alarm not sounding	Volume set to 'Off'; app notification issues; device malfunction.	Check volume settings in the app or on the device; ensure app notifications are enabled; contact support if issue persists.

SPECIFICATIONS

Feature	Detail
Model	TH017
Brand	Retekess
Power Source	4 AA batteries (included) or Type-C cable
Number of Compartments	21
Connectivity	Wi-Fi Wireless
Special Features	Automatic dispensing, Built-in alarm, Mobile connectivity, UV-resistant lid, Double mechanical lock, Power-off memory
Dimensions	32.2 x 32.1 x 6.4 cm
Weight	1.36 kg (Product), 1365 Grams (with packaging)
Country of Origin	China

WARRANTY AND SUPPORT

Retekess products are manufactured to high-quality standards. For warranty information, please refer to the documentation included with your purchase or visit the official Retekess website. If you encounter any issues or require technical assistance, please contact Retekess customer support through their official channels. Please have your model number (TH017) and purchase details ready when contacting support.

PRODUCT VIDEOS

No relevant official seller videos were found for this product based on the provided data.