

Rm-re

Generic 5-in-1 Anti-Aging Beauty Device User Manual

Model: Rm-re

1. INTRODUCTION

Thank you for choosing the Generic 5-in-1 Anti-Aging Beauty Device. This advanced device combines Radio Frequency (RF), Electrical Muscle Stimulation (EMS), LED Light Therapy, Sonic Vibration, and Warming Therapy to support your skincare routine. Please read this manual thoroughly before first use to ensure safe and effective operation.

2. SAFETY INFORMATION

- Do not use on broken, irritated, or infected skin.
- Consult a physician before use if you have any medical conditions, including but not limited to heart conditions, epilepsy, cancer, or if you are pregnant.
- Keep out of reach of children.
- Do not immerse the device in water. Clean with a damp cloth only.
- Discontinue use if you experience discomfort or adverse reactions.
- Use only the provided charging cable.

3. PRODUCT OVERVIEW

The device features a sleek dark gray design with an ergonomic grip for comfortable handling. It integrates multiple technologies for comprehensive skin care.

Key Features:

- **RF Radio Frequency:** Utilizes 13-pole RF technology for skin tightening and lifting.
- **EMS Microcurrent:** Delivers gentle electrical pulses to stimulate facial muscles.
- **LED Light Therapy:** Features 9 colored LED beads for various skin benefits.
- **Sonic Vibration:** Enhances product absorption and promotes circulation.
- **Warming Therapy:** Provides a soothing warmth during treatment.

- **Adjustable Intensity:** 5 intensity levels for customized treatment.

Included Components:

- Microcurrent Facial Device
- Large Care Head
- Eye Care Head
- Lip Care Head
- USB Charging Cable
- User Manual (this document)

5 major beauty features
Revive beautiful power



Figure 3.1: The 5-in-1 Anti-Aging Beauty Device. This image displays the device's main body and the treatment head, highlighting its compact and ergonomic design.

Hardcore anti-aging

Look for 5 Weight Technology

Multipole radio frequency

Hyperpulse electroporation

EMS Microcurrent Skin Exercise

LED red and blue color light
double color light phototherapy

Pulse vibration



Figure 3.2: Overview of the device's five core technologies: Multipolar Radio Frequency, Hyperpulse Electroporation, EMS Microcurrent Skin Exercise, LED Red and Blue Color Light Phototherapy, and Pulse Vibration.

4. SETUP

4.1 Charging the Device

1. Connect the provided USB charging cable to the charging port on the device.
2. Plug the other end of the USB cable into a compatible USB power adapter (not included) or a computer USB port.
3. The charging indicator lights will illuminate to show charging status. A full charge takes approximately 2 hours.
4. Once fully charged, disconnect the device from the power source.

4.2 Initial Preparation

- Ensure your skin is clean and dry before use.
- Apply a thin layer of conductive gel, serum, toner, lotion, or essence to the treatment area. This is essential for optimal performance of RF and EMS functions.

5. OPERATING INSTRUCTIONS

5.1 Powering On/Off and Mode Selection

1. Press and hold the Power button for 2 seconds to turn the device ON.
2. Press the Power button repeatedly to cycle through the available modes. The illuminated indicators on the device will show the active mode.
3. To turn OFF the device, press and hold the Power button for 2 seconds.

5.2 Intensity Adjustment

- The device has 5 intensity levels.
- Use the Intensity Adjustment buttons (usually marked '+' and '-') to increase or decrease the treatment intensity.
- Start with the lowest intensity and gradually increase to a comfortable level.

5.3 Treatment Modes

Brand new experience, stunning upgrade



**13 pole RF lead, wide coverage with 9
colored LED beads illuminating the metal body,
matte texture, skin friendly and smooth**

Figure 5.1: The 13-pole RF lead and 9 colored LED beads on the device's treatment head, designed for enhanced skin contact and light therapy.

3D photo electric heating technology reveals the mystery of youthful skin

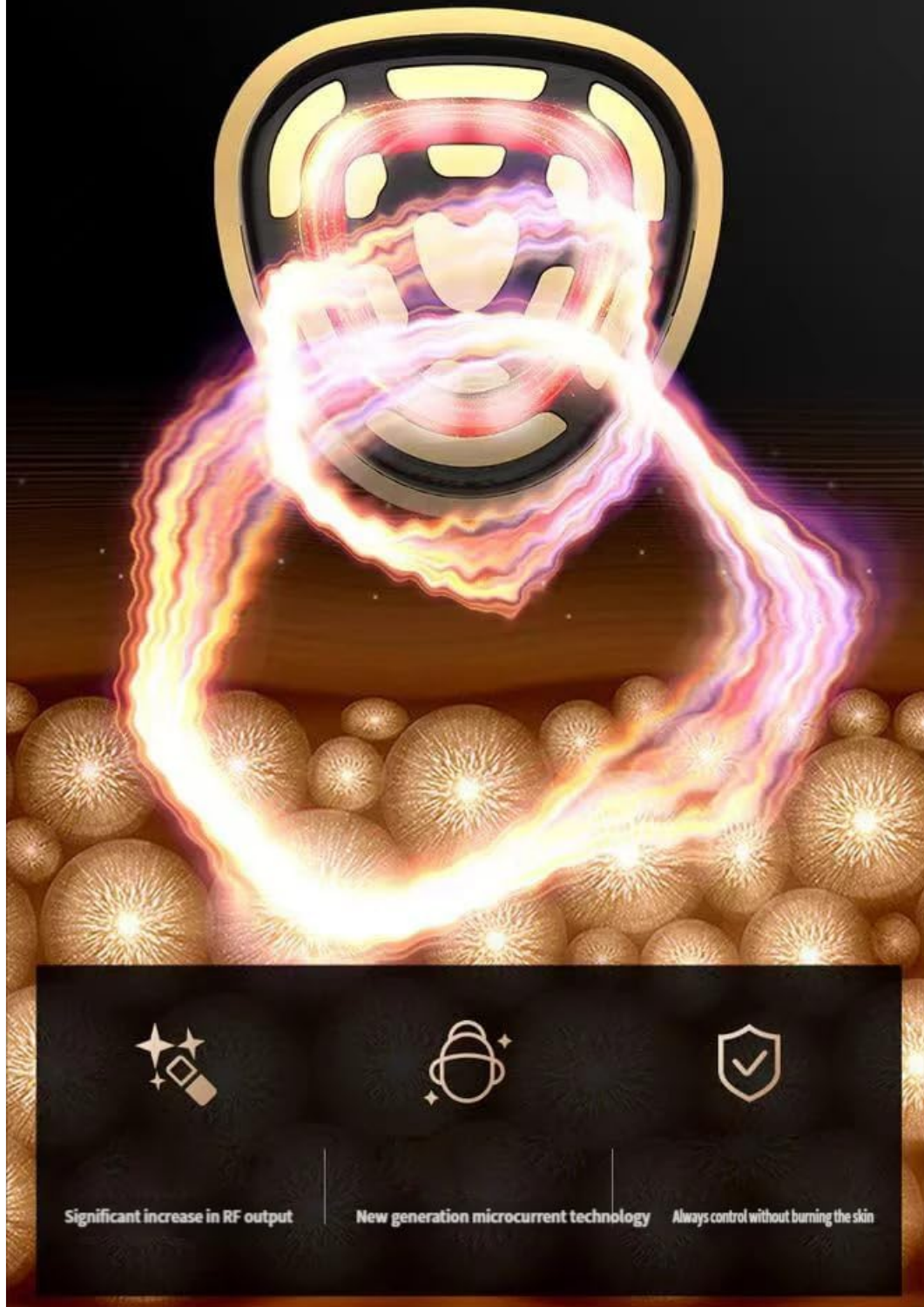


Figure 5.2: Visual representation of the 3D photo electric heating technology, showing its effect on skin for rejuvenation.

RM Rapid Care Model:

- Combines Pulse Current, Ultra-Pulse Electroporation, and Vibration.
- Recommended for quick skin care, 1-2 times a week, 10 minutes per session.
- Use with toner, lotion, and essence for optimal results.

Multi-effect combination care mode easily satisfies daily skin care



Quick skin care 1-2 times
a week, 10 minutes/time

RM Rapid Care Model

Pulse current + ultra-pulse electroporation +
radio frequency wave + LED + pulse
vibration

Use with toner, lotion, and Jinghua



受试者感觉 Subject feels 受试者感觉

Figure 5.3: Application of the RM Rapid Care Model. This image shows a user applying the device to their face, illustrating the recommended technique for the rapid care mode.

RE Nursing Model:

- Combines EMS Microcurrent and Pulse Vibration.
- Recommended for high micro care for facial expression muscles, once a day, 10 minutes per session.
- Use with toner, lotion, and essence.

High micro care for facial expression muscles

Once a day, 10 minutes per session

RE nursing model

EMS microcurrent+LED+pulse vibration

Use with toner, lotion and essence

Clearly pale and tight | Subject feedback



98%
Subject's feeling
Tight to



98%
Subject's feeling
Danwen

97%
Subject's feeling
Expansion bomb

Figure 5.4: Application of the RE Nursing Model. This image shows a user applying the device to their face, illustrating the recommended technique for the nursing mode.

5.4 LED Light Therapy

The device incorporates LED light therapy with different colors for targeted skin benefits.



Figure 5.5: LED Color Light Skincare. This diagram illustrates the scientific principle of lighting skincare, showing red light,

blue light, and red-blue alternating light modes and their respective wavelengths.

- **Red Light:** Stays bright for a long time. Typically associated with collagen production and anti-aging effects.
- **Blue Light:** Long lasting blue light. Often used for its calming properties and to address skin concerns.
- **Red-Blue Alternating Light:** Combines the benefits of both red and blue light.

6. MAINTENANCE

6.1 Cleaning the Device

1. Ensure the device is turned OFF and disconnected from the charger before cleaning.
2. Wipe the treatment head and body with a soft, damp cloth.
3. Do not use abrasive cleaners, alcohol, or harsh chemicals.
4. Allow the device to air dry completely before storing.

6.2 Storage

- Store the device in a cool, dry place, away from direct sunlight and extreme temperatures.
- Keep the device in its original packaging or a protective pouch to prevent damage.

7. TROUBLESHOOTING

- **Device not turning on:** Ensure the device is fully charged. If the issue persists, contact customer support.
- **No sensation during use:** Ensure sufficient conductive gel/serum is applied to the skin. Increase the intensity level gradually.
- **Charging indicator not working:** Check the USB cable and power adapter connection. Try a different USB port or adapter.
- **Skin redness or irritation:** Reduce intensity or frequency of use. Ensure proper hydration of the skin. If irritation persists, discontinue use and consult a dermatologist.

8. SPECIFICATIONS

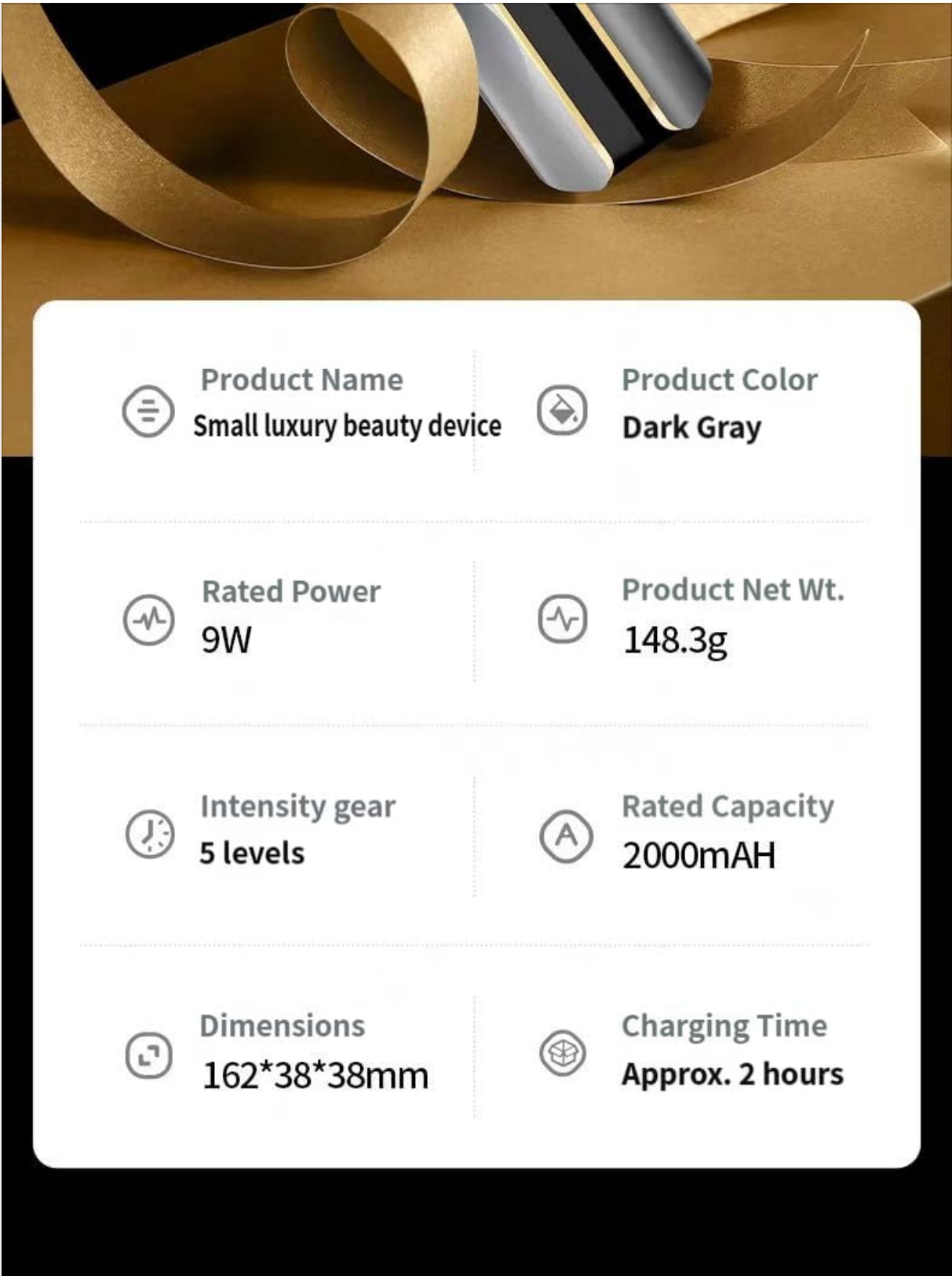


Figure 8.1: Summary of Product Specifications.

Feature	Specification
Product Name	5-in-1 Anti-Aging Beauty Device
Model Number	Rm-re
Product Dimensions	6.38 x 1.5 x 1.5 inches (162 x 38 x 38 mm)
Item Weight	0.5 Pounds (148.3g)

Color	Dark Gray
Rated Power	9W
Battery Capacity	2000mAh
Charging Time	Approx. 2 hours
Intensity Levels	5 levels
Material Type	ABS Plastic, Silicone, Metal Electrodes

9. WARRANTY AND SUPPORT

9.1 Warranty Information

This product comes with a standard manufacturer's warranty against defects in materials and workmanship. Please refer to your purchase documentation for specific warranty terms and duration. The warranty does not cover damage caused by misuse, accidents, unauthorized repairs, or normal wear and tear.

9.2 Customer Support

If you have any questions, concerns, or require assistance with your device, please contact our customer support team. Contact information can typically be found on the product packaging or the retailer's website where the product was purchased.