

COBEEK B0HC2D8122

COBEEK Ultrasonic Retainer Cleaner Machine (Model B0HC2D8122) - Instruction Manual

INTRODUCTION

This manual provides detailed instructions for the safe and effective use of your COBEEK Ultrasonic Retainer Cleaner Machine. This device utilizes advanced 52kHz ultrasonic technology and 254nm UV-C lights to provide thorough cleaning for various dental appliances and small items. Please read this manual carefully before operation and retain it for future reference.

SAFETY INFORMATION

- Do not operate the machine without water in the tank.
- Always use cold water for cleaning.
- The internal basket is removable and can be attached for convenience.
- It is not recommended to use this machine to clean soft silicone items.
- Ensure the power cord is securely connected and kept away from water.
- Keep out of reach of children.

PRODUCT FEATURES

- **Advanced 52kHz Ultrasonic Cleaning:** High-frequency ultrasonic waves generate micro-bubbles for deep cleaning in narrow gaps and hard-to-reach surfaces without damage.
- **Upgraded 6 UV-C Lights:** Equipped with 6 upgraded 254nm UV-C lights for 360° deep cleaning, odor removal, and surface stain reduction.
- **Dual Practical Cleaning Modes:** Features a 5-minute quick cleaning mode for daily maintenance and a 10-minute powerful cleaning mode for deep cleaning, both with UV light assistance.
- **250ML Large Capacity & Portable:** The 304 stainless steel tank accommodates retainers, dentures, mouth guards, toothbrush heads, and jewelry. Compact design for home or travel use.
- **Whisper-quiet Operation:** Operates below 45dB with non-slip rubber feet for stable and quiet use.
- **Top 3-Hole Design:** Convenient for cleaning toothbrushes and makeup brushes, also aids in heat dissipation and ventilation.

- **Digital Countdown Timer:** Provides clear indication of cleaning progress.

Professional Ultrasonic Cleaner

52kHz High Frequency Vibration with Pulse Cavitation Effect
to Removing Particles and Dirt from Hard-to-reach Areas



52000Hz
High Vibration



Powerful
cleaning



50 μ m
Cavitation Effect



Figure 1: Illustration of 52kHz high-frequency ultrasonic vibration and cavitation effect for deep cleaning.

254nm U-V Light Kill 99.99% of Dirt

UVC+UVA Double Ultraviolet Light Combination,
360° Deep Cleaning the Items



Figure 2: The cleaner's 6 UV-C lights providing 360° deep cleaning.

Newly Upgraded Design

New Features and Increased Capacity Make it a Good Substitute for Traditional Ultrasonic Cleaners

Product Compare			
Ultrasonic Frequency	✔ 52000Hz	42000Hz	45000Hz
UV Lights	✔ 6 UV Lights	NO	3 UV Lights
Display Screen	✔ YES	NO	NO
Holes for Cleaning	✔ 3 Holes	NO	2 Holes
Tank Capacity	✔ 250ML	150ML	200ML

Figure 3: Feature comparison of the COBEEK cleaner with other models.

PACKAGE CONTENTS

Verify that all items are present in the package:

- COBEEK Ultrasonic Retainer Cleaner Unit
- Power Adapter
- Removable Cleaning Basket
- Tweezers
- Cleaning Brush
- User Manual

SETUP

Follow these steps to set up your ultrasonic cleaner:

1. **Add Water:** Open the lid and add an appropriate amount of cold water to the stainless steel tank. Do not exceed the MAX fill line.
2. **Place Items:** Carefully place the items to be cleaned (e.g., retainer, denture, jewelry) into the removable basket, then lower the basket into the water. Ensure items are fully submerged.
3. **Connect Power:** Close the lid and plug the power adapter into the unit, then into a standard electrical outlet.
4. **Select Mode:** The digital display will show the default cleaning time. Press the power button to cycle through the cleaning modes.

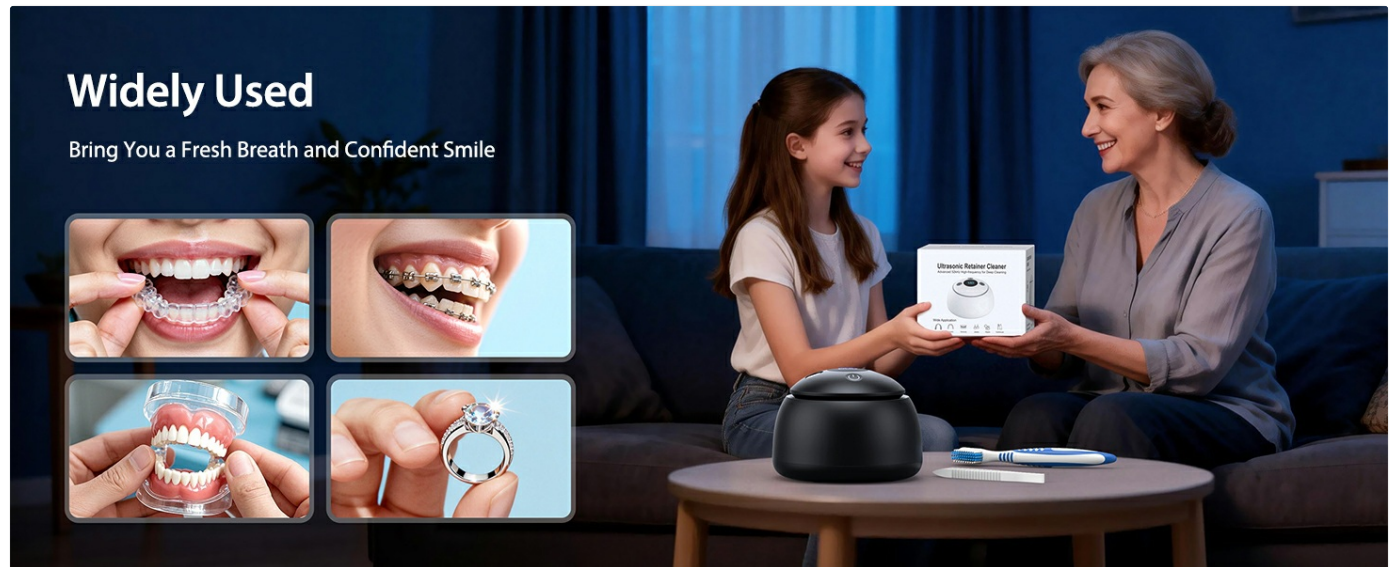


Figure 4: Visual guide for the 4 easy steps to use the retainer cleaner.

OPERATING INSTRUCTIONS

The COBEEK Ultrasonic Retainer Cleaner features two cleaning modes for optimal results:

- **5-Minute Quick Clean (300s):** Ideal for daily maintenance and light cleaning. This mode combines ultrasonic cleaning with UV-C light.
- **10-Minute Deep Clean (600s):** Suitable for thorough cleaning of heavily soiled items. This mode also combines ultrasonic cleaning with UV-C light.

To Operate:

1. After setup, press the power button once to activate the 5-minute quick clean mode. The digital display will show '300'.
2. Press the power button a second time to switch to the 10-minute deep clean mode. The digital display will show '600'.
3. The cleaning cycle will begin automatically. The digital display will count down the remaining time.
4. Once the cycle is complete, the machine will automatically shut off.
5. To manually stop a cycle, press and hold the power button.
6. To switch the UV light on or off during a cycle, long press the power button.

250ML Large Capacity & Wide Application

Ideal for Dental Supplies, Jewelry, Pacifier, Makeup Brush and More



Figure 5: Digital display and control button for selecting cleaning modes.

MAINTENANCE

Proper maintenance ensures the longevity and performance of your cleaner:

- **Water Replacement:** Always replace the cold water after each cleaning cycle to maintain hygiene and cleaning effectiveness.
- **Machine Rest:** Allow the machine to rest for at least 15 minutes after running two consecutive cleaning cycles to prevent overheating and prolong its lifespan.
- **Unit Cleaning:** Wipe the exterior of the unit with a soft, damp cloth. Do not submerge the main unit in water.
- **Tank Cleaning:** The stainless steel tank can be rinsed with water. For stubborn residue, use a mild detergent and the provided cleaning brush.
- **Storage:** Store the cleaner in a cool, dry place when not in use.

TROUBLESHOOTING

If you encounter any issues, please refer to the following common solutions:

- **Unit Not Turning On:** Ensure the power adapter is securely plugged into both the unit and a working electrical outlet. Check if the outlet has power.
- **No Ultrasonic Action:** Verify that there is water in the tank and that the unit is powered on and a cleaning mode is selected.
- **Poor Cleaning Results:** Ensure items are fully submerged. For heavily soiled items, use the 10-minute deep clean mode or pre-rinse items to remove large debris. Consider adding a small amount of specialized cleaning solution (not included) if recommended for your items.
- **Unusual Noise:** A slight buzzing sound is normal during ultrasonic operation. If the noise is excessively loud or unusual, ensure the unit is on a stable, flat surface and that no foreign objects are in the tank.

If the problem persists, please contact COBEEK customer support for assistance.

SPECIFICATIONS

Model	B0HC2D8122
Ultrasonic Frequency	52kHz
UV-C Wavelength	254nm
Tank Capacity	250ML
Noise Level	Below 45dB
Cleaning Modes	5-minute (Quick), 10-minute (Deep)
Tank Material	304 Stainless Steel
Package Dimensions	8 x 8 x 8 inches
Weight	13.76 ounces
Manufacturer	COBEEK

WARRANTY AND SUPPORT

COBEEK is committed to providing reliable products and excellent customer service.

- **2 Years Warranty:** Your COBEEK Ultrasonic Retainer Cleaner Machine is covered by a 2-year warranty from the date of purchase.
- **24 Hours Customer Support:** For any inquiries or issues, our customer support team is available 24 hours a day.
- **Lifetime Technical Support:** We offer lifetime technical support for your product.

Please refer to the contact information provided with your product or visit the official COBEEK website for support.



Figure 6: COBEEK's commitment to customer service and support.

OFFICIAL PRODUCT VIDEOS

COBEEK Ultrasonic Retainer Cleaner Machine Overview

Your browser does not support the video tag.

This video provides a general overview of the COBEEK Ultrasonic Retainer Cleaner Machine, demonstrating its features and basic operation. It highlights the ultrasonic cleaning process and the UV-C light function.

ASK A QUESTION ABOUT THIS MANUAL

Ask about setup, troubleshooting, compatibility, parts, safety, or missing instructions. Manuals+ will review the question and use this page's manual context to help answer it.

Name

Anonymous

Email

you@example.com

After a successful submission, we securely hash the supplied account or form email before sending it to Microsoft Advertising for conversion attribution. [Privacy details](#).

Question

Example: How do I reset this device or fix this error?

Details

Model number, symptoms, what you tried, or the section of the manual you are using.

Submit question