



# CINCOM CM-052G Percussion Massage Gun Operating Instruction

[Home](#) » [CINCOM](#) » CINCOM CM-052G Percussion Massage Gun Operating Instruction 

## Contents

- [1 CINCOM CM-052G Percussion Massage Gun](#)
- [2 Safety Precautions](#)
- [3 Feature](#)
- [4 Product Diagram](#)
- [5 Operating](#)
- [6 Charging Instruction](#)
- [7 Package Included](#)
- [8 Specifications](#)
- [9 FREQUENTLY ASKED QUESTIONS](#)
- [10 VIDEO – PRODUCT OVERVIEW](#)
- [11 References](#)
- [12 Related Posts](#)



**CINCOM CM-052G Percussion Massage Gun**



## Instruction

Thanks for purchasing our product. Please read the instructions carefully before use and keep it well for further reference. If you have any questions during use, please feel free to contact us through Amazon or mail us at [service@cincomhealth.com](mailto:service@cincomhealth.com).

## Safety Precautions

### Warnings & cautions


1. Those with any of the following conditions or persons who are receiving medical treatment should consult the doctor before using the machine:
  - Using pacemakers or other medical devices that are susceptible to electrical interference;
  - Suffering from malignant tumors;
  - Suffering from cardiac diseases;
  - Having serious peripheral neuropathy dysfunction or sensory disturbance caused by diabetes;
  - Being unsuitable to do the massage because of traumas on the body;
2. Keep it out of the reach of infants, children, and people without the ability to use it independently.
3. Do not immerse in water. Keep away from liquids or heat sources. Keep ventilation ports free from dust and debris.
4. Do not remove screws or attempt to disassemble.
5. Do not drop or misuse the Percussion Massage Gun.
6. Do not use the machine continuously for more than one hour, stop using it until you have a rest for 30 minutes.


7. Unplug the charger from the device after charging or before use.
8. it can only be used on soft tissues of the body, which will not cause pain or discomfort. Do not use it on any hard or bony area of the header body.
9. Fingers, hair, or any other part of the body should be kept away from the pole and back of the machine to avoid clipping.
  - Do not touch this part when the massage head is not installed and the device is on
  - Do not touch this part when the massage head is installed and the device is on
10. Abrasions may occur during use no matter how much intensity you choose. Always check the area for treatment and stop at the first sign of pain or discomfort.

## Feature

1. Powerful Brushless High-Torque Motor with Quiet Glide Technology
2. Rechargeable Lithium-Ion Battery { 4-6 Hours of Use Per Charge}
3. Multi-Speed Settings (4 Speed) for Optimal Pressure
4. Cordless and Lightweight (1 Lb)
5. One-button-control
6. Interchangeable Massage Heads (x4)

## Product Diagram

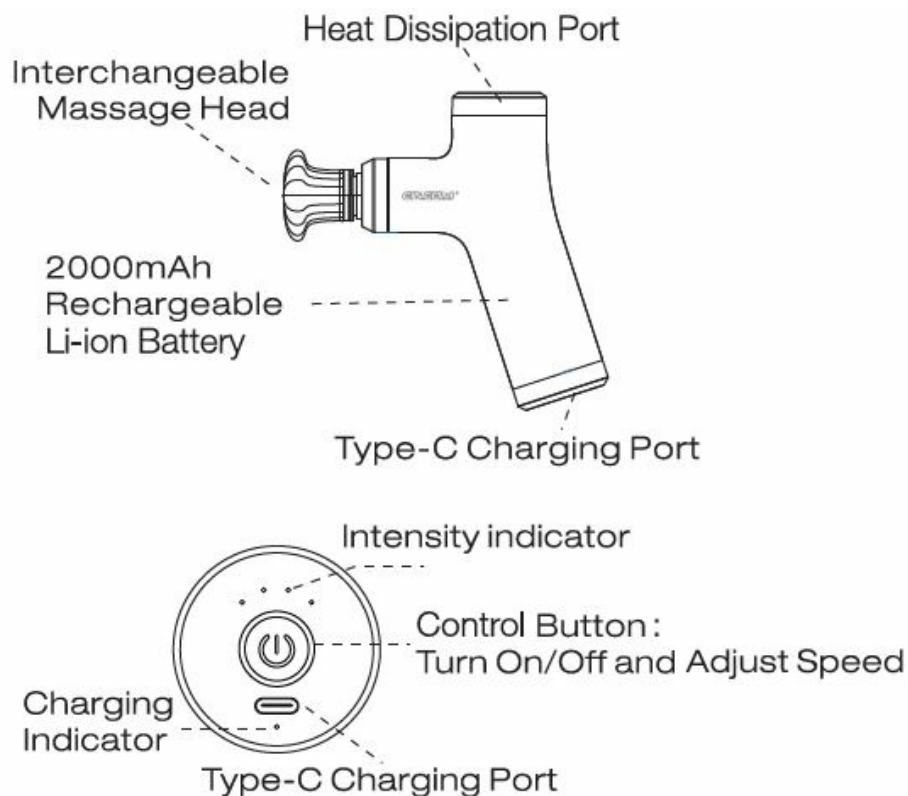


-  **No.1: Big Ball Head**  
Suitable for back, buttocks, etc.



- **No. 2: Flat Head**

Suitable for relaxing and muscle shaping.

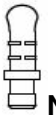

- **No. 3: Bullet Head**

Suitable for meridians, palms, soles, etc.



- **No. 4: Y-Shaped Head**

Suitable for neck, spine, Achilles tendon, etc

## Operating

1. Press the power button and hold for three seconds to turn on with level one. The blue lights indicate the speed level.
2. Adjust vibration speed via the power button. Maximum speed is level 4, minimum is level 1.
3. Press the power button and hold for three seconds will turn off the device.
4. Fully insert the applicator by pressing In. Remove It by pulling out,

**Notes:** The device automatically turns off when used for 15 minutes. Press the power button again to continue using it.

## Charging Instruction

- For first charging, please make sure the machine is fully recharged, which is good for the battery
- If the light of the power button is flashing, it means that the machine is out of power. Please stop using it and charge it.
- To recharge, connect the charging adapter to the device's charging port.
- When it is recharging, the indicator of the power button is red and keeps flashing: and when It Is fully recharged, the Indicator light changes from red to blue. It can't be used during the charging.
- The average running time is 4+ hours depending on the speed level and pressure applied during use.

## Maintenance, Gleaning &Storage

- Use a slightly damp towel to clean the device and dry it with a soft cloth.
- For storage or travel, Turn off the bottom power button and store it in the supplied portable box.

## Package Included

- 1 x Massage Gun
- 1 x USB-C Charging Cable
- 1 x Operating Instruction
- 4 x Massage Heads

## Specifications

Model	CM- 052G
Name	Mini Percussion Muscle Massager
AC/DC Adapter	100-240V 50/60Hz, DC Output: 5V 2A
Battery Capacity	2000Mah
Working Hour	4-6H
Charging Time	1500-2500 Rpm
No-load Speed	1100-2500 Rpm
Speed Levels	4

## Contact Us

We provide 2-year warranty and lifelong after-sale service for all products in or store, if you have any problems while using, feel free to contact us via email.

- Email: [service@cincomhealth.com](mailto:service@cincomhealth.com).

If you want to follow up our new products and latest promotional information, please visit our official website and Facebook account.



## Cincom

- [www.cincomhealth.com](http://www.cincomhealth.com)
- Email: [service@cincomhealth.com](mailto:service@cincomhealth.com)

## SHENZHEN CINCOM E-COMMERCE CO.,LTD

- **Add:** B806,11 lh Building , Gongyuandadi Community, Longx.Liang Road, Longcheng Street, Longgang District, Shenzhen, China
- **Web:** [www.cincomhealth.com](http://www.cincomhealth.com)
- **E-mail:** [service@cincomhealth.com](mailto:service@cincomhealth.com)

Made in China

## FREQUENTLY ASKED QUESTIONS

How does the CINCOM CM-052G Percussion Massage Gun provide deep muscle relief?

The massage gun offers 1500-2500 rpm percussion massage to reach deeper tissues, stimulate myofascial release, and reduce tension in muscles for quick pain relief and improved range of motion.

How many vibration speeds are available on the CINCOM CM-052G Percussion Massage Gun?

There are four vibration speeds available for a customizable massage experience.

What massage heads are included with the CINCOM CM-052G Massage Gun, and how do they enhance the massage experience?

Four massage heads are included, offering various massage experiences and allowing precise targeting of muscle groups, such as neck, shoulders, back, arms, and legs.

How portable is the CINCOM CM-052G Percussion Massage Gun, and what is its weight?

The massage gun is mini-sized, weighing only 1 lb, making it portable and easy to carry in a pocket for use at home, office, gym, or while traveling.

What is the power source of the CINCOM CM-052G Percussion Massage Gun?

The massage gun is battery-powered for convenience.

How is the CINCOM CM-052G Massage Gun recharged, and what is the battery life?

It features a modern USB C port for convenient and fast charging. A single charge lasts 4-6 hours.

How is the CINCOM CM-052G Massage Gun easy, quiet, and safe to use?

It has a single-button operation, with a press-and-hold to turn on, one short press to change speed, and automatic shutoff after 15 minutes for safety. The motor runs quietly, and the massage heads have rubber collars for reliability.

What is the sleek design of the CINCOM CM-052G Percussion Massage Gun like?

The mini massage gun has a sleek design, built with durable and skin-friendly aluminum material, providing a smooth and soft grip.

What is included in the packaging of the CINCOM CM-052G Massage Gun?

All accessories are packed in a gift box, including a charging cable.

What is the material used in the construction of the CINCOM CM-052G Massage Gun?

The massage gun is made of aluminum.

What are the dimensions of the CINCOM CM-052G Massage Gun?

The product dimensions are 7.87 x 7.42 x 2.36 inches.

What is the item weight of the CINCOM CM-052G Massage Gun?

The item weight is 16 ounces.

What is the model number of the CINCOM CM-052G Massage Gun?

The model number is CM-052G.

How long does a single charge of the CINCOM CM-052G Massage Gun last?

A single charge lasts 4-6 hours.

What is the safety feature regarding the massage heads of the CINCOM CM-052G Massage Gun?

The massage heads have rubber collars that grip the gun with no chance of coming loose, providing more reliability.

## VIDEO – PRODUCT OVERVIEW



[Download the PDF Link: CINCOM CM-052G Percussion-Massage-Gun Operating Instruction](#)  
[Percussion-Massage-Gun-Operating-Instruction.mp4](#)

## References

- [User Manual](#)