



COMPEX CX201FB01 Wireless Compression Boots Quick Start Guide

[Home](#) » [Complex](#) » COMPEX CX201FB01 Wireless Compression Boots Quick Start Guide 

Contents

- 1 COMPEX CX201FB01 Wireless Compression Boots
- 2 DESCRIPTION
- 3 SPECIFICATIONS
- 4 WHAT'S IN THE BOX
- 5 PRODUCT OVERVIEW
- 6 FEATURES
- 7 DIMENSIONS
- 8 INTENDED USE
- 9 CONTRAINDICATIONS
- 10 PRECAUTIONS
- 11 INSTRUCTIONS
- 12 FREQUENTLY ASKED QUESTIONS
- 13 VIDEO – PRODUCT OVERVIEW
- 14 References
- 15 Related Posts



COMPEX CX201FB01 Wireless Compression Boots



DESCRIPTION

Feel relief aid anywhere, anytime with Compex Ayre compression boots. Designed for elite and recreational athletes on the go, trainers, and coaches. Compex Ayre uses a dynamic pulsing compression up the limbs in waves for the duration of its use, which helps promote quick healing.

- **USE TO HELP TREAT / PREVENT:** Sore or aching muscles, pain and swelling of the muscles.
- **BEST FOR:** Aiding in the relief of sore or aching muscles from exercise or recreational activities helping you to be ready to go at full speed the next day.

SPECIFICATIONS

- **Brand:** Compex
- **Color:** Black
- **Special Feature:** Rechargeable
- **Product Dimensions:** 32 x 29 x 32 inches
- **Item model number:** CX201FB01
- **Power Source:** Battery Powered
- **Item Weight:** 3.49 Kilograms

WHAT'S IN THE BOX

- Wireless Compression Boots
- Quick Start Guide

PRODUCT OVERVIEW



Ayre™ Boots

- A 420D Nylon Boots
- B Hook and Loop
- C Zipper
- D C1, C2, C3, C4 (Chambers)



Ayre™ Control Unit


1. Power On/Off button
2. Lithium Ion Battery
3. Global Power Plug
4. SET button



FEATURES

- Ultra-portable, no tethers to a power outlet during compression
- 4 Separate compression chambers can be inflated separately or simultaneously
- 2 modes: Inflate the four chambers one by one or all four chambers at the same time
- Lithium-ion rechargeable battery for up to 3 hours of recovery
- High-quality exterior and interior linings make cleaning easy and help keep unwanted odors in check
- Integrated compression controls remove the bulk of a separate, wired unit
- Powerful compression: The pressure can be adjusted by increments of 10 from 0mmHg up to a max of 120mmHg

DIMENSIONS

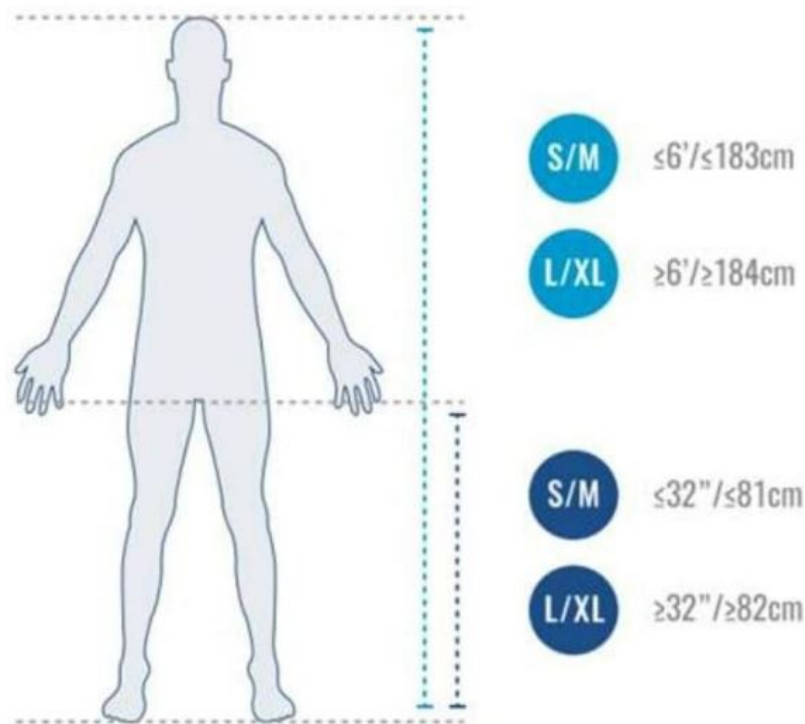

INTENDED USE

The Compex® Ayre™ is used to aid:

- Sore or aching muscles from exercise or recreational activities help you to be ready to go at full speed the next day.
- Sore or aching muscles, pain and swelling of the muscles, vein problems, diabetics blood flow.

This product is a non-medical device and is not used in medical conditions or medical environments.

CONTRAINDICATIONS

This device must not be used for the following conditions:

- Persons with suspected, active, or untreated: deep vein thrombosis, ischemic vascular disease, severe arteriosclerosis, pulmonary edema, severe congestive heart failure, thrombophlebitis or an active infection.

- On the legs cuffs would interfere with the following conditions: vein ligation, gangrene, dermatitis, open wounds, a recent skin graft, massive edema or extreme deformity of the leg.
- Not for use on any neuropathy. Do not use on extremities that are insensitive to pain.
- Not for use where increased venous or lymphatic return is undesirable.

WARNINGS

- Do not attempt to repair the device.
- Do not attempt to open or remove covers.
- Do not remove the pump unit from the cuff.
- Do not attempt to modify or change the device.
- NEVER attempt any service while the device is in use.
- Do not operate in a wet environment.
- Do not immerse in any liquid for any reason. For cleaning and disinfecting instructions refer to the “Clear and Disinfecting section”.
- Do not place the device in an autoclave for any reason.
- Not suitable for use in the presence of a flammable anesthetic mixture with air or with oxygen or nitrous oxide.
- If exposed to temperatures below 100(500F) allow the device to warm up to room temperature.
- Do not subject the device to extreme shocks, such as dropping the pump.
- Portable and mobile Radio Frequency Communication Equipment can be affected by this device.

PRECAUTIONS

- Stop using the device if swelling, skin irritation or any other unpleasant or painful sensation occurs, and consult a physician.
- If pulsation or throbbing occurs, loosen the cuffs immediately as the cuffs may be wrapped too tightly.
- For people with diabetes or vascular disease that require frequent skin assessment; Consult a physician.
- Consult a physician for people who use warming devices in combination with cuffs requiring frequent assessment as skin irritation may occur.
- For people positioned in the supine lithotomy position (with or without cuffs) for an extended period and require special attention to avoid extremity compartment syndrome; consult; a physician.

INSTRUCTIONS

READ THE COMPLEX MANUAL THOROUGHLY (available at compex.com)

- Put the boots on your legs
- Secure them tightly with the hook loop and zipper
- Hold the On/Off button for 1 second to turn on the device.
- The device will begin to operate in the default function/mode F1 which will inflate each of the four chambers to the default pressure setting. To change the mode and pressure settings, please read the full User Manual.
- Hold the On/Off button for 1 second to turn off the device.

CHARGING THE CONTROL UNIT

Make certain the charge pin is aligned with the global power plug charging receiver. The battery icon will flash during charging and will hold when fully charged. The device may be used while it is charging.

PRODUCT SUPPORT

- **Phone:** 1-877-266-7398 (877-COMPEX8)
- **Email:** service@compex.com

COPYRIGHT © 2020 by DJO, LLC Individual results may vary. Neither DJO, LLC nor any of its subsidiaries dispense medical advice. The contents of this document do not constitute medical advice. Rather, please consult your healthcare professional for information on the courses of treatment, if any, which may be appropriate for you.

DJO, LLC A DJO Global Company 1430 Decision Street Vista, CA 92081-8553 U.S.A. CX202IF02 [COMPEX.COM](https://www.compex.com)

FREQUENTLY ASKED QUESTIONS

What is the primary purpose of the COMPEX CX201FB01 Wireless Compression Boots?

The boots are used to aid sore or aching muscles from exercise or recreational activities, helping to prepare the user for the next day.

Is the COMPEX CX201FB01 a medical device?

No, it is mentioned as a non-medical device and is not used in medical conditions or medical environments.

What are the contraindications for using the COMPEX CX201FB01?

The device must not be used for conditions such as deep vein thrombosis, ischemic vascular disease, severe arteriosclerosis, pulmonary edema, severe congestive heart failure, thrombophlebitis, or active infection.

Where should the device not be used on the legs according to contraindications?

The device should not be used where cuffs would interfere with conditions like vein ligation, gangrene, dermatitis, open wounds, recent skin graft, massive edema, or extreme deformity of the leg.

What is the warning related to attempting to repair or modify the device?

Users should not attempt to repair, open, remove covers, modify, or change the device, and any service attempts should not be made while the device is in use.

Can the COMPEX CX201FB01 be used in a wet environment?

No, it should not be operated in a wet environment, and immersion in any liquid is not allowed.

What is the precaution related to using the device in the presence of flammable substances?

The device is not suitable for use in the presence of a flammable anesthetic mixture with air or with oxygen or nitrous oxide.

What should users do if skin irritation or unpleasant sensations occur?

Users should stop using the device and consult a physician if swelling, skin irritation, or any other unpleasant or painful sensation occurs.

How is the charging process described for the control unit of the COMPEX CX201FB01?

Users should align the charge pin with the global power plug charging receiver, and the battery icon will flash during charging, holding when fully charged. The device can be used while charging.

What is the power source for the COMPEX CX201FB01?

The power source is battery-powered.

What is the weight of the COMPEX CX201FB01?

The weight is 3.49 kilograms.

What special feature does the COMPEX CX201FB01 have?

The special feature is rechargeable.

What are the product dimensions of the COMPEX CX201FB01?

The product dimensions are 32 x 29 x 32 inches.

What is the item model number of the COMPEX CX201FB01?

The item model number is CX201FB01.

VIDEO – PRODUCT OVERVIEW



[Download the PDF link: COMPEX CX201FB01 COMPEX CX201FB01 Wireless Compression Boots Quick Start Guide](#)
[Wireless-Compression-Boots-Quick-Start-Guide.mp4](#)

REFERENCE: [COMPEX CX201FB01 Wireless Compression Boots Quick Start Guide-Device.Report](#)

References

- [User Manual](#)

[Manuals+](#), [Privacy Policy](#)