

# DAYTECH LC05 Caregiver Pager System User Manual

[Home](#) » [DAYTECH](#) » DAYTECH LC05 Caregiver Pager System User Manual 


## DAYTECH LC05 Caregiver Pager System



### Contents

- [1 Product Overview](#)
- [2 Instructions](#)
- [3 Documents / Resources](#)
- [4 Related Posts](#)

## Product Overview

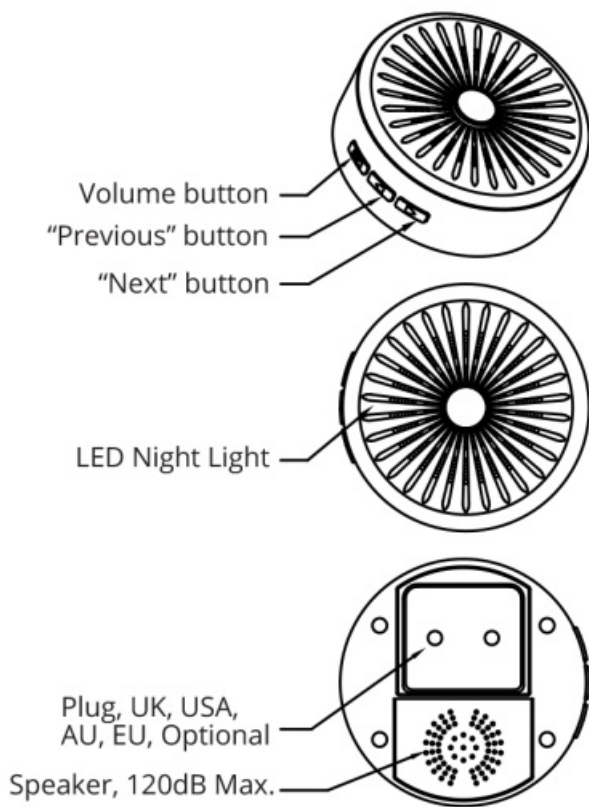

## Instructions

### 1. Volume adjustment

Each time you press the volume button on the receiver, the volume will increase by one level and the LED will flash. After the volume is adjusted to the maximum volume, press it again. The next level is the minimum volume, that is the 0dB mode, with only the LED light flashing, and a total of 5 levels loop.

### 2. Change the ringtone/learn to match code

1. Press the ringtone button (previous/next) on the receiver to tune to the user's favorite ringtone
2. Press and hold the volume button on the receiver for about 5 seconds, until the receiver's LED light flashes and a "ding" beep sounds, and it enters the learning mode (The learning mode time is only 8 seconds, after which it will automatically exit)
3. In the learning mode, quickly press the button on the transmitter, the receiver will emit two beeps of "ding-ding", and the setting is successful
4. Press the button on the transmitter again, and the ringing ringtone of the receiver is the ringtone set by the user
5. This method is suitable for changing ringtones and adding transmitters at the same time; multiple transmitters can correspond to the same ringtone or different ringtones, set to the same Ringtone Please press multiple transmitter buttons at the same time in the learning mode to set a different ringtone. After completing the pairing of the first transmitter, adjust it to the different ringtones, complete the pairing of other transmitters one by one
6. If the code matching operation is not completed within the learning time, please refer to the above steps and perform the operation again.

### 3. Code clearing/reset


Long press the receiver's next ringtone button until the receiver's LED light flashes and a "ding" beep sounds to complete the code clearing; The code between the transmitter and the receiver will be cleared. Press the

button of any transmitter again and it will automatically pair with the receiver. Other transmitters need the user adds the pairing by himself; the ringtone will be restored to the default “dingdong”.

#### 4. Lighting Mode

Press and hold the middle button for about 5 seconds. The LED is always on and enter the lighting mode; press any button on the receiver again and the night light flashes for 3 seconds to turn off the LED light, or press the transmitter button and the night light flashes for 3 seconds. Turn off the LED light.

## Documents / Resources

	<a href="#">DAYTECH LC05 Caregiver Pager System</a> [pdf] User Manual LC05 Caregiver Pager System, LC05, Caregiver Pager System, Pager System
---	--

[Manuals+.](#)