



ENCHEN Mint 3 Portable Water Flosser User Manual

[Home](#) » [ENCHEN](#) » ENCHEN Mint 3 Portable Water Flosser User Manual 

ENCHEN

**User Manual
for Portable Water Flosser Mint 3**




Please carefully read this User Manual before using the product.

Contents


- [1 Product Introduction](#)
- [2 Use Method](#)
- [3 Important Notes](#)
- [4 Care and Maintenance](#)
- [5 Warranty Instructions](#)
- [6 Instructions of Product Fault](#)
- [7 Basic Parameters of the Product](#)
- [8 Documents / Resources](#)
- [9 Related Posts](#)

Product Introduction

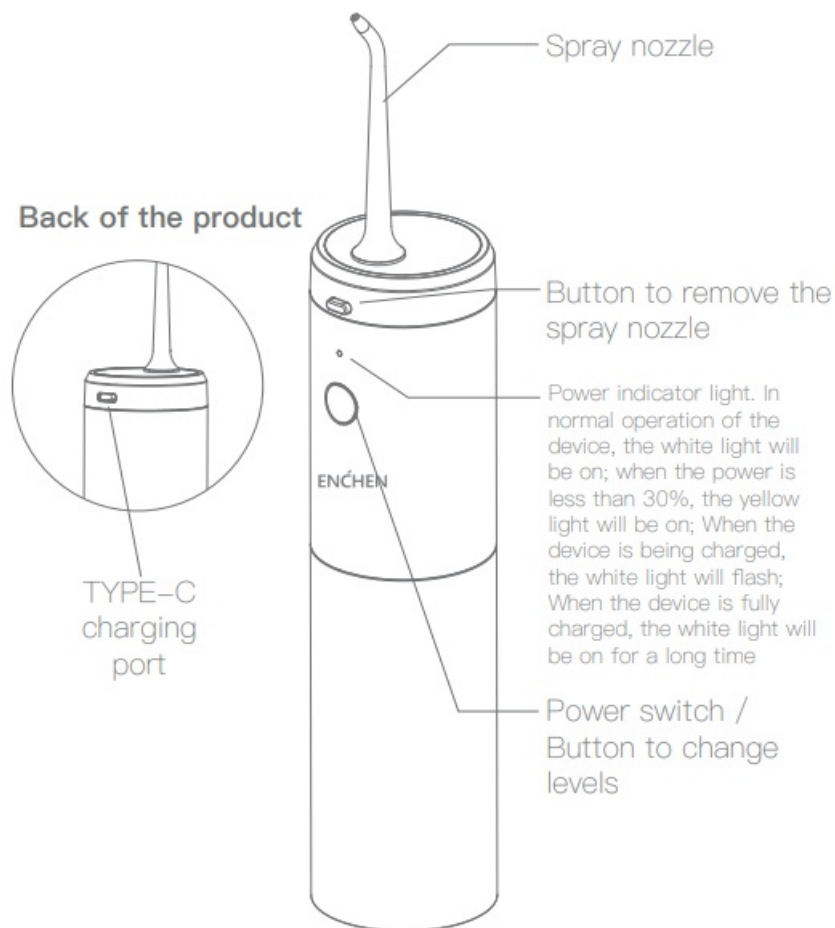


Use Method

1. Insert the spray nozzle into the nozzle jack and insert it firmly until it makes a clicking sound. Turn the water tank counterclockwise to remove it and fill it with water. The running tap water is recommended. If the water temperature is too low, use warm water of less than 40 ° C. Please do NOT add salt, toothpaste or chemicals into the tank. After the water is filled into the tank, turn the water tank clockwise into the device body and do NOT tilt or invert the tank during the process to avoid the water in the tank being poured out.

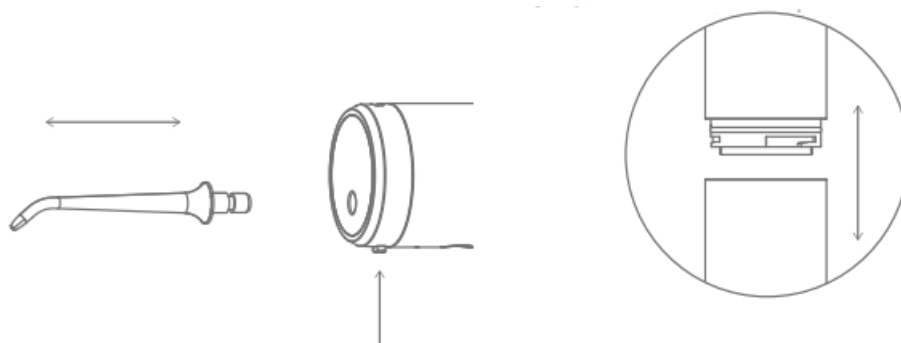

2. When using the product, lean forward slightly, open your mouth slightly, and let the water flow out from the mouth naturally or spit it out. Try your best to NOT stand upright nor sit to rinse your teeth (to avoid water splashing on your body or face). In order to achieve the best use effect, please start cleaning from the molars (back teeth), make cleaning at the direction toward the front teeth, clean teeth one by one along the gingival line, and pay attention to cleaning the tooth pit and fissure, wisdom teeth and other hidden parts. As the spray nozzle moves along the gingival line, pause briefly between two teeth to make sure the gum area and the part between those two teeth are handled to be clean.

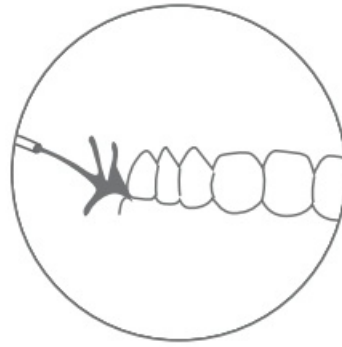

3. Quickly press the power switch to switch among the first level, the second level, the third level and off.

The first level (soft level):

Press the power switch, the device operates at soft level by default, and the operation indicator light is white.

The second level (standard level):

Within 3 seconds of its operation at the first level, press the power switch again, the device switches to this level, and the operation indicator light is white.

The third level (strong level):


Within 3 seconds of its operation at the second level, press the power switch again, the device switches to this level, and the operation indicator light is white.

Notes:

1. After 3 seconds of its operation at any level, press the power switch, and then the device stops operating.
2. In order to better maintain the service life of battery, the product built stops the protection function after 3 minutes of continuous operation.



4. Let the product operate for 2 or 3 seconds to drain the water after using it. Open the dustproof plug, firmly insert one end of the charging line into the charging port of the product, firmly insert the other end into the port of power adapter, and firmly insert the power adapter into the power socket.

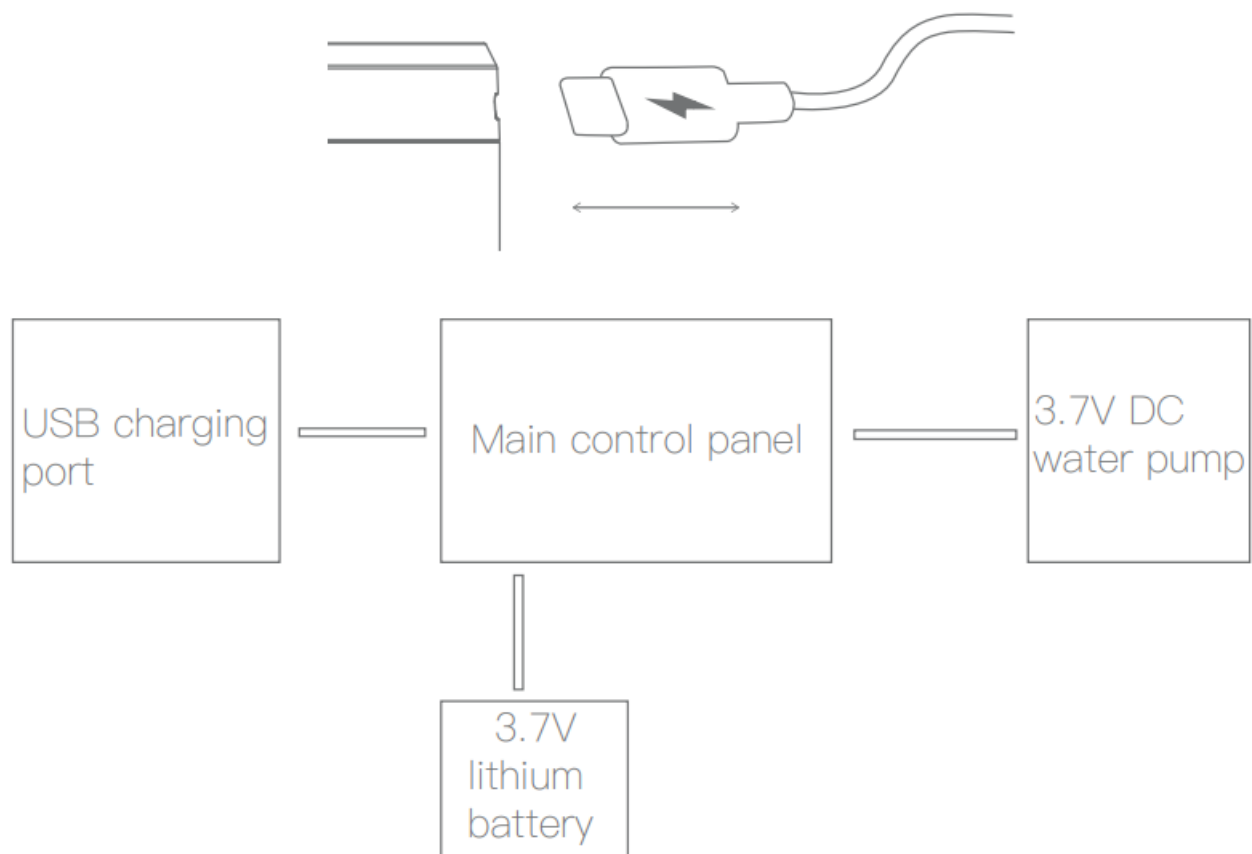

Important Notes

In order to reduce failure, electric shock, injury, fire, death, or damage to equipment or property, please follow the following safety precautions.

1. Do NOT disassemble the product or change the parts of the product without permission, so as not to damage the product or cause electric shock. The warranty will be invalid after unpermitted disassembly.
2. Do NOT put the water of more than 40 ° C into the container to clean the mouth, so as to avoid scalding.
3. If the product falls into the water when it is being charged, please unplug the power adapter from the power socket first, and then pull out the product to avoid the danger of electric shock
4. Do NOT put the product into fire or heat it, do NOT charge it at a place with high temperature or humidity, do NOT place the product at that place for a long time, and avoid direct sunlight. Do NOT pull out the power plug and charging plug with wet hands during charging.
5. The product is not equipped with a power adapter. When charging the product, please be sure to use a 5V 1A DC output adapter that meets the national standard.
Before using the power adapter, please check whether the voltage labeled meets local voltage requirements.
6. If the product is not used for a long time, please unplug the power plug and charging plug to avoid fire or other dangers.
7. Keep the product clean, please use running water or warm water to clean the mouth, and let the product operate for 2 or 3 seconds to drain the water after using it.
8. There is no replaceable or repairable part. Please return the product to the manufacturer, if any repair is required.
9. Please do NOT use the product, if it is not installed with a spray nozzle.
10. The product is only suitable for people over 6 years old. Children should be guided or accompanied by adults.
11. Please charge the product at room temperature.

12. When the product is charged at low temperature (below 15 ° C), its service time may be shortened; When the temperature is too low (below 0 ° C), charging may not be possible.
13. Use running tap water or warm water, and press the switch button after the spray nozzle is close to the tooth.
14. The built-in rechargeable battery may contain pollutants. When you discard the dental cleaner, please must remove the battery and handle it at an officially designated recycling center, or send the dental cleaner to an officially designated recycling center.
15. If the battery cannot be removed, you can send the dental cleaner to the maintenance center of Lvzhipin, the maintenance personnel will help you to remove the battery and handle it in environmentally-friendly way.
16. Do NOT remove the rechargeable battery of the dental cleaner at will. Only when discarding the dental cleaner can you remove the rechargeable battery.

Name and Content of Harmful Substances in the Product

Name of componen	Harmful substances or elements					
	Plumbum Pb	Mercury Hg	Cadmium Cd	Valentine chromium Cr(VI)	Polybrominated biphenyls PBB	Polybrominated diphenyl ethers PBDE
Control assembly	X	0	0	0	0	0
Power pack	X	0	0	0	0	0
pails p	0	0	0	0	0	0
Pipe component	0	0		0	0	0
Switch	0	0	0	0	0	0
Water pump	0	0	0	0	0	0
Transmission gear	0	0	0	0	0	0
Spray nozzle	0	0	0	0	0	0
Charging line	0	0	0	0	0	0
Screw	0	0	0	0	0	0

This form is prepared in accordance with regulations of SJ/T 11364.

○ indicates that the content of the harmful substance in all homogeneous materials of the component is below the limit specified in GB/T26572.

× indicates that the content of the harmful substance in at least one of the homogeneous materials of the component is above the limit specified in GB/T26572.

(The product conforms to requirements of RoHS environmental protection. At present, there is no mature technology to replace or reduce the plumbum content in electronic ceramics, optical glass, steel and copper alloy in the world.)

The number in the label represents that the normal environment-friendly use period of the product is 10 years in normal usage.

Care and Maintenance

Use clear water or neutral detergents to clean the product. Do NOT use corrosive or abrasive detergents,

because those will damage the product. And do NOT wash the product with hot water exceeding 50°C.
Do NOT immerse the product in water to clean it.
After each use, remove the spray nozzle and then clean the product.
The whole device can be rinsed under running water and then be stored.
Scrub the outer surface of the product with a soft and damp cloth to remove dirt, and wipe the relevant parts dry with a clean and dry cloth after each use. In addition, dental experts recommend changing the spray nozzle once a month.

Warranty Instructions

The three guarantees of product after-sales service are strictly in accordance with the Law of the People’s Republic of China on the Protection of the Rights and Interests of Consumers and the Law of the People’s Republic of China on Product Quality, and the service contents are as follows:

Warranty Policy

During the period of validity for three guarantees, you can enjoy return, replacement and repair services for free in accordance with this regulation.

Non-warranty Policy

Damage caused by failure to use, maintain or keep the product according to the instructions.
Damage caused by improper operation (such as cleaning the product in water), abnormal power supply, transportation and other accidents.
Damage caused by maintenances or modifications not approved by our company
Damage caused by other force majeure factors (such as lightning, earthquake, fire, flood, etc.).
Damage not caused by domestic or personal use (such as industrial or commercial use).
Aging, wear and so on caused by abnormal use of product.

Warranty Card of the Product

Certificate of Quality
Inspected

Product Name	Water Flosser	Warranty Period
Product Model		Within one year since the date of purchase
Date of Purchase	YY MM DD	Till YY MM DD
Customer’s Address		

Instructions of Product Fault

Fault	Cause Analysis	Solution
Not work	Just purchase the product	Connect it with power line to charge it for 3 hours till it is fully charged
	leave it unused for a few months	Consult customer service
Only work for a few minutes after each charge	Power is too low	Connect it with power line to charge it for 3 hours till it is fully charged
	Battery is at the end of service period	Consult customer service
Low water pressure	Dirt is stuck in the spray nozzle	Replace the spray nozzle
	There is dirt in the water pump	Consult customer service
No water comes out	The water in the tank is used up	Fill the tank with water
	There is dirt in the water pump.	Consult customer service

※If you have any other questions, please contact the customer service center

Basic Parameters of the Product

Product Name: Water Flosser

Product Model: Mint 3


Rated Voltage: 3.7V 

Rated Power: 4W



Shanghai Qiaopusi Smart Sci-Tech Development Co., Ltd.
Add: Room 103, 3th Bld, No.320 Donglan Rd, Minghang District,
Shanghai, China, 201102

Documents / Resources

	<p>ENCHEN Mint 3 Portable Water Flosser [pdf] User Manual Mint 3 Portable Water Flosser, Mint 3, Portable Water Flosser, Water Flosser, Flosser</p>
---	---