



DA FoCK-HPns EyePulse EMS Device User Guide

[Home](#) » [EyePulse](#) » DA FoCK-HPns EyePulse EMS Device User Guide 

DA FoCK-HPns EyePulse EMS Device



Contents

- 1 What is EMS or NMES technology?
- 2 What are the benefits of EMS or NMES?
- 3 DO NOT use EMS and our devices if you have:
- 4 How it is composed:
- 5 How to Use Our Device:
- 6 Command Center Specifications:
- 7 Technical Parameter
- 8 Latest Considerations
- 9 Documents / Resources
- 10 Related Posts

What is EMS or NMES technology?

Electrical Muscle Stimulation (EMS), also known as Neuro Muscular Electrical Stimulation (NMES), is the use of gentle electrical pulses to induce multiple muscle contractions.

Low-current electrical impulses stimulate the nerves first...

which then transmit the message to the surrounding muscle fibers, causing them to expand and contract.

This expansion and contraction is exactly what your body does when you exercise. When you use an EMS device on your face, it is as if each

muscle group is doing little crunches.

The result is a toned appearance with tighter, smoother skin.

One of our favorite aspects is that the benefits are immediate, with no recovery time.

What are the benefits of EMS or NMES?



Toning Facial Muscles



Firming Facial Skin



Stimulate collagen production



Improve Blood Circulation



Reduce Dark Circles



Deflating Eye Bags

DO NOT use EMS and our devices if you have:



If you're pregnant



Cancer Pathologies



Other Serious Medical Condition (check with your doctor)

How it is composed:

Our Device is composed of 3 basic Element similar to the functioning of our body:

1. The Command Center – the pulsating **“Heart”**
2. The Cables – the **“Blood Vessels”** of the tool
3. The StickyPods – the **“Muscles”** of the device



How to Use Our Device:

IMPORTANT: Before Using the device memorize the Steps and READ them 3 times.



Perform these 7 Steps perfectly to HAVE a Safe and Effective treatment.



The Command Center produces a small electric current, which passes through the Cables to the Sticky Pads.

The micro-current is transmitted to the skin through StickyPads.

Here are the steps to take:


1. Charge the Command Center for 4 hours the first time.

2. Then connect the cables, and to the cables connect 4 Sticky Pods.
3. Apply all 4 StickyPads to your face and lie down on the couch or bed in your home.
4. Turn the device on, and put mode 1, the lowest.
5. If you do not feel contractions turn the mode up every 30 seconds.
6. Hold the device for 15 minutes and then turn the device off.
7. Remove the StickyPads, apply Electro Gel to your face, and close the StickyPods in the sachet with the Protective Film.
8. Recharge the device and prepare it for the next treatment.

Command Center Specifications:

It has 8 Modes and 19 Levels here's how the buttons work:

1. Press + **button** to turn on.
2. Press **Button M** to change modes.
3. Press 2 times + **button** to increase the level.
4. Press 2 times – **button** to decrease the level.
5. Press – **button** to turn off.

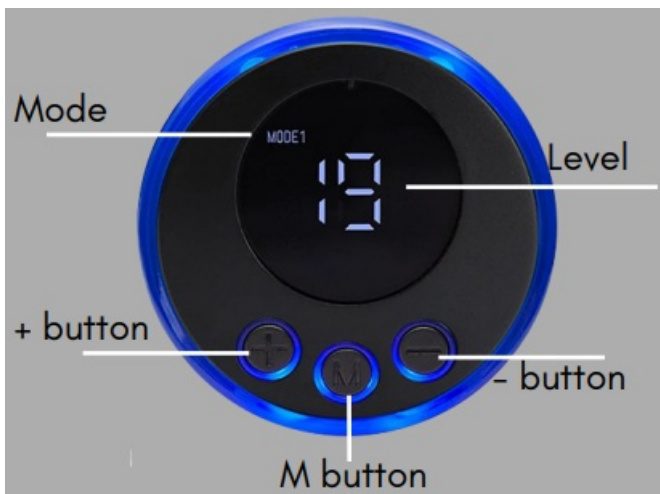

1. Can be used for about 160 minutes after one hour of charging. The higher the power, the higher the power consumption.
2. Cannot be used while charging, can only be used after charging.
3. The first time charge it for 4 HOURS.
4. After accessing it, it automatically turns off after 15 minutes.
5. Please set a lower intensity when using it for the first time, and adjust it after adapting to the low intensity.

Technical Parameter

Voltage:	3.7 V
Frequency:	1-100 Hz
Size	62x62x16 cm
Weight	31 g
Outout	35mA maximum
Battery capacity	200mAH
Working hours	15 minutes
Mode	8
Level	19

Issues that need attention

1. While charging, do not disconnect the host or turn it on, it may break.
2. When using 2 clips, they must be installed correctly, and adhere perfectly to the skin.

Latest Considerations

We thank you for relying on us.

Please note that Sticky Pads can be used for up to 6 treatments.

That' s why we offer you a 20% Discount on your next Sticky Pads order.


Just go here: paerly.uk/stickypads

and use the code **STICKY20**.

You can purchase packs of 6,12 and 24 pairs to continue using EyePulse.

A Greeting

Documents / Resources

	<p>EyePulse DAFoCK-HPns EyePulse EMS Device [pdf] User Guide</p> <p>DAFoCK-HPns EyePulse EMS Device, DAFoCK-HPns, EyePulse EMS Device, EMS Device, D evice</p>
-------------------------------------------------------------------------------------	------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------