

[Skip to content](#)

[Manuals+](#)

User Manuals Simplified.



HOMEDICS SP-180J-EU2 Cordless Double-Barrel Rechargeable Body Massager User Manual

[Home](#) » [Homedics](#) » HOMEDICS SP-180J-EU2
Cordless Double-Barrel Rechargeable Body Massager
User Manual

HOMEDICS®

LEADERS IN BEAUTY & WELLBEING



CORDLESS

DOUBLE-BARREL
RECHARGEABLE BODY MASSAGER
3 YEAR GUARANTEE
SP-180J-EU2

Contents [hide](#)

[1 PRODUCT FEATURES](#)

[2 INSTRUCTIONS FOR USE:](#)

[3 Cleaning](#)

[4 Storage](#)

[5 Documents / Resources](#)

[6 Related Posts](#)

PRODUCT FEATURES

1. Power – massage on/off
2. Charging port
3. Detachable straps



Fig. 1



Fig. 2



Fig. 3



Fig. 4



Fig. 5

INSTRUCTIONS FOR USE:

1. Your unit should arrive with a full charge. When you need to charge the Cleaning unit, plug the adapter into the jack on the unit, and plug the other end into a 100-240V mains outlet. The power button will illuminate red when charging and will change to green once fully charged. The unit should be charged after 5 hours of charge time. A full charge will last up to 2 hours.
Unplug the appliance and allow to cool before cleaning. Clean only with a soft, slightly damp sponge. Never allow water or any other liquids to come into contact with the appliance. Do not immerse in any liquid to clean.
Never use abrasive cleaners, brushes, glass/furniture polish, paint thinner, etc to clean.
Note: The product should only be charged using the adaptor supplied (SAW06C-050-1000GB). The adaptor must be removed from the socket when not in use. Adaptor output charging voltage of 5Vdc and 1A must not be exceeded. Storage Place the appliance in its bag or in a safe, dry, cool place. Avoid contact with sharp edges or pointed objects which might cut or puncture the fabric surface. To avoid breakage, DO NOT wrap the power cord around the appliance. DO NOT hang the unit by the cord.
2. This massager is versatile and can be used for the neck, shoulder, back, legs, arms, and feet (Fig. 1-3). For use on the neck, shoulders, or back, attach the adjustable straps to the ends of the unit (Fig. 4), and use the straps to hold the massager in the desired area. The double-barrel design allows the massager to roll up and down your body relaxing muscles and alleviating pain as it rolls.
3. To activate the massage action, short press the power button (Fig. 5), and the vibration waves will start on the lowest setting. Press 2 seconds to medium intensity and press 2 seconds again to experience the highest intensity. To turn the unit off, short press and turn it off.
4. When not in use, store your Cordless Double-Barrel Massager in its convenient drawstring storage bag.

Cleaning

Unplug the appliance and allow it to cool before cleaning. Clean only with a soft, slightly damp sponge. Never allow water or any other liquids to come into contact with the appliance. Do not immerse in any liquid to clean. Never use abrasive cleaners, brushes, glass/furniture polish, paint thinner, etc to clean.

Storage

Place the appliance in its bag or in a safe, dry, cool place. Avoid contact with sharp edges or pointed objects which might cut or puncture the fabric surface. To avoid breakage, DO NOT wrap the power cord around the appliance. DO NOT hang the unit by the cord.



FKA Brands Ltd

Manufacturer & UK Importer: FKA Brands Ltd, Somerhill Business Park,
Tonbridge, Kent TN11 0GP, UK

EU Importer: FKA Brands Ltd, 29 Earlsfort Terrace, Dublin 2, Ireland

Customer Support: +44(0) 1732 378557 |


support@homedics.co.uk

IB-SP180JEU2-0521-02



Register your product today at www.homedics.co.uk/product-registration

Documents / Resources

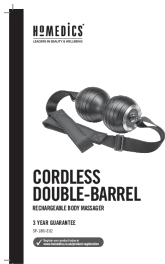

[HOMEDICS SP-180J-EU2 Cordless Double-Barrel Rechargeable Body Massager](#) [pdf] User Manual

SP-180J-EU2 Cordless Double-Barrel Rechargeable Body Massager, SP-180J-EU2, Cordless Double-Barrel Rechargeable Body Massager, Double-Barrel Rechargeable Body Massager, Rechargeable Body Massager, Body Massager, Massager

[Manuals+](#),

- [home](#)
- [privacy](#)