

**Hyperice  
Go 2 Percussion  
Massage Device**



# Hyperice Go 2 Percussion Massage Device Instruction Manual

[Home](#) » [Hyperice](#) » Hyperice Go 2 Percussion Massage Device Instruction Manual 

## Contents

- [1 Hyperice Go 2 Percussion Massage Device](#)
- [2 IMPORTANT SAFETY INSTRUCTIONS – ORIGINAL INSTRUCTIONS](#)
- [3 Operating](#)
- [4 Cleaning and storage](#)
- [5 Specifications](#)
- [6 Warranty](#)
- [7 Documents / Resources](#)
  - [7.1 References](#)

# Hyperice

**Hyperice Go 2 Percussion Massage Device**



TO REDUCE RISKS OF ELECTRIC SHOCK, FIRE, AND PERSONAL INJURY, OR PROPERTY DAMAGE, THIS DEVICE MUST BE USED IN ACCORDANCE WITH THE FOLLOWING WARNINGS, CAUTIONS, AND SAFETY INSTRUCTIONS.

## **IMPORTANT SAFETY INSTRUCTIONS – ORIGINAL INSTRUCTIONS**

Read all instructions before using the Hypervolt Go 2.

### **DANGER**

To reduce the risk of electric shock:

- Always unplug this device from the electrical outlet immediately after using and before cleaning.
- Do not reach for a device that has fallen into the water. Unplug immediately.
- Do not use it while bathing or in a shower.
- Do not place or store the device where it can fall or be pulled into a tub or sink. Do not place in or drop into water or other liquid.

### **WARNING**

To reduce the risk of burns, fire, electric shock, or injury to persons:

- A device should never be left unattended when plugged in. Unplug from outlet when not in use, and before putting on or taking off parts.
- Do not operate under a blanket or pillow. Excessive heating can occur and cause fire, electric shock, or injury to persons.
- This device can be used by children aged 8 years and above and persons with reduced physical, sensory, or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved.

- Use this device only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.
- Never operate this device if it has a damaged cord or plug, if it is not working properly if it has been dropped or damaged, or if it has gotten wet. Return the device to a service center for examination and repair.
- Do not carry this device by supply cord or use cord as a handle.
- Keep the cord away from heated surfaces.
- Never operate the device with the air openings blocked. Keep the air openings free of lint, hair, or any substance which may impair the flow of air.
- Never drop or insert any object into any opening.
- Do not use outdoors.
- Do not operate where aerosol (spray) devices are being used or where oxygen is being administered.
- Never operate on a soft surface such as a bed or couch where the air openings may be blocked.
- Do not use the device in close proximity to loose clothing or jewelry.
- Keep long hair away from the device while in use.
- Cleaning and user maintenance shall not be done by children without supervision.
- Do not remove screws or attempt to disassemble.
- Do not operate continuously for more than one hour. After one hour of use, allow the device to rest for 30 minutes before re-using.
- Unplug the device after charging or prior to use.
- Use the device only according to the Operating Instructions.
- When using the Hypervolt Go 2 use it only on dry, clean surfaces of the body by lightly pressing and moving across the body for approximately sixty(60) seconds per region. The Hypervolt Go 2 should only be used on soft tissue areas of the body and should not be used on any hard or bony areas of the body, including the head. Please immediately stop use of the Hypervolt Go 2 upon experiencing any pain or discomfort, apart from light muscle soreness. Bruising can occur from use of the Hypervolt Go 2, regardless of the pressure setting, and if bruising does occur you should not use the Hypervolt Go 2 on the bruised area until the bruise is fully healed. Use of the Hypervolt Go 2 should also be avoided on all bruises, contusions, rashes or irritated or injured areas of the skin until they are fully healed. When using the Hypervolt Go 2 be sure to keep fingers, toes, hair, and other body parts away from the back of the attachment head to avoid pinching or hair entanglement.

The Hypervolt Go 2 is intended to be charged with the Hyperice USB-C cord provided or by a 5V USB. Do not charge the Hypervolt Go 2 overnight or leave the Hypervolt Go 2 unattended when charging or in use. Use only with Listed/Certified ITE Power Supply.

Please do not use the Hypervolt Go 2, or any percussion device without first obtaining approval from your physician if any of the following apply:

- Pregnancy, diabetes with complications such as neuropathy or retinal damage, wear of pacemakers, recent surgery or injury, epilepsy or migraines, herniated disks, spondylolisthesis, spondylolysis, or spondylosis, recent joint replacements or IUD, metal pins or plates or any concerns about your physical health. Frail individuals and children should be accompanied by an adult when using any percussion or vibration device. These contraindications do not mean that you are not able to use a percussion or vibration device but we advise that before doing so you consult your physician. Ongoing research is being done on the effects of percussion massage for specific medical disorders that may result in the shortening of the list of contraindications as shown above.

The Hypervolt Go 2 contains a battery that must be complied with safely at an appropriate e-waste disposal or recycling facility.

### **CAUTION**

The battery used in this device may present a risk of fire or chemical burn if mistreated. Do not disassemble, heat above 100° C, or incinerate. Dispose of the used battery promptly. Keep away from children. Do not disassemble and do not dispose of in fire.

### **SAVE THESE INSTRUCTIONS**

The Hypervolt Go 2 is a handheld percussion massage device that delivers targeted pulses of pressure to care for muscles, relieve tension, provide a relaxing massage, accelerate warmup and recovery, and help maintain flexibility and range of motion.

## **Operating**

Press down the power/speed button on the back of the device. Hold for two seconds until the battery level indicator (LED band) illuminates, then release. Your Hypervolt Go 2 is now ready for use. Choose your desired speed (1-3) by pressing the power/speed button once per speed. The white lights on the back will indicate the speed setting. To power off, simply hold the power/speed button down for two seconds.

### **Changing head attachments**

Only change attachments after the speed setting is turned off and the head attachment stops moving. To change head attachments, pull the existing head attachment straight out to remove and insert the desired head directly into the opening while pushing firmly.


### **Charging**

Fully charge the battery for up to four hours before first use. To charge, connect the USB-C cord to the charging port on the bottom of the handle and plug the cord into a USB-A adapter. Then plug the adapter into the wall. The LED light band will pulsate to show the battery level and indicate active charging. The LED light band colors correspond to the charge level ranging from red (low) to green (fully charged). Full charge is indicated when the LED band remains green and fully illuminated. The battery may be recharged at any time and at any battery level. It is not recommended to fully deplete the battery to the red LED level. The average runtime is 2+ hours depending on speed level and pressure applied during use. It is recommended to turn the device OFF while it is not in use or during storage and travel.

## **Cleaning and storage**

Make sure the power is OFF and the battery charger is not attached. Use a damp, clean cloth and gently wipe your Hypervolt Go 2. Store in a clean, cool place, out of direct sunlight, when not in use.

## **Specifications**

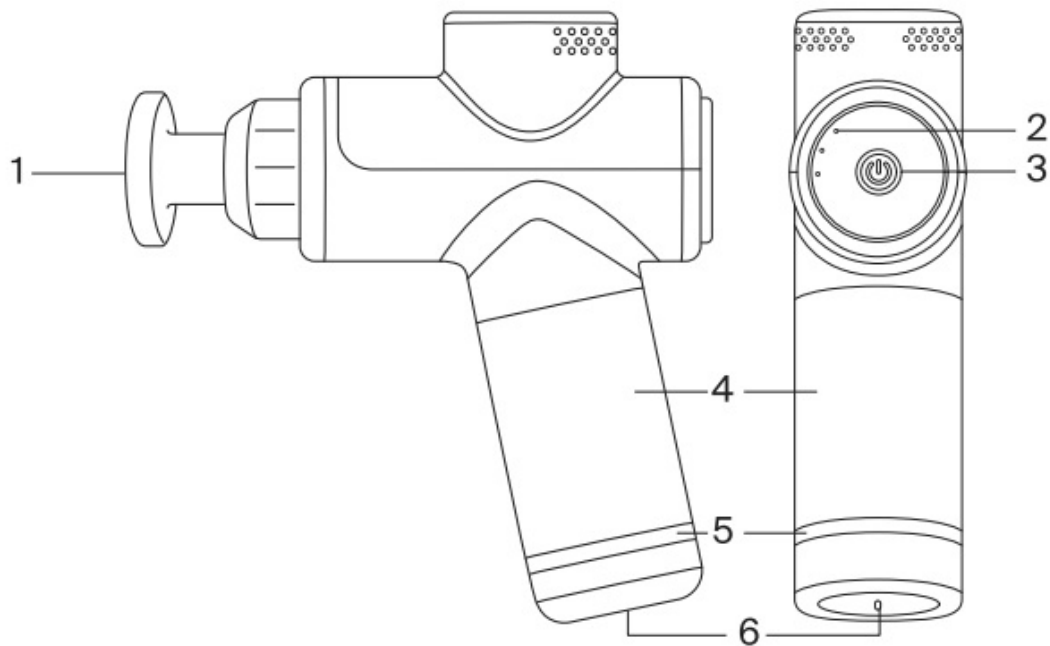

1. Interchangeable head attachments (2)
2. Speed indicators (3)
3. ON / OFF power button and speed setting button
4. Handle
5. Battery level indicator (LED band)
6. Charging port

5V = - 2.5A

1. Operating temperature: +32° F to + 113° F [0° C to +40° C]
2. Storage temperature: -4° F to +158° F [-20° C to +60° C]
3. Rechargeable lithium-ion battery 2500mAh
4. Frequency: Level 1 – 2000/33 Hz, Level 2 – 2400/40 Hz, Level 3 – 2800/47 Hz
5. Weight: 1.5 lbs/0.68 kg

## Warranty

This product is covered by a limited warranty from Hyperice. Please visit [hyperice.com/warranty](https://hyperice.com/warranty) to review the warranty in your country

## FOR MEXICO ONLY

- This device is not intended for use by persons (including children) whose physical, sensory, or mental capabilities are different or reduced, or who lack experience or knowledge unless such persons have been given supervision or training in the operation of the device by a person responsible for their safety.
- Children should be supervised to ensure that they do not use the devices as toys.
- The device must only be powered at the safety extra-low voltage indicated on the device. : The battery must be removed from the device before disposal. the device must be disconnected from the charger when the battery is removed. The battery is removed for safe disposal.


- Operation of this equipment is subject to the following two conditions: (1) this equipment or device may not cause harmful interference, and (2) this equipment or device must accept any interference, including interference that may cause undesired operation.

：Hypervolt Go 2			：55200 001-00			
Equipment name			Type designation (Type)			
Unit	Restricted substances and its chemical symbols					
	(Pb)	(Hg)	(Cd)	(Cr <sup>+6</sup> )	(PBB)	(PBDE)
	○	○	○	○	○	○
	—	○	○	○	○	○
	○	○	○	○	○	○
	○	○	○	○	○	○
	○	○	○	○	○	○
Note 1: “○” indicates that the percentage content of the restricted substances does not exceed the percentage of reference value of presences.						
Note 2: The “—” indicates that the restricted substance corresponds to the exemption.						

- [hyperice.com](https://www.hyperice.com)

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation. Hyperice, Inc., 525 Technology Drive, Suite 100, Irvine, California 92618 USA Copyright 2021 © Hyperice, Inc. All rights reserved. 170-00336 Rev 4

## Documents / Resources

	<a href="#">Hyperice Go 2 Percussion Massage Device</a> [pdf] Instruction Manual Go 2, Go 2 Percussion Massage Device, Percussion Massage Device, Massage Device, Device
---	---

## References

- [User Manual](#)

### [Manuals+](#). [Privacy Policy](#)

This website is an independent publication and is neither affiliated with nor endorsed by any of the trademark owners. The "Bluetooth®" word mark and logos are registered trademarks owned by Bluetooth SIG, Inc. The "Wi-Fi®" word mark and logos are registered trademarks owned by the Wi-Fi Alliance. Any use of these marks on this website does not imply any affiliation with or endorsement.