



# levenhuk LABZZ DM200 LCD Digital Microscope User Manual

[Home](#) » [levenhuk](#) » [levenhuk LABZZ DM200 LCD Digital Microscope User Manual](#) 

## levenhuk LABZZ DM200 LCD Digital Microscope User Manual




### Contents


- [1 Product Overview](#)
- [2 Getting started](#)
- [3 Use with a computer](#)
- [4 Specifications](#)
- [5 Care and maintenance](#)
- [6 Battery safety instructions](#)
- [7 Levenhuk International Lifetime Warranty](#)
- [8 Levenhuk prepared slides sets](#)
- [9 Documents / Resources](#)
  - [9.1 References](#)
- [10 Related Posts](#)

### Product Overview


## Levenhuk LabZZ DM200 LCD

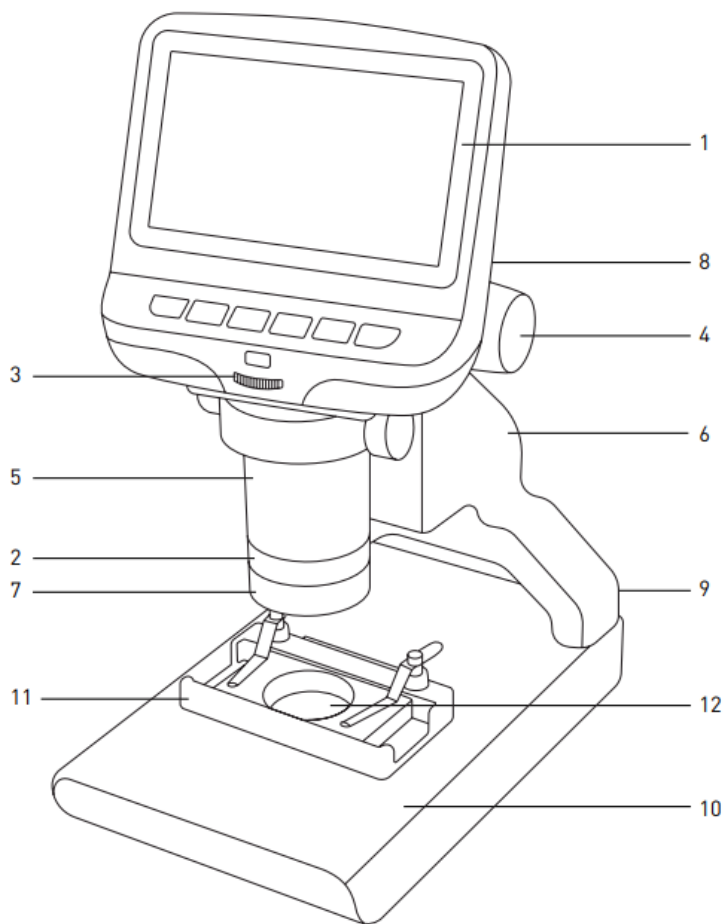

1. LCD screen
2. Objective
3. Brightness adjustment knob
4. Coarse focusing knob
5. Optical zoom ring
6. Pole
7. Upper illumination
8. LCD screen locking screw (not shown)
9. Power connector (not shown)
10. Stage
11. Movable block with specimen holders
12. Lower illumination

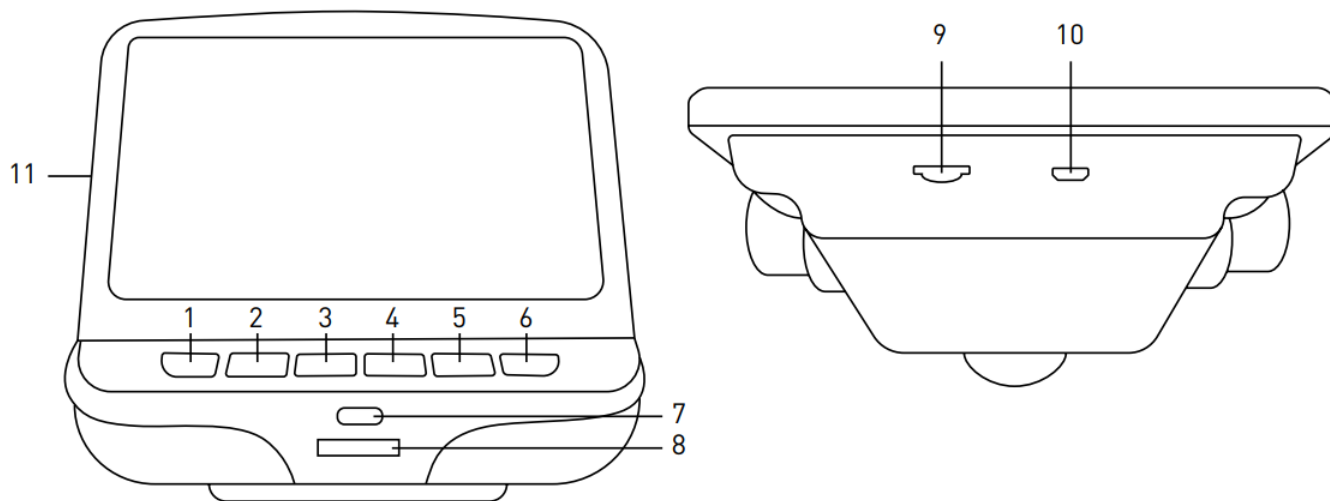


1. Power on/off
2. Menu
3. Up
4. Down
5. OK
6. Capture
7. Infrared sensor
8. Brightness adjustment wheel
9. SD slot
10. MiniUSB: output to PC
11. Reset (not shown)

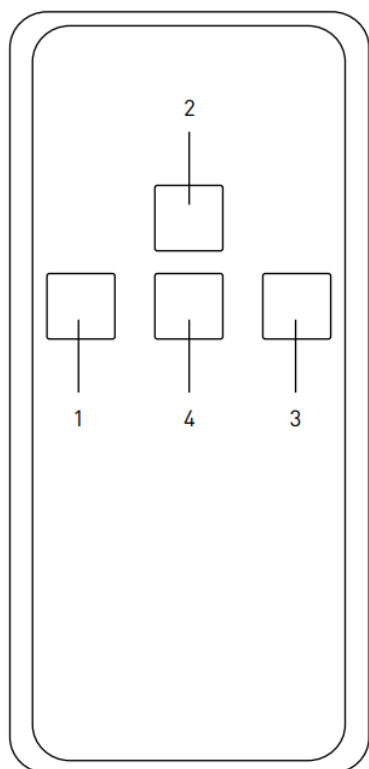

1. Capture

2. Zoom up
3. OK
4. Zoom down

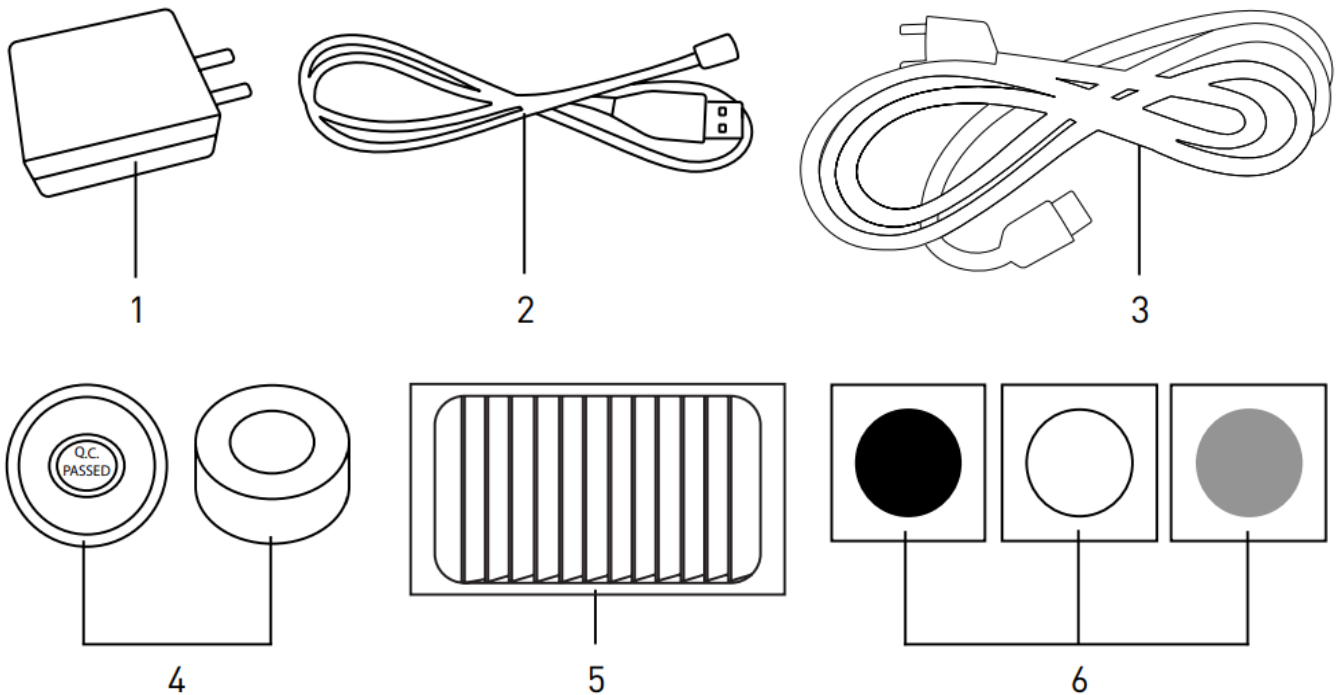

1. DC adapter
2. USB cable
3. Switch cable
4. Lens covers
5. 12 prepared slides
6. Black, white and uniform light plates

Please refer to fig. 1 (device), 2 (LCD screen), 3 (control panel), and 4 (accessories) and view all the details of the device.

## Getting started

- Unpack the microscope carefully and place it on a flat surface.
- Plug the power cable into the corresponding connectors on the base and on the LCD screen, the backlight and screen will turn on automatically.
- Press and hold the “Power on/off” button for 2 seconds on the screen to turn the screen on and off.
- Place the micropreparation on the movable block and fix it with the specimen holders. Place the movable block under the objective. The kit includes a set of ready-made micropreparations for study. Remove the lens cover labeled “Q.C. passed” from the bottom of the second lens cover.
- Adjust the brightness of the backlight and adjust the clarity of the image with the coarse focusing knob.

A microSD card (purchased separately) is required to save captured images and recorded videos. Turn the microscope off and insert the microSD card into the corresponding slot on the microscope body until it locks in place with a click. Do not apply excessive force when inserting the microSD card. Try rotating the card, if it does


not lock in place with ease.

If the LCD screen freezes, look for the “Reset” button on the back of the LCD screen. Use a narrow object, such as a straightened paperclip, to press the button. This will reset all of the settings to the factory defaults and restart the LCD screen. Saved photos or videos will not be deleted.

### Stage inset glass

- Frosted glass stage inset may be used for observations of transparent specimens or thin slices of objects. Place it on the stage and secure it in place with the locking screw. Lower illumination source should be used during observations with this glass stage inset.
- Additional black and white glass slides are supplied with the microscope. If the specimen is white or bright, then for better image contrast, use a black glass slide with only the top illumination on. Otherwise, use a white glass slide.

### Application menu

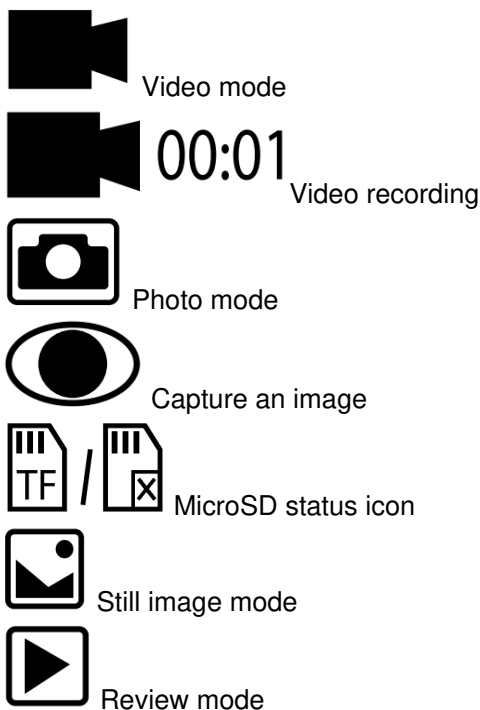

### Digital magnification

During observations, press the “Up” button to enlarge the image and the “Down” button to reduce the image.

### Capturing an image

1. Turn the microscope on and place a specimen on the stage.
2. Set the objective to the desired height using the coarse focusing knob.
3. Adjust the illumination brightness, so that the specimen is evenly lit.
4. Adjust the magnification and sharpness of the image by rotating the coarse focus knob and optical zoom ring.
5. To take a photo, press the “Capture” button.

## Recording a video

1. Press the “Video” button to switch to video mode. The icon in the upper left corner of the screen will change in confirmation.
2. Press the “OK” button to start recording. Note: do not hold down the “OK” button; instead, briefly press and release it.
3. Press the “OK” button again to stop recording.
4. Press the “Capture” button to take a picture while the microscope is recording a video. If you take a photo this way, the image resolution may not be as high as the video resolution – the image quality will be lower.

## View photos and videos

1. Insert the microSD card into the SD slot.
2. Press the “Video” button twice to turn on the viewing mode.
3. Press the “Right” and “Left” buttons to view the pictures and videos.  
Press the “OK” button to start playing the video.
4. Use the “Menu” button to provide additional options when viewing a still image or video. Press and hold the “Menu” button until options appear.

## Setup

Press the “Menu” button to go to the settings menu of the LCD screen. Use the control buttons to select the required parameters. Press the “Menu” button again to save the selected parameters. Do not turn off the power of the microscope while saving the settings. Reboot the microscope after applying the settings.

## Use with a computer

AMCAP, ViewPlayCap, WebCam Monitor, and Microscope Measure programs can be download from the official Levenhuk website.

1. Install the AMCAP (Windows 7/8), ViewPlayCap (Windows 10) and WebCam Monitor (Mac OS) software on the computer. The PC camera icon will appear on your desktop after the installation of the software is complete.
2. Use the “Power on/off” button to turn off the LCD screen.
3. Use the USB cable for connecting the LCD screen to the computer.
4. Click the PC camera icon on your computer’s desktop to start viewing.

## Specifications

|                   |                                 |
|-------------------|---------------------------------|
| <b>Product ID</b> | <b>76827</b>                    |
| <b>Model</b>      | <b>Levenhuk LabZZ DM200 LCD</b> |

|                               |                                |
|-------------------------------|--------------------------------|
| Digital magnification         | 17–220x                        |
| Optical magnification         | 17–55x                         |
| Optics material               | optical glass                  |
| Body                          | plastic                        |
| Stage                         | 120x180mm                      |
| Movable block                 | 83x43mm, with specimen holders |
| Focus                         | coarse, 93mm                   |
| Illumination                  | upper and lower, LED           |
| Rotatable color LCD screen    | 4.3"                           |
| Megapixels (sensor/image)     | 1/12                           |
| Max. resolution (image/video) | 4032x3024px/1920x1080px        |
| Image/video format            | *.jpg/ *.avi                   |
| Sensor                        | 1/4                            |
| Pixel size                    | 3µm                            |

|                             |   |
|-----------------------------|---|
| Sensitivity                 | highly sensitive sensor   |
| Spectral range              | 410–1100nm  |
| Manual settings             | white balance, exposure control   |
| Output (connectors)         | microSD, miniUSB  |
| Power supply (microscope)   | 110–220V; 5V, 1A via USB cable or power adapter; built-in 1800 mAh battery (working time: 2.5 hours, charging time: 10 hours) |
| Power supply (IR remote)    | 1pc CR2025 batteries  |
| Operating temperature range | –20...+70 °C (–4...+158 °F)   |

Levenhuk reserves the right to modify or discontinue any product without prior notice.

## Care and maintenance

- **Never, under any circumstances, look directly at the Sun, another bright source of light or at a laser through this device, as this may cause PERMANENT RETINAL DAMAGE and may lead to BLINDNESS.**
- Take necessary precautions when using the device with children or others who have not read or who do not fully understand these instructions.
- After unpacking your microscope and before using it for the first time check for integrity and durability of every component and connection.
- Do not try to disassemble the device on your own for any reason. For repairs and cleaning of any kind, please contact your local specialized service center.
- Protect the device from sudden impact and excessive mechanical force. Do not apply excessive pressure when adjusting focus. Do not overtighten the locking screws.
- Do not touch the optical surfaces with your fingers. To clean the device exterior, use only special cleaning wipes and special optics cleaning tools from Levenhuk. Do not use any corrosive or acetone-based fluids to clean the optics.
- Abrasive particles, such as sand, should not be wiped off lenses, but instead blown off or brushed away with a soft brush.
- Do not use the device for lengthy periods of time, or leave it unattended in direct sunlight. Keep the device



away from water and high humidity.

- Be careful during your observations, always replace the dust cover after you are finished with observations to protect the device from dust and stains.
- If you are not using your microscope for extended periods of time, store the objective lenses and eyepieces separately from the microscope.
- Store the device in a dry, cool place away from hazardous acids and other chemicals, away from heaters, open fire and other sources of high temperatures.
- When using the microscope, try not to use it near flammable materials or substances (benzene, paper, cardboard, plastic, etc.), as the base may heat up during use, and might become a fire hazard.
- Always unplug the microscope from a power source before opening the base or changing the illumination lamp. Regardless of the lamp type (halogen or incandescent), give it some time to cool down before trying to change it, and always change it to a lamp of the same type.
- Always use the power supply with the proper voltage, i.e. indicated in the specifications of your new microscope. Plugging the instrument into a different power outlet may damage the electric circuitry of the microscope, burn out the lamp, or even cause a short circuit.
- **Seek medical advice immediately if a small part or a battery is swallowed.**

## **Battery safety instructions**

- Always purchase the correct size and grade of battery most suitable for the intended use.
- Always replace the whole set of batteries at one time; taking care not to mix old and new ones, or batteries of different types.
- Clean the battery contacts and also those of the device prior to battery installation.
- Make sure the batteries are installed correctly with regard to polarity (+ and –).
- Remove batteries from equipment that is not to be used for an extended period of time.
- Remove used batteries promptly.
- Never attempt to recharge primary batteries as this may cause leakage, fire, or explosion.
- Never short-circuit batteries as this may lead to high temperatures, leakage, or explosion.
- Never heat batteries in order to revive them.
- Do not disassemble batteries.
- Remember to switch off devices after use.
- Keep batteries out of the reach of children, to avoid risk of ingestion, suffocation, or poisoning.
- Utilize used batteries as prescribed by your country laws

## **Levenhuk International Lifetime Warranty**

All Levenhuk telescopes, microscopes, binoculars, and other optical products, except for their accessories, carry a lifetime warranty against defects in materials and workmanship. A lifetime warranty is a guarantee on the lifetime of the product on the market. All Levenhuk accessories are warranted to be free of defects in materials and workmanship for six months from the purchase date. The warranty entitles you to the free repair or replacement of the Levenhuk product in any country where a Levenhuk office is located if all the warranty conditions are met.

For further details, please visit: [www.levenhuk.com/warranty](http://www.levenhuk.com/warranty)

If warranty problems arise, or if you need assistance in using your product, contact the local Levenhuk branch.

## **Levenhuk prepared slides sets**

The Levenhuk slides sets include thoroughly prepared specimens of various biomaterials as well as blank slides and cover glasses for you to make your own specimens. Prepared slides, blank slides and cover slips can be used with any microscope model.



#### Levenhuk Worldwide:

USA: [www.levenhuk.com](http://www.levenhuk.com)  
Azerbaijan: [www.levenhuk.az](http://www.levenhuk.az)  
Bulgaria: [www.levenhuk.bg](http://www.levenhuk.bg)  
Canada: [www.levenhuk.ca](http://www.levenhuk.ca)  
Czech Rep.: [www.levenhuk.cz](http://www.levenhuk.cz)  
Estonia: [www.levenhuk.ee](http://www.levenhuk.ee)  
Finland: [fi.levenhuk.com](http://fi.levenhuk.com)  
Germany: [www.levenhuk.de](http://www.levenhuk.de)  
Hungary: [www.levenhuk](http://www.levenhuk)  
Iceland: [www.levenhuk.is](http://www.levenhuk.is)  
India: [www.levenhuk.in](http://www.levenhuk.in)  
Italy: [www.levenhuk.it](http://www.levenhuk.it)  
Latvia: [www.levenhuk.lv](http://www.levenhuk.lv)  
Lithuania: [www.levenhuk.lt](http://www.levenhuk.lt)  
Netherlands: [www.levenhuk.nl](http://www.levenhuk.nl)  
Poland: [www.levenhuk.pl](http://www.levenhuk.pl)  
Romania: [www.levenhuk.ro](http://www.levenhuk.ro)  
Russia: [www.levenhuk.ru](http://www.levenhuk.ru)  
Slovakia: [www.levenhuk.sk](http://www.levenhuk.sk)  
Spain: [www.levenhuk.es](http://www.levenhuk.es)  
Turkey: [www.levenhuk.tr](http://www.levenhuk.tr)  
Ukraine: [www.levenhuk.ua](http://www.levenhuk.ua)  
United Kingdom: [www.levenhuk.uk](http://www.levenhuk.uk)  
EU: [www.levenhuk.eu](http://www.levenhuk.eu)

#### [levenhuk.com](http://www.levenhuk.com)

Levenhuk Inc. (USA): 928 E 124th Ave. Ste D, Tampa, FL 33612, USA, +1 813 468-3001,

[contact\\_us@levenhuk.com](mailto:contact_us@levenhuk.com)

Levenhuk Optics s.r.o. (Europe): V Chotejně 700/7, 102 00 Prague 102, Czech Republic, +420 737-004-919,

[sales-info@levenhuk.cz](mailto:sales-info@levenhuk.cz)


Levenhuk® is registered trademark of Levenhuk, Inc.

© 2006—2022 Levenhuk, Inc. All rights reserved.

20220124



## Documents / Resources

|   |  |
|---|--|
|  | <p><a href="#">Levenhuk LABZZ DM200 LCD Digital Microscope</a> [pdf] User Manual<br/>LABZZ, DM200 LCD Digital Microscope, DM200, LABZZ DM200, LCD Digital Microscope, LABZZ DM200 LCD Digital Microscope, Digital Microscope, Microscope</p> |
|---|--|

## References

- [Levenhuk's Official Finland Website!](#)
- [Levenhuk optical instruments store | Levenhuk - best optical equipment](#)
- [Официален уебсайт на Levenhuk в България](#)
- [Доживотна гаранция на Levenhuk – Официален уебсайт на Levenhuk в България](#)
- [Levenhuk optical instruments store | Levenhuk - best optical equipment](#)
- [Levenhuk lifetime warranty details | Levenhuk - best optical equipment](#)
- [Oficiální webové stránky Levenhuk pro Českou republiku](#)
- [Doživotní záruka společnosti Levenhuk – Oficiální webové stránky Levenhuk pro Českou republiku](#)
- [Die offizielle Website von Levenhuk in Deutschland](#)
- [Levenhuk Lebenslange Garantie – Die offizielle Website von Levenhuk in Deutschland](#)
- [Compre desde el sitio web oficial de Levenhuk en España](#)
- [Garantía internacional de por vida Levenhuk – Compre desde el sitio web oficial de Levenhuk en España](#)
- [Levenhuk's official website in Europe](#)
- [Levenhuk Lifetime Warranty – Levenhuk's official website in USA](#)
- [A Levenhuk élettartamra szóló szavatossága – A Levenhuk hivatalos magyarországi weboldala](#)
- [Levenhuk's official website in Europe](#)
- [Sito web ufficiale di Levenhuk in Italia](#)
- [Skyhunters Europe](#)
- [SkyHunters](#)
- [Levenhuk's official website in Europe](#)
- [Oficjalna witryna internetowa Levenhuk w Polsce](#)
- [Gwarancja bezterminowa Levenhuk – Oficjalna witryna internetowa Levenhuk w Polsce](#)
- [Официальный сайт Левенгук](#)
- [Поддержка - Гарантийное обслуживание Левенгук - Levenhuk Russia](#)
- [Optický přístroj Levenhuk - online prodejna optických přístrojů](#)
- [Интернет-магазин оптических приборов](#)

-  [Levenhuk's official website in Europe](#)

Manuals+.