


Snefe

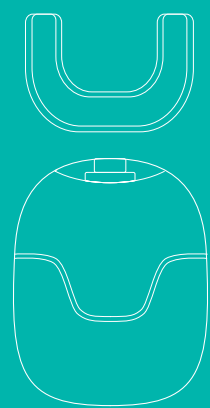

Smart U- toothbrush

Product Features

1. 360° Wrapped mouth.
2. IPX7 Waterproof.
3. Patent (u-shaped toothbrush Head).
4. 45° Scientific bass brushing method.

Toothbrush Specifications:

Product name : Snefe 360 ultrasonic electric toothbrush.
Product weight: 78g
Rated voltage: 5V
Rated power: 5W
Rated power: <1A
Charging time: 60min


List of Package:

Toothbrush host	x1
Toothbrush head	x1
USB-C Cable	x1
Instruction Manual	x1

01

How To Use

1. Please rinse up the brush head before use.
2. Apply appropriate amount of toothpaste on the four corners of the brush head.
3. Put the toothbrush into your mouth and bite gently, Press the button, choose the mode suitable you and do a brushing. (When brushing your teeth, swing the toothbrush from side to side to touch the teeth behind is recommended to perform two to three brush with the smart time).

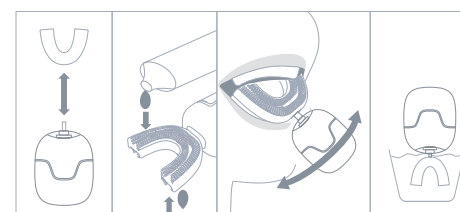

Note: the product pictures in this manual are for reference only, subject to physical objects.

02

Various Frequency Vibration

Snefe ultrasonic toothbrush has 4 working modes, including 3 different vibration frequencies and lights.
#1 mode: 5000 times/min, blue light, automatic shutdown in 55 seconds.
#2 mode: 8000 times/min, blue light, automatic shutdown in 55 seconds.
#3 mode: 15000 times/min, red light, automatic shutdown in 120 seconds.
#4 Cold light whitening mode: blue light no vibration, 5 minutes suggest to be used with whitening gel.

Tips:

Turn On: single click
Switch mode: single click
Close: extended press to turn off

DIY:
You can modify the Mouthpiece to achieve the appropriate length and width to fit your dental arch.



03

About Charging

Start Charging: Snefe toothbrush built in blue breathing light, it will be show in charging status.
Charging completed: The blue light will keep on , stop breathing battery indicator will keep in flash when in charging.

04

Certificate

Quality inspector: _____
Examination time: _____

Executive standard: GB 4706.1/G B4706.59

This product contains rechargeable batteries.

The treatment of batteries shall meet the requirements of local government.



05

Customer Service

If you have any suggestions or product problems when using , please send a message to our mailbox, and we will respond you and provide the best solution within 24 hours.

Email: lomom_@hotmail.com

06

07