




**MEDIC
THERAPEUTICS**

SHIATSU FOOT MASSAGER

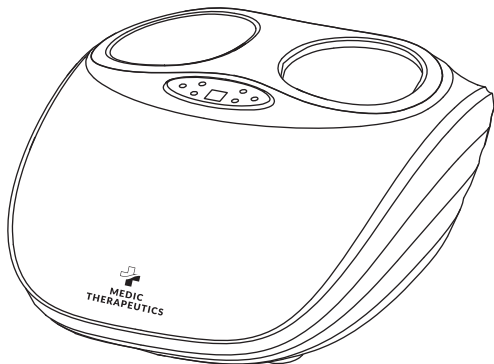

BASIC PRECAUTIONS:

1. Only use with parts included by the manufacturer.
2. Contact manufacturer for repair or servicing.
3. Do not use in damp or high-heat environments.
4. Always unplug after use.
5. Use only for 20-minute time increments.
6. If you start to feel unwell or uncomfortable during use, stop using immediately.
7. In instances of excessive temperature or strange order, discontinue use and unplug immediately.
8. Do not use the power cord to drag or move the device.
9. Device not suitable for children or physical compromised or who have paralysis.
10. Do not insert any foreign objects into the device.
11. Use with caution if you are heat sensitive in any way.

ADDITIONAL PRECAUTIONS:

1. This device is not to be used to treat or cure any medical conditions.
2. Do not use on pets.
3. Always consult a physician prior to use if any of the following pertain to you:
 - Pregnancy
 - Elderly
 - Skin rashes, infectious diseases, redness or scalds and burns
 - Heart Disease
 - Malignant Tumor
 - Osteoporosis
 - High or abnormal blood pressure
 - Cardiac Pacemaker or any other such medical device
 - Ongoing medical treatment


STORAGE AND CLEANING:

1. Prior to use, ensure the device has been unplugged and switched off.
Then plug into the outlet and switch on.
2. Use a clean, dry cloth to clean the surface of the device.
3. Do not use a corrosive or harsh cleaning agent or detergent on the device.
4. After use, always unplug and store in a cool, dry place.
5. To ensure overall safety or damage to the device,
do not use in a damp environment or in direct sunlight.


PRODUCT SPECIFICATIONS:

• Size of product itself:	360*400*210MM
• Net weight of product itself:	3.25KG
• Single package size:	375*225*435MM
• Net weight of product including package:	4.04KG
• Materials:	Plastic
• Certifications (FDA, UL, ETL) - with UL plug	
• List of all the parts and pcs included:	1x Shiatsu Foot Massager 1x Manual 1x Plug

PRODUCT DIAGRAM:

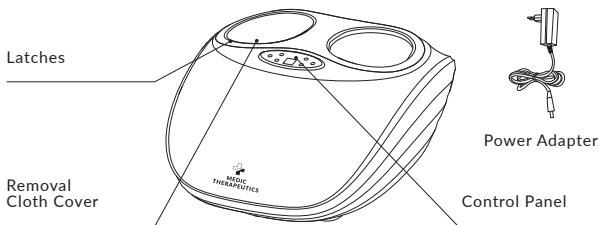

CONTROL PANEL DIAGRAM:

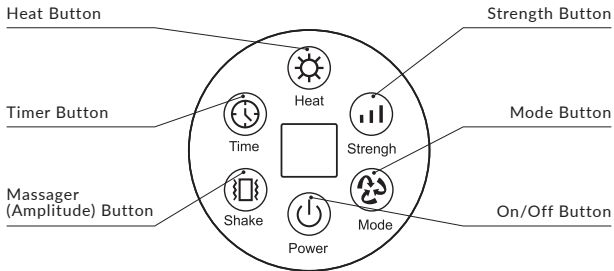

USAGE AND SETTINGS:

1. Power Button turns the device on or off. The default time for the massager is 15 minutes.
2. Use the Mode Button to switch between 3 types of operating modes.
P1 Setting: The strength for both the sole and pressure are adjustable, and the heat can be turned off or on.
P2 Setting: Only the pressure massage with adjustable strength is on. Heat can be on or off.
P3: Only the sole massager is on. Heat can be off or on.
3. The default of the Strength Button is set at midrange. To adjust, press the button to switch between 3 levels of pressure: L1 (Low), L2 (Midrange), L3 (High)

4. The default time is set for 15minutes. Use the Timer Button to switch between 10, 15, and 20minute time cycles.
5. The heating automatically switches on, indicated by the red LED light. To turn the heat off, press the Heat Button and the LED light will also switch off.
6. The default massager (amplitude) is set at midrange. To adjust, press the Massager (Amplitude) Button to switch between E1 (low) and E2 (high) massaging levels.

TROUBLESHOOTING:

PROBLEM:	CAUSE:	SOLUTION:
Unable to turn on.	1. Power adapter is not connected. 2. The plug is not plugged securely into the wall outlet. 3. The massage time exceeds the set time.	1. Re-Connect Power Adapter 2. Re-Plug into outlet. 3. Switch the power off and on again.
Suddenly stops during use.	The massage time is too long.	Turn off for one solid hour and try again.
Faint sound while in use.	1. General motor sound. 2. The sound of the rotator mechanism rotating.	Common sounds, perfectly normal while in use.

LIMITED LIFETIME WARRANTY

Your Medic Therapeutics Shiatsu Foot Massager is backed by a limited lifetime manufacturer's warranty. Medic Therapeutics will repair or replace your device at any time should it fail due to a defect in material or workmanship, subject to the certain limitations.

This limited warranty does not cover any damage that results from unauthorized or improper use, service, or repair. Further, it does not cover damage caused by accident, impact, negligence, or normal wear and tear. Should you discover your Medic Therapeutics Shiatsu Foot Massager is not functioning properly, please send your device to our repair center for evaluation. If your product cannot be repaired or serviced, we will reserve the right to change it for a similar or newer model.

Please note that a flat rate of \$35.00 will be charged to cover a service evaluation fee and return shipping of your device. All warranty claims must be accompanied by a copy of your proof of purchase from an authorized retailer. Please send your device, proof of purchase, and a check or money order in the amount of \$35.00 made out to Medic Therapeutics to:

Address:
Medic Therapeutics Service Center
3069 Taft Street
Hollywood, FL 33009

Contact:
warranty@medictherapeutics.com