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4

HEAD EYE MASSAGER

(Air pressure+ Vibration + Hot compress
+Built-in music+ Voice broadcast)

INSTRUCTION MANUAL

Please read this manual carefully before using the product This manual contains service instructions, please keep it in a safe place.

Basic Product Information

Products


Intelligent head eye massager adopts microcomputer chip control technology, through the air pressure/vibration/hot compress effect on the head and eye area, it helps the head and eye muscles to relax, and brings you an all-round comfortable experience.

This product is suitable for all kinds of people who use computers, cell phones and other electronic

Product Parameters

[Product Features] Head massage [Battery Capacity] 1000mAh
[Massage Mode] Five modes adjustable [Working Time] 15 minutes
[Rated Voltage] Dc 5V [Rated Current]1A
[Packing List] Host *1 USB cable *1 Instruction manual *1
[Product Features] Air pressure Hot compress Vibration massage Built-in music
Voicebroadcast

Product mix

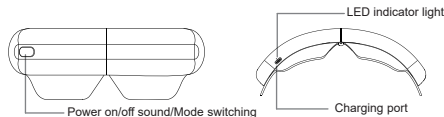(Note: The figure above is a schematic diagram, product appearance shall be subject to the actual)

1

Preparation and charging before use

Pre-use Preparation

Aim the head massager at your forehead, put on the machine, adjust the elastic band to the most comfortable size, and adjust the machine up and down to the most comfortable position.

Please remove glasses, contact lenses, hair clips, hair accessories, etc. before wearing to avoid discomfort during use. Do not use this product with wet hands. Do not fall or step on this product.

Charging

The USB plug of the USB charging cable is inserted into the USB socket of the power adapter (5V-1A), then the power adapter (5V-1A) into the power outlet, and finally insert the plug of the USB charging cable into the charging port of the product.

Warm reminder

During charging, the indicator light flashes, and the indicator light is always on when fully charged.

During charging, the product can not be turned on and work, after charging is completed, please disconnect the power adapter.

If the product is charging while working, the product will stop working automatically.

Usage

Turn the product on or off

Press the "⏻" key for 2 seconds to turn on the machine. After the machine is turned on, the indicator light is lit and the machine starts to work, default is automatic mode. The working time is 15 minutes, playing music.

While the product is working, keep pressing the "⏻" button for 2 seconds to turn off the machine or the machine will turn off automatically after 15 minutes of use.

Switching Mode

Press the "⏻" key for less than one second to switch modes. Operate once to switch to the second mode; operate twice to switch to the third mode; and so on.

Music control

Press the "🔊" key twice in a row for less than - seconds to turn the music off or on.

2

Mode introduction

Automatic mode = Air pressure massage + hot compress massage + vibration massage + music

Decompression mode = Air pressure massage + hot compress massage + music

Vitality mode = Air pressure massage + vibratory massage + music

Comfort mode = Hot compress massage + music

Sleep mode = Air pressure massage + music

(Note: Because the product version is constantly updated, the actual mode may differ from the manual)

Transportation and storage

Storage environment

Room temperature, no corrosive gases, well ventilated, dry room.

Transport

This product has simple shock-proof facilities, suitable for air, railway, road and ship transport, should avoid rain, exposure, inversion and collision and other bad conditions.

Cleaning and storage

Please use a dry soft cloth to clean the product (please turn off the machine when cleaning, do not use benzene or diluent and other corrosive cleaning agents, otherwise it will cause fading, surface aging, structural damage and other problems).

Please put the product in a well-ventilated, dry room, and cover it with a protective cover when it is not used for a long time. When the storage time exceeds 3 months, please open the box for inspection and charge for more than half an hour.

Caveat

1. Do not use when there are symptoms of eye disease. People who have performed eye surgery, suffered from glaucoma, cataract, retinal detachment and other major eye diseases are prohibited from using.

2. Do not wash this product directly with water. Please wipe it with cotton soaked in alcohol or water.

3. This product has a hot surface, please pay attention to the use of people who are not sensitive to heat. Minors, the disabled, the elderly, etc., should be used under adult supervision.

3

4. For sensitive skin, atopic dermatitis, cosmetic materials placed on the site and such special groups, please use this product under the guidance of a doctor.

5. Do not disassemble this product by yourself, This product contains batteries, do not put dry high temperature (above 50° C) or into the fire, otherwise there will be explosion risk.

6. This product is a personal health care massage product, can not replace medicine or medical equipment.

7. This product is not a toy and should not be used by children alone unless parents or guardians have set it up in advance, or have given children adequate guidance on its use.

8. The use of this product will block the line of sight, do not use in walking, driving and other behaviors.

9. After using this product, the eyes need to return from the relaxed state to the tense state, and there will be short-term blurred vision, which is a normal phenomenon. Please close your eyes for 1-3 minutes, and your vision will gradually become clear.

10. The storage and use temperature of this product is -10°c ~ 50°c, it is prohibited to store and use in extreme environment (low dry temperature negative 10°C or more than 50°C).

11. This product should be charged at least once every 3 months when it is not used for a long time.

12. Due to product hardware and software upgrade, product functions, please refer to the actual.

Common problems and solutions

Problems encountered	Possible cause	Treatment method
Automatic shutdown unable to turn on	Low battery	Please charge
The charging indicator light is of	The UsB port is improperly connected to the power supply	Plug in the usB charging cable/power cable

If the above processing is carried out and there is an abnormal phenomenon, please turn off the power and contact the dealer you purchased or the company.

Please do not disassemble or repair the product without authorization, and the company will not be responsible for any product damage, personal injury or other losses.

4