MASSAGE GUN / PISTOLET DE MASSAGE SMGB A1

(GB) (IE)

MASSAGE GUN

Operating instructions

(FR) (BE)

PISTOLET DE MASSAGE

Mode d'emploi

(CZ)

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(SK)

MASÁŽNY PRÍSTROJ

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(DK)

MASSAGEPISTOL

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(HU)

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(DE) (AT) (CH)

MASSAGE GUN

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(NL) (BE)

MASSAGE GUN

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(PL)

MASSAGE GUN - PISTOLET DO MASAŻU

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(ES)

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(IT)

PISTOLA MASSAGGIANTE

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
(SI)

MASAŽNA PIŠTOLA

Navodila za uporabo

IAN 391896_2201

(DE) (AT) (BE)
(NL) (CZ) (PL)
(SK) (DK) (HU)

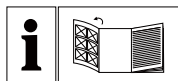GB IE

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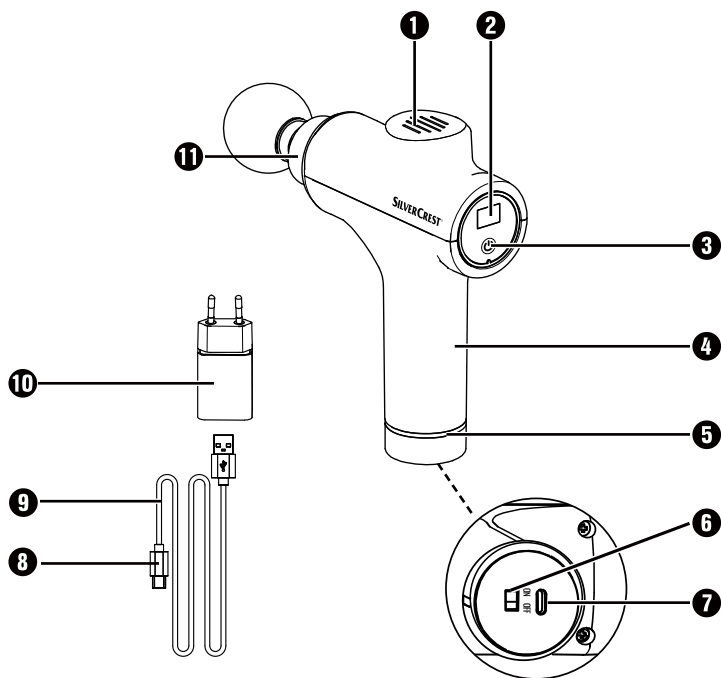
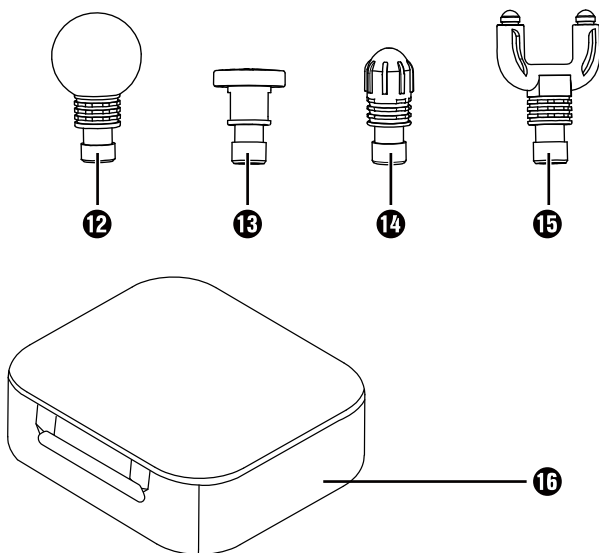
A**B**

Contents

Introduction	2
Intended use	2
Warning symbols used	2
Notes on trademarks	3
Safety information	3
Safety guidelines for rechargeable batteries	6
Package contents	7
Disposal of the packaging	8
Parts description	8
Before first use	9
Charging the appliance	9
Using the appliance	10
Selecting/fitting attachments	10
Examples of use	11
Switching the appliance on/off and selecting intensity	12
Intensity levels	13
Cleaning and care	13
Storage	14
Troubleshooting	14
Disposal	14
Technical specifications	15
Kompernass Handels GmbH warranty	16
Service	18
Importer	18

Introduction

Congratulations on the purchase of your new appliance.

You have selected a high-quality product. The operating instructions are part of this product. They contain important information about safety, usage and disposal. Before using the product, please familiarise yourself with all operating and safety instructions. Use the product only as described and for the specified range of applications. Please also pass these operating instructions on to any future owner.

Intended use

This appliance and its attachments are intended exclusively for massaging the muscles of the neck, shoulders, arms, back (excluding the spine), thighs, calves and buttocks. Only use the attachments intended for this purpose. Never use the appliance directly on bones. The appliance is not intended for use on the face, head, neck or genital area. Do not use the appliance on your stomach or chest or directly on your heart. Never use the appliance to massage babies, small children, people who are incapacitated or animals.

This appliance is intended solely for personal, private use. It is not to be used for medical or commercial applications and is not a substitute for medical treatment. It is not intended for use in medical/therapeutic or commercial fields. Any other use or use beyond that specified is considered improper use.

The manufacturer accepts no responsibility for damage caused by failure to observe these instructions, improper use or repairs, unauthorised modifications or the use of unapproved replacement parts. The operator bears sole liability.

Warning symbols used

The following warnings are used in these operating instructions (where applicable):

DANGER

A warning notice at this hazard level indicates an imminently hazardous situation.

If the hazardous situation is not avoided, it could result in death or serious physical injury.

- Follow the instructions in this warning to prevent a fatality or serious personal injury.

WARNING

A warning at this hazard level indicates a potentially hazardous situation.

Failure to avoid this hazardous situation could result in injuries.

- Follow the instructions in this warning notice to prevent injuries.

CAUTION

A warning at this hazard level indicates a risk of property damage.

Failure to avoid this situation could result in property damage.

- ▶ Follow the instructions on this warning label to prevent property damage.

NOTE

- ▶ A note provides additional information which will assist you in using the appliance.

Notes on trademarks

- USB® is a registered trademark of USB Implementers Forum, Inc.
All other names and products may be trademarks or registered trademarks of their respective owner

Safety information

To ensure safe operation of the appliance, follow the safety guidelines set out below.

RISK OF ELECTROCUTION!

- ▶ Charge the appliance only with the mains adapter supplied.
- ▶ Connect the mains adapter only to correctly installed and earthed mains sockets. Ensure that the power supply corresponds with the information on the type plate of the appliance/mains adapter.
- ▶ Make sure that the power cable does not become wet or moist. Route the cable in such a way that it cannot be crushed or damaged.
- ▶ Keep the power cable away from hot surfaces.
- ▶ Do not perform any repairs to the appliance. All repairs must be performed by our Customer Service Unit or by a qualified specialist technician.
- ▶ Never open the appliance housing. There are no parts inside the appliance which require user maintenance. In addition, this will invalidate the warranty.

RISK OF ELECTROCUTION!

- ▶ Switch off the appliance after each use and before cleaning it.
- ▶ Use the appliance in dry interior rooms only; never use it outside.
- ▶ Never immerse the appliance in water or other liquids!
- ▶ Do not use the appliance while taking a bath or shower. Do not store or keep the appliance in places from where it can fall into the bath tub, sink or other containers filled with water or other liquids.
- ▶ Never reach for the appliance if it has fallen into water.
- ▶ Do not use the appliance if the power cable is damaged, the mains adapter is defective, the appliance is not working properly, if the attachments show signs of damage, if it has fallen into water or is otherwise damaged.

Do not use the mains adapter if it has a damaged plug or cable. If the plug contacts of the mains adapter are defective or bent, the mains adapter cannot be repaired and must be scrapped.

WARNING! RISK OF INJURY!

- ▶ This appliance may be used by children aged 8 years and above and by persons with limited physical, sensory or mental capabilities or lack of experience and knowledge, provided that they are under supervision or have been told how to use the appliance safely and are aware of the potential risks.
- ▶ Cleaning and user maintenance tasks may not be carried out by children unless they are supervised.
- ▶ Children must not be allowed to play with the appliance.
- ▶ Never leave children unattended with the packaging material. Keep children away from the packaging and the appliance. There is a risk of suffocation!

⚠ WARNING! RISK OF INJURY!

- ▶ People with pacemakers, artificial joints or electronic implants should consult their doctor before using the appliance.
- ▶ Do not use the appliance around your heart, face or genital area.
- ▶ Do not use the appliance if you suffer from one or more of the following diseases or medical conditions: circulatory disorders, varicose veins, open wounds, bruises, cracked skin, phlebitis, thrombosis, diabetes or slipped disc.
- ▶ Never use the appliance on children, animals, during pregnancy, on your head, throat or other parts of the body that are exhibiting swelling, burns, inflammation, wounds or other sensitive areas.
- ▶ Do not use the appliance under blankets or pillows. Never use the appliance in the vicinity of petrol or other inflammable materials.
- ▶ Do not use the appliance while sleeping or after taking substances that may reduce your ability to react quickly (e.g. painkillers or alcohol).
- ▶ Never use the appliance directly on bones, e.g. on your spine, shins, the tops of your feet or similar places on your body.
- ▶ Use the appliance only on dry and clean parts of your body.
- ▶ Stop using the appliance immediately if you find the massage unpleasant or if you experience pain.
- ▶ Uninterrupted application of the appliance should not last any longer than 15 minutes. After use, switch off the appliance for at least 15 minutes and allow it to cool down.
- ▶ Do not insert your fingers or any other objects into the opening for the attachments. Keep long hair away from moving parts.
- ▶ Use the appliance only for its intended purpose.

CAUTION! PROPERTY DAMAGE!

Use the mains adapter/the appliance indoors only.

- ▶ Use the appliance only with the accessories supplied (mains adapter and attachments). Make sure that the attachments can always move freely.
- ▶ Protect the appliance and attachments from shocks, dust, chemicals or electromagnetic fields. Keep the appliance away from heat sources (radiators, ovens, etc.) and do not expose it to extreme fluctuations of temperature.
- ▶ Do not use any aggressive cleaning agents or solvents. They could damage the surface.
- ▶ Keep the ventilation slits free of dust and do not cover them.

NOTE

- ▶ No user action is required to switch the product between 50 and 60 Hz. The product switches automatically to either 50 or 60 Hz as required.

Safety guidelines for rechargeable batteries

The product is equipped or supplied with a lithium-ion battery. The battery in this appliance cannot be replaced.

RISK OF EXPLOSION!

- ▶ Do not open or short-circuit the battery. This could cause the battery to overheat and explode. There is a risk of fire!
- ▶ Never throw the battery into a fire or water.
- ▶ Do not subject the battery to mechanical stress.
- ▶ Charge the battery only with the supplied mains adapter (type GQ12-050200-ZG)
- ▶ Never leave the appliance unattended during charging.

⚠ WARNING! RISK OF INJURY DUE TO BATTERY LEAKAGE!

- ▶ High temperatures can affect the battery. Avoid direct sunlight and proximity to heat sources such as radiators.
- ▶ Under extreme conditions, the battery cells may leak. Avoid contact with skin, eyes and mucous membranes. In the event of accidental contact, rinse the affected areas with plenty of clean water and contact a doctor.
- ▶ Leaking or damaged batteries can cause chemical burns to the skin. Always wear protective gloves if you have to touch the batteries.

Package contents

The appliance is supplied with the following components as standard:

- Massage gun
 - 4 attachments (1 ball, 1 spiral, 1 U-shaped and 1 flat head)
 - Storage bag
 - Mains adapter with USB charging cable
 - Operating instructions
- 1) Remove all parts of the appliance and the operating instructions from the box.
 - 2) Remove all packaging materials and any stickers from the appliance.
 - 3) Remove the protective film from the display.

NOTE

- ▶ Check contents for completeness and damage immediately after unpacking. Contact Customer Service if necessary.

Disposal of the packaging

The packaging materials have been selected for their environmental friendliness and ease of disposal, and are therefore recyclable.

Dispose of packaging materials that are no longer needed in accordance with applicable local regulations.

Dispose of the packaging in an environmentally friendly manner. Note the labelling on the packaging and separate the packaging material components for disposal, if necessary. The packaging material is labelled with abbreviations (a) and numbers (b) with the following meanings:

1 - 7: plastics, 20 - 22: paper and cardboard, 80 - 98: composites.

Parts description

Figure A (see the fold-out page):

- ❶ Ventilation slits
- ❷ Display/intensity display
- ❸ Button  (turn device on/off / set intensity level)
- ❹ Handle
- ❺ Ring display
- ❻ Main switch ON/OFF
- ❼ Connection socket
- ❽ USB plug (type C)
- ❾ USB charging cable (USB type-A to USB type-C)
- ❿ Mains adapter
- ⓫ Socket for attachments

Figure B:

- ❿ Ball head
- ⓫ Flat head
- ⓬ Spiral head
- ⓭ U-shaped head
- ⓮ Storage bag

Before first use

Before you can use the appliance, it needs to be charged for at least 3 ½ hours.

Charging the appliance

CAUTION

- Use only the supplied mains adapter 10 to charge the battery.
- Do not use the appliance while it is charging.

- 1) Make sure that the appliance is switched off. The main switch 6 must be set to OFF.
- 2) Insert the USB plug 8 into the USB connection socket 7 on the underside of the handle 4 (fig. 1):

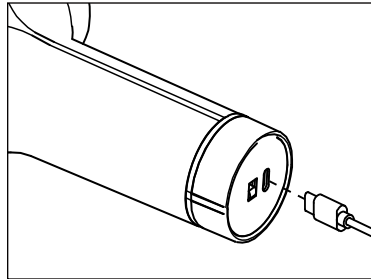
Fig. 1

- 3) Connect the USB charging cable 9 to the mains adapter 10 and plug the mains adapter 10 into a suitable mains socket.
- ✓ When the ring display 5 lights up green, the appliance is fully charged.

The battery level is displayed as follows:

Ring display 5	Meaning
Lights up red	Appliance switched on, battery charge less than 20 %.
Flashes red	Appliance is charging, battery level less than 30 %.
Flashes green	Appliance is charging, battery level at least 30 %.
Lights up green	Appliance switched on, battery fully charged.

- 4) Disconnect the appliance from the power supply and unplug the USB charging cable 9 from the appliance.

A full battery charge lasts up to 4 hours of operation at a medium intensity level. The battery can be recharged at any time and at any battery level.

CAUTION

- The appliance has a deep discharge protection device. When the battery charge level is too low, the display **2** shows the message **E !**. Charge the appliance immediately to avoid damage to the battery.

Using the appliance

You can use this appliance to give yourself or others a soothing and effective trigger point massage. Trigger point massages can relax or stimulate the muscles to relieve tension, pain and fatigue. The massage gun provides a powerful, intensive massage for the thighs, calves, shoulders, arms and the back region (excluding the spine). Depending on the extent of the tension in the part of the body to be treated, a maximum treatment time of approx. 15 minutes is recommended.

A massage should always feel pleasant. A slight reddening of the treated areas and even slight bruising are a normal side effect, as the massage specifically promotes blood circulation. However, if you experience excessive skin reactions, discontinue use and consult a doctor.

Selecting/fitting attachments

- 1) Select the appropriate attachment for the desired treatment and body part as follows:

Attachment	Suitable for:
 12	Uniform and gentle massage to relax smaller muscle groups, e.g. upper arms, shoulders, waist and legs.
 13	Large areas all over the body, such as buttocks, thighs and back (excluding spine).
 14	Targeted and selective massage of deep muscle tissue, e.g. soles of the feet, palms of the hands, muscles of the shoulder blades.
 15	Uniform, 2-point stimulation of deep muscles in the lower back (except spine), neck, Achilles tendon.

- 2) Insert the desired attachment into the socket ⑪ and press it as far as it will go (Fig. 2):

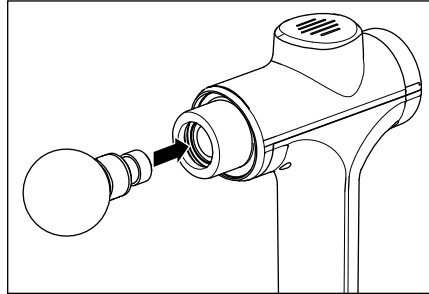
Fig. 2

- 3) To change the attachment, carefully pull it out of the socket ⑪ horizontally, turning it slightly if necessary.

Examples of use

NOTE

- Do not apply the attachments directly to a bone, only to muscle tissue.

Examples of application on arms (fig. 3), thighs (fig. 4), sole of the foot (fig. 5) and Achilles tendon (fig. 6).

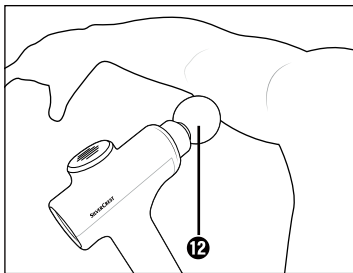
Fig. 3

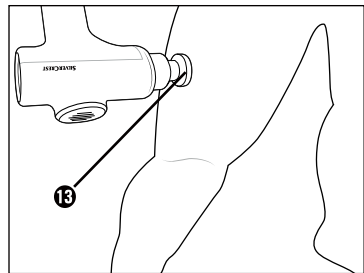
Fig. 4

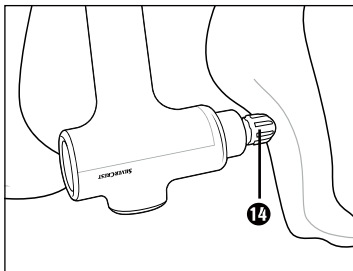
Fig. 5

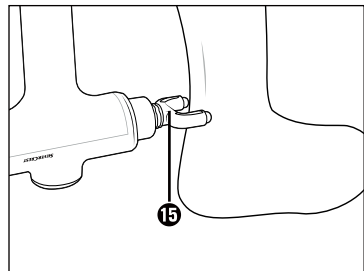
Fig. 6

Switching the appliance on/off and selecting intensity

- 1) Switch the appliance on by moving the main switch **6** on the underside of the handle **4** to the ON position.
 - ✓ The ring display **5** lights up green.
- 2) Press the button **1** **3** to start the appliance.
 - ✓ The display **2** shows the intensity level **1**.
The appliance starts on intensity level **1**.
- 3) To increase the intensity level, press the button **1** **3** briefly.
 - ✓ Each press increases the intensity by one level. This is shown on the display **2**.

NOTE

- After reaching the highest intensity level **05**, if you press the button **1** **3** again, the appliance stops working and switches into standby mode.

- 4) Massage the desired body part by gently moving the attachment back and forth over the muscle area.

NOTE

- Do not apply too much pressure, otherwise the appliance will automatically switch off to protect the motor from a possible overload. If this safety device has been triggered, **E3** will appear on the display **2**.
You need to restart the appliance.

- 5) After you have finished using the appliance, switch it off by pressing and holding the button **1** **3** for around 2 seconds and then move the main switch **6** to the OFF position.

NOTE

- The appliance is equipped with an automatic switch-off to prevent over-heating. The appliance switches off automatically after 15 minutes of continuous operation. Allow the appliance to cool down for at least 15 minutes before using it again.

Intensity levels

The intensity of the massage varies depending on the level selected, as follows:

Level	Intensity
00	Off
01	1800 rpm – 30 Hz
02	2100 rpm – 35 Hz
03	2400 rpm – 40 Hz
04	2700 rpm – 45 Hz
05	3000 rpm – 50 Hz

Cleaning and care

DANGER! ELECTRIC SHOCK!

- ▶ Do not clean the appliance while it is charging. Switch off the appliance before cleaning.
- ▶ Ensure that no moisture can penetrate the device during cleaning.

ATTENTION! PROPERTY DAMAGE!

- ▶ Do not use abrasive or aggressive cleaning materials. These can damage the surface of the appliance!
- 1) Switch the appliance off by moving the main switch **6** on the underside of the handle **4** to the OFF position.
 - 2) Make sure that the appliance is not connected to the mains adapter **10**.
 - 3) Wipe off the appliance using a damp cloth. Do not use any solvents or cleansers which can damage plastic. If necessary, add a little mild detergent to the cloth and wipe off afterwards with fresh water. Make sure that no water or moisture can get into the appliance!
 - 4) To hygienically clean the attachments, wash them with water and a mild detergent and then rinse them with clear water.
 - 5) Allow everything to dry thoroughly before storing the appliance or putting it back into operation.

Storage


Store the cleaned appliance in the storage bag 16 in a dust-free and dry location out of direct sunlight. Charge the battery at regular intervals to ensure a long battery life.

Troubleshooting

Fault	Cause	Solution
The appliance is not working.	The appliance is switched off.	Set the main switch 6 to ON.
	The appliance is defective.	Contact the Customer Service department.
E1 appears on the display 2.	The deep discharge protection has been triggered. Battery is too weak.	Charge the appliance immediately.
E3 appears on the display 2.	The motor overload protection has been triggered.	Restart the appliance. Do not press too hard while massaging.

If the malfunction cannot be corrected using one of the suggestions above, or if you have found some other kind of malfunction, contact Customer Service.

Disposal

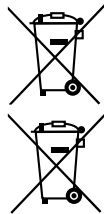The adjacent symbol of a crossed-out dustbin means that this appliance is subject to Directive 2012/19/EU. This directive states that this appliance may not be disposed of in the normal household waste at the end of its useful life, but must be taken to specially set-up collection locations, recycling depots or disposal companies.

The disposal is free of charge for the user. Protect the environment and dispose of this appliance properly.

If it is possible to do so without destroying the old appliance, remove the old batteries or rechargeable batteries before returning the appliance for disposal and take them to a separate collection point. In the case of permanently installed rechargeable batteries, you must indicate during disposal that the appliance contains a battery. The integrated rechargeable battery in this appliance cannot be removed for disposal.

Your local community or municipal authorities can provide information on how to dispose of the worn-out product.

The product and its packaging are recyclable. They are subject to an extended manufacturer responsibility and will be collected separately.

Technical specifications

Mains adapter 10

Manufacturers	DongGuan City GangQi Electronic Co.,Ltd. Privately Operated Industrial District, ShiYong, 532600 Heng Li Town, Dong Guan, Guangdong, P.R.CHINA Commercial registration number: 914419007929512928 Importer: KOMPERNASS HANDELS GMBH, BURGSTRASSE 21, 44867 BOCHUM, DEUTSCHLAND, Registergericht AG Bochum, Registernummer: HRB 4598
Model	GQ12-050200-ZG
Input voltage	100~240 V ~ (AC)
Input AC frequency	50/60 Hz
Output voltage	5.0 V $\equiv$ (DC)
Output current	2.0 A
Output power	10.0 W
Average operating efficiency	80.1%
No-load power consumption	0.06 W
Input current	max. 0.4 A
Efficiency class 6	
Polarity of the USB Port	P1  P4 
Safety transformer, short-circuit proof	
Switched-mode power supply	
Protection class	II /  (double insulation)
Nominal ambient temperature (ta)	45 °C
Protection type	IP20: Protection against ingress of solid objects with a diameter of more than 12.5 mm

Appliance	
Input voltage	5.0 V $\equiv$ (DC)
Input current	2.0 A
Protection class	III /  (Protection by low voltage)
Integrated battery (Li-Ion)	2500 mAh / 7.4 V $\equiv$ (DC) / 18.5 Wh
	detachable supply unit

Kompernass Handels GmbH warranty

Dear Customer,

This appliance has a 3-year warranty valid from the date of purchase. If this product has any faults, you, the buyer, have certain statutory rights. Your statutory rights are not restricted in any way by the warranty described below.

Warranty conditions

The warranty period starts on the date of purchase. Please keep your receipt in a safe place. This will be required as proof of purchase.

If any material or manufacturing fault occurs within three years of the date of purchase of the product, we will either repair or replace the product for you or refund the purchase price (at our discretion). This warranty service requires that you present the defective appliance and the proof of purchase (receipt) within the three-year warranty period, along with a brief written description of the fault and of when it occurred.

If the defect is covered by the warranty, your product will either be repaired or replaced by us. The repair or replacement of a product does not signify the beginning of a new warranty period.

Warranty period and statutory claims for defects

The warranty period is not prolonged by repairs effected under the warranty. This also applies to replaced and repaired components. Any damage and defects present at the time of purchase must be reported immediately after unpacking. Repairs carried out after expiry of the warranty period shall be subject to a fee.

Scope of the warranty

This appliance has been manufactured in accordance with strict quality guidelines and inspected meticulously prior to delivery.

The warranty covers material faults or production faults. The warranty does not extend to product parts subject to normal wear and tear or to fragile parts which could be considered as consumable parts such as switches, batteries or parts made of glass.

The warranty does not apply if the product has been damaged, improperly used or improperly maintained. The directions in the operating instructions for the product regarding proper use of the product are to be strictly followed. Uses and actions that are discouraged in the operating instructions or which are warned against must be avoided.

This product is intended solely for private use and not for commercial purposes. The warranty shall be deemed void in cases of misuse or improper handling, use of force and modifications / repairs which have not been carried out by one of our authorised Service centres.

Warranty claim procedure

To ensure quick processing of your case, please observe the following instructions:

- Please have the till receipt and the item number (IAN) 391896_2201 available as proof of purchase.
- You will find the item number on the type plate on the product, an engraving on the product, on the front page of the operating instructions (below left) or on the sticker on the rear or bottom of the product.
- If functional or other defects occur, please contact the service department listed either by **telephone** or by **e-mail**.
- You can return a defective product to us free of charge to the service address that will be provided to you. Ensure that you enclose the proof of purchase (till receipt) and information about what the defect is and when it occurred.

You can download these instructions along with many other manuals, product videos and installation software at www.lidl-service.com.

This QR code will take you directly to the Lidl service page (www.lidl-service.com) where you can open your operating instructions by entering the item number (IAN) 391896_2201.

Service

GB Service Great Britain
Tel.: 0800 404 7657
E-Mail: kompernass@lidl.co.uk

IE Service Ireland
Tel.: 1800 101010
E-Mail: kompernass@lidl.ie

IAN 391896_2201

Importer

Please note that the following address is not the service address. Please use the service address provided in the operating instructions.

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