*Bob
AND
Brad.*

User Manual

Q2

Massage Gun

Model: QL/MINI.Q2-D

Maintenance

The battery in this device is not designed to be removed or exchanged by users. If any problem with the power supply occurs, please contact us directly.

Do not try to open the device, repair it or fix it by taking the device apart yourself. Tampering with the device's mechanics in any way will result in an automatic expiration of the product's warranty.

Warnings and Safety Cautions

Please be aware of the following to ensure a safe use of the product.

- The battery cannot be removed from the product. Do not disassemble the product or the battery.
- Please check the massage gun and the charging cable before using them. If you notice damages, please report it to us or the retailer and do not use the product.
- Do not hang the product by the power cable.
- This product is not a medical product and is not intended to be used as a treatment for medical conditions.
- This product is only intended for healthy adults, seniors, children, or people with limited cognitive capacity should not use this product or can only use it under supervision.
- Please consult your physician in case you have concerns over your health or experience pain or discomfort while using the product.
- Avoid use the product on your body parts with little muscle coverage, like your head, shins, and back of your hands and feet.
- Avoid use the product on or around your sensitive body parts like neck, eyes, teeth, genitals, or breast implants.
- Do not use the product in any area where there is a risk of bone fracture.
- Use this product gently. Using this product with too much force or too high a speed might cause bruising. Check the massaged area frequently and stop using the product as soon as you experience any pain or discomfort.
- Be careful that your fingers, hair, or other body parts might get caught between the massage head and the shaft of the massage gun.

- Do not use the product in water or get it wet. Use the product only on dry and clean skin.
- Do not put the product in fire or heat it up in any way.
Do not charge, use, or store the product near heat sources or under extremely hot conditions.
- If the battery leaks or emits an odor, keep it away from fire.
- Do not use the product if the battery:
 - emits an odor;
 - gives off excessive heat;
 - shows signs of deformation, discoloration, or any of abnormalities.
- If the battery is being charged, stop charging immediately.

***THIS PRODUCT IS NOT A MEDICAL DEVICE.**

FCC Statement

This device complies with Part-15 of the FCC Rules.

Operation is subject to the following two conditions:

- (1) This device may not cause harmful interference.
- (2) This device must accept any interference received, including interference that may cause undesired operation.

FCC Statement NOTE

This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in residential installation. This equipment generates, uses, and can radiate radio frequency energy, and if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause any harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help with the supplied fastener.

DO NOT hang the unit by the power cord.

Disclaimer

Please contact your physician before using the Iliad and Iliad QIL/MINI.Q2-D if you have any physical health concerns or have any of the following conditions:

- ⓐ Metal implants like pins or plates in your body.
- ⓑ Recent joint replacements
- ⓒ Pregnancy
- ⓓ Intrauterine device
- ⓔ Diabetes
- ⓕ Epilepsy
- ⓖ Migraine
- ⓗ Hemistated disks
- ⓘ Spondylolisthesis (a vertebra slipping forward)

★THIS PRODUCT IS NOT A MEDICAL PRODUCT


IC Warning

This device contains licence-exempt transmitter(s)/receiver(s) that comply with Innovation, Science and Economic Development Canada's licence-exempt RSS(s). Operation is subject to the following two conditions:

- (1) This device may not cause interference.
- (2) This device must accept any interference, including interference that may cause undesired operation of the device.

Main Structure and Fittings

Removable and Replaceable Massage Head

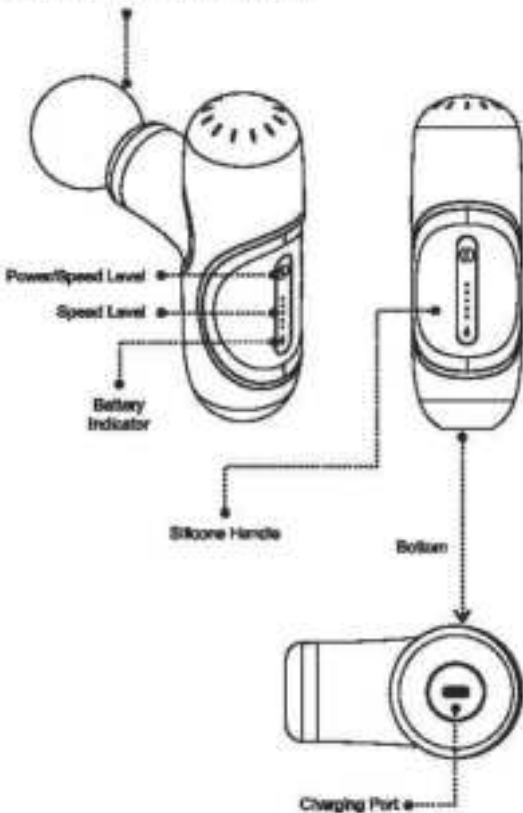**Ball Head**

For large muscle groups like neck, glutes, abdomin muscles and joint exercises.

**Air Cushion Head**

For relaxation of sensitive body parts.

**Ballistic Head**

For calves, lower back, trigger points, and small muscle areas like feet and wrists.

**U-Shaped Head**

For spine, neck, and shoulder. For sensitive muscles and soft tissues.

**Flat Head**

For all body parts.

**Rubber Blower**

Each massage head is pre-installed with a rubber blower.


If it is damaged, please replace it with the provided spare one.

Operation Steps

1. When the device is off, insert the desired massage head with pressure into the massage gun. To move the massage head, pull it out when the device is stopped.

Do not use the product when it is being charged.

Note: Position the U-shaped head for installation as shown in the figure below (top view) to prevent the risk of the massage head falling out while the massage gun is being used. When the buckle is inserted into the slot and a click is felt, the installation of the U-shaped head is in place.

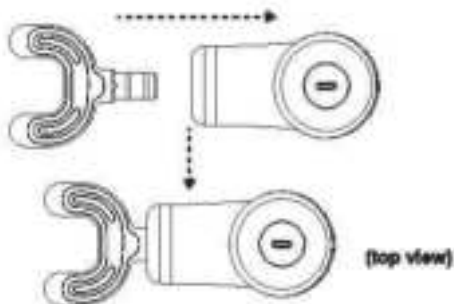2. Turning the device ON/OFF: Press and hold the " ① " power button for 2 seconds.

3. When the message gun is on, press the power button " ① " to choose a suitable message level and the level is indicated by the speed level indicator above and the power button.

To pause the message gun, press the power button for 2 seconds.

To resume the use, short press the power button.

Note

A. After 10 minutes of operation, the device automatically turns off.

B. LED light indicates the corresponding strike frequency level, the lowest at level 1 and the highest at level 5.

C. When the device is standby mode (0 Level), it will automatically shut down without operation for 20 seconds.

Battery Display

When the message gun is being used:

- If the battery indicator is blue, the device battery level is above 50%;
- If the battery indicator is red, the device battery level is between 30% and 20%;
- If the battery indicator flashes red, the device battery level is below 20%. Please charge your message gun.

When the message gun is being charged, the battery indicator is green and it turns off once the device is fully charged. If the green battery indicator turns off soon after the charging starts, it means the device battery is full and does not require to be charged.

Charge Massage Gun

1. Please charge the product fully before the first use. It takes about 3 hours for a full charge.

2. Connect the product to a power source (5V/2A) with the provided type C cable to charge it. To charge it from a power socket, an adapter is needed for the type C Cable (a regular phone adapter is enough). It will take longer time to charge the product if the power is below 5V/2A.

Use with 5V=2A
charger output

5V=2A

3. When a massage gun is being charged, the battery indicator on the handle flickers regularly, indicating the power level.

4. Once the battery indicator stops flickering and remains on, the massage gun is fully charged.

5. The battery can be recharged at your preference, please make sure to turn off the massage gun before charging it.

6. The massage gun does not work while being charged.

Technical Details

QL/MINI.Q2-D Mini Massage Gun

Speed: 1800 - 3000 r/min

Rated Power: 25W

Dimensions: 14.6x6.6x4.8cm/ 5.7" x3.3" x1.8"

Battery:

Rated Voltage: 11.1V=

Battery Capacity: 2500mAh

Charging Time: 210 minutes

Troubleshooting

PROBLEM	CAUSE	SOLUTION
Massage gun doesn't turn on.	Power is off.	Long press the Power Button to turn on the massage gun.
	Battery is too low.	Charge the massage gun.
	Massage gun shuts-off after 10 minutes.	Long press the Power Button to turn on the massage gun for another 10 minutes.
Massage gun shuts off automatically.	The massage gun's overload protection will switch off the massage gun automatically when too hot.	Wait for the massage gun to cool down before pressing it back on.
	Massage gun shuts off after 10 Minutes.	Long press the Power Button to turn on the massage gun for another 10 minutes.
Massage gun fails to run.	The massage head is inserted too tight.	Plug out the massage head a little, then power it on.
If you encounter any other problems, contact: support@gobanfitness.com		

User Guide for Gob and brad QU/MINI.Q2-D massage gun

Do not use the product when it is being charged

Guarantee

This product comes with a one-year warranty from the date of purchase. Warranty will be voided due to any of the following:

- Damaged due to improper usage or storage.
- Damaged by attempting to repair or disassemble the device.
- Period of warranty exceeded.

UK Importer

UK	REP
----	-----

EVATOST CONSULTING LTD
Office 101 32 Threadneedle Street, London,
United Kingdom, EC2R 8AY
contact@evatost.com

Unlike the Bob and Brad massage gun before, during, or after your workout for maximum benefits. This massage gun can help accelerate muscle tissue responsiveness by increasing hydration, range of motion and function to tight muscles.

Indication	Action	Speed	Pressure	Duration	Muscle Position	Movement
Warm Up / Activation	Activate and prime muscles prior to exercise. Reduce risk of injury.	Level 3	Light to Moderate	30 sec or less	Relaxed. Move along the muscle	Quick, high speed over area.
Mild Workout / Myofascial Mobility	Help muscle recovery and reactivation.	Level 4	Heavy (especially for large muscle groups)	8 - 12 sec	Target specific area.	Localized
Recovery	Relax muscles and aid recovery. Down-regulate or turn off muscles.	Level 1 or 2	Moderate to Heavy	90 sec or more	Move along the muscle.	Slow movement along the muscle
Pain Modulation	Reduction in pain to increase movement.	Level 1	Light	10 to 30 sec	Random. Localized area.	Over particular area.
Body Awareness	Improve movement or mechanics of a joint. Activation of a muscle	All Levels	Light to Moderate	10 to 30 sec	Relaxed or with movement	Over specific target area or movement over area.

CHEST

2 Minutes

ATTACHMENT: BALL HEAD
PROTOCOL: Sweep from A to B on each side for 1 minute
DURATION: 2 minutes

ABS

2 Minutes

ATTACHMENT: BALL HEAD
PROTOCOL: Sweep from A to B on each side for 1 minute
DURATION: 2 minutes

UPPER BACK

4 Minutes

ATTACHMENT: BALL HEAD
STEP 1: Upper back
PROTOCOL: Sweep from A to B on each side for 1 minute
DURATION: 2 minutes

UPPER BACK

4 Minutes

STEP 2: Lower back
PROTOCOL: Sweep from A to B on each side for 1 minute
DURATION: 2 minutes

SHOULDERS

2 Minutes

ATTACHMENT: BALL HEAD
PROTOCOL: Sweep from A to B on each side for 1 minute
DURATION: 2 minutes

BICEPS

2 Minutes

ATTACHMENT: BALL HEAD
 PROTOCOL: Sweep from A to B on each side for 1 minute
 DURATION: 2 minutes

TRICEPS

4 Minutes

ATTACHMENT: BALL HEAD
 PROTOCOL: Sweep from A to B on each side for 2 minutes
 DURATION: 4 minutes

FOREARM

2 Minutes

ATTACHMENT: BALL HEAD
 PROTOCOL: Sweep from A to B on each side for 1 minute
 DURATION: 2 minutes

HANDS

1 Minute

ATTACHMENT: BALL HEAD
 PROTOCOL: Sweep from A to B on each side for 30 seconds
 DURATION: 1 minutes

HIP FLEXORS

2 Minutes

ATTACHMENT: BALL HEAD
 PROTOCOL: Sweep from A to B on each side for 1 minute
 DURATION: 2 minutes

GLUTES

4 Minutes


ATTACHMENT: BALL HEAD
 PROTOCOL: Sweep from A to B
 on each side for 2 minutes
 DURATION: 4 minutes

HAMSTRINGS

2 Minutes

ATTACHMENT: BALL HEAD
 PROTOCOL: Sweep from A to B
 on each side for 1 minute
 DURATION: 2 minutes

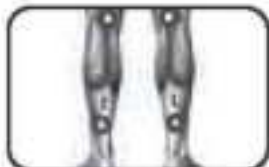CALVES

2 Minutes

ATTACHMENT: BALL HEAD
 PROTOCOL: Sweep from A to B
 on each side for 1 minute
 DURATION: 2 minutes

SHINS

1 Minute

ATTACHMENT: BALL HEAD
 PROTOCOL: Sweep from A to B
 on each side for 30 seconds
 DURATION: 1 minute

FEET

4 Minutes

ATTACHMENT: BALL HEAD
 PROTOCOL: Sweep from A to B
 on each side for 2 minutes
 DURATION: 4 minutes

WEEE Instructions

Correct Disposal of this product

This marking indicates that this product should not be disposed with other household wastes throughout the EU. To prevent possible harm to the environment or human health from uncontrolled waste disposal, recycle it responsibly to promote the sustainable reuse of material resources. To return your used device, please use the return and collection systems or contact the retailer where the product was purchased. They can take this product for environmental safe recycling.

Battery Instructions

- ❶ Do not submerge the product into water or get it wet.
- ❷ Do not charge the product near an open flame or in any extremely hot conditions.
- ❸ Do not use or store the product near heat sources like furnaces or heaters.
- ❹ If the battery leaks or emits an odor, remove it immediately from close proximity to an open flame.
- ❺ Do not put the product into a fire or heat it up in any way.
- ❻ Do not disassemble the battery in any way.
- ❼ Do not use the product if, while the battery is in use or being charged, it emits an odor, gives off excessive heat, shows signs of deformation, discoloration, or any other abnormalities. If the battery is in use or being charged, remove it from the apparatus or charger immediately.

RoHS Declaration

Our Product is compliant with RoHS 2.0 (2011/65/EU Directive) regarding restriction of hazardous substances, the raw materials and components of our product meet the following toxic and hazardous substance limits.

Cadmium (Cd): < 100 ppm	Lead (Pb): < 1000 ppm
Mercury (Hg): < 1000 ppm	Hazardous Substances (SVHC): < 1000 ppm
Polychlorinated Biphenyls (PCB): < 1000 ppm	Polychlorinated Biphenyl Esters (PCBE): < 1000 ppm
Hexo-Chlorocyclohexane (HCHC): < 1000 ppm	Hexo-Chloro-cyclohexane (HCHC): < 1000 ppm
Chloro-biphenyls (CBP): < 1000 ppm	Chlorobiphenyl esters (CBPE): < 1000 ppm

A Dedicated Team Just For You

Every Bob and Brad product is backed by a 1-year warranty and an award-winning customer service team. If there are any quality problems with the product you ordered, we will try our utmost to help you get the problems resolved, as well as offer a replacement, return, or refund service under the warranty. Simply email us the problem, and we'll take care of it.

Troubleshoot

Help you solve any issues with your device.

Product Knowledge

Answer questions about your device and best practices.

Initiate Returns

Provide instructions for replacements, returns and refunds.

Place New Orders

Help with placing orders online or via authorized resellers.

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EVATOST CONSULTING LTD
Office 101 32 Threadneedle Street, London,
United Kingdom, EC2R 8AY
contact@evstool.com

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GLUTES

4 Minutes


ATTACHMENT: BALL HEAD
PROTOCOL: Sweep from A to B
on each side for 2 minutes
DURATION: 4 minutes

HAMSTRINGS

2 Minutes

ATTACHMENT: BALL HEAD
PROTOCOL: Sweep from A to B
on each side for 1 minute
DURATION: 2 minutes

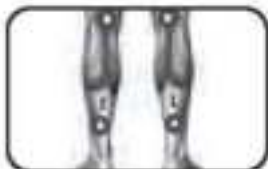CALVES

2 Minutes

ATTACHMENT: BALL HEAD
PROTOCOL: Sweep from A to B
on each side for 1 minute
DURATION: 2 minutes

SHINS

1 Minute

ATTACHMENT: BALL HEAD
PROTOCOL: Sweep from A to B
on each side for 30 seconds
DURATION: 1 minutes

FEET

4 Minutes

ATTACHMENT: BALL HEAD
PROTOCOL: Sweep from A to B
on each side for 2 minutes
DURATION: 4 minutes

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EC Declaration of Conformity

We (UK Importer is responsible for this declaration)

EVATOST CONSULTING LTD

(Company name)

Office 101 32 Threadneedle Street, London, United Kingdom, EC2R 8AY

(Company address)

declare under our sole responsibility that the product

Tradename: BOB AND BRAD

Equipment: Massage Gun

Model No. : QL/MINI.Q2-D

ASIN: B06MBYSFC7/B08XPGTZLL/B08B6KQD96X

to which this declaration relates is in conformity with the essential requirement set out in the Council Directive on the Approximation of the laws of the Member States relating to LVD Directive(2014/35/EU) & EMC Directive(2014/53/EU) & RoHS(2002/95/EC) product is responsible to affix CE marking, the following standards were applied:

EN 60335-1:2012+A11:2014+A13:2017

EN 60335-2-32:2003+A1:2008+A2:2015

EN 62233:2008

EN 55014-1:2017

EN 55014-2:2015

EN 61000-3-2:2014

EN 61000-3-3:2013

Full Name: Taozdul

Position: Manager

Signature:

Date: 2021.5.6

EC Declaration of Conformity

We (EU Importer is responsible for this declaration)

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(Company address)

declare under our sole responsibility that the product

Tradename: BOB AND BRAD

Equipment: Massage Gun

Model No. : QL/MINLO2-D

ABIN: B08M8YSFC7/B08XPGTZLL/B08BGKD96X

to which this declaration relates is in conformity with the essential
requirement set out in the Council Directive on the Approximation of the
laws of the Member States relating to LVD Directive(2014/35/EU) &
EMC Directive(2014/53/EU) &RoHS(2002/95/EC) product is
responsible to affix CE marking, the following standards were applied:

EN 60335-1:2012+A11:2014+A13:2017

EN 60335-2-32:2003+A1:2008+A2:2015

EN 62233:2008

EN 55014-1:2017

EN 55014-2:2015

EN 61000-3-2:2014

EN 61000-3-3:2013

Full Name: Tacodru Position: Manager

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EN 61000-3-2:2014

EN 61000-3-3:2013

Full Name: Tacórrul Position: Leiter

Signature: _____ Date: 2021.5.6