



**Affordable. Reliable. Home Improvement.**

**COMMODE**  
**MODEL:FBL760639**

# VEVOR

Affordable. Reliable. Home Improvement.

## COMMODE

**MODEL:FBL760639**



This is the original instruction, please read all manual instructions carefully before operating. VEVOR reserves a clear interpretation of our user manual. The appearance of the product shall be subject to the product you received. Please forgive us that we won't inform you again if there are any technology or software updates on our product.

## **WARNING:**

Read the manual completely to avoid harm, product damage, or warranty avoidance. Seek professional help if you encounter discomfort, pain or worsening symptoms. If necessary, use supervision in order to assemble and use this product correctly and safely. Do not exceed the weight limit.

Regularly check for secure assembly and contact support if in doubt. Ensure the product's legs are even and locked before use, and discontinue usage if there are signs of wear, damage, or missing parts. The product is to be used on flat, stable surfaces. Do not disassemble if not necessary, as it might compromise the longevity of the product. Do not modify the product unless specified in the manual.

## **Physiotherapist**

1. Before first use, adjust the height and the width of the rail to your toilet.
2. Back up to the toilet with your assistive device (if available) until the back of your legs is touching the toilet.
3. Reach back for the rail, hold firmly, and slowly lower yourself into a seated position on the toilet.
4. To return to a standing position, place your hands on the rails firmly and push yourself up into a standing position.
5. Once you are standing, gradually let go of the rails, and only then, grab the handle of your assistive device (if available).

## **IMPORTANT**









1. Don't exceed the maximum weight capacity.
2. Use on stable surfaces only.
3. Regularly inspect for signs of wear, damage, or loose parts.
4. Clean regularly using a mild detergent and warm water, and rinse off any soap residue. After cleaning, dry thoroughly.
5. Before use, confirm with your healthcare provider that the product meets your specific needs. The following guidelines are not a substitute for medical advice.

## **SAVE THESE INSTRUCTION**

## PART LIST

When unpacking this product, check to make sure the parts listed below are included and carefully inspect them for any damage that may have occurred

during transit. If any part is missing or damaged, do not attempt to assemble or use the product.


NO	PICTURE	NAME	QTY	NO	PICTURE	NAME	QTY
①		Main frame	1	⑤		Leg B	1
②		Backrest	1	⑥		Suction cup leg A	1
③		Handrail tube	1	⑦		Suction cup leg B	1
④		Leg A	1	⑧		Bucket	1

## PRODUCT PAREMETERS

<b>Model</b>	<b>FBL760639</b>
<b>Main Materials</b>	Aluminum
<b>Seat Height</b>	480-570mm/18.9"-22.4"
<b>Seat Rated Load</b>	227kg/500lbs

## INSTALLATION AND OPERATION INSTRUCTIONS

### 1. Installation Procedures

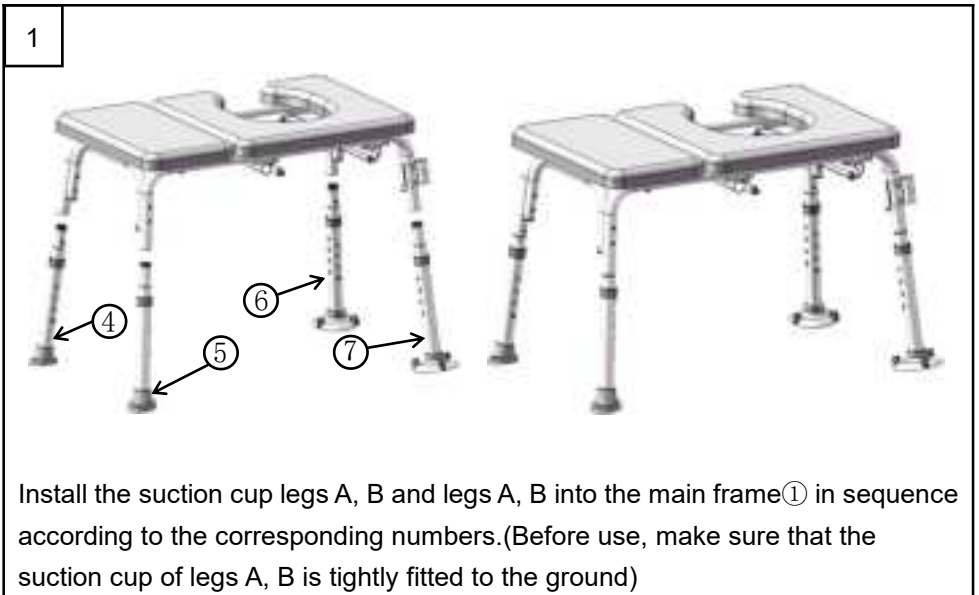


2



Press the buttons on the handrail assembly and insert them into the 8-shaped tube of the main frame until the buttons pop out.

3

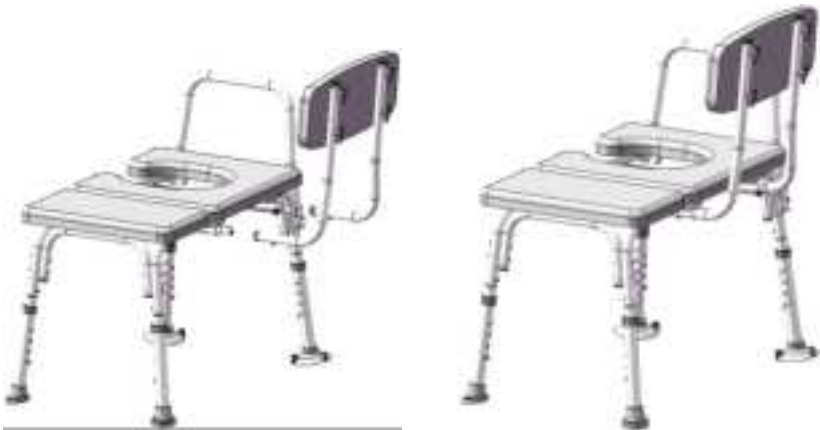

Press the buttons on the backrest assembly and insert it into the vertical tube under the main frame until the buttons pop out.

4



Place the bucket into the seat tube from the front of the main frame.  
As shown in the figure, the product assembly is completed.

## 2.Instruction

1



Press the button and pull out the leg  
to adjust the height (8 levels of height adjustment)

## **Zhongshan Bliss Medical Instrument Co., Ltd**

**Address:** 3-6F & 7F-1, No. 63, Tongxing East Road, Yumin Community  
Xiaolan, Zhongshan, Guangdong, 528414, CHINA

<b>UK</b>	<b>REP</b>
-----------	------------

YH CONSULTING LIMITED. C/O YH Consulting  
Limited Office 147, Centurion House, London  
Road, Staines-upon-Thames, Surrey, TW18 4AX

<b>EC</b>	<b>REP</b>
-----------	------------

E-CrossStu GmbH  
Mainzer Landstr.69,  
60329 Frankfurt am Main.

