




1.If you feel uncomfortable when using this product, please stop using it • Product Accessories immediately and consulting a doctor. 2. When using, please do not make metal substances (such as belts, bracelets,  $watches, necklaces, etc.) \, in \, contact \, with \, the \, pad.$ 3.It is recommended that the same part of the using time not more than 30minutes, please do not paste the patch on the body for a long time.  $\ensuremath{\mathsf{4.When}}$  using, please keep the skin dry, and ensure there is no scar on the position that needs to be massged . 5. Please do not use this product during intense exercise. 6.Do not put this product in high temperature, high humidity and dusty 7. Do not disassemble, repair or modify the product without permission. 8.Do not dispose of the product and waste batteries. 9.If you need to move or remove the product, please be sure to turn off the 10. The elderly must use the product in the company of others. 11.Recommended age:16 years old. 12. May appear phenomenon of rapid heartbeat in the first use, please remove the pad; wait for the heartbeat normal and then use again, if still appear this phenomenon, please stop using. Massage Patch X2 13. The product package does not contain a power adapter, charging cable of Controller X 2 this product is TYPE-C. 14. Charging box and the controller is not waterproof, please do not wash, if


Charging Box X1

the surface is dirty, please use alcohol cotton to wipe it.

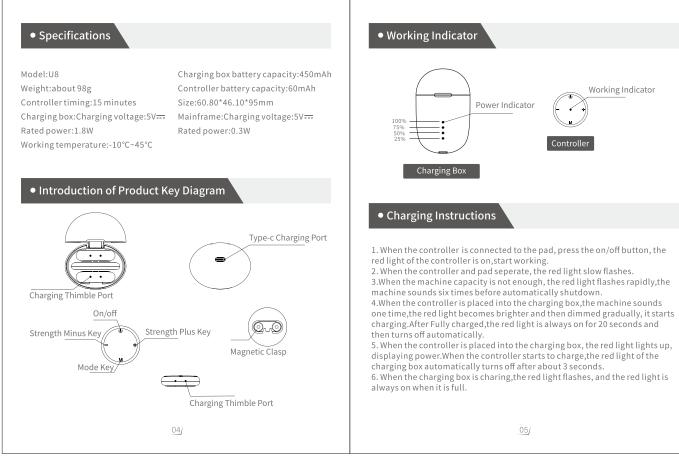


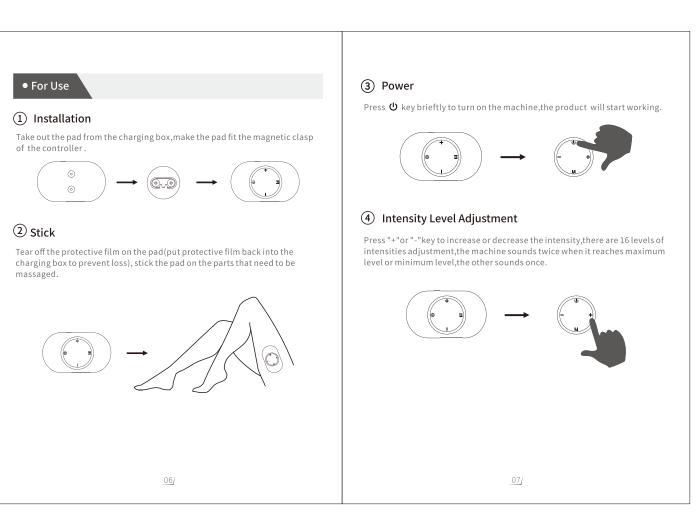


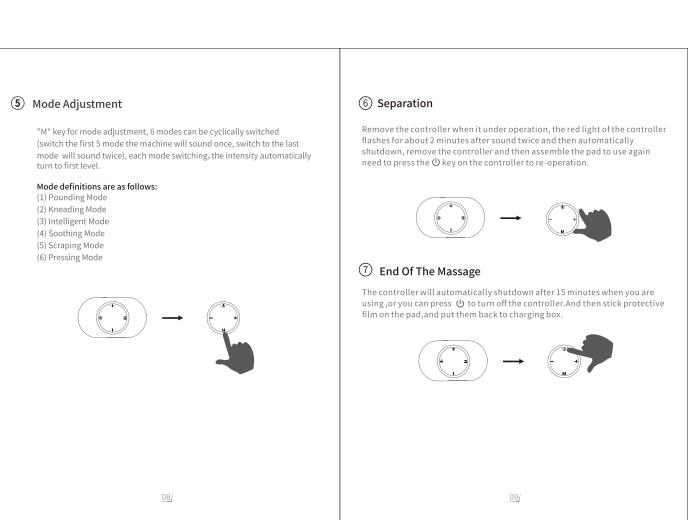

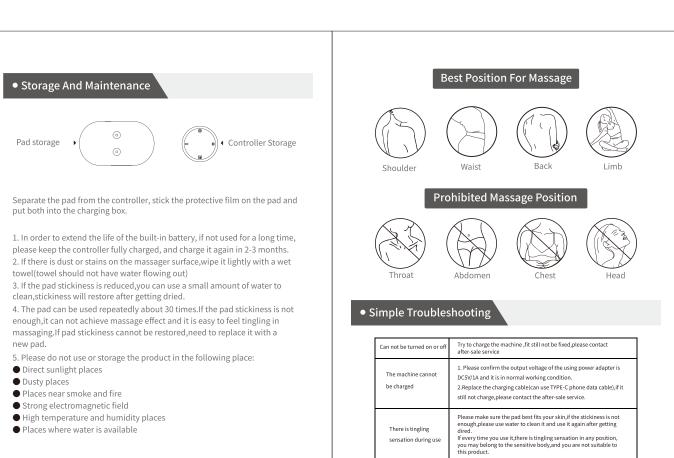

11/

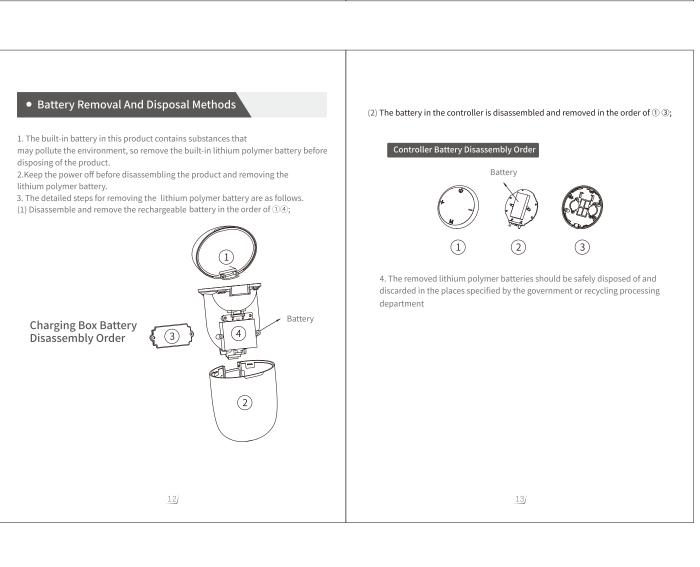

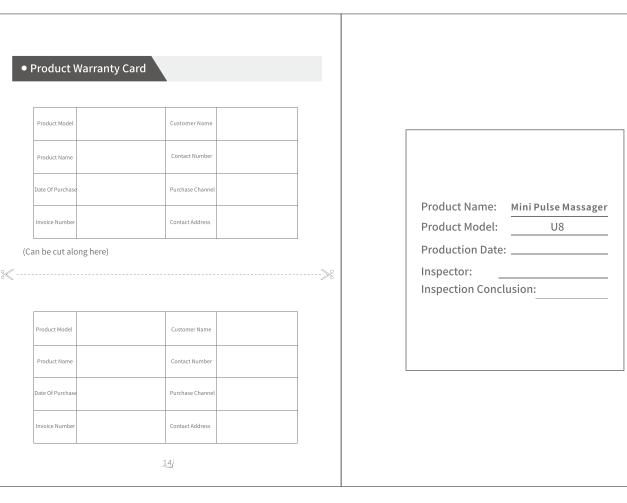

PDF