




FLOSSY®

## USER MANUAL

Portable Oral Irrigator

Model Name: Flossy

INSTRUCTION (EN) >>

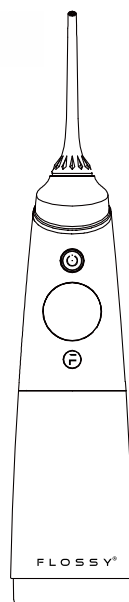

Thankyou for choosing Flossy Oral Irrigator

If you have any problems or questions, please feel free to contact us by [info@owlyvision.com](mailto:info@owlyvision.com) . We will try our best to resolve it.

Thank you for your time and wish you the best.

Note: Please choose the mode which is suitable for you firstly, then put the nozzle into your mouth and turn on the unit. For the first-time user and people who have sensitive teeth, it's suggested to choose the "soft" mode or "DIY" mode at first.

The logo features a stylized 'F' inside a circle, with the word 'FLOSSY' in a sans-serif font below it, followed by a registered trademark symbol (®).

FLOSSY®

## Contents

EN

Contents .....	01
Dangerous .....	02
Product Accessories .....	04
Operating Procedures .....	05
Instruction for the modes .....	06
Operating Step .....	07
Recommended Techniques .....	08
Disassembly Example .....	09
Battery Maintenance and Handling .....	09
Problems and solutions .....	10
Comment Problems and solutions .....	11
After-sales Services .....	12


## READ ALL INSTRUCTIONS BEFORE

### DANGERS:


#### To reduce the risk of electrocution:

- ◆ Do not handle charger with wet hands.
- ◆ Do not immerse in water or other liquid.
- ◆ Check the charger cord for damage before use.
- ◆ Keep charger and cord away from heated surfaces.
- ◆ Use this product only as indicated in these instructions or as recommended by your dental professional.
- ◆ Do not reach for a product that has fallen into water. Unplug immediately.
- ◆ Do not plug this device into a voltage system which is different from specific voltage range.
- ◆ Do not use if you have trauma on your tongue or in your mouth.
- ◆ If you are in dental treatment, you should consult your dentist before using this instrument.
- ◆ Do not direct water flow under the tongue, into the ear, nose or other delicate organs.
- ◆ It's normal phenomenon that caused gum bleeding for the first time use. It is recommended to start from a lower pressure level.
- ◆ Never be left unattended when plugged in. Close supervision is necessary when it used by, on, or near children or invalids.

**Product Description**

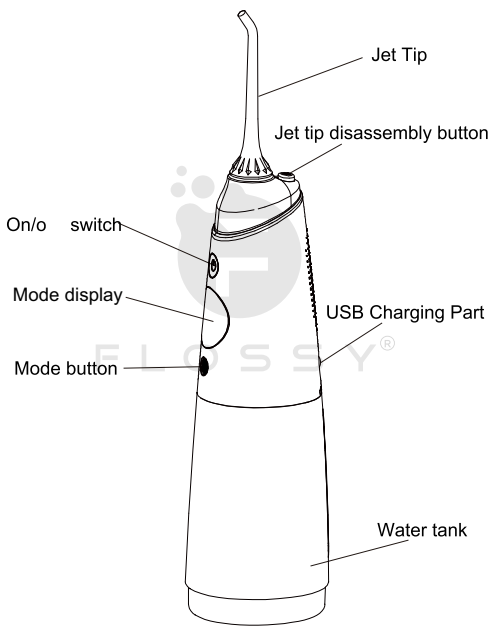

## Product Accessories

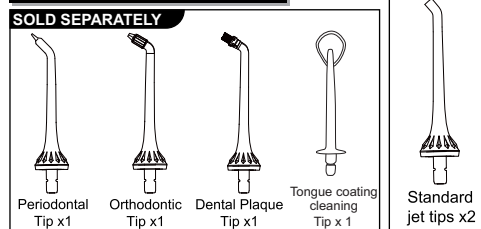

**Standard Tips:** For teeth daily cleaning, ultra-fine pulsed water flow, it massages gums and get all the food particles out of your teeth.

**Periodontal Tip:** Deeply rinsing into the interdental fissure and subgingival site. Relieve periodontitis, gingivitis, teeth shrinking, reduce the thickness of pocket.

**Orthodontic Tip:** Suitable for the person with the false tooth, wearing orthotics, dental braces, teeth socket. It removes hard-to-reach plaque around orthodontic brackets and wires.


**Dental Plaque Tip:** For implants, dental bridge, false tooth crowns. It has 4 thin tufts of bristles to access stubborn plaque around implants, crowns, bridges and other areas.



Charging Cable x1

Replacement Nozzle / Accessories can be purchased from  
**OWLYVISION on Bol.com**

## Instruction for the modes

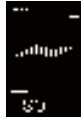

Soft: for people who have sensitive teeth and first-time-user. Pressure range: 20-50 PSI .



Medium: for daily cleaning. Pressure range: 50-80 PSI



Strong: for deep cleaning. Pressure range: 70-110 PSI



Pulse: massage for your teeth and gum. Pressure range: 70-100PSI




DIY Mode: Customize the water pressure / speed as you need. Long press the "Mode" button to enter the DIY mode, water pressure will change from min to max with pattern in OLED display. Press the "Mode" button again to choose the water pressure as you need. Pressure range: 8-110 PSI.

## Operating Procedures


### Usage instruction

Please ensure the machine is fully charged battery before first time usage, do not use it while charging.

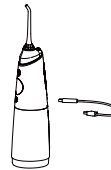


### Installation and Disassembly of Jet Tip

Press the release button and insert the nozzle into spray head hold until you hear a "click" ring.

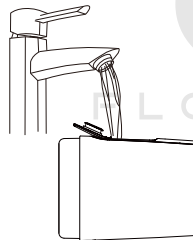

Open the water tank lid and fill with the water in it. The choose the mode which suitable for you.

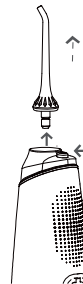

After using, press the On/ Off button to turn off the unit, press the nozzle button and take out the nozzle.




## Operating Step

- ◆ Press the mode button continuously, the mode switches among the sequence: SOFT-MEDIUM-STRONG-PULSE-DIY; This device has a memory function. So in the off state, you can also select the mode suitable for the user and then turn it on again.
- ◆ Irrespective of the mode, after operating for 120 seconds the product will automatically switch off.
- ◆ When power is insufficient, the indicated light will flash until power off. During charging, this product cannot be used.
- ◆ The charger is suitable for: AC 100-240VAC, 50/60Hz, output 5V  $\equiv$  1A
- ◆ Keep the device vertical for the best effect, rotate the jet tip dial in order to adjust the direction of the water flow.
- ◆ It's normal phenomenon that caused gum bleeding for the first time use. It is recommended to start from a lower pressure level.




## Recommended Techniques

Holding the device at downward, aim the jet tip to the teeth by rotating the tip base before switching on the device.

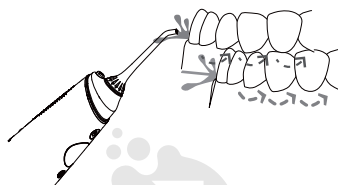


When the water flow and the gums are at a 90 degree angle, slightly close your mouth in order to avoid the water spraying out. In order to achieve the best results. Please start from molar area to the incisors. Slide the water flow along the gums, as well as stay a while between teeth, until the area around the gums and the spaces between the teeth are fully cleaned.

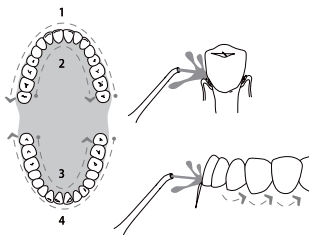


## Disassembly Example

### Following Usage

Please switch off the power, as well as empty the water tank.

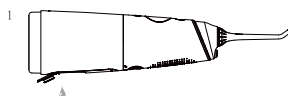

Using hands as indicated in Diagram 2 use your left hand to hold the product body, and your right hand to hold the water tank. At the same time rotate downwards in an anti-clockwise direction, in order to remove the water tank.

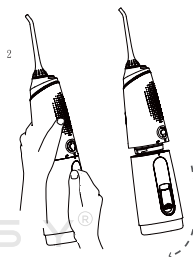

## Battery Maintenance And Handling

### Battery Maintenance

Prior to initial usage the device must be fully charged, for approximately 4 hours. If to be stored for a long-period of time, please fully charge the device in advance. In order to maximize the usage life of the battery, please promptly recharge when power is insufficient.

## Battery Handling


The device includes a lithium battery which cannot be changed. At the end of the battery life, please dispose in a local recycling centre.

## Comment Problems and solutions

Malfunction	Cause	Handling method
Water leakage between the jet tip and screw part.	1.The jet tip cannot be installed properly. 2.The seal ring socket of the jet tip is damaged.	1.Insert the jet tip into the screw part until you hear a clicking sound. 2.Contact the retailer for repair and install the jet tip properly
Inadequate pressure	Seal ring is loose or damaged	Replace the seal ring
Inadequate pressure	Low voltage usage	Use after charging
Device cannot be switched on	Low voltage usage	Use after charging
Device is not spraying water	Filtration net is clogged	Wet and clean the Filtration net
Water tank sides are leaking water	Air entry hole is producing water	Standard phenomenon high air pressure principle

If the above methods cannot resolve these problems, please consult with the retailer or proceed to the designated after-sales repair site for repairs.

Any question please contact us at

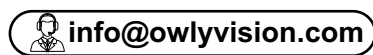

Thank you.

