Xiaomi Water Flosser 2

User Manual · 01

Manual de usuario de Xiaomi

Water Flosser 2 · 17

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Water Flosser 2 · 33

Xiaomi Munddusche 2

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Flosser 2 · 65

Xiaomi Waterflosser 2

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Xiaomi 2 · 97

Manual do utilizador do Irrigador

oral Xiaomi 2 · 113

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příručka · 177

Read this manual carefully before use, and retain it for future reference.

Safety Instructions

WARNINGS:

- Do not use this product unless you understand the operating procedures and safety instructions to prevent injuries resulting from improper operation.
- Users with periodontal disease or who have undergone oral surgery in the past six months should use this product under the guidance of a dentist.
- This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.
- Children should be supervised to ensure that they do not play with the appliance.
- Unplug the power adapter and charging cable before use, cleaning or maintenance.
- Do not expose the appliance or battery to excessive temperatures.
- This product contains a built-in rechargeable battery. Do not throw it into fire, heat or charge, use, or place it in a high-temperature environment, otherwise it will cause overheating, fire or explosion.
- Do not plug and unplug the power adapter and charging cable with wet hands, as it may result in electric shock.
- If this product is not in use for an extended period of time, please unplug the power adapter and charging cable, otherwise it may lead to electric shock or electrical fire caused by insulation deterioration.
- If this product is not used for an extended period of time, please fully charge it before storage and subsequently maintain a full charge every six months.
- This product contains a non-replaceable battery.
- Be aware of the risk of terminals of the battery-operated appliance or battery being short circuited by metal objects.

- This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved.
- Cleaning and user maintenance shall not be made by children without supervision.
- Do not let children play with it.
- This appliance contains a battery that is non-replaceable.
- Use only with <A520-050100U-EU1> power supply unit.


⚠ Cautions:

- In case of oral or lingual wounds, please seek medical advice before using this product to avoid unnecessary injury.
- Be sure to empty the water tank and allow it to idle for a few seconds after use to drain any remaining water, thus preventing the growth of bacteria, unpleasant odors, and malfunctions.

⚠ Note:

- Illustrations of product and accessories in the user manual are for reference purposes only. Actual product and functions may vary due to product enhancements.

Product Overview

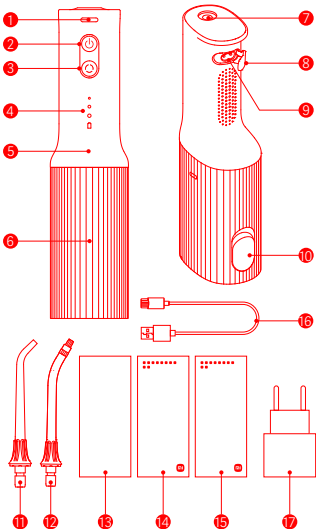- ① Tip Release Button
- ② Power Button
- ③ Mode Button
- ④ Indicators
 - Gentle mode
 - Standard mode
 - Powerful mode
 - Battery level
- ⑤ Flosser Body
- ⑥ Water Tank

- ⑦ Tip Hole
- ⑧ Rubber Plug
- ⑨ Charging Port
- ⑩ Water Tank Cover
- ⑪ Standard Tip
- ⑫ Orthodontic Tip
- ⑬ Quick Start Guide
- ⑭ User Manual
- ⑮ Warranty Notice
- ⑯ Charging Cable
- ⑰ Power Adapter

How to Use

Cautions:

- The product undergoes testing before leaving the factory. Please rest assured that it is normal for some of the water for testing to remain inside, and feel free to use it.
- When using the product for the first time, it is common for the gums to experience itching or bleeding, which can be attributed to the gums not being accustomed to the stimulation from the flosser. Typically, this condition will resolve within a period of 1–2 weeks.
- For first-time users, it is recommended to start with the gentle mode and use the standard tip. Maintain this routine for a week before switching to another mode. If any discomfort persists for over 2 weeks, stop using this product and seek advice from your dentist.

Charging

- The product should be fully charged before its first use.
- The provided charging cable is a specialized accessory designed for the product. It requires approximately 4 hours to achieve a full charge, please keep it properly.
- Before charging, ensure that both the product and your hands are dry. Insert the charging cable into the charging port, connect the power adapter, and activate the power supply to initiate charging.

Note:


- The product should be charged within the temperature range of 0°C to 35°C at room temperature. Charging is not possible when the temperature falls below 0°C.

Battery level indicator descriptions

Flosser status	Battery level	Battery level indicator
Operating	40%-100%	Keeping white
	20%-40%	Keeping orange
	10%-20%	Quickly blinking orange
	< 10%	Quickly blinking orange (unable to turn on)
Charging	100%	Keeping white
	40%-100%	Blinking white
	< 40%	Blinking orange

Filling the water tank

- Hold the flosser in a horizontal position, open the water tank cover, and fill it with an appropriate amount of water, then close the water tank cover.

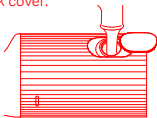
⚠ Cautions:

- It is recommended to use warm water not exceeding 40°C or running tap water.
- The flosser should not be filled with toothpaste, salt, solid mouthwash tablets, mouthwash containing isopropyl myristate, or other chemicals that may clog or damage the flosser.

Flosser tips

- Standard tip: Featuring a very thin water jet, this standard tip eliminates food debris lodged deep within interdental spaces and fissures, thereby preventing tooth decay and dental calculus. It is suitable for daily oral cleaning as well as for individuals with deep cleaning needs.
- Orthodontic tip: This tip is capable of cleaning intricate areas, including those with dental braces, which is suitable for individuals undergoing orthodontic treatment.

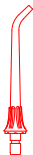Standard Tip

Orthodontic Tip

Note:


- If the top of the tip is deformed, please replace it with a new one immediately to avoid potential product failure. It is recommended to replace the tip every 3–6 months.

Assembly & disassembly of the flosser tip

- Insert the tip vertically into the tip hole, and it is installed in place when you hear a "click" sound.
- Detach the tip by simultaneously pressing the tip release button and pulling the tip upwards.

⚠ Caution:

Do not press the tip release button while in use, as it may result in ejection of the tip and potential accidental injury.

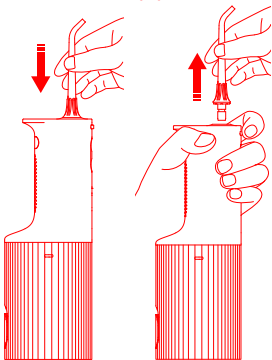Turning on/off

- Press the power button  to turn on/off the flosser.

Notes:


- The flosser is provided with a 2-minute timing function. The flosser will automatically turn off after a 2-minute continuous operation. To continue using the flosser, simply press the power button again.
- The standard mode allows for 50 seconds of uninterrupted cleaning with a full tank of water, while the duration of usage may vary for different modes.

Switching modes

- Press the mode button  to cycle through gentle, standard, and powerful modes. The indicator corresponding to the mode on the flosser body will light up once the mode is enabled.

Note:

- The flosser has a mode memory function and can remember the mode used before it was turned off. Once turned back on, the flosser will automatically restore the mode.

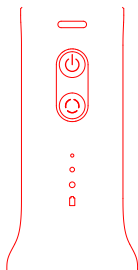
- Gentle mode:
Gentle water pressure for delicate cleaning of your teeth.
- Standard mode:
Standard water pressure for effective cleaning.
- Powerful mode:
Powerful water pressure for deep cleaning.

How to use the flosser

1. Hold the flosser in a vertical position and lean over the wash basin. Press the mode button to select your desired mode, then insert the tip into your mouth, targeting the interdental space or the gum line. Press the power button to turn on the flosser. Keep your lips closed to avoid water splashing, while maintaining a small gap to facilitate the smooth flow of water from your mouth into the wash basin.
2. Keep the water jet to maintain a 90° angle with the gum line. It is recommended to initiate the cleaning from the posterior teeth and gradually guide the tip along the gum line towards the front teeth, pausing briefly at the interdental space.
3. After each use, it is recommended to empty the water tank and allow the flosser to idle for a few seconds in order to drain the remaining water inside the flosser. Open the water tank cover for air drying or pull the water tank downward and place both the flosser body and water tank in a well-ventilated area. Once the inside of the water tank is dry, push it back up into position and close the water tank cover.

Cautions:


- Be sure to turn on the flosser with the tip in the mouth.
- Do not rinse directly against the gums.
- Do not apply excessive pressure on the teeth or gums with the tip, as this may trigger the motor jam protection under the powerful mode. When the protection is activated, the flosser will automatically pause and all mode indicators will blink white for 10 seconds before turning off. To restore normal functionality, simply press the power button to restart the flosser.

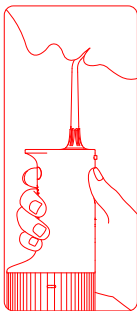Care & Maintenance

Daily cleaning

1. After usage, turn on the flosser and allow it to idle for a few seconds to facilitate the drainage of any residual water inside the flosser.
2. Hold the flosser body with one hand, while simultaneously pulling down the water tank with the other hand.
3. Rinse the water tank with water or clean it using a damp cloth or brush.
4. After cleaning, wipe the flosser body and water tank dry and store in a cool location for drying.

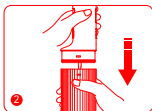Notes:

- The water temperature for cleaning should not exceed 40°C.
- Please do not use corrosive liquid, thinner, alcohol, etc. to clean the flosser.
- The flosser body and water tank can be wiped clean with a damp cloth or rinsed with water, ensuring thorough drying afterwards.
- The interior of the water tank can be cleaned using a brush in case of any scale formation.
- Please ensure that both the interior of the water tank and the exterior surface of the flosser body are dried before reinstalling the water tank into the flosser body.
- The flosser should be cleaned and stored after drying if it is not used for an extended period of time.

Regular maintenance

- Fill the water tank with warm water that does not exceed a temperature of 40°C. Add one tablespoon of white vinegar and turn on the flosser until the water tank is completely emptied. Once the water tank has been emptied, continue operating the flosser for approximately one minute until it automatically turns off. At this point, refill the water tank with warm water at a temperature not higher than 40°C and turn on the flosser again until the water tank is emptied again.

Note: To ensure optimal performance of the product, it is recommended to perform regular maintenance every 1–3 months.

Environmental Protection



When the product reaches the end of its service life, please remove, dispose of and recycle the built-in battery according to local laws and regulations of the country or region where the product is used. If the battery cannot be removed, dispose of it in strict accordance with the regulations on classified domestic waste disposal in your area.

⚠ WARNINGS:



- Do not expose the battery to fire to prevent it from exploding or releasing toxic and hazardous substances.
- Do not short-circuit the battery to prevent it from combusting and causing a fire.
- Do not disassemble the product unless you want to dispose of it, otherwise, it will be damaged.
- The battery must be removed from the appliance before it is scrapped.
- The appliance must be disconnected from the supply mains when removing the battery.
- The battery is to be disposed of safely.

Removing the rechargeable battery



1. Remove the water tank, unplug the water pipe, carefully pry open the 3 screw plugs at the bottom of the flosser, unscrew the 3 screws in order, and then remove the bottom cover of the flosser.

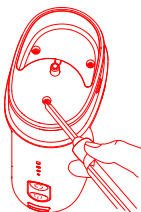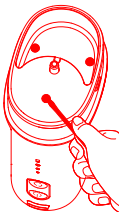
2. After removing the bottom cover, proceed to unscrew the 2 screws at the bottom of the flosser body and carefully pry open the top cover of the flosser.

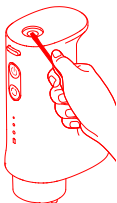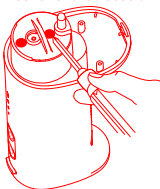3. Unscrew the 3 screws on the top of the flosser, remove the tip release button, pull out the charging port rubber plug, and pry out the inner mechanism.

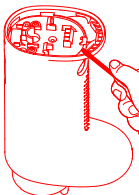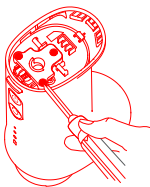4. Cut the 3 connecting wires (one by one) between the battery and the main board of the flosser body, and remove the battery.

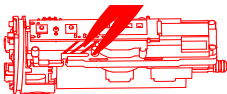⚠ Warnings:

- Please exercise caution when removing the battery to prevent any potential short-circuiting between the positive and negative terminals of the battery.
- The product shouldn't be connected to power once the battery has been removed.

Troubleshooting

Errors	Possible Causes	Solutions
The battery level indicator does not light up while charging	Incorrect power supply connection	Connect the power supply correctly
Too short operating time after charging	Too short charging time	Charge the flosser
The flosser cannot be turned on	Low battery level	Charge the flosser
Too low water pressure	Tip clogged	Replace the tip
	Tip's outlet damaged	Replace the tip
	Water pipe clogged	Clean the bottom of the water pipe

Specifications

Name: Water Flosser

Model: MEO705

Rated Voltage/Power: 3.7 V $\approx$ 12 W

Rated Input: 5 V $\approx$ 1 A

IP Rating: IPX7

Charging Time: Approx. 4 hours

Battery Type: KAN 18650 2200

Power Consumption Off mode: $\leq$ 0.3 W

Time to Enter the Condition Off mode: $\leq$ 20 min

Power Adapter

Input: 100–240 V~, 50/60 Hz, 0.2 A			Output: 5.0 V $\overline{\rule{0.5pt}{0.4pt}}\rule{0.5pt}{0.4pt}$ 1.0 A, 5.0 W		
Information Published	Value and Precision			Unit	
Manufacturer's Name	Jiangxi Jian Aohai Technology Co., Ltd.			—	
Model Identifier	A520-050100U-EU1			—	
Input Voltage	100–240			V~	
Input AC Frequency	50/60			Hz	
Output Voltage	5.0			V $\overline{\rule{0.5pt}{0.4pt}}\rule{0.5pt}{0.4pt}$	
Output Current	1.0			A	
Output Power	5.0			W	
Average Active Efficiency	73.77			%	
Efficiency at Low Load (10%)	64.588			%	
No-load Power Consumption	0.075			W	

Europe-EU Declaration of Conformity

We, Xiaomi Communications Co., Ltd., hereby declare that this equipment is in compliance with the applicable Directives and European Norms, and amendments. The full text of the EU declaration of conformity is available at the following internet address:

<http://www.mi.com/global/service/support/declaration.html>

For detailed e-manual, please go to
www.mi.com/global/support/user-guide

WEEE Disposal and Recycling Information

All products bearing this symbol are waste electrical and electronic equipment (WEEE as in directive 2012/19/EU) which should not be mixed with unsorted household waste. Instead, you should protect human health and the environment by handing over your waste equipment to a designated collection point for the recycling of waste electrical and electronic equipment, appointed by the government or local authorities. Correct disposal and recycling will help prevent potential negative consequences to the environment and human health. Please contact the installer or local authorities for more information about the location as well as terms and conditions of such collection points.

Manufacturer: Xiaomi Communications Co., Ltd.
Address: #019, 9th Floor, Building 6, 33 Xi'erqi Middle Road, Haidian District, Beijing, 100085, China

For further information, please go to
www.mi.com

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Výrobce: Xiaomi Communications Co., Ltd.
Xiaomi Technology Netherlands B.V.
Xiaomi Technology UK Limited