Non-Contact Infrared Thermometer

Model JXB-311 / Ref No 09-313


These brief instructions allow you to quickly begin using your Infrared Thermometer for Body Temperature Measurements. Please refer to the Instruction Manual for more detailed instructions covering all features of the thermometer, including Object/Surface and Room Temperature measurements.

- The infrared technology used in this thermometer is safe and is not harmful in any way. However, we recommend avoiding looking into the lens during use.

The Infrared Thermometer offers a selection of settings for person, object/surface and room measurements.

To toggle between modes:

- Press and release the **POWER / SCAN** button to turn the thermometer on. The unit will run a self-test and all symbols on the display will momentarily appear.
- Press and release the **MODE** button to scroll through and select person, object/surface and room temperatures.

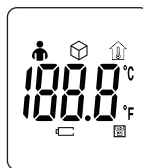Self-test

Person Mode

To measure body temperature using your Infrared Thermometer, follow these steps:

- Press and release the **POWER / SCAN** button to turn the thermometer on. The unit will run a self-test and all symbols on the display will momentarily appear, Fig. 1.
- The display will show dashed lines, the selected mode, and the temperature mode (°F/°C). The thermometer is ready for measurement, Fig. 2.
- Press and release the **POWER / SCAN** button while aiming the infrared lens at the individual's forehead, holding the thermometer approximately 0.5" – 1.0" away.
- When the thermometer is in the correct position, it will slightly vibrate, indicating the positioning sensor has detected the ideal reading location and the measurement has been taken.
- The results will be displayed, Fig. 3.
- If the results are at or above 100.4°F (38°C), the device will vibrate 3 times, indicating a fever may be present.
- After approximately 30-seconds of non-use, the thermometer will automatically shut off.

Fig. 1

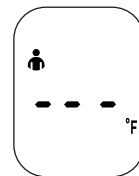
Fig. 2

Fig. 3

**IF YOU HAVE ANY QUESTIONS OR ARE EXPERIENCING DIFFICULTIES,
PLEASE CALL THE TOLL-FREE CUSTOMER SERVICE HELP LINE FOR FURTHER
ASSISTANCE AT 1-866-326-1313.**

Distributed by:
Veridian Healthcare
1175 Lakeside Drive
Gurnee, IL 60031
www.veridianhealthcare.com
Made in China
#93-2162 07/20
©2020 Veridian Healthcare, LLC

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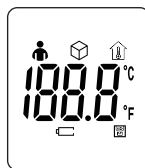
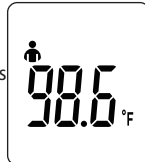
Distributed by:
Veridian Healthcare
1175 Lakeside Drive
Gurnee, IL 60031
www.veridianhealthcare.com
Hecho en China
#93-2162 07/20
©2020 Veridian Healthcare, LLC