

The Coslus logo is a teal-colored circle with the word "COSLUS" in white, uppercase, sans-serif font. The letter "O" is stylized with a white cross-like shape inside it.

# COSLUS

ACTIVATE YOUR  
**24-MONTH**  
WARRANTY FOR FREE



\* Please contact customer service if you're unable to provide the order number

[www.coslus.com/warranty](https://www.coslus.com/warranty)



Win Free Accessories



Covered by 24-Month Warranty



Friendly Responsive Customer Service

**Expire time:** 7 days after the delivery date

 +1 855 594 3950 (9am-5pm Mon. - Fri. (EST))

 [support@coslus.com](mailto:support@coslus.com)


**FLIP →**

# QUICK GUIDE

Scan QR code to watch the operation video. >>




To minimize environmental impact, we recommend using your phone adapter (5V 1A/2A). Please contact us if you need a new one.

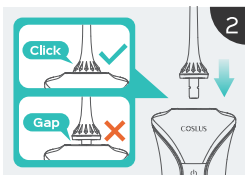


## Install nozzle

Note: Don't press nozzle release button during using.

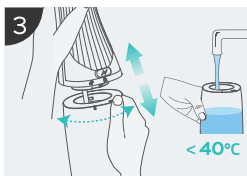


## Fill water

Note: It is recommended to fill warm water for sensitive gums.

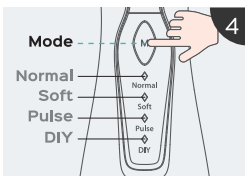


## Adjust mode

Note: Press mode button to select mode.

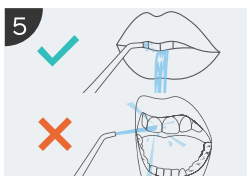


## Put nozzle tip into mouth

Note: Aim the nozzle closely against the gum line, open your mouth slightly.

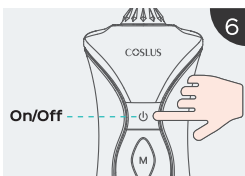

## Start irrigating

Note: Short press On/Off button to turn it on.

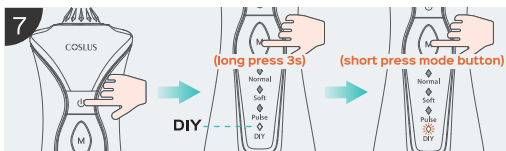

## DIY Mode Instruction: including 8 levels

Power on → adjust mode button to DIY mode → long press 3s until light flash → short press mode button to select water pressure → lock DIY mode (light flash and then off)