

COSLUS

ACTIVATE YOUR **24-MONTH** WARRANTY FOR FREE



* Please contact customer service if you're unable to provide the order number

www.coslus.com/warranty



Win Free Accessories



Covered by 24-Month Warranty



Friendly Responsive Customer Service

Expire time: 7 days after the delivery date

☎ +1 855 594 3950 (9am-5pm Mon. - Fri. (EST))

✉ support@coslus.com




FLIP →

QUICK GUIDE

Scan QR code to watch the operation video. >>




To minimize environmental impact, we recommend using your phone adapter (5V 1A/2A). Please contact us if you need a new one.

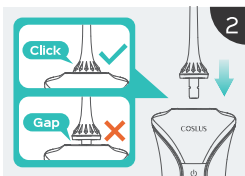


Install nozzle

Note: Don't press nozzle release button during using.

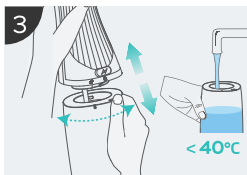


Fill water

Note: It is recommended to fill warm water for sensitive gums.

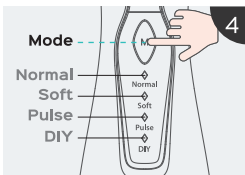


Adjust mode

Note: Press mode button to select mode.

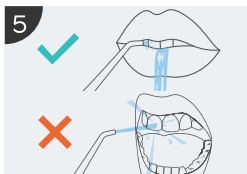


Put nozzle tip into mouth

Note: Aim the nozzle closely against the gum line, open your mouth slightly.

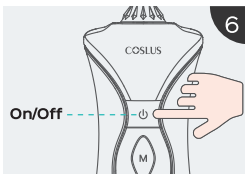

Start irrigating

Note: Short press On/Off button to turn it on.

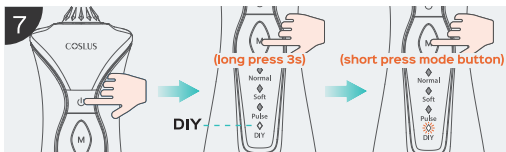

DIY Mode Instruction: including 8 levels

Power on → adjust mode button to DIY mode → long press 3s until light flash → short press mode button to select water pressure → lock DIY mode (light flash and then off)