
*Bob
AND
Brad.*

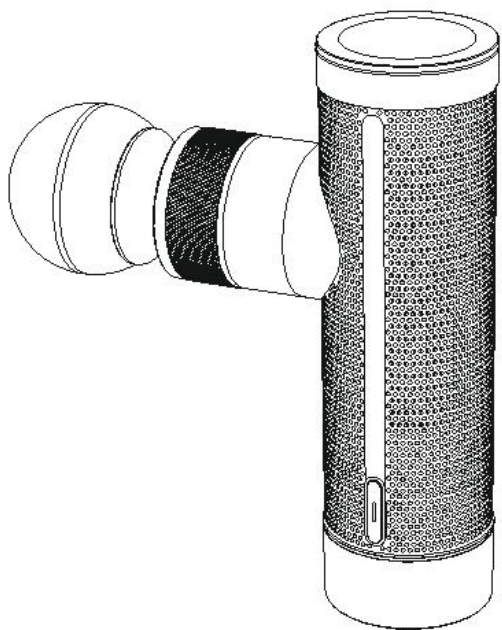User Manual

A7

Massage Gun

Modell: T1035A

Do Not Replace Head During Use.

To avoid any risk of injury or damage, it is strictly forbidden to change the massage head while using the massage gun. Always make sure the device is turned off and unplugged before attempting to remove or insert the massage head. Once the device is safely powered down and unplugged, firmly grasp the inserted head and pull it out to remove it.

To attach the new head, insert its stem into the unit and press firmly until it is securely in place.

Maintenance

The battery in this device is not designed to be removed or exchanged by users. If any problem with the power supply occurs, please contact us directly.

Do not try to open the device, repair it or fix it by taking the device apart yourself. Tampering with the device's mechanics in any way will result in an automatic expiration of the product's warranty.

Warnings and Safety Cautions

Please be aware of the following to ensure a safe use of the product.

- The battery cannot be removed from the product. Do not disassemble the product or the battery.
- Please check the massage gun and the charging cable before using them.
If you notice damages, please report it to us or the retailer and do not use the product.
- Do not hang the product by the power cable.
This product is not a medical product and is not intended to be used as a treatment for medical conditions.
- This product is only intended for healthy adults. Seniors, children, or people with limited cognitive capacity should not use this product or can only use it under supervision.
- Please consult your physician in case you have concerns over your health or experience pain or discomfort while using the product.
- Avoid use the product on your body parts with little muscle coverage, like your head, shins, and back of your hands and feet.
- Avoid use the product on or around your sensitive body parts like neck, eyes, teeth, genitals, or breast implants.
- Do not use the product in any area where there is a risk of bone fracture.
Use this product gently.
- Using this product with too much force or too high a speed might cause bruising.
- Check the massaged area frequently and stop using the product as soon as you experience any pain or discomfort.
- Be careful that your fingers, hair, or other body parts might get caught between the massage head and the shaft of the massage gun.
- Do not use the product in water or get it wet. Use the product only on dry and clean skin.
- Do not put the product in fire or heat it up in any way.
- Do not charge, use, or store the product near heat sources or under extremely hot conditions.

- DE • If the battery leaks or emits an odor, keep it away from fire.
- FR • Do not use the product if the battery:
- emits an odor;
 - gives off excessive heat;
- ES - shows signs of deformation, discoloration, or any of abnormalities.
- IT • If the battery is being charged, stop charging immediately.
- **THIS PRODCUT IS NOT A MEDICAL DEVICE.**

Fcc Statement

This device complies with part 15 of the FCC Rules.

Operation is subject to the following two conditions:

- (1) This device may not cause harmful interference.
- (2) This device must accept any interference received, including interference that may cause undesired operation.

Fcc Statement NOTE

This equipment has been tested and found to comply with the limits for a class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in residential installation. This equipment generates, uses, and can radiate radio frequency energy, and if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause any harmful interference to radio or television receptions, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help with the supplied fastener.
- **DO NOT hang the unit by the power cord.**

Disclaimer

Please contact your physician before using the Bob and Brad T1035A if you have any physical health concerns or have any of the following conditions:

- Metal implants like pins or plates in your body.
- Recent joint replacements
- Pregnancy
- Intrauterine device
- Diabetes
- Epilepsy
- Migraine
- Herniated disks
- Spondylolisthesis (a vertebra slipping forward)


● THIS PRODUCT IS NOT A MEDICAL PRODUCT

IC Warning

This device contains licence-exempt transmitter(s)/receiver(s) that comply with Innovation, Science and Economic Development Canada's licence-exempt RSS(s) . Operation is subject to the following two conditions:

- (1) This device may not cause interference.
- (2) This device must accept any interference, including interference that may cause undesired operation of the device.

Main Structure and Fittings

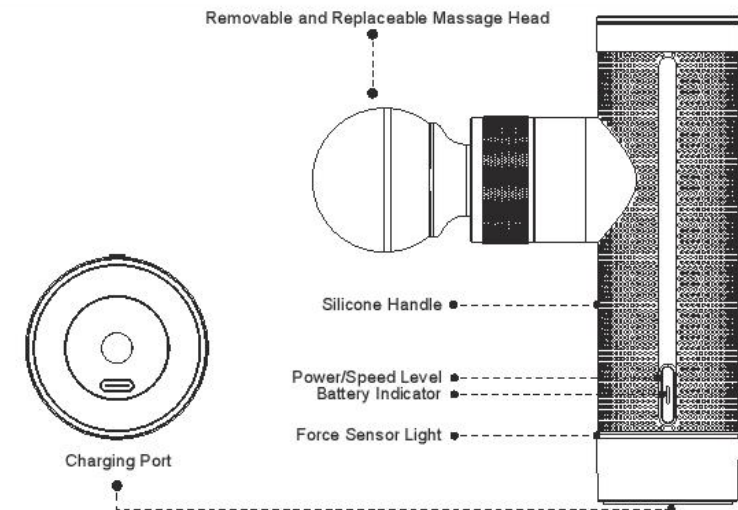Aluminium Alloy Flat Head
For all body parts.

Leg Massage Head
For simulating hand massage with even intensity on curved areas such as arms and thighs.

Ball Head
For large muscle groups like quads, glutes, extreme muscles and joint soreness.

Air Cushion Head
For relaxation of sensitive body parts.

Fan Shaped Silicone Head
Suitable for massaging shoulder girdles, shoulder blades.

DE Operation Steps

FR 1.Align the fitting with the front nozzle mounting port and push to install so that it fits the bearing completely.

ES To remove, hold the fitting, place your fingers on either side of the plastic connector and pull off the fitting.

IT Position the ball head for installation as shown in the figure below (top view).

Do not use the product when it is being charged.

(top view)

2.Turning the device ON/OFF: Press and hold the power button for 0.6 seconds.

3.When the massage gun is on, press the power button "⏻" to choose a suitable massage level and the level is indicated by the speed level indicator above and the power button.

To pause the massage gun, press the power button for 0.6 seconds.

To resume the use, short press the power button.

Note

A.After 10 minutes of operation, the device automatically turns off.

B.LED light indicates the corresponding strike frequency level, the lowest at level 1 and the highest at level 5.

C.When the device is standby mode (0 Level), it will automatically shut down without an operation for 20 seconds.

Battery Display

When the massage gun is being used:

- If the battery indicator is green, the device battery level is above 50%;
- If the battery indicator is orange, the device battery level is above 21%. Advice charge your massage gun.
- If the battery indicator is red, the device battery level is below 2%.Please charge your massage gun promptly.
- If the battery indicator flashes red, the device battery level is below 1%. Please charge your massage gun immediately.

Charge Massage Gun

1.Please charge the product fully before the first use. It takes about 4.5 hours for a full charge.

2.Connect the product to a power source (5V/2A) with the provided type C cable to charge it. To charge it from a power socket, an adapter is needed for the type C Cable

(a regular phone adapter is enough) .It will take longer time to charge the product if the power is below 5V/2A

Use with 5V $\overline{\text{---}}$ 2A
charger output

5V $\overline{\text{---}}$ 2A

3. Once the battery indicator stops flickering and remains on, the massage gun is fully charged.

4. The battery can be recharged at your preference. please make sure to turn off the massage gun before charging it.

5. The massage gun does not work while being charged.

Technical Details

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Speed: 1600 - 3000r/min

Rated Power: 25W

Dimensions:

9x5x16.5cm/ 3.5 " X1.9 " X6.5 "

Battery:

Rated Voltage: 10.8V $\overline{\text{---}}$

Battery Capacity: 2550mAh

Charging Time: About 4.5 Hours

Charge Input:5V $\overline{\text{---}}$ 2A

Troubleshooting:

PROBLEM	CAUSE	SOLUTION
Massage gun doesn't turn on.	power is off.	Long press the power Button to turn on the massage gun.
	Battery is too low.	Charge the massage gun.
	Massage gun shuts-off after 10 minutes.	Long press the power Button to turn on the massage gun for another 10 minutes.
Massage gun shuts off automatically.	The massage gun's overheat protection will switch off the massage gun automatically when too hot.	Wait for the massage gun to cool down before powering it back on.
	Massage gun shuts off after 10 Minutes.	Long press the power Button to turn on the massage gun for another 10 minutes.
Massage gun fails to run.	The massage head is to inserted too tight.	Plug out the massage head a little, then power it on.
If you encounter any other problems, contact: support@bobandbrad.com		

Guarantee

This product comes with a one-year warranty from the date of purchase. Warranty will be voided due to any of the following:

- Damaged due to improper usage or storage.
- Damaged by attempting to repair or disassemble the device.
- Period of warranty exceeded.

Utilize the Bob and Brad massage gun before, during, or after your workout for maximum benefits. This massage gun can help accelerate muscle tissue responsiveness by increasing hydration, range of motion, and function to tight muscles.


Indication	Action	Speed	Pressure	Duration	Muscle Position	Movement
Warm Up/ Activation	Activate and prime muscles prior to exercise. Reduce risk of injury.	Level 3	Light to Moderate	30 sec or less	Relaxed. Move along the muscle.	Quick, high speed over area.
Mid Workout / Myofascial Mobility	Help muscle recovery and reactivation.	Level 4	Heavy (especially for large muscle groups)	8-12 sec	Target specific area.	Localized
Recovery	Relax muscles and aid recovery.	Level 1 or 2	Moderate to Heavy	90 sec or more	Move along the muscle.	Slow movement along the muscle
Pain Modulation	Reduction in pain to increase movement.	Level 3	Light	10 to 30 sec	Relaxed. Localised spot.	Over painful areas.
Body Awareness	Improve movement or mechanics of a joint. Activation of a muscle. Draw attention to a muscle	All Levels	Light to Moderate	10 to 20 sec	Relaxed or with movement	Over specific target area or movement over area.

CHEST

2 Minutes

ATTACHMENT: **BALL HEAD**
PROTOCOL: Sweep from A to B on each side for 1 minute
DURATION: 2 minutes

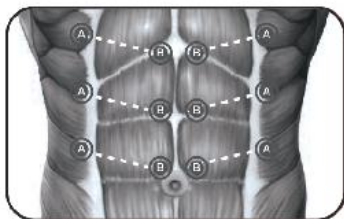ABS

2 Minutes


ATTACHMENT: **BALL HEAD**
PROTOCOL: Sweep from A to B on each side for 1 minute
DURATION: 2 minutes

UPPER BACK

2 Minutes

ATTACHMENT: **BALL HEAD**
STEP 1: **Upper back**
PROTOCOL: Sweep from A to B on each side for 1 minute
DURATION: 2 minutes

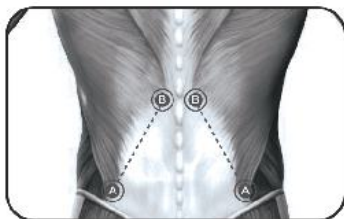LOWER BACK

2 Minutes


ATTACHMENT: **BALL HEAD**
STEP 2: **Lower back**
PROTOCOL: Sweep from A to B on each side for 1 minute
DURATION: 2 minutes

SHOULDERS

2 Minutes


ATTACHMENT: **BALL HEAD**
PROTOCOL: Sweep from A to B on each side for 1 minute
DURATION: 2 minutes

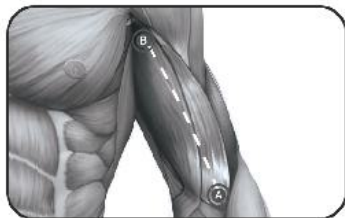BICEPS

2 Minutes


ATTACHMENT: **BALL HEAD**
PROTOCOL: Sweep from A to B
on each side for 1 minute
DURATION: 2 minutes

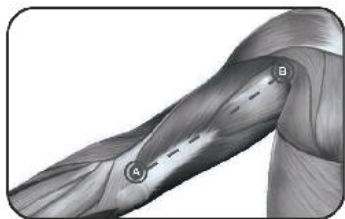TRICEPS

4 Minutes

ATTACHMENT: **BALL HEAD**
PROTOCOL: Sweep from A to B
on each side for 2 minutes
DURATION: 4 minutes

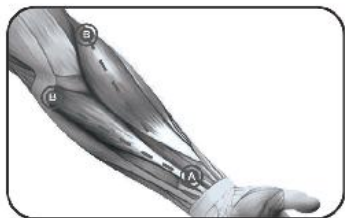FOREARM

2 Minutes

ATTACHMENT: **BALL HEAD**
PROTOCOL: Sweep from A to B
on each side for 1 minute
DURATION: 2 minutes

HANDS

1 Minute


ATTACHMENT: **BALL HEAD**
PROTOCOL: Sweep from A to B
on each side for 30 seconds
DURATION: 1 minute

HIP FLEXORS

2 Minutes


ATTACHMENT: **BALL HEAD**
PROTOCOL: Sweep from A to B
on each side for 1 minute
DURATION: 2 minutes

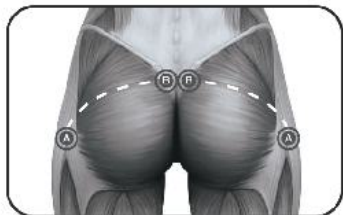GLUTES

4 Minutes

ATTACHMENT: **BALL HEAD**
PROTOCOL: Sweep from **A** to **B** on each side for 2 minutes
DURATION: 4 minutes

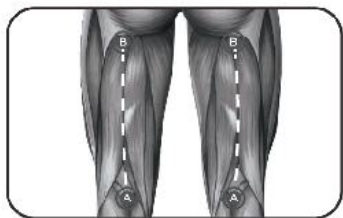HAMSTRINGS

2 Minutes

ATTACHMENT: **BALL HEAD**
PROTOCOL: Sweep from **A** to **B** on each side for 1 minute
DURATION: 2 minutes

CALVES

2 Minutes

ATTACHMENT: **BALL HEAD**
PROTOCOL: Sweep from **A** to **B** on each side for 1 minute
DURATION: 2 minutes

SHINS

1 Minute

ATTACHMENT: **BALL HEAD**
PROTOCOL: Sweep from **A** to **B** on each side for 30 seconds
DURATION: 1 minute

FEET

4 Minutes

ATTACHMENT: **BALL HEAD**
PROTOCOL: Sweep from **A** to **B** on each side for 2 minutes
DURATION: 4 minutes

WEEE Instructions

Correct Disposal of this product

This marking indicates that this product should not be disposed with other household wastes throughout the EU. To prevent possible harm to the environment or human health from uncontrolled waste disposal, recycle it responsibly to promote the sustainable reuse of material resources. To return your used device, please use the return and collection systems or contact the retailer where the product was purchased. They can take this product for environmental safe recycling.

Battery Instructions

- 1 Do not submerge the product into water or get it wet.
- 2 Do not charge the product near an open flame or in any extremely hot conditions.
- 3 Do not use or store the product near heat sources like furnaces or heaters.
- 4 If the battery leaks or emits an odor, remove it immediately from close proximity to an open flame.
- 5 Do not put the product into a fire or heat it up in any way.
- 6 Do not disassemble the battery in any way.
- 7 Do not use the product if, while the battery is in use or charging, it emits an odor; gives off excessive heat; shows signs of deformation, discoloration, or any other abnormalities. If the battery is in use or being charged, disconnect the battery from the power source.
- 8 Charge for about three hours before its first use to a full charge. The device can be connected or disconnected.
- 9 Connect the output end of the connection cable to the bottom of the massage instrument's DC charging port. Then, connect the input end of the connection cable to a mobile power adapter or cell phone charger (5V 2A; if less than 2A, the charging time will be longer).
- 10 The power indicator light will flash regularly, indicating that the battery is being charged.

RoHS Declaration

Our product is compliant with RoHS 2.0 (2011/65/EU Directive) regarding restriction of hazardous substances and raw materials, and components of our product meet the following toxic and hazardous substance limits.

Cadmium (Cd): < 100 ppm	Lead (Pb): < 1000 ppm
Mercury (Hg): < 1000 ppm	Hexavalent Chromium: (Cr VI) < 1000 ppm
Polybrominated Biphenyls (PBB): < 1000 ppm	Polybrominated Diphenyl Ethers (PBDE): < 1000 ppm
Bis(2-Ethylhexyl) phthalate (DEHP): < 1000 ppm	Benzyl butyl phthalate (BBP): < 1000 ppm
Dibutyl phthalate (DBP): < 1000 ppm	Diisobutyl phthalate (DIBP): < 1000 ppm

A Dedicated Team Just For You

Every Bob and Brad product is backed by a 1-year warranty and an award-winning customer service team. If there are any quality problems with the product you ordered, we will try our utmost to help you get the problems resolved, as well as offer a replacement, return, or refund service under the warranty. Simply email us the problem, and we'll take care of it.

Troubleshoot

Help you solve any issues with your device.

Product Knowledge

Answer questions about your device and best practices.

Initiate Returns

Provide instructions for replacements, returns and refunds.

Place New Orders

Help with placing orders online and/or authorized resellers.

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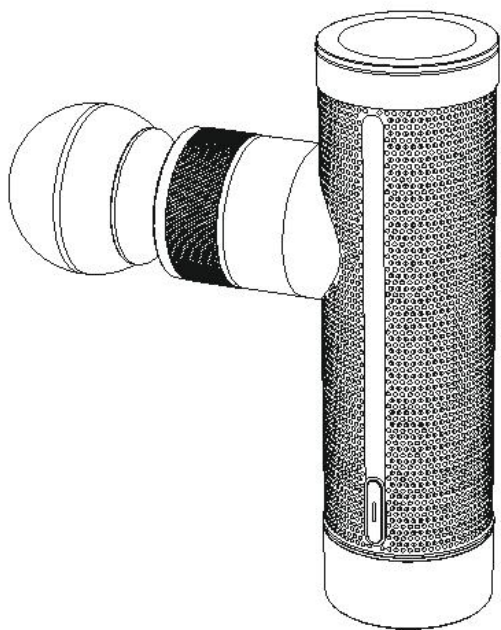
Tension nominale : 10.8V=

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Troubleshooting:

PROBLEM	CAUSE	SOLUTION
Massage gun doesn't turn on.	power is off.	Long press the power Button to turn on the massage gun.
	Battery is too low.	Charge the massage gun.
	Massage gun shuts-off after 10 minutes.	Long press the power Button to turn on the massage gun for another 10 minutes.
Massage gun shuts off automatically.	The massage gun's overheat protection will switch off the massage gun automatically when too hot.	Wait for the massage gun to cool down before powering it back on.
	Massage gun shuts off after 10 Minutes.	Long press the power Button to turn on the massage gun for another 10 minutes.
Massage gun fails to run.	The massage head is inserted too tight.	Plug out the massage head a little, then power it on.
If you encounter any other problems, contact: support@bobandbrad.com		

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UK Importer:

KOVAASSOCIATES LTD

72a Kingston Road, Portsmouth, Hampshire, PO2 7PA, United Kingdom

KOVA.Compliance@outlook.com

0044 (0) 7542988998

Manufacturer:

Shenzhen Laige Technology Co.Ltd.

Room 301, 3/F, Building 2, Yunli Intelligent Park, Bantian Street, Longgang District, Shenzhen, Guangdong, 518129, China

Support@bobandbrad.com

EC Declaration of Conformity

We(UK Importer is responsible for this declaration

KOVA ASSOCIATES LTD

(Company name)

72a Kingston Road, Portsmouth, Hampshire, PO2 7PA, United Kingdom

(Company address)

KOVA.Compliance@outlook.com

(Email Address)

declare under our sole responsibility that the product

Tradename:BOB AND BRAD

Equipment: Massage Gun

ModelNo.:T1035A

ASIN: B0D41NJDC2/B0D41R1NTH

to which this declaration relates is in conformity with the essential requirement set out in the Council Directive on the Approximation of the laws of the Member States relating to LVD Directive(2014/35/EU) & EMC Directive(2014/30/EU) &RoHS(2002/95/EC) product is responsible to affix CE marking, the following standards were applied:

EN 60335-1:2012+A11:2014+A13:2017

EN 60335-2-32:2003+A1:2008+A2:2015

EN 62233:2008

EN 55014-1:2017

EN 55014-2:2015

EN 61000-3-2:2014

EN 61000-3-3:2013

Full Name: Taozirui Position: Manager

Signature: Taozirui Date: 2021.5.6

EC Declaration of Conformity

We (EU Importer is responsible for this declaration)

XDH Tech

(Company name)

2 Rue Coysevox Bureau 3, Lyon, 69001, France

(Company address)

xdh.tech@outlook.com

(Email Address)

declare under our sole responsibility that the product

Tradename:BOB AND BRAD

Equipment: Massage Gun

ModelNo.:T1035A

ASIN: B0D41NJDC2/B0D41R1NTH

to which this declaration relates is in conformity with the essential requirement set out in the Council Directive on the Approximation of the laws of the Member States relating to LVD Directive(2014/35/EU) & EMC Directive(2014/30/EU) &RoHS(2002/95/EC) product is responsible to affix CE marking, the following standards were applied:

EN60335-1:2012+A11:2014+A13:2017

EN60335-2-32:2003+A1:2008+A2:2015

EN 62233:2008

EN 55014-1:2017

EN 55014-2:2015

EN 61000-3-2:2014

EN61000-33:2013

Full Name: Taozirui Position: Manager

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