

# OWNER'S MANUAL

## Massage gun
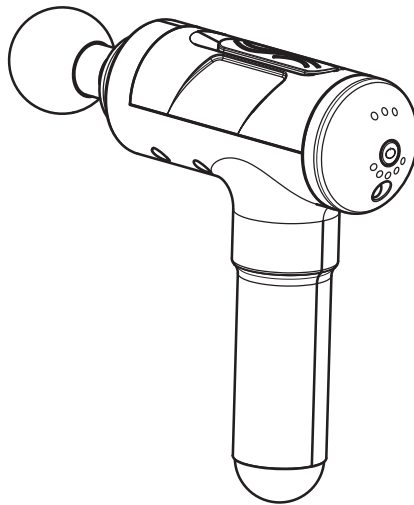
# MASSAGEPISTOOL

## SPECIFICATIES

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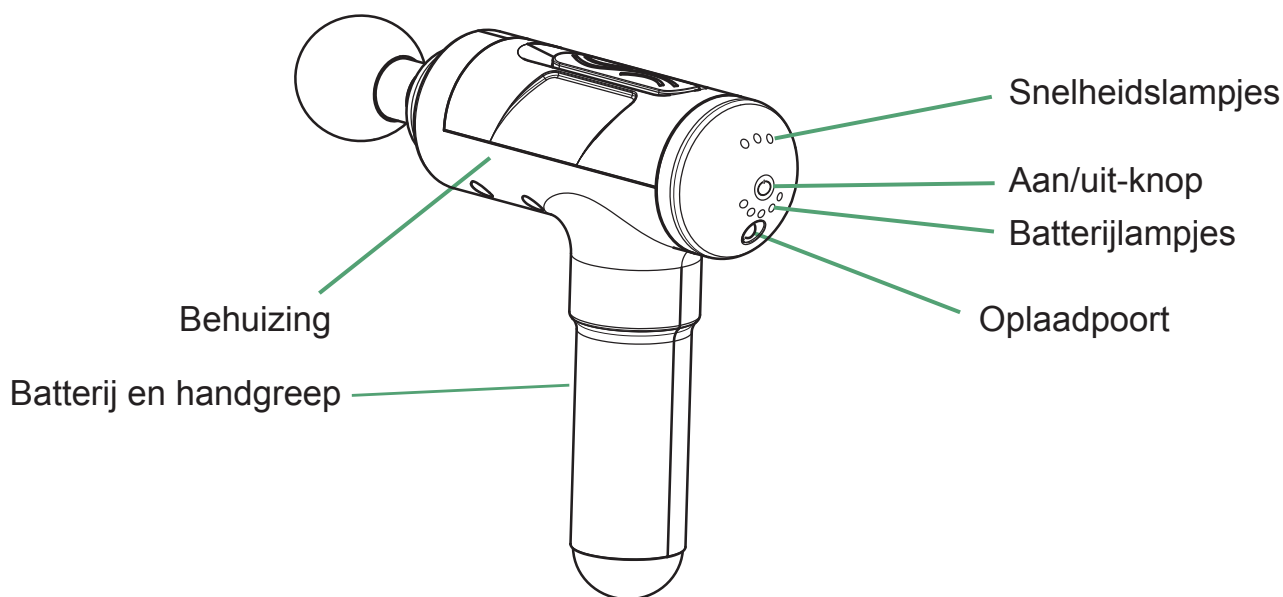
# MASSAGE GUN

## SPECIFICATIONS

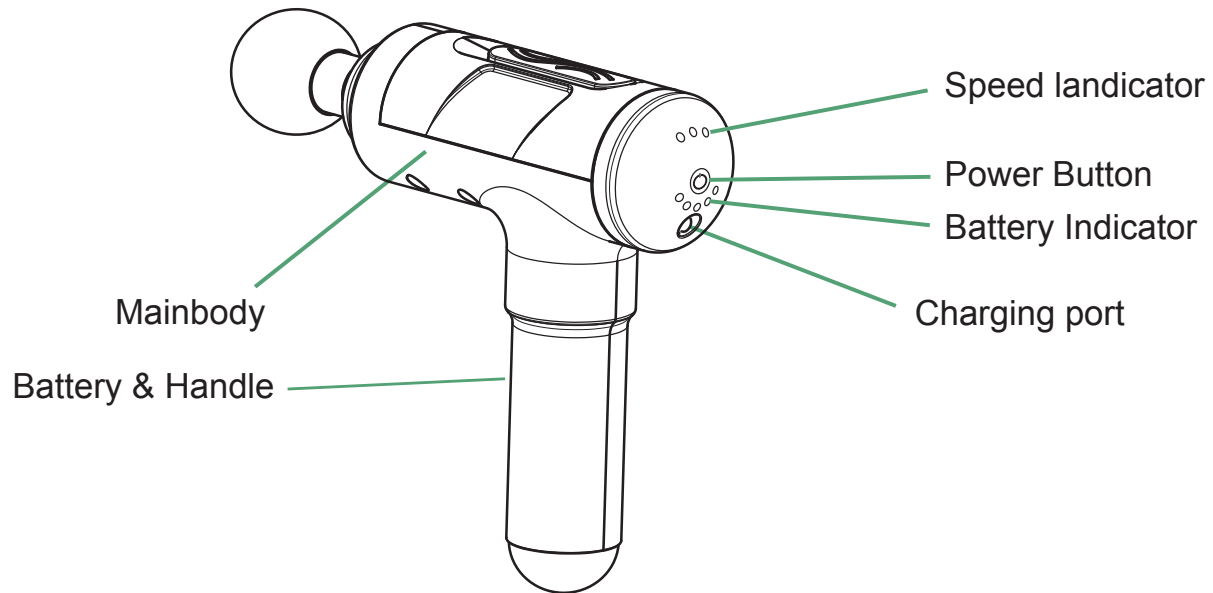1. Powerful High-Torque Motor with Quiet Glide Technology
2. Built-In Rechargeable Lithium-Ion Battery
3. Multi-Speed Setting for Optimal Pressure
4. Cordless and Lightweight including charger
5. Easy Speed Adjusting

### **Vibration Speed:**

Level 1- Level 3

### **Vibration Frequency:**

Level 1-33.3Hz, Level 2-40.0Hz, Level 3-53.3Hz

### **Charge:**

110-240V 50/60Hz 12V 1A

# USAGE INSTRUCTIONS

PLEASE READ PRIOR TO USE

## CHARGING

- Fully charge battery for up to six hours before first use either attached or detached from the device
- To charge, connect the DC end of the supplied adapter to the battery's charging port and plug AC side cable into wall outlet
- One battery light means 20% level.
- The massager should be charged when there is only one battery light left
- The battery lights stop changing and stay on means it is full charged
- It is not recommended to fully deplete the battery during using
- The average run time is 2+ hours depending on speed level and pressure applied during use

## HOW TO USE

- Press the Power button and hold it until the 5 battery Lights up
- Press the Power button once to turn device on to level one. twice to turn device to level two and three times for level three. pressing the power button again will turn off
- Power Off: Press and hold until the button lights off

## MAINTENANCE, CLEANING STORAGE

- Use a slightly damp towel to clean device and a dry with a soft cloth
- For storage or travel, remove battery and store in supplied carton box

# DISCLAIMER

## Warnings & Caution

- Do not immerse in water. Keep away from liquids or heat sources. Keep ventilation ports free from dust and debris
- Do not remove screws or attempt to disassemble
- Do not operate continuously for more than one hour. allow device to rest for 30 minutes before re-using
- Unplug the charger from device battery after charging or prior to use
- For adult use only do not use if injured. Consult your doctor before using this product

## SAFETY INSTRUCTIONS

DANGER: TO REDUCE RISKS OF ELECTRIC SHOCK, FIRE, AND PERSONAL INJURY, THIS PRODUCT MUST BE USED IN ACCORDANCE WITH THE FOLLOWING INSTRUCTIONS

- For adult use only
- Use only on a dry clean surface of the body over clothing by lightly pressing and moving across the body for approximately 60 seconds per region
- Use the Massage Gun only on the body's soft tissue as desired without producing pain or discomfort. Do not use on head or any hard or bony area of the body
- Use only the applicators that provide the best desired effect
- Bruising can occur regardless of control setting or pressure applied Check treated areas frequently and immediately stop at first sign of pain or discomfort
- Keep fingers hair or any other body part away from the shaft and backside of the applicator as pinching may occur
- Do not place any objects into the ventilation holes of the Massage Gun
- Do not immerse in water or allow water to enter the ventilation holes of the Massage Gun
- Do not drop or misuse the Massage Gun
- Only recharge with supplied 12 VDC Massage Gun charger
- Carefully examine device and battery before each use
- Do not tamper or alter Massage Gun in any way
- Never leave the Massage Gun operating or charging unattended

# DISCLAIMER


PLEASE DO NOT USE THE MASSAGE GUN OR ANY VIBRATION DEVICE WITHOUT FIRST OBTAINING APPROVAL FROM YOUR DOCTOR IF ANY OF THE FOLLOWING APPLY

Pregnancy diabetes with complications such as neuropathy or retinal damage, wear of pace-makers, recent surgery epilepsy or migraines, herniated disks, spondylolisthesis, spondylolysis, or spondylosis, recent joint replacements or IUDs, metal pins or plates or any concerns about your physical health frail individuals and children should be accompanied by an adult when using any vibration device. these contra indications do not mean that you are not able to use a vibration or exercise device but we advise you to consult a doctor first

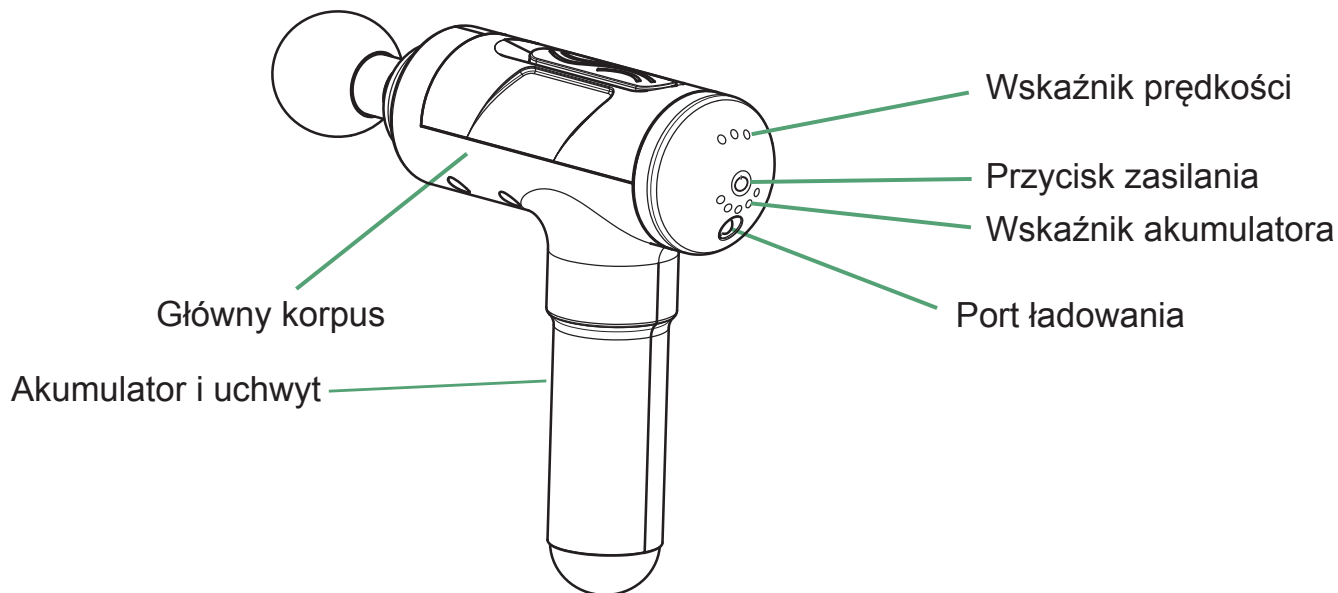Ongoing research is being done to the effects of vibration exercise for specific medical disorders. This will likely result in the shortening of the list of contra indications as shown above. Practical experience has shown that integration of vibration exercise into a treatment plan is advisory in a number of cases. this must be done on the advice of and in accompaniment with a doctor specialist or physiotherapist

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