COSLUS

Dual stream Technology

Oral Irrigator
Munddusche

Model

C20

USER MANUAL
Benutzerhandbuch

www.coslus.com
Tel: +1(855) 594-3950
(8am-5pm Mon. - Fri. (EST))
E-mail: support@coslus.com
Made in China

CONFORMS TO STD.
UL60335-1&60335-2-52

01

What's in the box ?

EN	① Main unit x 1 ② Standard jet tip x 2 ③ Orthodontic tip x 1 ④ Periodontal tip x 1	⑤ Tongue scraper x 1 ⑥ User manual x 1 ⑦ USB cable x 1 (adapter is not included)
DE	① Hauptgerät x 1 ② Standard-Düsen spitze x 2 ③ Kieferorthopädische Spitze x 1 ④ Parodontale Spitze x 1	⑤ Zungenschaber x 1 ⑥ Bedienungsanleitung x 1 ⑦ USB-Kabel x 1 (Adapter ist nicht enthalten)
FR	① Unité principale x 1 ② Embout de jet standard x 2 ③ Embout orthodontique x 1 ④ Embout parodontal x 1	⑤ Racleur de langue x 1 ⑥ Manuel d'utilisation x 1 ⑦ Câble USB x 1 (adaptateur non inclus)
ES	① Unidad principal x 1 ② Boquillas estándar x 2 ③ Boquilla de ortodoncia x 1 ④ Boquilla periodontal x 1	⑤ Raspador (lingua) x 1 ⑥ Manual de usuario x 1 ⑦ Cable USB x 1 (adaptador no incluido)
IT	① Unità principale x 1 ② Punta jet standard x 2 ③ Punta ortodontica x 1 ④ Punta parodontale x 1	⑤ Raschia lingua x 1 ⑥ Manuale d'uso x 1 ⑦ Cavo USB x 1 (adattatore non incluso)

03

 IMPORTANT SAFETY CAUTIONS

EN All information in this part is provided to instruct you to use the appliance in a correct and safe way so that you can avoid injuries to yourself and others as well as damage to property.

- The appliance can be used by children aged 8 years and above. Person with reduced physical, sensory or mental capabilities or lack of experience and knowledge are not allowed to use unless they have guardian supervision or instruction concerning safe using.
- Children shall not play with the appliance. Cleaning and maintenance shall not be made by children without supervision.
- The appliance is only to be used with the power supply unit provided with the appliance.
- Never alter the appliance. Also, do not try to disassemble or repair it.
- The battery must be removed from the appliance before it is scrapped.
- Unplug the appliance from power supply when removing the battery.
- The battery should be disposed of safely.
- Regarding how to remove the battery, please refer to disposal instruction.
- Do not fill the tank with water warmer than 40°C.
- If the external flexible cable or cord is damaged, please ask the manufacturer, service agent or qualified professional person for help to avoid hazard.
- Keep the cable or cord away from heated surfaces.
- Please do not use it if you lost the nozzle.
- Do not submerge the device in water.
- Do not use for any purpose other than oral cleaning.
- Keep nozzle away from small children. The nozzle is small enough to be swallowed and can be lodged in the throat.
- Certain types of mouthwash may damage the appliance by causing the housing or water reservoir to crack. Therefore, do not fill water tank with mouthwash.
- People who is suffering from tooth or gum pain may not be able to use the appliance and should consult their dentist before use.
- Do not press the nozzle too firmly against your teeth or gums, or it may cause injury.

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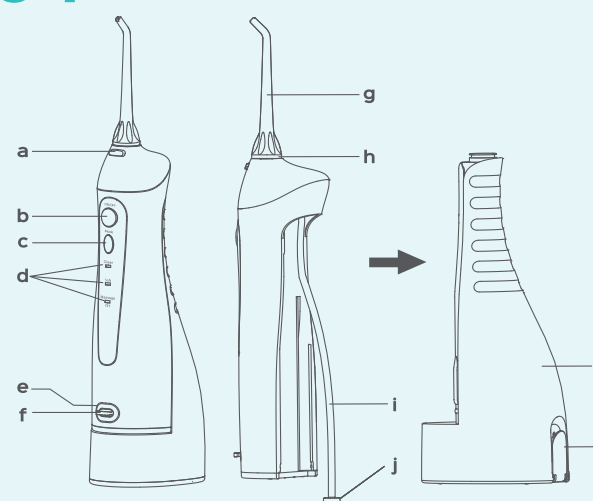
EN

How to Charge (1)

1. Take charging base cover down, securely plug USB cable into charging port of the appliance. (Make sure no liquid remaining inside the port before plugging).

2. Connect the charging adaptor securely with the power outlet.

3. "CH" indicator blinks during the appliance is charging and it will turn off after it's fully charged.

4. Cover charging base again after charging is finished.

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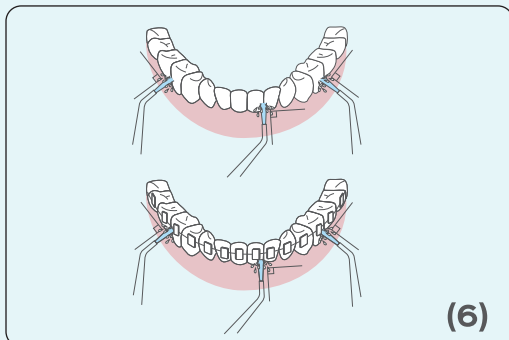
EN 6. Press "On/Off" button to start irrigating. (6)
1. Direct water flow at 90° to the teeth and gums
2. Direct water flow slowly along the teeth
3. Direct water flow along the gum line with nozzle tip at 90°
4. Clean between the teeth and around any braces, crowns or bridgework thoroughly

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
• After Using (8)

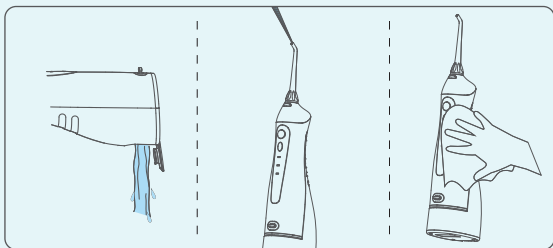
EN **Note:** To prevent bacteria breeding, make sure there is no remaining water inside. If you're not going to use the appliance for a long time, be sure to wipe and dry it before storing.

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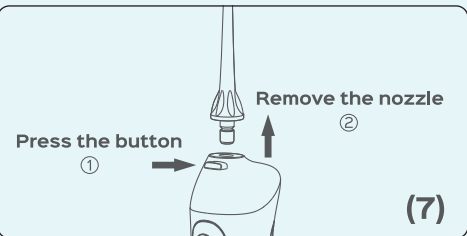
EN 7. Press "On/Off" button to turn it off, then press nozzle release button to replace other nozzle you need. (7)

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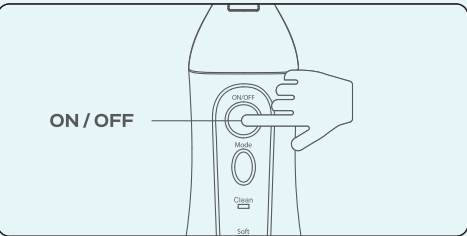EN 8. After irrigating is finished, press "On/Off" button again to turn off. (7)

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EN 1. Open water tank cap and dispose of water inside the tank.
2. Press "On/Off" button to turn on and drain water inside the appliance.
3. Press "On/Off" button again to turn it off, and close water tank cap.
4. Wipe the appliance with dry cloth.

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
CLEANING TIP INSTRUCTION

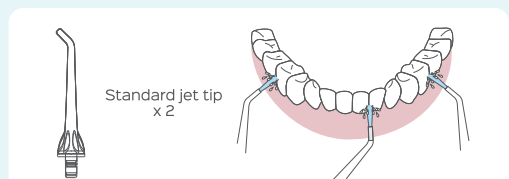
EN **Standard Jet Tip:** daily clean teeth, gum line, fissure, remove food residue, stain, plaque, etc.

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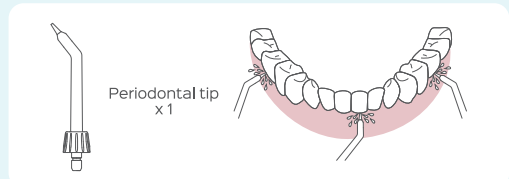EN **Periodontal Tip:** clean gingival sulcus, especially suitable for those who are suffering from periodontal diseases. When cleaning periodontal pocket, please use the lowest water pressure. Direct water flow along the gum line with nozzle tip at 45° and turn it on.

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
EN **Note:** It is recommended to replace cleaning tip every 6 months for better cleaning experience.

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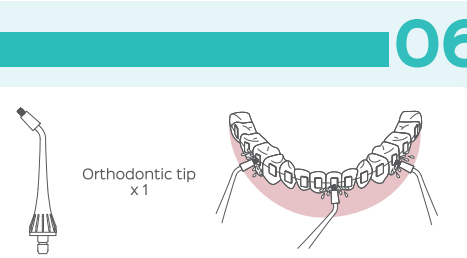
EN **Orthodontic Tip:** clean bridge and teeth with braces, suitable for those who wear braces or have dental implant.

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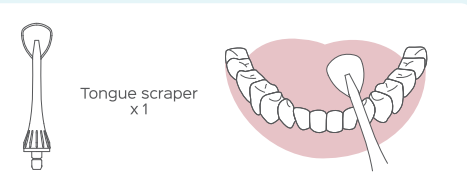EN **Tongue Scraper:** clean tongue surface, remove bad breath, food residue and thick tongue fur. When cleaning tongue, please use the soft mode. Lean forward to put scraper on your tongue and turn it on, then move gently from back to front. You can adjust water pressure according to your need.

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MAINTENANCE

EN Use only water or neutral cleaning agent to clean the appliance. Never use caustic or abrasive agents (e.g., vinegar or lime scale removers), since they may cause damage to the appliance. Do not use water over 40°C. Keep the appliance away from high temperature and avoid direct sunlight.

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EN • **Main unit**
1. Wipe the dirt by cloth with soapy water.
2. Do not submerge the main unit in water to clean.
3. Dry after washing with water.
• **Water tank cap**
Wipe off any grime and dirt with a cloth.
• **Nozzle**
1. Wash with water and wipe with soft cloth.
2. Replace the nozzle every 6 months.
• **Water suction tube**
1. Wash with water and wipe with soft cloth.
2. Do not bend, pull or twist the water suction tube.
• **Tank**
1. Wash with water and wipe with soft cloth.
2. Dry the tank when not using for a long period.

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08 TROUBLESHOOTING EN


Problem	Cause	Solution
Water does not come out	Negative pressure caused by not using for a long time.	Solution 1: Full fill the tank with water, invert the device and shake it to eliminate negative pressure.
	You've just purchased the appliance or it has been left unused for a period.	Solution 2: Install the nozzle tip on the main unit, suck out the air from nozzle tip with mouth to eliminate negative pressure
It does not work	You've just purchased the appliance or it has been left unused for more than 3 months.	Charge the appliance for at least 4 hours.
It only works for a few minutes even after charging	The battery has reached the end of its service life.	Contact service center for support.
	Charging time was insufficient.	Charge the unit for at least 4 hours.
Water pressure is very low	There is dirt stuck in the nozzle.	Replace the nozzle.
	The nozzle is deformed.	Clean the filter.

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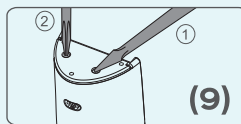EN 1. Remove 2pcs rubber feet on the bottom of the appliance. (9)
2. Use screwdriver counter-clockwise to loosen bottom cap. (9)
3. Slide water tank down off the main unit and unscrew two screws beside nozzle suction tube. (10)
4. Open top cover and put inner part out of the main unit, then unscrew all screws of the inner part. (11)
5. Cut off battery wire and pull out batteries from the inner. (12)

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Warm Tip: If you have any question about recycling, please contact your local waste management office.

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09 DISPOSAL

EN The rechargeable battery inside your appliance is easy to remove for recycling. At the end of its service life and prior to disposal, please remove the rechargeable battery following the instructions below.

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