
*Bob
AND
Brad.*

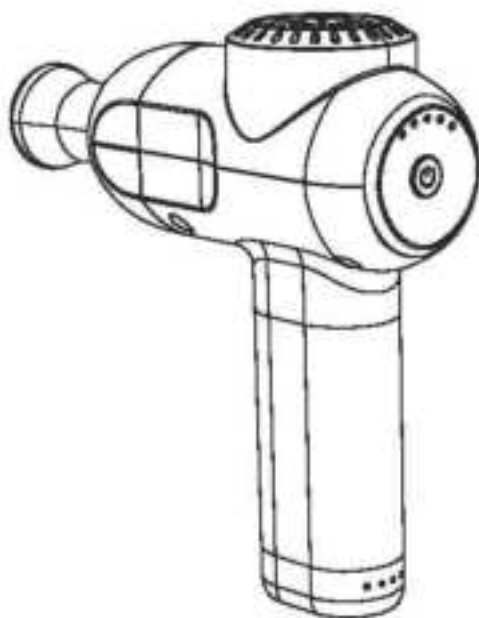User Manual

QL/DMS.C2-D
Massage Gun

Maintenance

The battery in this device is not designed to be removed or exchanged by users. If any problem with the power supply occurs, please contact us directly.

Do not try to open the device, repair it or fix it by taking the device apart yourself. Tampering with the device's mechanics in any way will result in an automatic expiration of the product's warranty.

Warnings and Safety Concerns

Before connecting the device to the power supply, please ensure that the supply voltage stated on the rating plate is compatible with the unit. Please charge the new batteries before the first use.

Do not submerge the device in water or any other liquid.

Do not use the device for more than 30 minutes at a time. Wait 5 minutes prior to using the device again, following a 30-minute session with the device.

Do not apply to your head or near your genitals.

Do not use the device to supplement or replace any medical treatment. Please consult your physician before using the device if you are experiencing unexplained pain, if you're receiving medical treatment, or using medical equipment.

FCC Statement

1. This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions:

- (1). This device may not cause harmful interference.
- (2). This device must accept any interference received, including interference that may cause undesired operation.

FCC Statement NOTE:

This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in residential installation. This equipment generates, uses, and can radiate radio frequency energy, and if not installed and used in accordance with the instructions, may cause harmful interference to radio communications.

However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause any harmful interference to radio or television receptions, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help with the supplied fastener.
- DO NOT hang the unit by the power cord.

DISCLAIMER

Please contact your physician before using the Bob and Brad QL/DMS.C2-D if you have any physical health concerns or have any of the following conditions:

- Metal implants like pins or plates in your body.
- Recent joint replacements
- Pregnancy
- Intrauterine device
- Diabetes
- Epilepsy
- Migraine
- Herniated disks
- Spondylolisthesis (a vertebra slipping forward)


★ THIS PRODUCT IS NOT A MEDICAL PRODUCT

IC Warning:

This device contains licence-exempt transmitter(s)/receiver(s) that comply with Innovation, Science and Economic Development Canada's licence-exempt RSS(s). Operation is subject to the following two conditions:

- (1) This device may not cause interference.
- (2) This device must accept any interference, including interference that may cause undesired operation of the device.

Main structure and the fittings

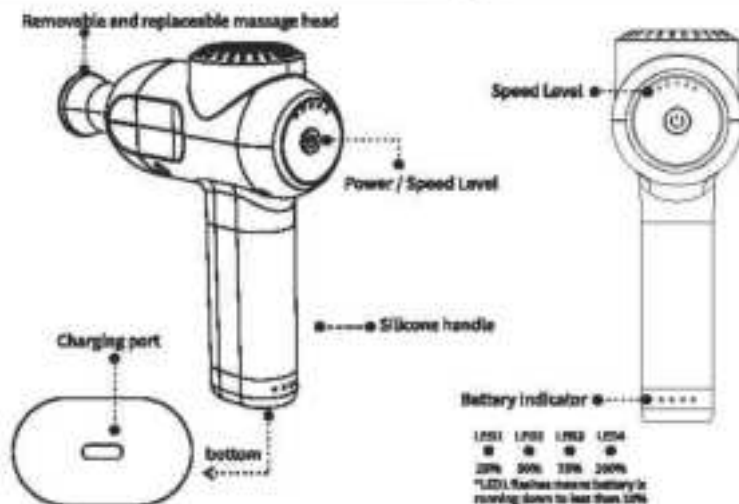Ball Head

For large muscle groups like quads, glutes, extreme muscles and for joint soreness.

U-Shaped Head

For spine, neck, and shoulder. For sensitive muscles and soft tissues.

Air Cushion Head

For relaxation of sensitive body parts.

Flat Head

For all body parts.

Bullet Head

For joints, deep tissues, trigger points, and small muscle areas like feet and wrists.

Rubber Sleeve

Each massage head is pre-installed with a rubber sleeve. If it is damaged, please replace it with the provided spare one.

Operation steps

1. When the device is off, insert the desired massage head with pressure into the message gun. To remove the massage head, pull it out when the device is stopped.

Do not use the product when it is being charged.

Note: Position the U-shaped head for installation as shown in the figure below (top view) to prevent the risk of the massage head falling out while the massage gun is being used. When the buckle is inserted into the slot and a click is felt, the installation of the U-shaped head is in place.

2. Turning the device ON/OFF: Press and hold the "⊙" power button for 2 seconds.
3. When the massage gun is on, press the power button "⊙" to choose a suitable massage level and the level is indicated by the speed level indicator above the power button.
To pause the massage gun, press the power button for 2 seconds.
To resume the use, short press the power button.

Note

- A. After 10 minutes of operation, the device automatically turns off.
- B. LED light indicates the corresponding strike frequency level, the lowest at level 1 and the highest at level 5.
- C. When the device is standby mode (0 level), it will automatically shut down without any operation for 20 seconds.

Battery display

LED1	LED2	LED3	LED4
●	●	●	●
25%	50%	75%	100%

* LED 1 flashing means the battery is running down to less than 10%

When the first indicator light is on, it indicates that 25% of the battery remains.

When the second indicator light is on, it indicates that 50% of the battery remains.

When the third indicator light is on, it indicates that 70% of the battery remains.

When the fourth indicator light is on, it indicates that 100% of the battery remains.

Note

1. When the first indicator light flashes, the battery is running down to less than 10% and a recharge is required.
2. It is not recommended that the battery power be below depleted below 25%.

Charge Massage Gun

1. Please charge the product fully before the first use. It takes about 3 hours for a full charge.
2. Connect the product to a power source (5V/2A) with the provided type C cable to charge it. To charge it from a power socket, an adapter is needed for the type C cable (a regular phone adapter is enough).
It will take longer time to charge the product if the power is below 5V/2A.

Use with 5V $\equiv$ 2A
charger output

5V $\equiv$ 2A

3. When a massage gun is being charged, the battery indicator on the handle flickers in blue regularly, indicating the power level.
4. Once the battery indicator stops flickering and remains on, the massage gun is fully charged.
5. The battery can be recharged at your preference. Please make sure to turn off the massage gun before charging it.
6. The massage gun does not work while being charged.

Technical details:

QL /DMS.C2-D Massage Gun

Speed: 2000 - 3200 r/min

Rated Power: 25W

Dimensions: 18.3x6.2x13.2cm/ 7.2" x2.4" x5.1"

Battery Information:

Rated Voltage: 7.4V---

Battery Capacity: 2500mAh

Charging Time: 210 minutes

Troubleshooting:

PROBLEM	CAUSE	SOLUTION
Message gun doesn't turn on.	Power is off.	Long press the Power Button to turn on the message gun.
	Battery is too low.	Charge the message gun.
	Message gun shuts-off after 10 minutes.	Long press the Power Button to turn on the message gun for another 10 minutes.
Message gun shuts off automatically.	The message gun's overheat protection will switch off the message gun automatically when too hot.	Wait for the message gun to cool down before powering it back on.
	Message gun shuts off after 10 Minutes.	Long press the Power Button to turn on the message gun for another 10 minutes.
Message gun fails to run.	The message head is too inserted too tight.	Plug out the message head a little, then power it on.
If you encounter any other problems, contact: support@bobandbrad.com		

Guarantee

This product comes with a one-year warranty from the date of purchase. Warranty will be voided due to any of the following:

- Damaged due to improper usage or storage.
- Damaged by attempting to repair or disassemble the device.
- Period of warranty exceeded.

Utilize the Bob and Brad massage gun before, during, or after your workout for maximum benefits. This massage gun can help accelerate muscle tissue responsiveness by increasing hydration, range of motion and function to tight muscles.

Indication	Action	Speed	Pressure	Duration	Muscle Position	Movement
Warm Up / Activation	Activate and prime muscles prior to exercise. Reduce risk of injury.	Level 3	Light to Moderate	30 sec or less	Relaxed. Move along the muscle.	Quick, high speed over area.
Mid Workout / Myofascial Mobility	Help muscle recovery and reactivation.	Level 4	Heavy (especially for large muscle groups)	8-12 sec	Target specific area.	Localised
Recovery	Relax muscles and aid recovery. Down-regulate or turn off muscles.	Level 1 or 2	Moderate to Heavy	90 sec or more	Move along the muscle.	Slow movement along the muscle
Pain Modulation	Reduction in pain to increase movement.	Level 1	Light	10 to 30 sec	Relaxed. Localised spot.	Over painful area.
Body Awareness	Improve movement or mechanics of a joint. Activation of a muscle. Draw attention to a muscle	All Levels	Light to Moderate	10 to 20 sec	Relaxed or with movement	Over specific target area or movement over area.

CHEST

2 Minutes

ATTACHMENT: BALL HEAD

PROTOCOL: Sweep from A to B on each side for 1 minute

DURATION: 2 minutes

ABS

2 Minutes


ATTACHMENT: BALL HEAD

PROTOCOL: Sweep from A to B on each side for 1 minute

DURATION: 2 minutes

UPPER BACK

4 Minutes

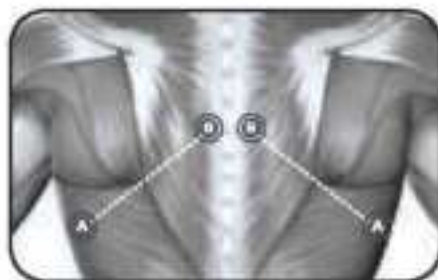ATTACHMENT: BALL HEAD

STEP 1: UPPER back

PROTOCOL: Sweep from A to B on each side for 1 minute

DURATION: 2 minutes

UPPER BACK

4 Minutes

STEP 2: Lower back

PROTOCOL: Sweep from A to B on each side for 1 minute

DURATION: 2 minutes

SHOULDERS

2 Minutes

ATTACHMENT: BALL HEAD

PROTOCOL: Sweep from A to B on each side for 1 minute

DURATION: 2 minutes

BICEPS

2 Minutes

ATTACHMENT: BALL HEAD

PROTOCOL: Sweep from A to B on each side for 1 minute

DURATION: 2 minutes

TPICEPS

4 Minutes

ATTACHMENT: BALL HEAD

PROTOCOL: Sweep from A to B on each side for 2 minute

DURATION: 4 minutes

FOREARM

2 Minutes

ATTACHMENT: BALL HEAD

PROTOCOL: Sweep from A to B on each side for 1 minute

DURATION: 2 minutes

HANDS

1 Minute

ATTACHMENT: BALL HEAD

PROTOCOL: Sweep from A to B on each side for 30 seconds

DURATION: 4 minute

HIP FLEXORS

2 Minutes

ATTACHMENT: BALL HEAD

PROTOCOL: Sweep from A to B on each side for 1 minute

DURATION: 2 minutes

SHINS

1 Minute

ATTACHMENT: BALL HEAD

PROTOCOL: Sweep from A to B on each side for 30 seconds

DURATION: 4 minute

GLUTES

4 Minutes

ATTACHMENT: FLAT HEAD

PROTOCOL: Sweep from A to B on each side for 2 minutes

DURATION: 4 minutes

HAMSTRINGS

2 Minutes

ATTACHMENT: BALL HEAD

PROTOCOL: Sweep from A to B on each side for 1 minute

DURATION: 2 minutes

CALVES

2 Minutes

ATTACHMENT: BALL HEAD

PROTOCOL: Sweep from A to B on each side for 1 minute

DURATION: 2 minutes

FEET

4 Minutes

ATTACHMENT: BALL HEAD

PROTOCOL: Sweep from A to B on each side for 2 minutes

DURATION: 4 minutes

WEEE Instructions

Correct Disposal of this product	
	This marking indicates that this product should not be disposed with other household wastes throughout the EU. To prevent possible harm to the environment or human health from uncontrolled waste disposal, recycle it responsibly to promote the sustainable reuse of material resources. To return your used device, please use the return and collection systems or contact the retailer where the product was purchased. They can take this product for environmental safe recycling.

Battery Instructions

- ⓘ Do not submerge the product into water or get it wet.
- ⓘ Do not charge the product near an open flame or in any extremely hot conditions.
- ⓘ Do not use or store the product near heat sources like furnaces or heaters.
- ⓘ If the battery leaks or emits an odor, remove it immediately from close proximity to an open flame.
- ⓘ Do not put the product into a fire or heat it up in any way.
- ⓘ Do not disassemble the battery in any way.
- ⓘ Do not use the product if, while the battery is in use or being charged, it emits an odor; gives off excessive heat; shows signs of deformation, discoloration, or any other abnormalities. If the battery is in use or being charged, remove it from the apparatus or charger immediately.

RoHS Declaration

Our Product is compliant with RoHS 2.0 (2011/65/EU Directive) regarding restriction of hazardous substances, the raw materials and components of our product meet the following toxic and hazardous substance limits.

Cadmium (Cd): < 100 ppm	Lead (Pb): < 1000 ppm
Mercury (Hg): < 1000 ppm	Hexavalent Chromium: (Cr VI) < 1000 ppm
Polybrominated Biphenyls (PBB): < 1000 ppm	Polybrominated Diphenyl Ethers (PBDE): < 1000 ppm
Bis(2-Ethylhexyl) phthalate (DEHP): < 1000 ppm	Benzyl butyl phthalate (BBP): < 1000 ppm
Dibutyl phthalate (DBP): < 1000 ppm	Dibenzyl phthalate (DBP): < 1000 ppm

A Dedicated Team Just For You

Every Bob and Brad product is backed by a 1-year warranty and an award-winning customer service team. If there are any quality problems with the product you ordered, we will try our utmost to help you get the problems resolved, as well as offer a replacement, return, or refund service under the warranty. Simply email us the problem, and we'll take care of it.

Troubleshoot

Help you solve any issues with your device.

Product Knowledge

Answer questions about your device and best practices.

Initiate Returns

Provide instructions for replacements, returns and refunds.

Place New Orders

Help with placing orders online and/or authorized resellers.

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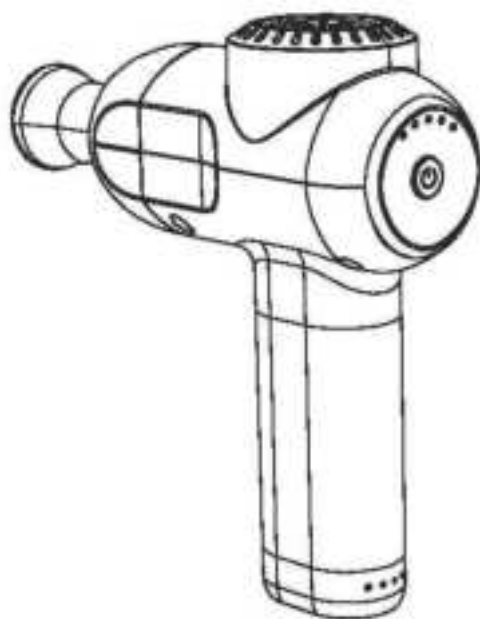
3001 Esperanza Crossing Austin, TX 78758

EC Declaration of Conformity

We (UK Importer is responsible for this declaration)

EVATOST CONSULTING LTD

(Company name)

Office 101 32 Threadneedle Street, London, United Kingdom, EC2R 8AY

(Company address)

declare under our sole responsibility that the product

Tradename: BOB AND BRAD

Equipment: Massage Gun

Model No. : QL/DMS.C2-D

ASIN: B08CKWVYMF

to which this declaration relates is in conformity with the essential requirement set out in the Council Directive on the Approximation of the laws of the Member States relating to LVD Directive(2014/35/EU) & EMC Directive(2014/30/EU) &RoHS(2002/95/EC) product is responsible to affix CE marking, the following standards were applied:

EN 60335-1:2012+A11:2014+A13:2017

EN 60335-2-32:2003+A1:2008+A2:2015

EN 62233:2008

EN 55014-1:2017

EN 55014-2:2015

EN 61000-3-2:2014

EN 61000-3-3:2013

Full Name: Taozirul

Position: Manager

Signature: Taozirul

Date: 2021.5.6

EC Declaration of Conformity

We (EU Importer is responsible for this declaration)

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(Company address)

declare under our sole responsibility that the product

Tradename: BOB AND BRAD

Equipment: Massage Gun

Model No. : QL/DMS.C2-D

ASIN: B08CKWVYMF

to which this declaration relates is in conformity with the essential requirement set out in the Council Directive on the Approximation of the laws of the Member States relating to LVD Directive(2014/35/EU) & EMC Directive(2014/30/EU) & RoHS(2002/95/EC) product is responsible to affix CE marking, the following standards were applied:

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