
Withings

Smart Body Analyzer WS-50

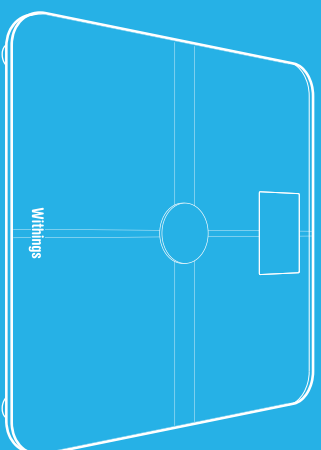Quick Installation Guide

FR Guide d'installation **CZ** Průvodce instalací
DE Installationsanleitung **IT** Guida all'installazione

NEED HELP

BESOIN D'AIDE | BRAUCHEN SIE HILFE | POTŘEBUJETE POMOC | TI SERVE AIUTO




→ <http://support.withings.com/ws50>

Thank you for choosing the Smart Body Analyzer

Merci d'avoir choisi le Smart Body Analyzer | Wie bedanken uns für ihre Wahl zum Smart Body Analyzer
Děkujeme vám, že jste si vybrali Smart Body Analyzer | Grazie per aver scelto lo Smart Body Analyzer

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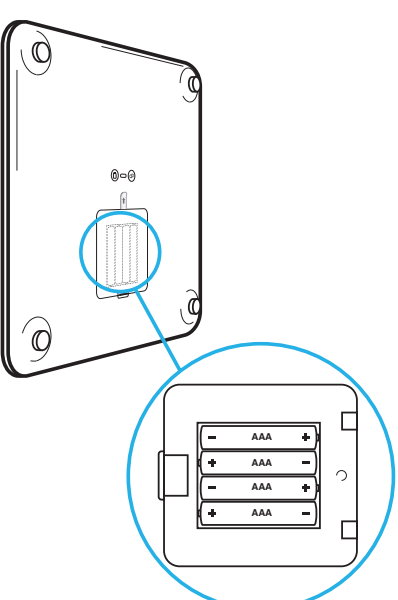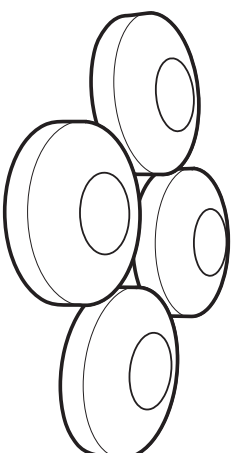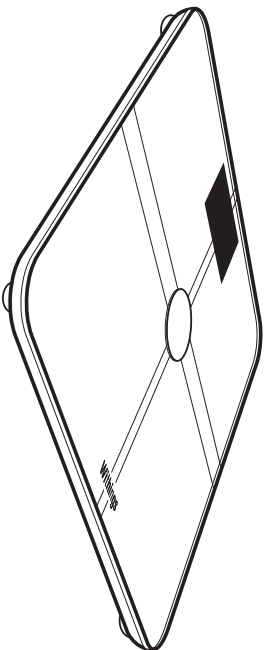
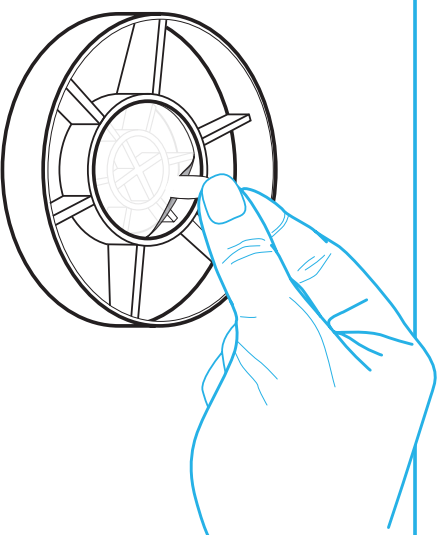
2

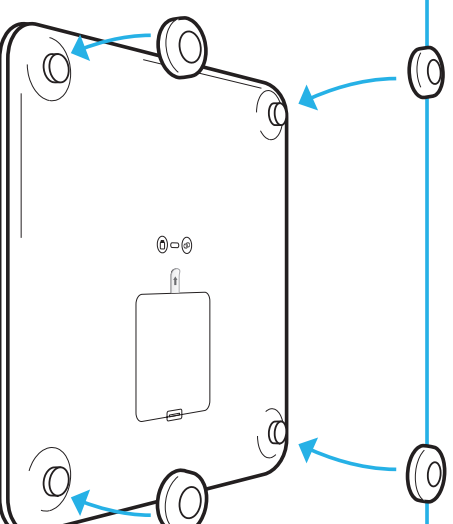
INSTALLATION


INSTALLATION | INSTALLATION
INSTALACE | INSTALLAZIONE

Before starting, you need:

FR Avant de commencer, vous devez avoir à disposition :
DE Bevor Sie beginnen, müssen Sie bereitlegen :
CZ Než začnete, potřebujete mít k dispozici :
IT Prima di iniziare, è necessario avere a disposizione:

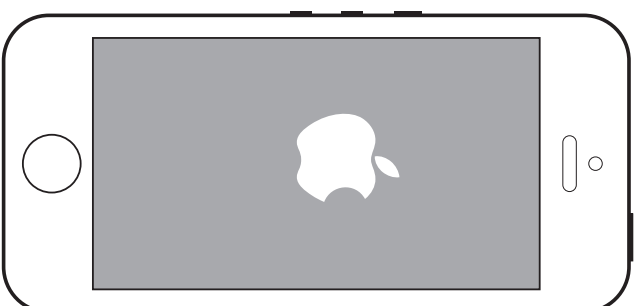iPhone,
iPad & iPod touch

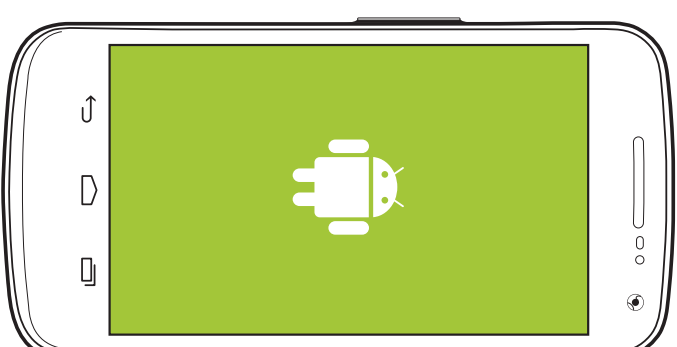Android


Installation using your smartphone / tablet only. Computer-free.

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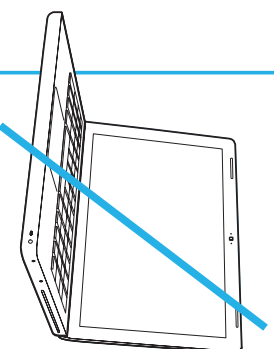
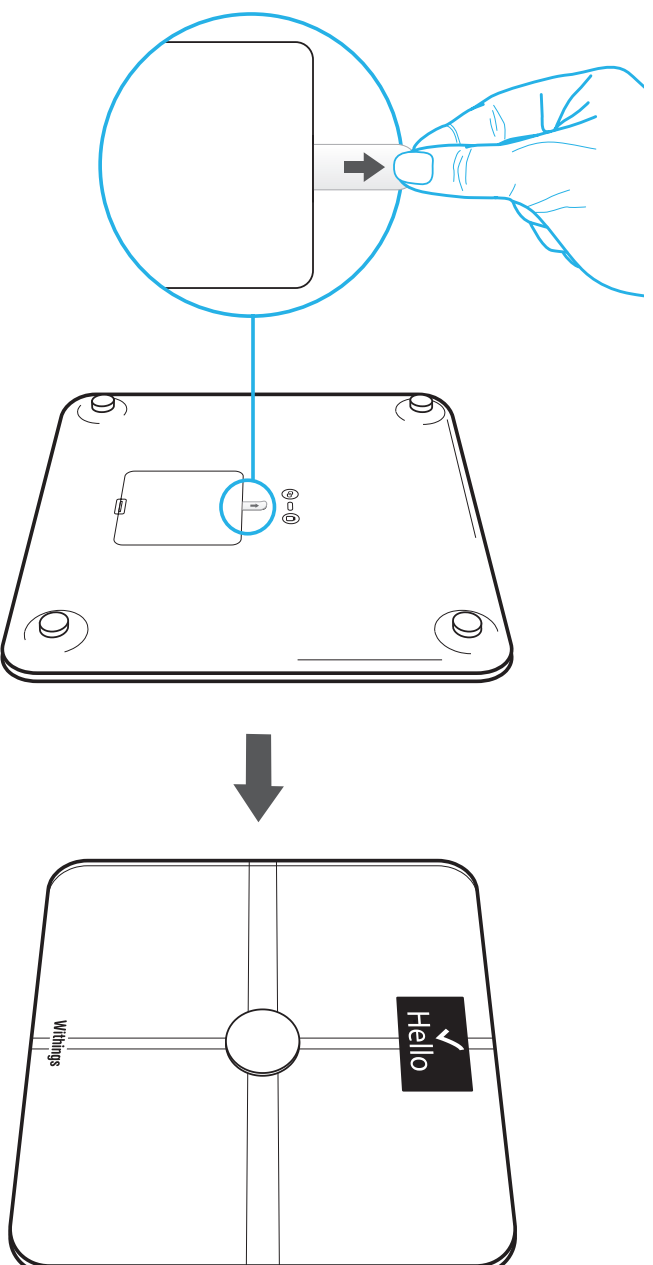
2

Press the pairing button for two seconds

FR Appuyez 2 secondes sur le bouton d'appairage
DE Drücken Sie 2 Sekunden lang die Pairing-Taste
CZ Stiskněte na 2 sekundy párovací tlačítko
IT Premere 2 secondi il tasto di accoppiamento

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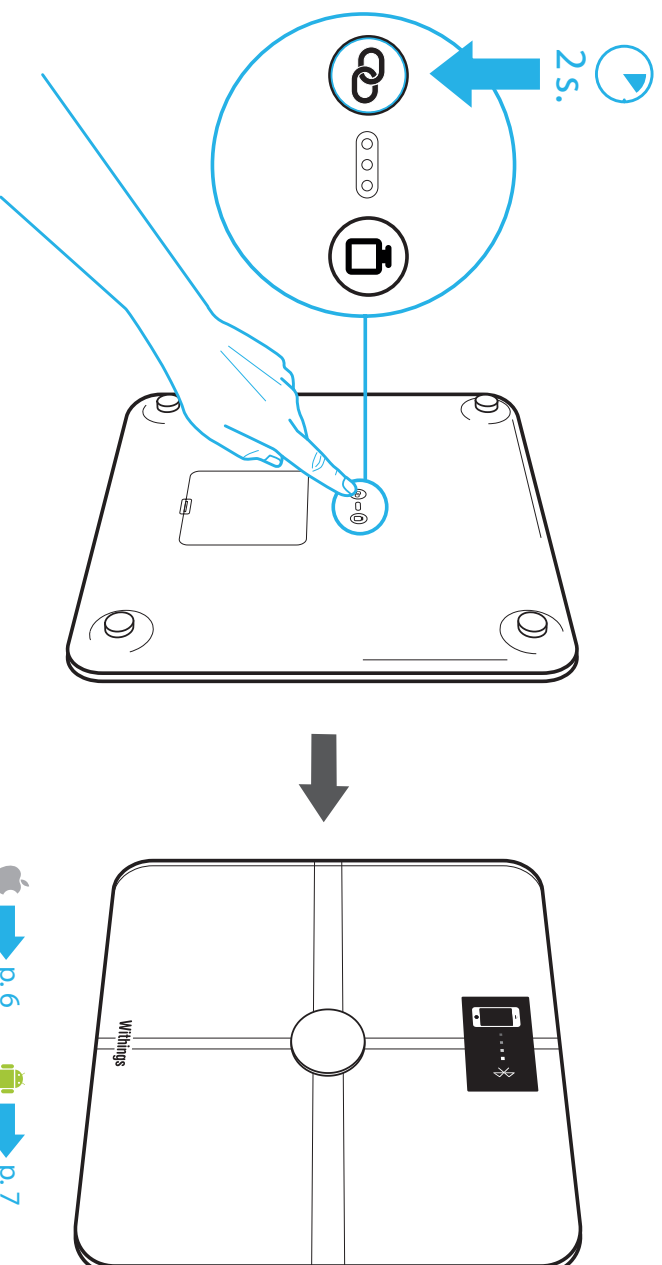IT Regolazioni / Bluetooth / ON → "Withings WS50"

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4


Install and launch the Withings app from the App Store then follow the instructions

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
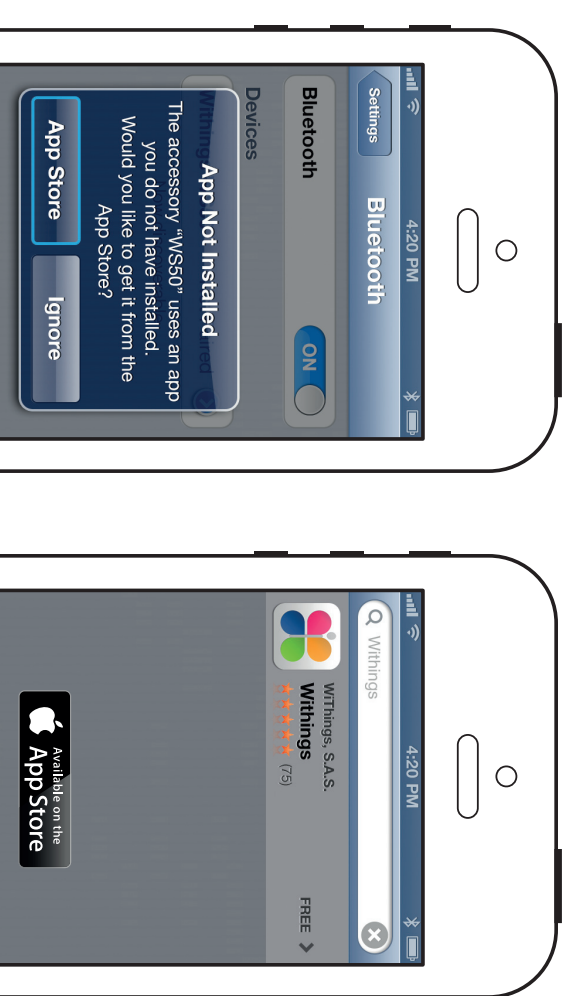
4

Install and launch the Withings app from the Google Play Store then follow the instructions

FR Installez et lancez l'appli Withings, puis suivez les instructions
DE Installieren und starten Sie die Withings-App und folgen Sie dann den Anweisungen
CZ Nainstalujte si z Google Play Store aplikaci Withings a postupujte podle pokynů
IT Installare e lanciare l'app Withings, poi seguire le istruzioni

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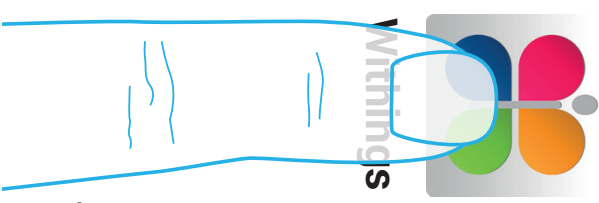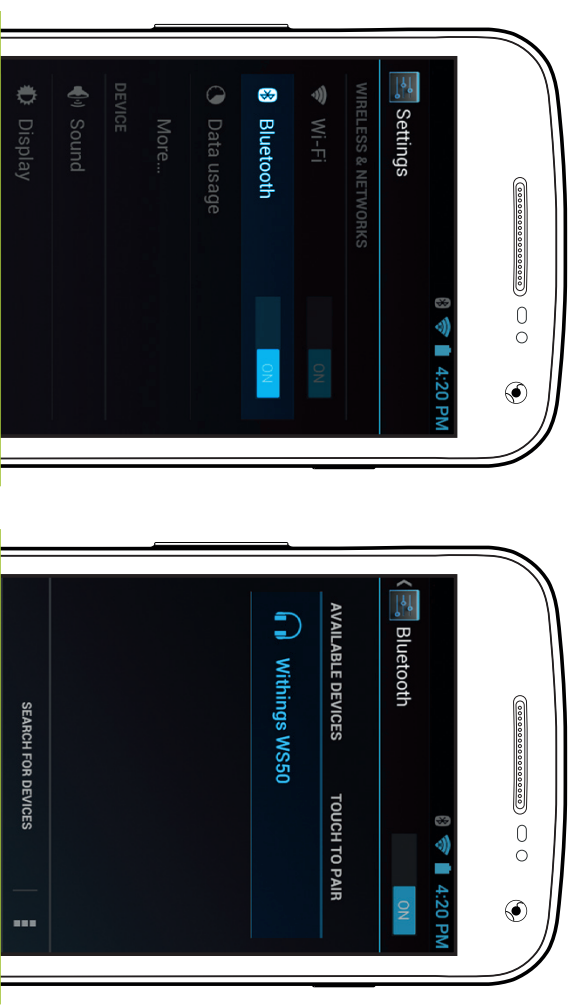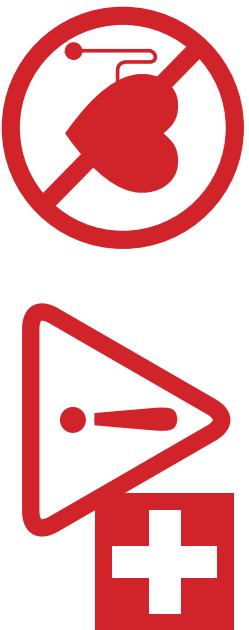
- This product is not suitable for people wearing a pacemaker or any other internal medical device.
- The Withings scale data must be considered purely as reference and not as a medical opinion.

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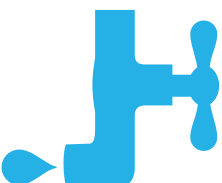
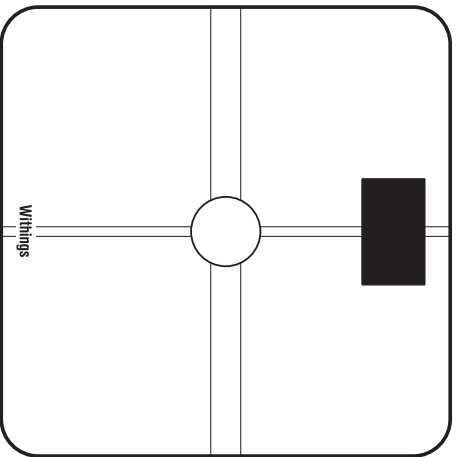
- Use only a cloth dampened with water to clean your scale.

The use of any detergents or cleaners may cause irreversible damage to the coating that protects the electrodes used to measure your body fat and heart rate.

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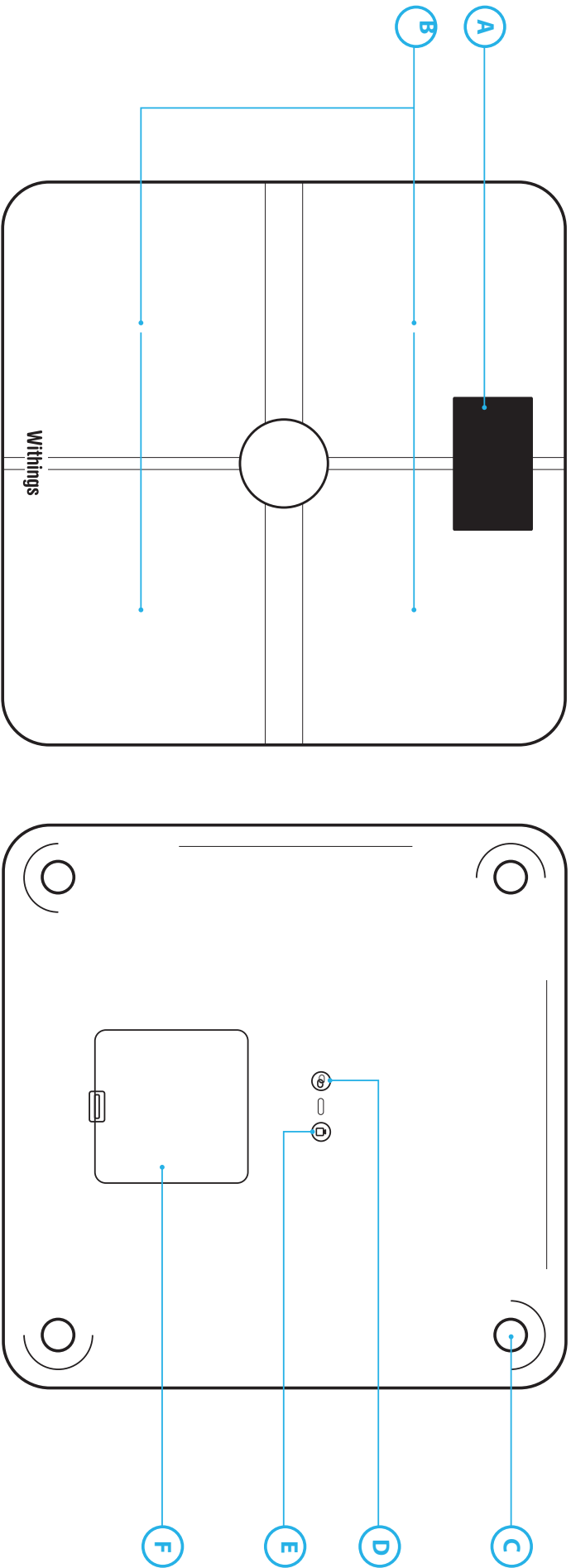
FCC ID: XNA-WS-50

This device complies with Part 15 of the FCC Rules.
Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) This device must accept any interference received, including interference that may cause undesired operation.

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Withings

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©2013 Withings. All right reserved. Všechna práva vyhrazena.

iPad, iPhone and iPod touch are trademarks of Apple inc.

Android is a trademark of Google Inc.

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iPad, iPhone a iPod touch jsou ochranné známky společnosti Apple Inc.
Android je ochranná známka společnosti Google Inc.