




# Smart Oral Irrigator User Manual

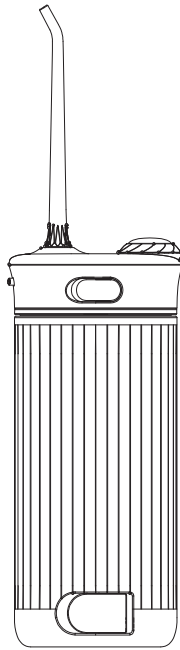

## Product Introduction

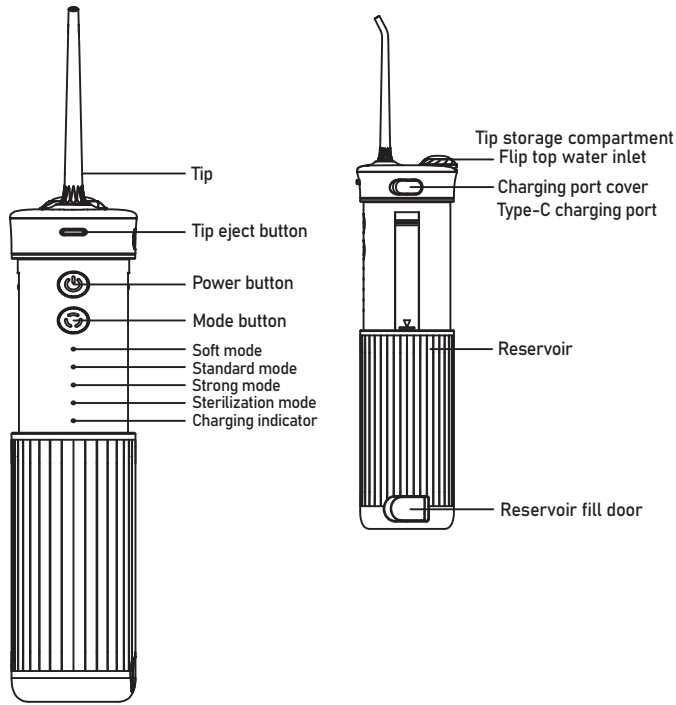

\*Replace Every 3 months  
\*\*Replace Every 6 months

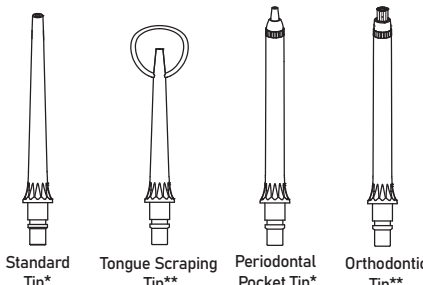

## Important Safeguards

### IMPORTANT SAFEGUARDS

When using electrical products, especially when children are present, basic safety precautions should always be followed, including the following:

#### READ ALL INSTRUCTIONS BEFORE USING. SAVE THESE INSTRUCTION.

##### DANGER:

To reduce the risk of electrocution:

- Do not place or store product where it can fall or be pulled into a tub or sink.
- Do not place in or drop into water or other liquid.
- Do not handle the user supplied USB wall charger with wet hands.
- Do not place or drop the user supplied USB wall charger into water or other liquid.
- Check the charging cord and user supplied USB wall charger for damage before the first use and during the life of the product.

## Basic Parameters

Main Material	ABS, PC	Charging	USB Type-C
Product size	140x64x45mm	Charging time	2~4 hours
Tip diameter	0.65mm	Waterproof rating	IPX 7
Rated power	10W	Rated input	5V-1A
Battery capacity	1100mAh	Water Reservoir	7.6oz/225ml
Water pressure range		50-120 PSI	
Pulse frequency		1800-2000 times/min	
Package Content		Unit*1, Tip*4, Charging cable*1, Manual*1	

## Warning

### WARNING

To reduce the risk of burns, electrocution, fire, or injury to persons:

- This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.

- Children should be supervised to ensure that they do not play with the appliance.
- Use this product only as indicated in these instructions or as recommended by your dental professional.
- Do not use this product if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into any liquid. Contact the seller.
- Keep the charging cord and user supplied USB wall charger away from heated surfaces.
- Never drop or insert any object into any opening or hose.
- Do not use outdoors or operate where aerosol (spray) products are being used or where oxygen is being administered.
- Do not expose the battery or appliance to excessive temperatures.
- A user supplied USB-A wall charger (power adapter) with the following specifications is required to charge this device. Rated Output: **5.0 VDC, 1.0 A**
- Do not plug this device into a voltage system that is different from the voltage system specified on the device or charger.
- Unplug USB wall charger when not in use.
- Do not direct water from the Water Flosser under the tongue, into the ear, nose or other delicate areas. This product is capable of producing pressures that may cause serious damage in these areas.
- If your physician advised you to receive antibiotic premedication before dental procedures, you should consult your dental professional before using this product or any other oral hygiene aid.
- Fill reservoir with warm water or other dental professional recommended solutions only.
- Do not use iodine, saline (salt solution), or water insoluble concentrated essential oils in this product. Use of these can reduce performance and will shorten the life of the product.
- Remove any oral jewelry prior to use of this product.
- Do not use if you have an open wound on your tongue or in your mouth.
- There are no consumer serviceable electrical parts in this product and it does not require maintenance.
- Remove charging cord from the unit and replace charging port cover prior to filling the reservoir.

## Before Get Started




### Charging Unit Prior to First Use

Charge unit for approximately 4 hours prior to first use. Plug the user supplied USB wall charger into a wall outlet and connect the USB charging cable to the unit (see images). Refer to the IMPORTANT SAFEGUARDS section for required USB wall charger specifications

### Ensure contacts on BOTH the cord and the unit are clean and dry, and the reservoir is empty and collapsed.


The LED light will emit a slow glow on/off when unit is charging. Once the unit is fully charged and the charger is connected the LED will remain on until the charger is disconnected. A full charge will last up to 4 weeks of typical use. When the batteries are in need of charging, the charge indicator light emit a rapid triple blink during and after use to indicate that you should recharge the unit. A complete recharge will take 2-4 hours. If you use your Water Flosser once a day or less, it is recommended you charge every 3-4 weeks. Put the protective charging port cover back into place before use. **The charging port cover should only be removed for charging.**

## Getting Started

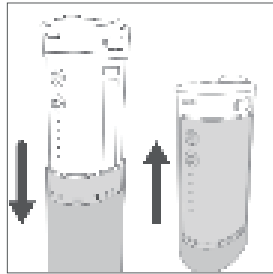


### 1.Filling the Reservoir

Remove the charging cable from the unit and replace charging port cover prior to filling reservoir and product use. Pull down on reservoir to open. Unit will "click" into place when fully extended. Lift flip top on reservoir and fill the reservoir with lukewarm water. Collapse by pushing reservoir up. Reservoir will "click" into place when fully closed.

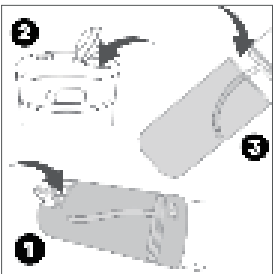

### Three ways to fill water:

1. Reservoir fill door
2. Flip top water inlet (tip storage compartment)
3. Remove the reservoir from the unit

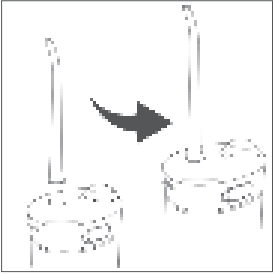


### 2.Inserting and Removing Tips

Insert tip at the top of the unit. Press firmly until the tip clicks into place. To remove the tip from the unit, press the tip eject button and pull tip from the unit.


#### CAUTION:

**Do not press tip eject button while unit is running.**

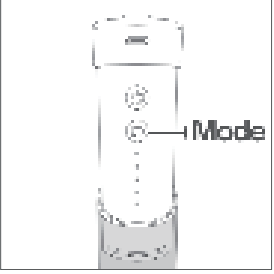


### 3.Adjusting the Pressure Setting


1. Three modes:
  - Soft, Standard, Strong pressure settings.
2. Unit will remember last setting when turned off.
3. Unit will cycle through the modes when the button is repeatedly pressed.

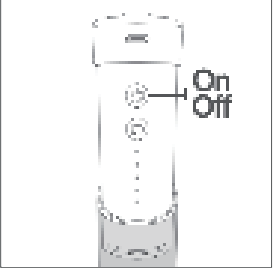


### 4.Press POWER button to start flossing. (automatically power off in 70s)

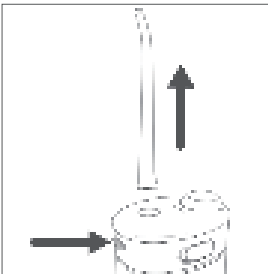


### 5.Hold the tip eject button to remove tip.

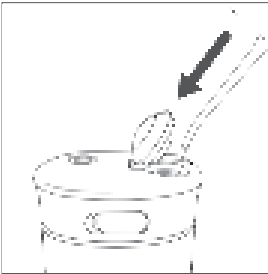

### 6.Put the tip into tip storage compartment.

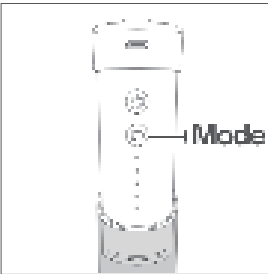


### 7.Press MODE button to Sterilize. (automatically turn off after 90s)

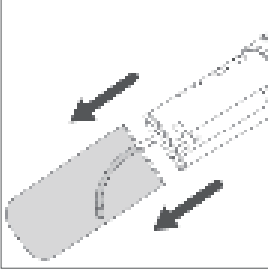

### 8.When Finished

Turn the unit OFF. Empty any liquid left in the reservoir. If desired, the reservoir may be removed from the power handle by simply pulling down on the reservoir to remove; Or just open the fill door on the reservoir.

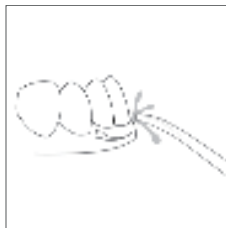

### Recommended Technique

Lean low over sink and place tip in mouth. Aim the tip toward gumline, and turn the unit ON. For best results, start with the back teeth and work toward the front teeth. Glide tip along gumline and pause briefly between teeth. Continue until you have cleaned the inside and outside of both the upper and lower teeth. Direct the jet stream at a 90-degree angle to your gumline. Slightly close lips to avoid splashing but allow water to flow freely from mouth into the sink. Keep unit upright during use for best results.

## Tip Usage




#### Periodontal Pocket Tip

The Periodontal Pocket tip is specifically designed to deliver water or anti-bacterial solutions deep into periodontal pockets. To use the Periodontal Pocket tip, set the unit to the soft mode. Place the tip against a tooth at a 45-degree angle and gently. Place the tip under the gumline, into the pocket.



#### Tongue Cleaner Tip

Set to soft mode, place tip in the center/middle of your tongue about half way back. Pull forward with light pressure. Increase pressure as you prefer.

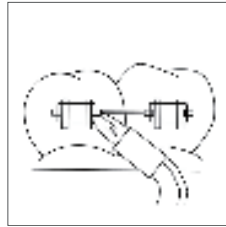

#### Orthodontic Tip

To use the orthodontic tip, gently glide tip along gumline, pausing briefly to lightly brush area between teeth and all around orthodontic bracket, before proceeding to the next tooth.

## Using Mouthwash and Other Solutions

The water flosser can be used to deliver mouthwash and antibacterial solutions. For ease of adding other solutions, remove the reservoir (see "When Finished" section) and add from opening on reservoir. If you are adding additional water you may replace the reservoir and fill reservoir fully with the fill door. After using any special solution, rinse unit to prevent clogging by partially filling reservoir with warm water and running unit with tip pointed into sink until unit is empty.

## Cleaning

### Never immerse unit in water.

Before cleaning, unplug from the electrical outlet. Clean the product when necessary by using a soft cloth and mild nonabrasive cleanser to wipe exterior. The water reservoir is removable for easy cleaning. It is recommended that the reservoir is removed and cleaned weekly.

### Removing Hard Water Deposits / General Cleaning

Hard water deposits may build up in your unit, depending on the mineral content of your water. If left unattended, it can hinder performance. Cleaning internal parts: add 1 tablespoon of white vinegar to a full reservoir of warm water. Point the handle and tip into sink. Turn unit ON and run until reservoir is empty. Rinse by repeating with a full reservoir of clean warm water. This process should be done every 1 to 3 months to ensure optimal performance.

## Battery Care, Disposal/Recycling/Removal

### Battery Care

If the Water Flosser is to be stored for an extended period of time (six months or more), be sure to fully charge it prior to storage. Do not expose the battery or appliance to excessive temperatures.

## Disposal/Recycling

This product contains a rechargeable lithium ion battery. The battery is not replaceable. At the end of product life, do not dispose of in household garbage or recycling bins. Please recycle at your nearest recycling center or take to hazardous waste collection point to comply with local and state regulations. Please contact your local waste administrative department for any questions. Note: The rechargeable battery is not replaceable. Remove only at end of product life. This process is not reversible and voids the product warranty.

### Battery Removal

To remove the rechargeable battery for recycling, please discharge the battery fully before removing. Observe basic safety precautions when you follow the procedure below. Be sure to protect your eyes, hands, fingers and surface on which you work. Do not disassemble this product unless you want to discard it, follow the basic safety precautions for operation, and disconnect the power supply before removing the battery. Make sure that the battery is exhausted before removing the battery.



- 1.Insert a cutter to pry off the top of the unit.
- 2.Use a small flat-head screwdriver to pry open the screw cork.
- 3.Use a philips screwdriver to loosen 4 screws inside the hole and 1 screw beside the water tube.
- 4.Slide battery out of chassis. Cut wires individually (avoid cutting positive and negative at the same time to prevent shorting).

## Troubleshooting

FAQ	Problem Causes	Solution
Press the on/off button, the motor does not work and the indicator light does not light up	The battery is exhausted and the device cannot be turned on	Must charge the product for 4 hours
The product has low water pressure	Inlet pipe filter is clogged	Clean the water inlet filter with a sharp object such as a toothpick
	Tip is clogged	Clean the tip or replace the tip
	Working in a low pressure mode	Press the mode button to an appropriate mode
Plug in the charging cable, the battery power indicator does not light up, and the product cannot be turned on	The battery is exhausted and cannot be turned on	Plug in the charging cable with good contact and charge for 2 hours, and check whether the function returns or not.
Plug in the charging cable, the battery power indicator does not light up, and the product can be turned on	Poor charging contact	Re-plug and connect the charging cable, or directly replace the adapter and charging cable and try charging again
The device is overheating and the lights is flashing	Do not use the Water Flosser for more than 5 times continuously.	To prevent overheating, please allow the device to cool down for at least 30 minutes before using it again.

If following the above method does not solve the problem, please contact our customer service for assistance.  
PS. When using the device for the first time or after a long period of non-use, you may experience the following phenomenon: After pressing the power button, you will need to wait for a moment before the tip can spray water. In this case, please follow the instructions in the illustration:



Add an appropriate amount of water to the tip installation port, reinstall the tip, and then turn on the device.

## 12 Months Warranty

12 months warranty against materials and manufacturing defects from the date of purchase.We will replace at its discretion the product, provided the product has not been abused, misused, altered or damaged after purchase, was used according to instructions, and was used only with accessories or consumable parts approved by the seller. This limited warranty excludes accessories or consumable parts such as tips, etc. Use with improper voltage will damage product and VOID warranty.