
HYUNDAI

Electric toothbrush USB charging
(User Manual)

HHA732301 / HHA732302



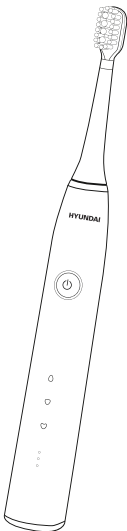
GENUINE PRODUCT OF
HYUNDAI CORPORATION

PRECAUTIONS

- Check the power cord regularly, and do not use the appliance in case of damage or failure.
- This appliance is not suitable for children under the age of 14.
- It is recommended that the appliance shall be placed out of reach of children.
- Do not place the handle in water or other liquids.
- Do not insert any item in the handle when the brush head is extracted.
- Do not touch the plug with a wet hand to avoid electric shock.
- Do not repair and modify the appliance by yourself, or use the accessories or chargers of other products.
- If you are in oral therapy, please consult your dentist before using the appliance.

PRODUCT DESCRIPTION

Electric toothbrush take care of your teeth and gums, cleaning them with outstanding effect.

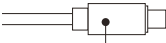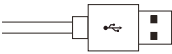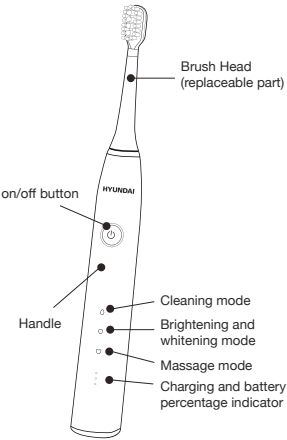Type C Cable

CHARGING

Please connect the A-end plug of the USB power cord to the adapter into the power socket, and then connect the handle to the Type C plug of the USB power cord to charge. The first charging should continue for 2 hours. This appliance can be used continuously for about 30 days when it is fully charged (calculated as 2 times/day, 2 minutes each time). During daily use, please keep the Type-C charging port of the handle clean and avoid toothpaste and other stains from entering it. This appliance is water-resistant for the use in the bathroom.

HOW TO USE THE APPLIANCE

Squeeze a proper amount of toothpaste onto the wetted bristles.

To avoid splashing the toothpaste, please place the toothbrush in the mouth before pressing it on.

According to the Bass technique, it is recommended that the brush head be aimed at the teeth obliquely at an angle of 45°, and the bristles be gently and slowly moved, staying on each tooth for 2-3 seconds


Brush your teeth for at least 2 minutes each time, and then press it off.


CLEANING ADVICE

After use, keep the appliance operating and place the brush head under running water for a few seconds to clean, then turn off the appliance. Clean the handle with a slightly damp cloth and do not use chemical cleaners

BRUSH HEAD REPLACEMENT

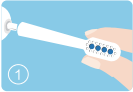
When replacing the brush head, please extract the old brush head firstly

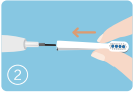Insert the new brush head in a proper direction

Insert the new brush head to the very bottom till it's fixed

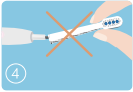Please do not insert it from an inclined direction

SWITCHING METHOD OF THE THREE MODES

1. Take out the Electric toothbrush

2. Press On/Off button (Default Mode 1)

3. Release On/Off button

4. Press On/Off button
Switch to Mode 2/3
(Cycle: Mode 1-Mode 2-Mode 3-Off-Mode 1...)

SMART POWER MANAGEMENT

- During charging: Three LED lights flicker and display in green.
- Fully charged: Three LED lights are bright green.
- Smart identification of battery percentage when operating:
Three LED lights are turned off one by one, indicating the current battery percentage
- Power shortage prompt function:
When the appliance is in power shortage (less than 20% of the battery), the charging LED light automatically flashes in orange, prompting the user to charge.
- Low battery protection function:
When the battery voltage is extremely low, in order to protect the battery from being over-discharged, the appliance will be automatically disabled.

SMART TIMING FUNCTION

Pause (0.5 seconds) after operating for 30 seconds to remind to change position and stop operating automatically in 2 minutes.

SMART MEMORY FUNCTION

The modes can be changed with the button within 5 seconds after turn on the appliance. It can be turned off after 5 seconds and remember the current working mode. The next time the machine is turned on, the working mode is automatically entered

SPECIAL TIPS

In the first few days of using the Electric toothbrush, there may be a slight bleeding in the gums, which usually stops after a few days .If the gum continues bleeding for 2 weeks please stop using it immediately and seek medical attention.

RECYCLING

1.Environment friendly use period

Remarks: All electric and electronic products sold in China must be marked with the Environment Friendly Use Period. The Environment Friendly Use Period only refers to the use period of the product under normal use conditions, temperature, and humidity, instead of the service life of consumables like the battery.

2. Declaration of hazardous substance

Note: Please do not dispose the appliance with general domestic waste, but hand it over to the recycling center designated by the government to achieve environmental protection.

Danger: There is a built-in rechargeable battery in the toothbrush handle. Do not heat it, throw it into fire, or place, charge or use it in a high-temperature environment, otherwise overheating, fire or explosion will be caused.

© 2023 The Investment Company Europe B.V.
HYUNDAI is a trademark of HYUNDAI Corporation holdings, Korea. Licensed to The Investment Company Europe B.V. All rights reserved.

Licensed by Hyundai Corporation Holdings, Korea.
Imported/Distributed by The Investment Company Europe B.V. exclusively developed for sales in the Netherlands, Belgium, and Luxembourg.

PATENT PENDING

The artwork and design of this package are protected by European copyright law and may not be reproduced, distributed, displayed, published or used for any purpose without prior written permission. It is not permitted to alter, or reproduce any trademarks or copyright notice from this package.

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Lizenziert von HYUNDAI Corporation Holdings, Korea