human touch®


humantouch.com
800.355.2762

©2017 Human Touch®, LLC. Patents pending and intellectual property rights enforced.
No medical claims are warranted or implied by the use of this product.
Human Touch is a trademark of Human Touch®, LLC.
iJoy-MA-NA-A3

Register your product **ONLINE** today
www.humantouch.com/register

iJOY®
massage anywhere

**A POWERFUL,
PORTABLE MASSAGE
ON THE GO.**

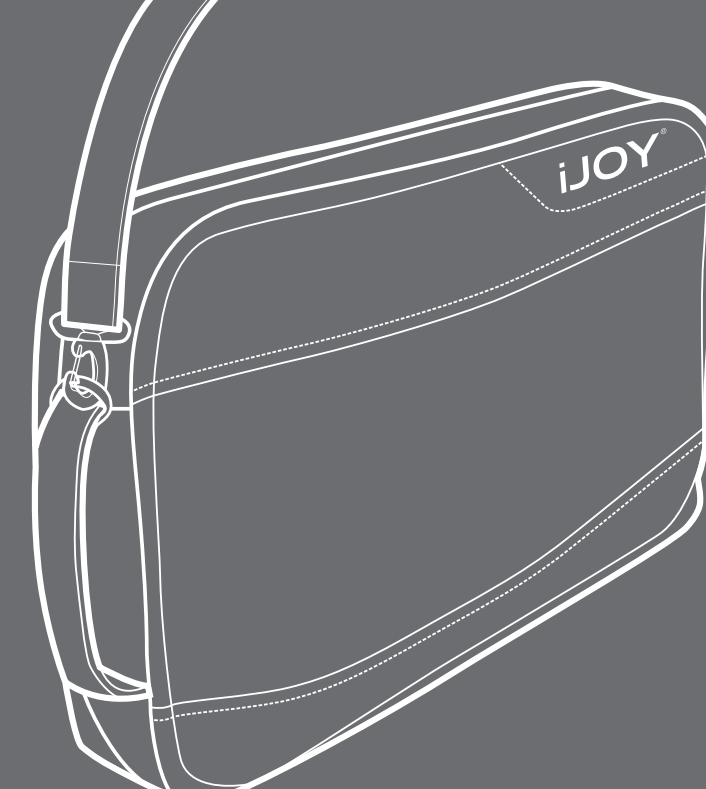 **human touch**
be your best. **feel your best.™**

use & care manual

Register your product **ONLINE** today
www.humantouch.com/register

ENG

USE & CARE MANUAL

FRC

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ESN

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For **more than 35 years**, Human Touch has established itself as an industry leader that is changing the lives of its customers with products and services that help them **feel, perform, and live their best**. Developed in cooperation with the medical community and industry experts, we deliver products that blend patented technology, craftsmanship and intelligent ergonomics with superior design to enable therapeutic relief for everyone.

Our industry-leading massage chairs, Perfect Chair® recliners and evolving set of innovative solutions share attributes in their ease-of-use, attention to performance and quality, and design to adapt to most any lifestyle. These attributes are highly valued by our customers and garner broad industry award and recognition. As a result, Human Touch empowers people to enjoy greater health and vitality in their professions, pursuits and everyday lives. And with a population increasingly aware of the benefits of wellness, Human Touch extends its brand leadership and available **solutions for the benefit of everyone**.

Proudly endorsed by:

welcome to **human touch**®

Your new Human Touch® iJoy® Massage Anywhere is the powerful, portable solution for a massage anytime and anywhere you need to recharge and soothe tired muscles. Its revolutionary FlexGlide® orbital massage technology delivers a stress-reducing experience like never before, offering a seamless glide to help you bask in the glow of ultimate relaxation.

The sleek, compact Massage Anywhere design moves from indoor to outdoor use with a rechargeable lithium-ion battery for two hours of cordless massage and an easy-to-carry handle and adjustable shoulder strap for maximum portability.

Experience awesome techniques used by massage and chiropractic professionals—all at the touch of a button. To ensure you don't overlook any of your product's features, and that you use it safely, read this manual carefully. Then start to enjoy the benefits of your Massage Anywhere—anywhere you go.

keep it safe

IMPORTANT SAFETY INSTRUCTIONS

When using an electrical appliance, basic precautions should always be followed, including the following:

READ ALL INSTRUCTIONS BEFORE USING THIS APPLIANCE.

 **DANGER:** To reduce the risk of electrical shock Always unplug this product immediately after using and before cleaning.

 **WARNING:** To reduce the risk of burns, fire, electric shock or injury to persons:

- Never use any power source other than a standard 120 volt wall outlet.
- Do not put hands or fingers near the massage mechanism or rollers while this product is running. Touching the mechanism can result in injury from the rollers squeezing together.
- Do not operate under blanket or pillow. Excessive heating can occur and cause fire, electric shock, or injury to persons.
- Do not use on an infant, invalid, sleeping, or unconscious person.
- Use this product only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.
- Never operate this product if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water. Return this product to a service center for examination and repair.
- If you feel any discomfort, discontinue use immediately.
- Never drop or insert any object into any opening.
- Do not operate where aerosol (spray) products are being used or where oxygen is being administered.

- Do not massage any area of the body that is swollen, inflamed or covered with skin eruptions.
- Do not use on unexplained pain.
- Do not use this product when standing on a damp floor, or when any part of the body is in contact with plumbing or any similar ground.
- Do not stand on appliance.
- Anyone with a medical condition, including those who are or may be pregnant, should consult a physician before using any massage product.
- Unplug the product before cleaning.
- Do not operate this product while driving.
- Comfortable clothing is recommended while using this product.
- Do not sit on the controller.
- When using this product after it has been in storage, check to see that the product is operating properly prior to use.
- Pay attention to any abnormal condition, such as smoke or odor. Should any abnormal condition occur, turn off the product immediately and discontinue use.
- If you experience any shock, stop using the product immediately.

Do not attempt to repair this product yourself.

- This product is designed for personal in-home use. Commercial use will void warranty.

SAVE THESE INSTRUCTIONS

This product is a rechargeable Li-ion battery that has been certified to the UL2054 and CSA E62133 standards.

- Do not dismantle, open or shred secondary cells or batteries.
- Do not expose cells or batteries to heat or fire. Avoid storage in direct sunlight.
- Do not short-circuit a cell or a battery. Do not store cells or batteries haphazardly in a box or drawer where they may short-circuit each other or be short-circuited by other metal objects.
- Do not remove a cell or battery from its original packaging until required for use.
- Do not subject cells or batteries to mechanical shock.
- In the event of a cell leaking, do not allow the liquid to come in contact with the skin or eyes. If contact has been made, wash the affected area with copious amounts of water and seek medical advice.
- Do not use any charger other than that specifically provided for use with the equipment.
- Observe the plus (+) and minus (–) marks on the cell, battery and equipment and ensure correct use.
- Do not use any cell or battery which is not designed for use with the equipment.
- Do not mix cells of different manufacture, capacity, size or type within a device.
- Battery usage by children should be supervised.
- Seek medical advice immediately if a cell or a battery has been swallowed.
- Always purchase the battery recommended by the device manufacturer for the equipment.
- Keep cells and batteries clean and dry.
- Wipe the cell or battery terminals with a clean dry cloth if they become dirty.
- Secondary cells and batteries need to be charged before use. Always use the correct charger and refer to the manufacturer’s instructions or equipment manual for proper charging instructions.
- Do not leave a battery on prolonged charge when not in use.
- After extended periods of storage, it may be necessary to charge and discharge the cells or batteries several times to obtain maximum performance.
- Retain the original product literature for future reference.
- Use only the cell or battery in the application for which it was intended.
- When possible, remove the battery from the equipment when not in use.
- Dispose of properly.
- Use only the charger model QD08-160040-A to charge the batteries.

Note: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:


- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

FCC ADAPTOR STATEMENT This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.


 **WARNING:** Changes or modifications to this unit not expressly approved by the party responsible for compliance could void the user’s authority to operate the equipment.

get ready

Unpack the box

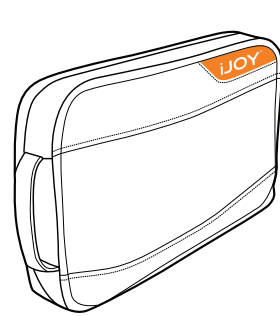
Massage Anywhere

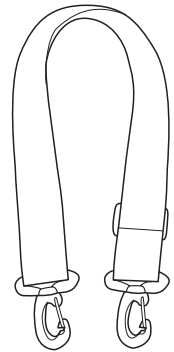Shoulder Strap

Battery Pack

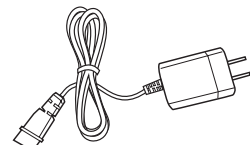Charger

1. Charge it


IMPORTANT: You must charge the battery for a minimum of eight hours before initial use.

- Connect the battery pack to the charger and plug the charger into a standard 120 volt wall outlet. The charging LED on the battery turns red. When the battery is fully charged (after about eight hours), the LED turns off.
- Connect the battery to the battery compartment cable, insert it fully into the compartment, then close the zipper to secure it.

IMPORTANT: See *Understanding Battery Life* on page 11 for additional battery information.

2. Power on



Turn the On/Off switch to On.

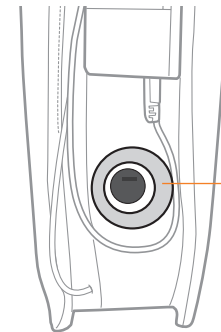On/Off switch

Grab and Go!



Convenient handles on each end make it easy to pick up your massager and go. Or, attach the adjustable shoulder strap and adjust the length to carry your massager on the side, like a briefcase.

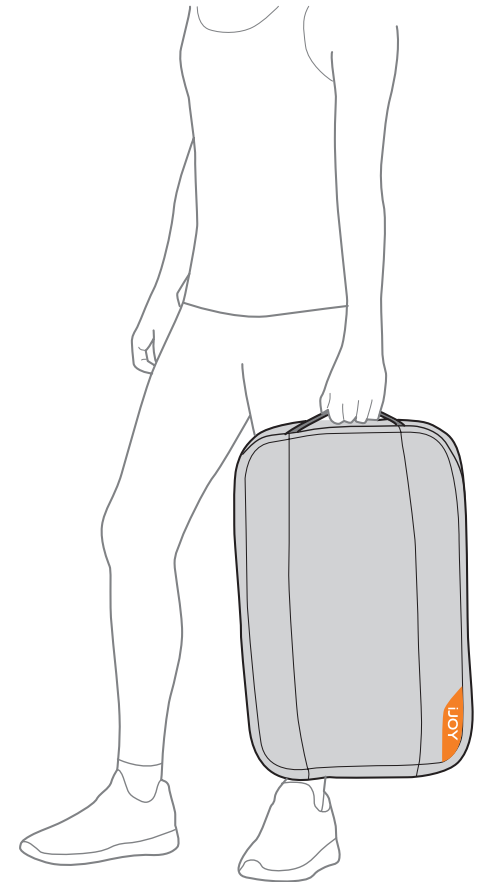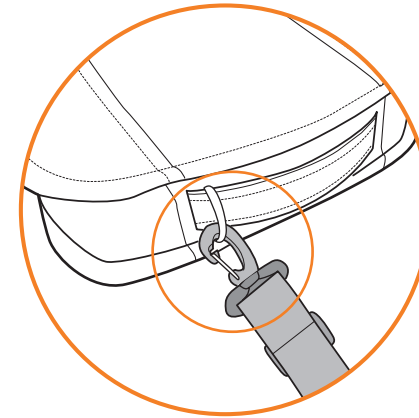Your Massage Anywhere comes fully assembled. Simply remove and charge the battery, reinstall it, and you're set to go!

choose a position

- A** Place the massager vertical against a wall and lean back for a penetrating full back massage.
- B** Place the massager horizontal against a wall, select Kneading or Compression, then lean back for a deep-kneading lower back massage.
- C** Place the massager on the ground to use on your thighs or calves.

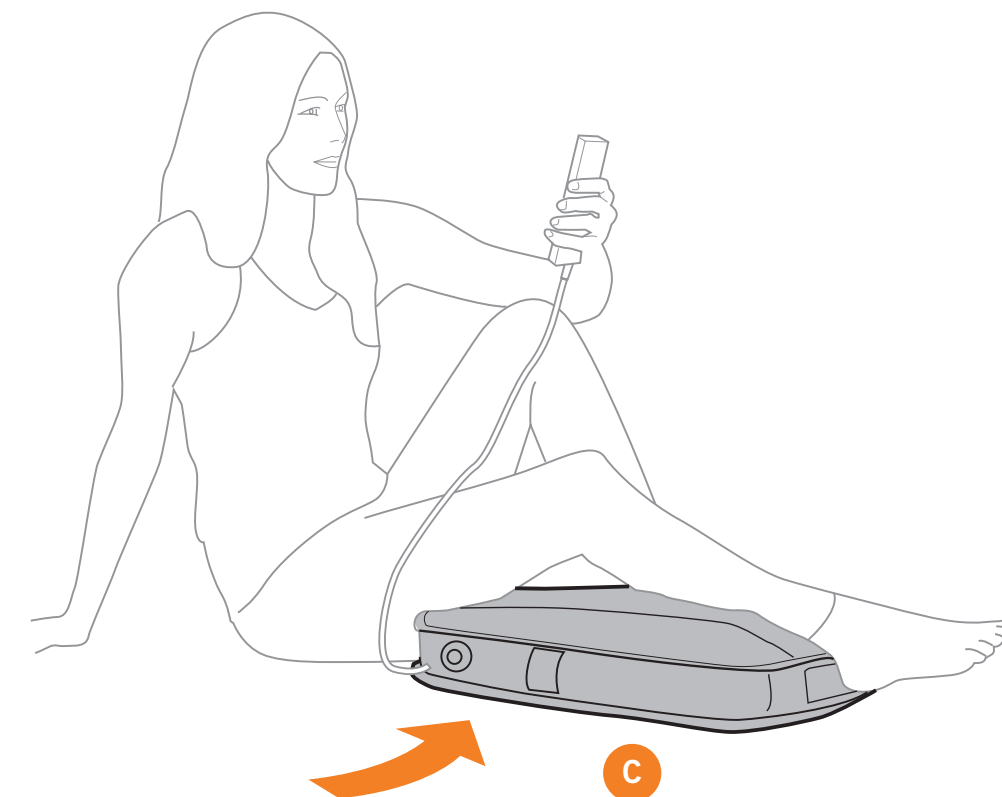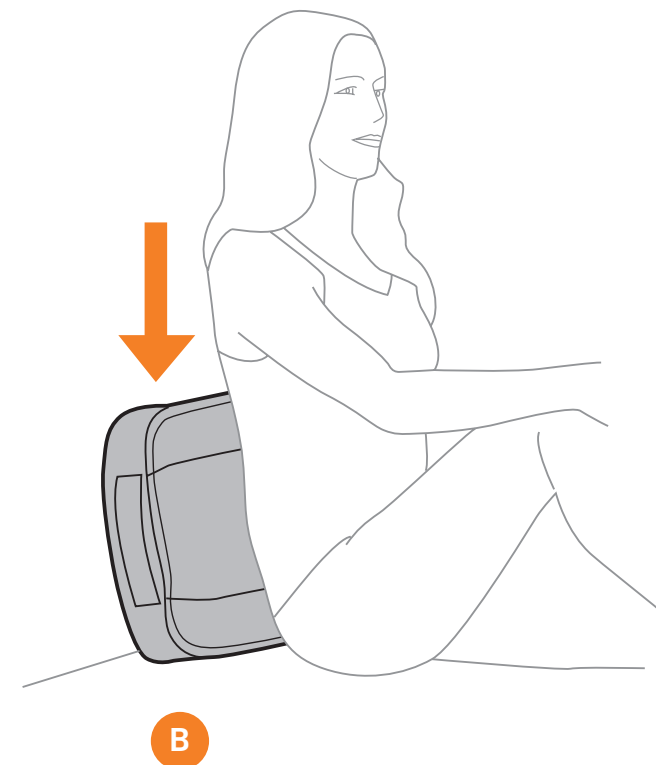
Use your versatile **Massage Anywhere** in a variety of angles and positions to get the perfect massage.

take control

Stop — Immediately stops the massage.


Play — Press to toggle through three auto programs, which combine massage techniques in different ways for ultimate relief.

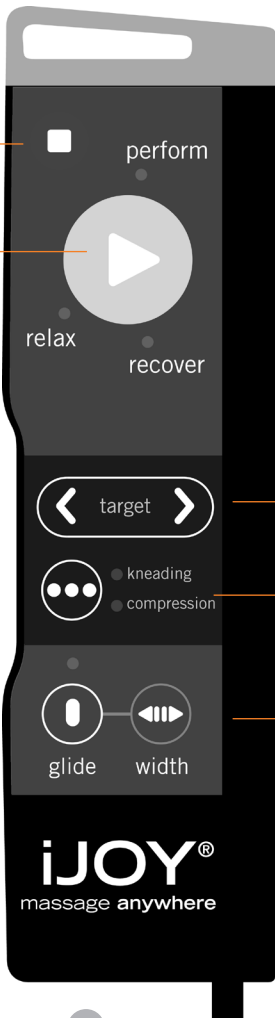Target — When enjoying a Kneading or Compression massage, adjust the massage position to just the right spot.

Width — Focus or diffuse massage pressure during a Glide massage. Select the width setting after selecting Glide.

 *Your Massage Anywhere automatically limits your massage to an ideal 15-minute interval. To restart, simply re-select your favorite massage setting.*

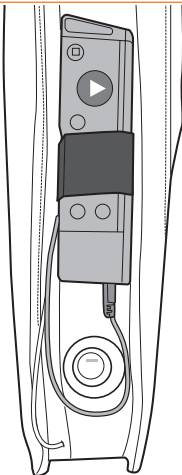Attaching the remote to the massager.
Slip the controller into the elastic loop until the loop snaps into the notch. Keep the cord tidy by looping and tucking it in behind the controller.

Using your controller while mounted on the massager
Just press the **Play** button to start the massage. Or, remove it from the elastic loop for hand-held control.

Now choose a massage technique to suit your mood. Your iJOY **Massage Anywhere** lets you customize your massage using a wide variety of massage techniques. Choose from two stationary massage techniques and a rolling massage or select one of the three auto massage programs and let the massager take control. No matter how you use it, your massager offers maximum massage flexibility and refreshing benefits.

Understanding Targeted Benefits

Targeted Benefit	What It Does	Benefits
Kneading 	Kneads your back in a small, circular motion on both sides of your spine simultaneously.	Improves deeper blood and lymph circulation and helps bring vital nutrients to the spinal area, also helps clear out the toxins from muscle and nerve tissue in the area.
Compression 	Kneads your back in a side-to-side motion on both sides of your spine simultaneously.	Relieves tension and soreness by lifting and stretching muscles. Improves circulation and helps bring vital nutrients to the spinal area.
Glide 	Rolls gently up and down your back near your spine.	Warms and loosens your muscles, preparing your back for deeper massage. Relieves muscle tension, improves circulation and relaxes the back.

care and maintenance

The product is specially designed and constructed to be maintenance free. No periodic lubrication or servicing is required. Your product should only be serviced by an authorized service center. Foreign objects should not be inserted between the massage nodes or the unit housing. Care should be taken not to overload the unit.

STORAGE:

Your product has been carefully designed to give you many years of trouble-free service. If not in use for a long period of time, we recommend that the battery be unplugged, and that the unit be placed in an environment free of dust and moisture. Do not store near heat or open flame. Do not leave this product exposed to direct sunlight for long periods of time. This may cause fading or damage.

CLEANING:

Clean with a slightly damp cloth. Never soak or submerge in water.

Understanding battery life

- After initial use, a fully drained battery must be charged for approximately twelve hours to provide maximum duration massage function.
- A fully charged battery will provide approximately two hours of massage under normal load (e.g., leaning against the massager); continuous heavier load use (e.g., laying on top of the massager) will reduce run-time.
- If the unit is fully charged and stored with the power switch set to ON, the battery will fully drain in approximately one week.
- If the unit is fully charged and stored with the power switch set to OFF, the battery will fully drain in approximately eight months.

Specifications*

CHARGER:

Class 2, Model No. QD08-160040-A

Input: AC 120V 60Hz 8.0W

Output: 16VDC 400mA

Timer: 15 Min.

Product Weight: 10.5 Lbs.

Maximum Load Weight: 200 Lbs

BATTERY:

Rechargeable Li-Ion 12.8V 3Ah 38.4Wh

SPEC-26650-3000mAh X4C 41FpR27/65

**Specifications are subject to change without prior notice.*

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