

# Fairywill



   @FairywillCare

 [www.fairywill.com](http://www.fairywill.com)

 [support@fairywill.com](mailto:support@fairywill.com)

## E10

Sonic Electric Toothbrush


Thank you for choosing the Fairywill E10 sonic toothbrush. Please read the instructions carefully before using the product, and keep the user manual in a safe place for future reference.

## 1. Features:

- A Waterproof and non-slip design.
- Replaceable brush heads with unique 3D shape bristles design and a bristle replacement reminder, to help keep your teeth whiter and healthier.
- The high-frequency vibration effectively removes plaque, giving your teeth and gums a deep cleaning, and preventing gingivitis.
- USB charging, travel friendly and compatible with global voltage.
- A low-power reminding function, when the voltage of the battery is low on power, the indicator light will keep flashing quickly or the toothbrush is unable to turn on during the operating process, as a reminder to charge the toothbrush.
- A 2 minute smart timer with interval pausing, that gives a reminder to change the brushing location every 30 seconds.
- There are 5 modes: White, Clean, Sensitive, Polish and Massage settings built into the toothbrush. There is also a mode memory function built into the brush, that means when you turn on the toothbrush, the default starting mode stays the same as the mode that was set at the time the toothbrush was last shutoff.

Model: E10  
Sonic Electric Toothbrush

## 2. Designs:

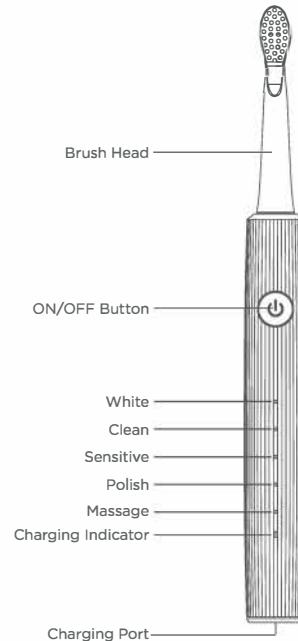

### 3. Operations:

#### ★ Using the toothbrush

- Slide the toothbrush head down to the handle and make sure it is locked in the correct position.
- Turn on the toothbrush by pressing the “Power” button.
- Shift between the 5 modes (“White, Clean, Sensitive, Polish and Massage”) by continuously pressing the “Power” button.
- Wet the bristles and apply an appropriate amount of toothpaste onto the brush.
- Place the bristles with toothpaste against the teeth at a slight angle towards the gum line.
- Gently move the brush head across the teeth on the front, back, and all chewing surfaces. The brush should be allowed to glide across the teeth with minimal pressure. Do not apply too much pressure to the teeth/gums, or apply the bristles to one area of the mouth for an extended period of time.
- The brush has a 2 minute smart timer as this is recommended by dentists as the optimal time for brushing the teeth and there is interval pausing after every 30 seconds as a reminder to move to another quadrant of the mouth.
- Move to the next section when you feel the pause and the toothbrush will automatically power off after 2 minutes.
- To avoid splashing toothpaste, only turn on the toothbrush after you have placed it against your teeth and only remove it from the mouth after it is turned off.

#### ★ Changing the toothbrush head

- Install: Slide the toothbrush head down toward the handle. Please make sure it is in the correct position.
- Detach: Firmly hold the base then, pull the toothbrush head out away from the handle.

#### ★ Charging the toothbrush

- If the toothbrush is in low battery condition when turning it on for the first time, please ensure it is fully charged before usage.
- If the toothbrush is not in use for more than 3 months, please fully charge it according to the charging instructions before using it.
- If the indicator light flashes quickly during the operation process, it indicates the toothbrush is in a low battery condition, please charge the toothbrush as soon as possible.
- When you are going to charge your toothbrush, please insert the USB cable into the charging port in the bottom of the toothbrush and ensure the brush is connected to the electric power safely.
- The charging indicator light will be white and flashes slowly during the charging process. When the battery is fully charged, the indicator will turn on for two minutes and not flash anymore, after 2 minutes, it will turn off automatically. After the toothbrush is fully charged, it is possible to use it for two minutes per session, twice a day for 25 days.
- There is a waterproof function built into the toothbrush, so it can be used in the bathroom safely.

#### ★ Modes Introduction

- The first mode: White  
In this mode, the white indicator light turns on in the “White” mode position, the swing of the brush is large and the vibration intensity is the strongest. It is suitable for users who have high requirements for cleaning effect and do not have sensitive teeth.
- The second mode: Clean  
In this mode, the white indicator light turns on in the “Clean” mode position, the swing of the brush is normal and the vibration intensity is in medium speed. It is suitable for normal oral conditions, to gain the strongest cleaning effect, and maintain healthy oral hygiene.

- The third mode: Sensitive  
In this mode, the white indicator light turns on in the “Sensitive” mode position, the brush swing is soft and the vibration intensity is the softest of the five modes, which is suitable for users with sensitive teeth.
- The fourth mode: Polish  
In this mode, the white indicator light turns on in the “Polish” mode position, the bristles’ swing and the vibration intensity changes rapidly. The simulated effect can be used for additional treatment of the incisor area which has a whitening effect.
- The fifth mode: Massage  
In this mode, the white indicator light turns on in the “Massage” mode position, the bristle swing and the vibration intensity changes slowly. This is for gentle stimulation, designed to keep gums firm and healthy.

### 4. Specifications:

Working temperature: 0-50°C

Working voltage: DC3.7V

Bristles material: Nylon

Input voltage: DC5.0V

Rated power: 2W

Waterproof rating: IPX7

### 5. Precautions:

- Do not disassemble, repair, or alter the product without authorization.
- People who have severe periodontal disease or oral surgery should consult a doctor before use.
- Toddlers, infants and people who are unable to operate should be supervised by a parent or adult.
- Do not soak the product in water for an extended period of time.
- Do not place the toothbrush in high temperature condition or expose it to the sun.
- Do not force the bristles hard against the teeth or gums.
- Not suitable for children under 12 years old.