

# SPACETOUCH

ADVANCED SKINCARE THERAPY

# JUPITER



[spacetouch.com](http://spacetouch.com)

# JUPITER USER GUIDE

Hello!

Thank you for purchasing Jupiter and welcome to the Spacetouch family! There is no better investment than investing in yourself! You are now part of a growing community of millions of people who use advanced LED skincare therapy to achieve flawless skin.

As you unbox your new device, take a few moments to review the manual and become familiar with its features before using your Jupiter LED device. The instructions provide important details on how to safely operate your device.

Spacetouch is committed to providing you with the best skincare experience possible. If you have any questions or concerns, our friendly and knowledgeable customer service team is here to help. You can reach us at [Service@Spacetouch.com](mailto:Service@Spacetouch.com).

# TABLE OF CONTENTS

---

What Is Jupiter?	03
How Does It Work?	03
Technical Specifications	05
Compliance	05
Package Contents	05
How To Use Jupiter	06
Contraindications, Warnings, and Precautions	07
Light Sensitivity Test	09
Limited Lifetime Warranty	09
Register Your Device	10

## **WHAT IS JUPITER?**

Jupiter is a home-use LED light device that combines all of the latest technologies in luxury skincare and medical devices in one easy-to-use device.

Designed by plastic surgeons and dermatologists, this product provides easy and accessible quality skin care treatments; all from the comfort of your home. Jupiter offers 2 modes of treatment- a blue and a red infrared wavelength that work to protect your skin for the long term. Both treatments can be used on your face and neck. Jupiter is rechargeable and travel-friendly, meaning that it provides powerful results wherever you are, allowing you to receive your LED light treatments from the comfort of your own home.

## **HOW DOES IT WORK?**

LED therapy is a non-invasive treatment that utilizes light at precise wavelengths to penetrate the skin's surface without causing any damage. The light stimulates natural chemical reactions in the skin cells, triggering the rejuvenation and healing processes. Developed by leading plastic surgeons and dermatologists, Jupiter is a home-use LED light device that effectively addresses a plethora of various skin

concerns via its 2 available modes of treatment:

- The blue light setting is designed to improve skin texture and reduce sebaceous hyperplasia (enlarged oil glands). It can help with reducing the appearance of actinic keratoses (sun spots), acne, and even acne scarring. The blue light treatment is accompanied by freezing technology. Via this technology, the surface of the device, or "the probe," will cool down to 7°C. This cooling effect tightens pores, slows blood circulation, eases inflammation, helps kill the P. bacteria on the skin's surface, and can be used in tandem with other products as part of an acne treatment plan.
- The red light setting is designed to reduce signs of aging by assisting in the reproduction of cells, collagen, and fiber and reaches the hypodermis level, the deepest and thickest layer of the skin. The red light treatment is accompanied by heating technology. During treatment, the surface of the device, or "the probe" will heat up to 40°C, providing topical heat to open up the pores and increase blood flow to the treatment area. This innovative treatment combines two cutting edge

anti-aging therapies - light and heat -  
for rapid and enduring effects to  
produce fast and long-lasting results.

## **TECHNICAL SPECIFICATIONS**

Jupiter weight: 230 Grams

Product measurement: 20\*6.5\*6.5C

Input power: 10W

Charging time: Approximately 3 hours

Battery type: Lithium 3.7V 900mAh

## **COMPLIANCE**

All Spacetouch LED devices are approved by the U.S. Food and Drug Administration (FDA) and comply with the international standard of basic safety for medical devices IEC 60601 and with the relevant collateral and particular standards for non-laser light source medical equipment for home use IEC 60825-1 and IEC 62471

Photobiological Safety of Lamps and Lamp Systems. Spacetouch devices also conform to the application of risk management to medical devices ISO 14971, and biological evaluation of medical devices in tests for invitro cytotoxicity ISO 10993-5 and irritation and skin sensitization ISO 10993-10.

## **PACKAGE CONTENTS**

Jupiter LED device

Charging dock

Charging cable

# HOW TO USE JUPITER

## TURNING THE DEVICE ON

To turn Jupiter on, hold the device's button down until the device's lights turn red. To switch to the blue mode, or to switch between modes give a short tap to the button.

To turn the device off, hold the button down until the lights turn off.

## USING THE DEVICE

You can use the device morning or night – whichever you prefer! Use the device on the face and neck areas, in circular movements. During the blue light treatment, you may hear the sound of an internal fan, do not be concerned this is expected as part of the freezing process. A single treatment lasts for 8 minutes and the device will turn off once the treatment is complete. You can do multiple treatments on different areas of the face. It is recommended to use Jupiter daily for the first 30 days, and then once a week.

Finish with your favorite moisturizer or serum for maximum results. To clean the device, wipe off the top with a damp cloth or alcohol.

## **CHARGING THE DEVICE**

Connect the charging dock to power via the cable that came in your box. Only use the original cable provided as other cables may damage the device. Now that your charging dock is set up, place Jupiter on the charging dock to charge the device. When your device is fully charged, the device will beep 3 times and the lights on the device will turn off, indicating that it is now ready for use. Charging the device can take up to 3 hours.

## **CONTRAINDICATIONS, WARNINGS, AND PRECAUTIONS**

### **CONTRAINDICATIONS**

DO NOT use the device if you are pregnant.

DO NOT use the device while showering, bathing, or swimming.

DO NOT drop the device or subject it to excessive force.

DO NOT continue to use the device if you develop a skin irritation; consult a doctor if skin irritation persists.

DO NOT use the device if you suffer from Lupus erythematosus, Albinism, photosensitive disorders, photosensitive eczema, or any medications that can cause photosensitivity unless you have consulted with a doctor. LED lights can cause severe skin reactions for the conditions listed.



DO NOT use the device if you suspect that it is damaged or needs repair. Instead, contact Spacetouch for a replacement.

## **WARNINGS**

DO NOT use the device without removing all skincare products and/or makeup prior to treatment. Some ingredients in these products may cause light sensitization or reduce the efficiency of the treatment.

DO NOT tamper with or attempt to modify the device in any way. Modifying any Spacetouch device will result in a voided warranty.

DO NOT continue using the device if it begins to overheat or feel uncomfortably warm on your skin, turn off it immediately and contact customer support.

DO NOT continue using the device if you encounter prolonged skin redness or discomfort that lasts for more than 2 hours after treatment. This could be due to light sensitivity or another condition. Please consult your doctor for advice.

## **PRECAUTIONS**

Before using the Jupiter device, we highly recommend that you conduct a light sensitization test. Stop using the device if you experience any adverse effects from the treatment and consult with your doctor if you have any

concerns about the product and/or treatment.

## **LIGHT SENSITIVITY TEST**

**Step 1:** Power the device and place it on your clean face.

**Step 2:** Use the device for one round of treatment- this should take 8 minutes.

The device will turn off once the treatment is complete. At any point, if you experience uncomfortable heat or discomfort, stop using the device immediately.

**Step 3:** Wait for 24 hours before using the device again. If you experience pain, discomfort, or redness that lasts for more than 2 hours, please consult a doctor for advice as you might have light sensitivity.


## **LIMITED LIFETIME WARRANTY**

The lifetime warranty for Spacetouch devices guarantees coverage for manufacturing defects. To activate the warranty, you must register your device on our website, **[www.spacetouch.com](http://www.spacetouch.com)**, within 7 days of receiving the order.

Failure to register the device correctly and timely will invalidate the warranty.

Please note that the warranty does not cover damages caused by intentional or unintended misuse, use of the device for purposes other than its intended use, or any damage beyond normal wear and tear. Certain restrictions may apply.

Register your device by  
scanning the QR code below

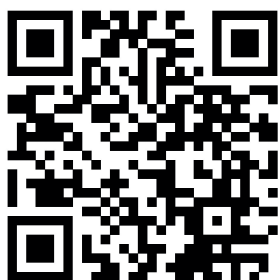



@spacetouch